**READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE
USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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Questions or comments?

Call our US-based customer service toll free at 1-800-536-0366.

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ENGLISH

Content

1. IMPORTANT SAFETY INSTRUCTIONS	2	6. Care, maintenance and disposal.....	4
2. Package Contents	3	7. Troubleshooting Guide	4
3. Before First Use	4	8. FCC Compliance Information.....	5
4. Parts and Controls.....	4	9. Specifications	5
5. To Use	4	10. Warranty	5

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

Dear Valued Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit www.registerbeurer.com or contact customer service at 1-800-536-0366.

Please read these instructions for use carefully and keep them for later use.

Make them accessible to other users and observe the information they contain.

With kind regards,
Your Beurer Team

1. IMPORTANT SAFETY INSTRUCTIONS

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

WARNING

WARNING indicates a hazardous situation which, if not avoided, could result in death or serious injury.

CAUTION

CAUTION indicates a hazardous situation which, if not avoided, may result in minor or moderate injury.

NOTICE

NOTICE addresses practices not related to personal injury, such as product and/or property damage.

**WARNING To reduce the risk of serious personal injury:**

- Keep batteries away from children. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately. Only adults should handle batteries.
- Swallowing batteries and/or battery fluid can be extremely dangerous. Keep the battery and scale out of the reach of children and disabled persons. Should any person swallow a battery and/or battery fluid, please call 911 immediately.
- Should battery fluid leak and come into contact with your eyes or skin, immediately rinse with plenty of clean water. Please call 911 immediately.
- Batteries should not be charged or reactivated by any other means.
- Batteries should not be taken apart, thrown in the fire or short circuited.

**CAUTION: To prevent personal injury and/or product/property damage:**

- Keep scale away from moisture, dust, chemicals, toiletries, liquid cosmetics, extreme hot and cold temperatures, and sources of heat (e.g., open flame, radiators).
- Prevent impact with hard objects.
- This product does not and is not intended to provide a medical diagnosis. For medical advice, contact a licensed physician. If you have an emergency, call 911 immediately.
- Consult your physician or healthcare provider before beginning a weight reduction or exercise program.
- Do not step on the edge of the scale; otherwise it may tip.
- Do not jump on the scale.
- Strong electromagnetic fields (e.g. cell phones) may impair the accuracy of the scale.
- The scale is not intended to be used for commercial purposes.
- Repairs may only be performed by an authorized service center. Contact Customer Service with any questions regarding operation, usage, or repair of this product.
- Packaging materials are a deadly hazard for children and can cause suffocation. Remove all packaging materials immediately and keep them away from children at all times.
- Do not drop the scale or drop any objects on it.
- Do not use the unit on slippery surfaces such as wet floors. There is a danger of slipping.
- Do not step on the scale if you weigh more than 400 lb (180 kg / 28 st).
- Do not use this product for any purpose other than weighing human beings.
- Place the scale on level flooring, otherwise it may tip.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- Do not mix different types of batteries together (e.g. Alkaline and Carbon-zinc) or old batteries with fresh ones. Always replace batteries as a simultaneous set.
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions.)
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.


2. Package Contents

- PS 25 Bathroom Scale
- 2x AAA 1.5 V batteries
- Instructions for use

3. Before First Use

1. Remove the Battery Compartment Cover located underneath the scale and insert two "AAA" size batteries (included) according to the polarity markings inside the compartment. Replace Battery Compartment Cover. When "Lo" appears in the LCD Display, the batteries need to be replaced.
2. Place the scale on a hard, level floor (no carpeting or rugs); a hard surface is required for accurate measurement.
3. Before first use or after moving the scale, reset the scale to true zero by performing the following steps:
 - a. Step onto the scale for one second, then quickly step off.
 - b. The LCD Display will show "0.0 Lb" (or kg/ st as selected).
 - c. After approx. 15 seconds, the LCD Display will shut off.
 - d. You can now use the scale normally, as described in the "To Use" section below.

4. Parts and Controls

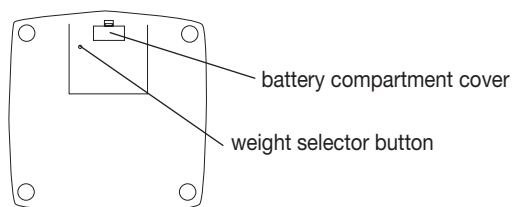5. To Use

1. Press the Weight Selector Button, located underneath the scale next to the Battery Compartment, to show "lb" (pounds), "kg" (kilograms), or "st" (stones) in the LCD Display. Results are shown in increments of 0.2 lb/ 0.1 kg/ 0 st : 1 lb).
2. Place the scale on a flat, level, hard surface to ensure accurate measurement.
3. Carefully step onto the scale.
4. Stand still and distribute your weight evenly on both feet.
5. After a few seconds, the scale will stabilize and your weight will flash on the LCD Display. Step off the scale; the LCD Display will shut off after about 15 seconds.

6. Care, maintenance and disposal

- Clean the scale with a damp cloth and a mild detergent. Do not immerse the scale in water or rinse it under running water.
- Store the scale in a cool, dry place when not in use; ensure no objects are placed on the scale during storage.

Disposal:

Observe the local regulations for material disposal. Dispose of the device in accordance to local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

7. Troubleshooting Guide

Error	Possible causes	Remedy
FFFF	Capacity of 400 lb (180 kg / 28 st) exceeded.	Do not exceed maximum weight.
Lo	Low Power.	Replace batteries.
No display	Batteries are depleted or inserted incorrectly.	Replace batteries and ensure correct polarity.

Scale shows unrealistic result	Scale placed on soft surface or movement took place during measurement.	• Wait until scale switches off automatically and repeat measurement. • Place the scale on hard surface (no carpeting or rugs). • Stand still during measurement.
--------------------------------	---	---

8. FCC Compliance Information

Changes or modifications not expressly approved by the party responsible for compliance could void the user’s authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

9. Specifications

Battery type	2 x AAA 1.5 V batteries
Measuring range	11 lb to 400 lb; 5 kg to 180 kg; 0 st: 11 lb to 28 st
Measurement Increments	0.2 lb; 0.1 kg; 0 st: 1lb
Transport / storage conditions	32 °F to 122 °F (0 °C to 50 °C), 30 % to 80 % relative humidity
Operating conditions	32 °F to 104 °F (0 °C to 40 °C), 30 % to 80 % relative humidity
Scale size	13.3 x 13.1 x 1.6 in (338 x 335 x 41 mm)
Net Weight of the scale	4.2 lb (1.9 kg)

Subject to technical modifications.

10. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer Scale, model PS 25, excluding the batteries, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Scale, model PS 25, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer’s recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period

the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: www.beurer.com

Distributed by:
Beurer North America LP
1 Oakwood Blvd. Suite 255
Hollywood, FL 33020 USA
www.beurer.com

Questions or comments?

Call our US-based customer service toll free at 1-800-536-0366 or contact info@beurer.com

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