READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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Questions or comments? Call our US-based customer service toll free at 1-800-536-0366 or contact info@beurer.com

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Contents

1. IMPORTANT SAFETY INSTRUCTIONS.....	2	6. Care, Maintenance and Disposal	4
2. Battery Handling Safety Precautions	3	7. Troubleshooting Guide	5
3. Parts and Controls.....	3	8. Technical Specifications.....	5
4. To set up.....	4	9. FCC Compliance Information.....	5
5. To Use	4	10. Warranty	6

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

1. IMPORTANT SAFETY INSTRUCTIONS

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
 WARNING	Indicates a hazardous situation, which, if not avoided, could result in death or serious injury.
 CAUTION	Indicates a hazardous situation, which, if not avoided, may result in minor/moderate injury.
NOTICE	Addresses practices not related to personal injury, such as product and/or property damage.

WARNING

- Never leave a child unattended while on the scale, even for a moment.
- Ensure the baby lies down in the tray with minimal movement. Excessive movement could result in a fall, potentially causing serious injury.
- Packaging materials are a deadly hazard for children and can cause suffocation. Remove all packaging materials immediately and keep them away from children at all times.
- This product is not a toy. Keep it out of reach of children.
- People with disabilities, such as impairments, activity limitations, and participation restrictions, or who are physically frail, should be assisted by another person when using this scale.
- This product does not and is not intended to provide a medical diagnosis. For medical advice, contact a licensed physician.
- Consult your physician or healthcare provider before beginning a weight reduction or exercise program.

CAUTION

- Do not weigh more than one person on the scale at one time.
- Do not step on the edge of the scale; otherwise it may tip.
- Do not jump on the scale.

- Do not step or stand on the scale when your body and/or feet are wet, for example after taking a bath or shower. There is a danger of slipping.
- Do not use the unit on slippery surfaces such as wet floors.
- If the glass platform is damaged, there could be risk of danger or injuries caused by broken glass.
- This scale has a capacity of 396 lbs (180 kg). Do not place more than this maximum weight limit on the scale.
- Ensure the equipment is placed on a flat, level and stable surface. Ensure that there is no risk of the equipment tipping over.

NOTICE

- Keep scale away from moisture, dust, chemicals, toiletries, liquid cosmetics, extreme hot and cold temperatures, and sources of heat (e.g., open flame, radiators).
- Prevent impact with hard objects.
- Do not place any object on the scale when it is not in use.
- Do not drop the scale or drop any objects on it.
- Strong electromagnetic fields (e.g. cell phones) may impair the accuracy of the scale.
- The scale is not intended to be used for commercial purposes.

2. Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- Do not mix different types of batteries together (e.g. Alkaline and Carbon-zinc carbon or rechargeable batteries) or old batteries with fresh ones. Always replace batteries as a simultaneous set.
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions.)
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.


3. Parts and Controls

Front - Figure 1

1. LCD/Hold Button
2. LCD Display
3. Power/Tare Button
4. Tray

Rear - Figure 2

5. Weight Selector Button
6. Battery Compartment
7. Tray tabs

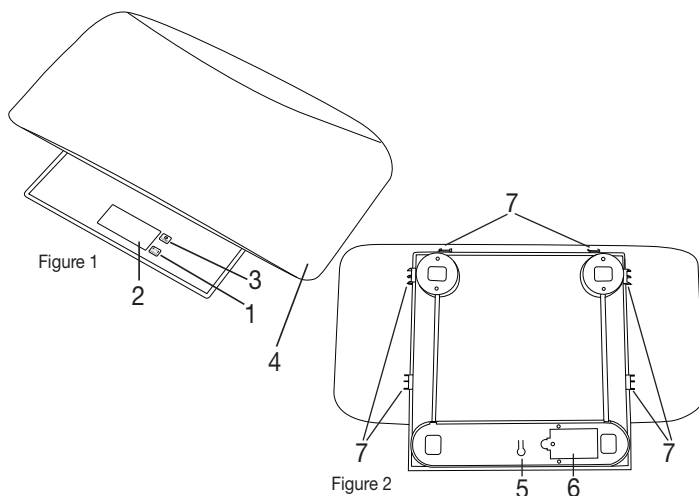
4. To set up

How to install the batteries

Remove the Battery Compartment Cover located underneath the scale by unscrewing the lock screw counterclockwise and lifting the cover away by the Tab. Insert three "AAA" (LR03) size batteries (included) according to the polarity markings inside the compartment. Replace Battery Compartment Cover. When "Lo" appears in the LCD Display, the batteries need to be replaced.

How to install the tray

Slide the Baby Tray onto the scale by aligning the slots on the underside of the tray with the sides of the scale and pushing the tray onto the scale until the tabs on the bottom of the tray contact the bottom edge of the scale.

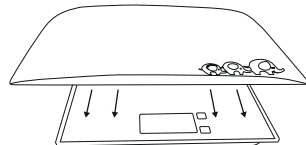Where to place the scale

Place the scale on a flat, level, stable surface to ensure accurate measurement.

5. To Use

1. Press the **Power/Tare** Button to activate the scale. When zeroes appear on the LCD Display, the scale is ready for use.

2. With the scale powered on, press the **Weight Selector** Button located underneath the scale next to the battery compartment, to select through "lb" (pounds), "oz" (ounces), "lb:oz" (pounds ounces), or "kg" (kilograms).

NOTE: The default is set to "lb:oz" (pounds ounces). The scale will remember your last selection as long as the batteries are not replaced.

3. Short press the **LCD/Hold** Button to set the desired display orientation.

4. Carefully place the baby onto the tray, taking extra care with the baby's head. If you do not want to place the baby directly on the weighing surface and prefer to use a towel or blanket, place the towel/blanket on the scale first. The scale features an automatic **TARE** function for weights ranging from 0 to 500 grams. Please note that the weight will remain visible for 5 seconds before resetting to zero. If you place anything over 500 grams on the scale, briefly press the **Power/Tare** button to execute the Tare function manually. "-----" will briefly appear on the LCD Display and the weight resets to zero.

5. This scale has both automatic and manual **HOLD** functions. The Auto Hold function displays the weight when it stabilizes enough to detect an exact value and will display for about 30 seconds.

NOTE: If the baby is moving around too much, the Auto Hold function will not activate. In this case, press the **LCD/Hold** Button when the baby is as calm as he/she will become. Within three seconds, the numbers stop flashing and the weight is displayed.

6. To shut the scale off, press and hold the **Power/Tare** Button. The scale will also shut off by itself after about 30 seconds.

To Weigh a Standing Person

1. Slide the tray off the scale (if attached).

2. Briefly press the **Power/Tare** Button to activate the scale.

3. Briefly press the **LCD/Hold** Button to orient the LCD Display so it is readable while you are standing on it with the LCD Display toward the front.

NOTE: **Power/Tare** Button will become inoperable when the scale is being used in this manner to prevent inadvertent activation from a user's feet if they are standing on top of the **LCD/Hold** Button.

4. When the weight values in the scale read all zeroes, stand motionless on the scale with your feet apart until the weight stabilizes and appears in the LCD Display.

5. The scale will shut off by itself after about 30 seconds.

6. To reactivate the scale, tap your foot on the scale.

NOTE: To reactivate the **Power/Tare** Button, reorient the LCD Display by briefly pressing the **LCD/Hold** Button.

6. Care, Maintenance and Disposal

Care and Maintenance

To keep your scale in the best condition and protect the unit from damage follow the directions below: The scale should occasionally be cleaned with a damp cloth. Do not use aggressive cleaning agents. The unit is not waterproof, do not immerse the unit in water and do not spill water over the unit. If your scale

accidentally gets wet immediately clean off the water.

Repairs may only be performed by an authorized service center. Contact Customer Service at 1-800-536-0366 or info@beurer.com with any questions regarding operation, usage, or repair of this product.

Store the unit in a safe and dry location. Do not expose the unit to direct sunlight, extreme hot or cold temperatures or humidity. Thoroughly dry any moisture off the unit before storing. Store the unit safely out of the reach of children.

Disposal

Dispose of the device in accordance with local regulations.

If you have any questions, please contact the local authorities responsible for waste disposal.

7. Troubleshooting Guide

Error	Possible Cause	Remedy
EEE	Capacity of 396 lb (180 kg) exceeded	Do not exceed maximum weight
Lo	Low Power	Replace batteries
No display	Batteries are depleted or inserted incorrectly	Replace batteries and ensure correct polarity
Scale shows unrealistic result	Scale placed on soft surface or movement took place during measurement	• Wait until scale switches off automatically and repeat measurement • Place the scale on hard surface (no carpeting or rugs) • Stand still during measurement

8. Technical Specifications

Battery type	3 x AAA batteries
Measuring range	1.1 lb to 396 lb (0.5 kg to 180 kg)
Resolution	0.01 lb (5 g) increments (0 lb to 44 lb / 0 kg to 20 kg) 0.1 lb (50 g) increments (44 lb to 396 lb / 20 kg to 180 kg)
Temperature Ranges	Operation: 41°F to 95°F (5°C to 40°C), relative humidity 20 - 85% Storage: 68°F to 140°F (20°C to 60°C), relative humidity 20 - 90%
Measuring units:	lb, oz, lb:oz, kg
Dimensions (with tray):	20.3 x 12.9 x 2.8 in (515 x 327 x 71 mm)
Dimensions (without tray):	11.9 x 11.9 x 1.1 in (302 x 302 x 27 mm)
Net weight of scale (with tray and batteries):	4.4 lb (2.01kg)

9. FCC Compliance Information

Multifunctional Baby Scale BY 82

Responsible Party – U.S. Contact Information:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

United States

1-800-536-0366

info@beurer.com

FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

10. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer Multifunctional Baby Scale, Model BY 82, excluding the batteries, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Baby Scale, Model BY 82, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement. For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn battery; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: www.shop-beurer.com

Distributed by:

Beurer North America LP

1 Oakwood Boulevard, Suite 255

Hollywood, FL 33020

United States

www.shop-beurer.com

Questions or Comments?

Call our US-based customer service toll free at 1-800-536-0366.

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