READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT

Keep this manual in a safe location for future reference

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 **Questions or Comments? Call our US-based customer service toll free at 1-800-536-0366.**

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Content

1. IMPORTANT SAFETY INSTRUCTIONS	2	7. Care, Maintenance and Disposal	4
2. Package Contents	3	Care and Maintenance	4
3. Parts and Controls.....	3	8. Troubleshooting Guide	4
4. Introduction	3	9. Technical Specifications	4
5. To Set Up.....	3	10. FCC Compliance Information.....	5
6. To Use	3	11. Warranty	6

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact Beurer customer service. Please see warranty for service contact.

1. IMPORTANT SAFETY INSTRUCTIONS

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ AND UNDERSTAND THIS ENTIRE MANUAL INCLUDING THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

WARNING

Indicates a hazardous situation which, if not avoided, could result in death or serious injury.

IP 21

Device in protection class 2.

Protected against foreign objects equal to or greater than 12.5 mm in diameter and against vertically falling drops of water.

For home use only.

Read the instructions before starting work and/or operating devices or machines.

WARNING:

- Before use, ensure that all packaging materials are removed and that there is no visible damage to the unit or accessories. When in doubt, do not use the unit and contact customer service.
- Do not immerse the unit into water and do not use it in wet rooms such as a bathroom, laundry room, or pool area.
- Do not use after taking painkillers or drugs or drinking alcohol.
- Keep children away from packaging materials (risk of suffocation).
- Do not touch the unit with wet hands while it is plugged in; do not allow any water to be sprayed onto the unit. Operate only when the unit is completely dry.
- Ensure that you only insert and remove the power adapter with dry hands and that you only press the On / OFF / Timer button with dry hands.
- Keep the power cord away from heated objects and flames.
- Do not disconnect the power adapter from the electrical outlet by pulling on the cord.


- Do not use the device in the presence of flammable anesthetic gas connections with air, oxygen or nitrogen oxide.
- People with retinal diseases as well as diabetics need to be examined by an ophthalmologist before using the daylight lamp.
- Do not use if you have an eye disease such as cataracts, glaucoma, diseases of the optic nerve in general and inflammation of the vitreous humor.
- If you have health concerns of any kind, contact your doctor before using this product.
- The light of the daylight lamp is only intended for use on the human body.
- Connect the unit only to the power voltage listed on the rating plate.
- Do not use it in the vicinity of people who are helpless or children under age 3.
- Do not use the device or adapter if it shows signs of damage or does not function properly. In these cases, contact Customer Service.
- Do not position the device so that it is difficult to unplug the power adapter.
- Strangulation resulting from baby or child entanglement in power adapter output cables is possible.
- Do not make any changes to the device without the manufacturer's permission.
- Only operate the device with the power adapter provided.

2. Package Contents

- 1 x Daylight Lamp
- 1 x AC-Adapter
- 1 x Storage Bag
- 1 x Stand (adjustable)
- 1 x Instruction Manual

3. Parts and Controls

- 1 LED Screen
- 2 Color temperature button
- 3 Timer display (in minutes)
- 4 ON / OFF / Timer button
- 5 Brightness button (dimmer)
- 6 Stand
- 7 AC-Adapter connector

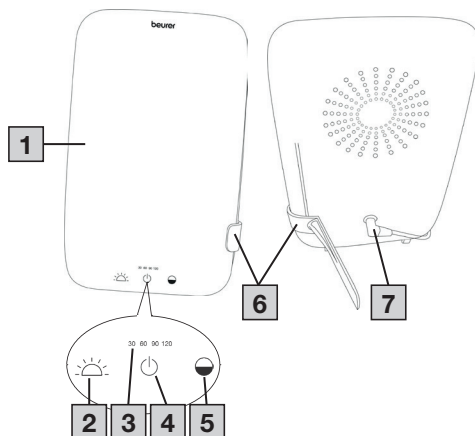4. Introduction

Why use a daylight lamp?

Natural spectrum daylight may help to reduce the effects associated with winter blues, jet lag, shift work and seasonal time changes and boost your sense of well-being.

5. To Set Up

1. Attach the Stand to the side of the lamp so that the clip in the Stand snaps into the groove on the rear side of the lamp. The Stand accommodates both vertical and horizontal lamp orientation.
2. Place the lamp 5-14 inches (13 - 35 cm) away from you on a flat, level, stable, hard surface.
3. Plug the round Power Cord plug into the AC Adapter Socket on the rear of the lamp and connect the adapter to a 110 V AC electrical outlet. Ensure the power cord is routed so that it cannot be tripped over.

6. To Use

1. Touch the **On / Off / Timer Button** to switch on the lamp; the **LED screen** will illuminate on its brightest setting and with cool white color temperature.
2. To adjust brightness and color temperature, touch the **Brightness** and **Color Temperature Buttons**, respectively. Brightness can be set to one of four intensity levels, while color temperature can be set to one of three options as listed in the table below.

Therapy	6500° Kelvin	Cool White
Active	4300° Kelvin	Neutral White
Relax	3000° Kelvin	Warm White

3. To adjust the timer, touch the **On / Off / Timer Button** again. Each touch will set one of four timer durations (30, 60, 90 or 120 minutes) as indicated by the corresponding LED scale. When the set timer duration elapses, the lamp will shut off. After use, the timer duration is saved and the lamp will start up again at the last setting. If no LEDs are illuminated, the timer is inactive.
4. Suggested usage duration depends on the distance; the closer you are to the light, the less time you should use it.

Lux	Distance	Duration
10,000	5 inches (13 cm)	½ hour per day
5,000	8 inches (20 cm)	1 hour per day
2,500	14 inches (35 cm)	2 hours per day

5. Stay close to the lamp with your eyes open since the effect depends on your eyes/retina receiving the light. While oriented toward the light, you can go about your normal activities, such as reading, writing, using your phone, etc.

NOTE: Do not stare directly at the light since this could overstimulate the retina.

6. Use the Daylight Lamp as many days as you want. However, treatment is most effective if used for at least seven consecutive days. The most effective treatment is between 6 AM and 8 PM.
7. Start with shorter periods of usage and increase over the course of a week.

NOTE: Eye aches and mild headaches may occur for the first few uses but will subside after continued usage, as your nervous system becomes accustomed to the new stimulation.

8. Touch and hold the **On / Off / Timer Button** for about two seconds to shut the lamp off manually. Unplug the **AC Adapter** from the electrical outlet when the lamp is not in use.

7. Care, Maintenance and Disposal

Care and Maintenance

- Before cleaning the unit, always switch it off, unplug it, and let it cool.
- Use a slightly damp cloth to clean the device. Do not use abrasive detergents.
- If you are not going to use the device for an extended period of time, disconnect it and store it in a dry place, out of the reach of children.
- Do not cover or store the unit while it is warm.
- Protect the unit from strong impacts.
- This product has no user-serviceable parts; return the product to customer service for any repair or replacement.

Disposal

Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

8. Troubleshooting Guide

Problem	Possible cause	Solution
Device does not light up.	On / Off / Timer Button switched off.	Switch on the On / Off / Timer Button .
	No electricity.	Connect the power cord properly.
	No electricity.	The power cord is faulty. Contact customer service.
	LEDs have reached the end of their service life or are faulty.	Contact customer service.

9. Technical Specifications

Model	TL 35
Type	GUS506

Dimensions (LxWxH)	approx. 6.1 x 0.8 x 9.8 in (15.5 x 2.1 x 25.0 cm)
Weight	approx. 18.3 oz (520 g)
Lighting elements	LED
Power	max. 15 Watts
Light intensity*	Therapy: 10,000 lux (at a distance of approx. 5.1 in (13 cm)) Active: 7,000 lux (at a distance of approx. 3.9 in (10 cm)) Relax: 3,000 lux (at a distance of approx. 7.9 in (20 cm))
Radiation	Output of radiation beyond the visible spectrum (infrared and UV) is so low that it is harmless to eyes and skin.
Operating conditions	32 °F to 95 °F (0 °C to 35 °C), 15–90 % relative humidity
Storage conditions	-4 °F to 140 °F (-20 °C to 60 °C), 15–90 % relative humidity
Product classification	External power supply, Protection class II, IP21
Color temperature	Therapy: 6,500 °K ± 400 Active: 4,300 °K ± 400 Relax: 3,000 °K ± 400
Limit short wavelength	Therapy 435-680 nm
Expected life of the device	10,000 hours
Power Adapter Model	MP12M-120150-AU (Mitra Power Solutions Co., Ltd.)
Power Adapter Input	100 –240 V ~50/60 Hz, 0.5 A Max
Power Adapter Output	12 V DC, 1.5A Max
Polarity	 Center Positive

Maximum output of optical radiation for the TL 35

Radiance Output	Risk group classified by IEC 62471	Max. value
E_{UVA} : Eye UV-A	Exempt Group	0.00007 W·m ⁻²
E_S : Actinic UV skin & eye	Exempt Group	0.00011 W·m ⁻²
E_{IR} : Infrared radiation hazard exposure limits for the eye	Exempt Group	0.00409 W·m ⁻²
L_{IR} : Retinal thermal (mild visual irritation)	Exempt Group	N/A
L_B : Blue light	Exempt Group	59.63 W·m ⁻² ·sr ⁻¹
L_R : Retinal thermal	Exempt Group	63.22 W·m ⁻² ·sr ⁻¹

Subject to technical modifications. The serial number is located on the device.

*Brightness: 10,000 Lux (this light intensity value is for information only).

10. FCC Compliance Information

Daylight Lamp TL 35

Responsible Party – U.S. Contact Information:

Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020
United States
1-800-536-0366

FCC Compliance Statement

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

1. This device may not cause harmful interference, and
2. This device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

11. Warranty

Limited Lifetime Warranty For Original Purchaser

Your Beurer daylight lamp, model TL 35, excluding the battery, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer daylight lamp, model TL 35, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at **1-800-536-0366** or at **info@beurer.com** to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you.

Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: www.shop-beurer.com

Distributed by:
Beurer North America LP
1 Oakwood Blvd. Suite 255
Hollywood, FL 33020, USA
www.shop-beurer.com

For Canada only:
Manufacturer:
Beurer GmbH, Söflinger Str. 218
89077 Ulm, Germany

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

Assembled in China
German Design and Engineering

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