DRUMM

[WARNING]

- DO NOT use if pregnant, have health or medical problems, are suffering from mental illness or have recently had surgery or metallic implants. If you have any doubts about using this device, please consult your local Doctor or Medical Professional before use.
- DO NOT place your fingers around the piston head whilst in use & avoid getting too close to hair or other sensitive areas that could become trapped in the mechanism.
- DO NOT use any power adapters or battery chargers other than those supplied or approved for use by the manufacturer.
- Disconnect the power adapter prior to use. Please ensure that the power switch is OFF when not in use.
- The device is fitted with a 15 min safety timer & will automatically shut-off after 15 minutes of continuous use.
- DO NOT get the device wet & avoid contact with water. Use a soft dry cloth to wipe & clean the surface of the device.


[WARNING]


- DO NOT block the air vents of the motor.
- DO NOT use on head, face, bones or genitals. This device is only suitable for muscle massage.
- DO NOT massage the same spot for longer than 3 minutes to avoid discomfort.
- DO NOT use for longer than 1 hour intervals. After 1 hour of continuous use, please rest for 30 minutes.
- NOT SUITABLE FOR CHILDREN. This device is for adult use only.
- DO NOT disassemble the device or attempt to remove the battery.
- Please store responsibly in a dry & ventilated area.
- Please save these instructions for future reference.

[EN]

[T-1]

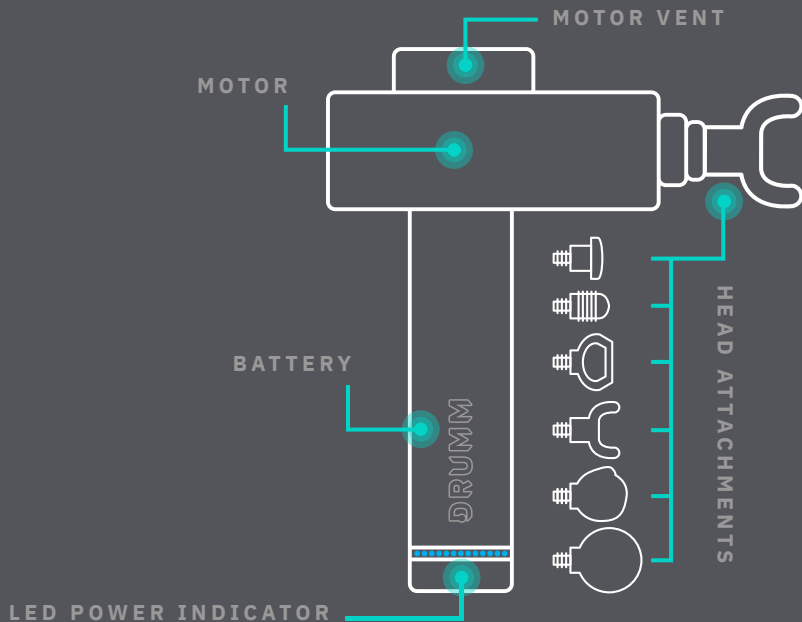[T-1]

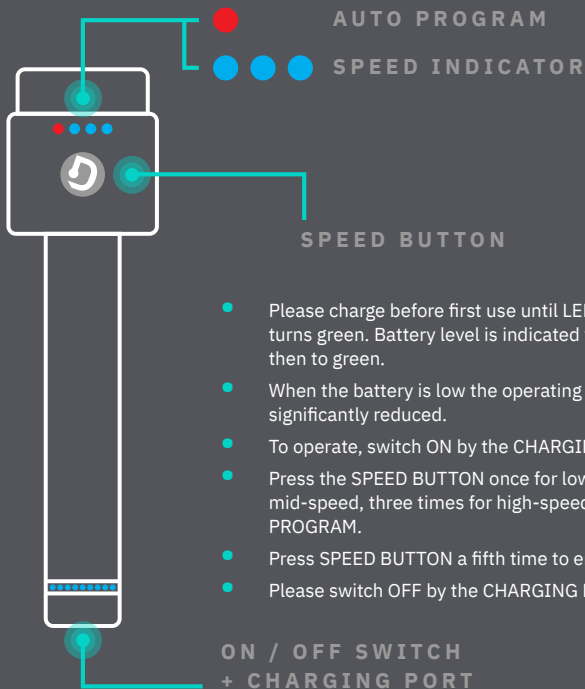- Please charge before first use until LED POWER INDICATOR turns green. Battery level is indicated from red, to blue & then to green.
- When the battery is low the operating intensity will be significantly reduced.
- To operate, switch ON by the CHARGING PORT.
- Press the SPEED BUTTON once for low-speed, twice for mid-speed, three times for high-speed & four times for AUTO PROGRAM.
- Press SPEED BUTTON a fifth time to enter standby.
- Please switch OFF by the CHARGING PORT when not in use.


[EN]

[HOW TO USE]

- Start up the device away from your body.
- To begin, make light contact & float across the muscular area to get used to the feeling of percussive therapy.
- Apply comfortable pressure for you & allow the device to do the work.
- Start slow & increase the speed & pressure only when you feel comfortable doing so. You should not be sore from using the device.
- We recommend using the device for 15 minute sessions at a time & only treating one specific area for no longer than 3 minutes to avoid discomfort.

[TECHNICAL DATA]

Size:	222mm x 160mm x 55mm
Frequency:	30 - 43 Hz / 1800 - 2600 Percussions Per Minute
Amplitude:	8mm
Battery:	1500mAh
Motor:	DC 16.8 V (Brushless)
Weight:	0.9kg
Noise:	50db (Distance of 30cm)

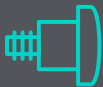
FLAT HEAD

Massage & relax all muscles.


BULLET HEAD

Target joints & trigger points.

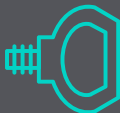
SPADE HEAD

Penetrate shoulder, neck
& pelvic muscle aches.

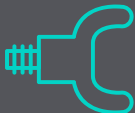
FORK HEAD

Grip spinal, forearm
& leg groups.

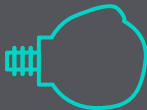
BOXING HEAD

Add more punch to
larger muscle groups.

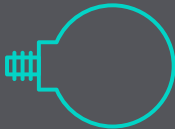BALL HEAD

Massage & relax
larger muscles.

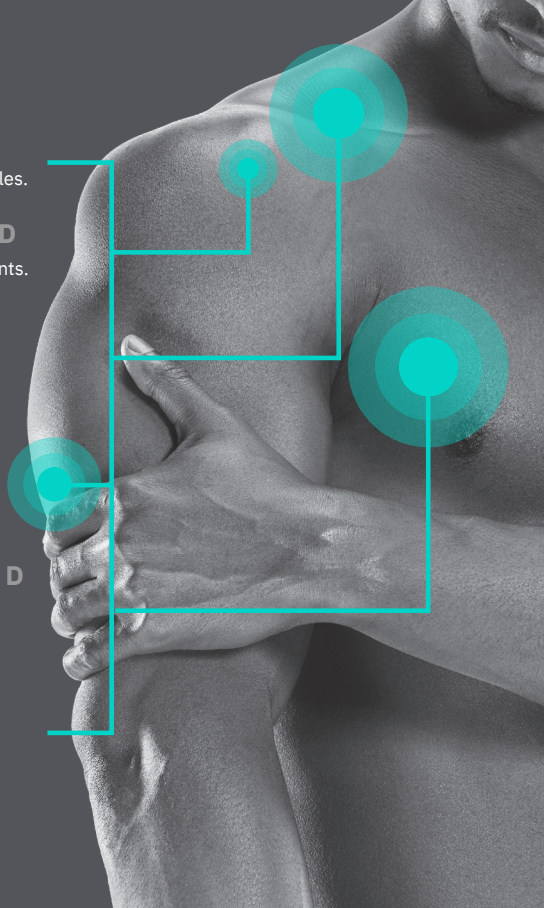[WARRANTY]

The Bluefin warranty covers any technical faults with your item. If you need to make a warranty claim contact us via email at support@bluefintrading.co.uk

Your warranty does not cover the following:

- General wear and marks from use.
- Water damage.
- Damage due to improper use.
- Chassis damage and scratching due to improper use or dropping.
- Electronic damage due to improper use or dropping.

Your warranty begins on receipt of delivery, technical faults are covered for 12 months.

If you have any questions regarding your warranty please email us at:
support@bluefintrading.co.uk

William Vaughan
CEO, Bluefin Trading Ltd.

BLUEFIN
— FITNESS —

[EU-DECLARATION OF CONFORMITY]

For the latest & up-to-date EU declaration of conformity for this product visit:
www.bluefinfitness.com/manuals/drummT1_DoC.pdf

Importer Address:

Bluefin Trading Ltd, Keelham Farm, Hebden Bridge, UK, HX7 8TG

www.bluefinfitness.com

Waste Electrical and Electronic Equipment

Do not dispose of this device with domestic waste. The components of the device must be separately recycled or disposed of according to local regulations. The system operator may also return any waste electrical and electronic equipment to the manufacturer at the end of their lifetime. Charges incurred for transport to the manufacturer will be at the expense of the owner of the waste electrical and electronic equipment.

Battery Directive

Do not dispose of any batteries with domestic waste. Please recycle or safely dispose of batteries in accordance with local regulations. Any batteries supplied by the manufacturer or distributor can be returned at the end of their lifetime for recycling and safe disposal. Charges incurred for transport to the manufacturer or distributor will be at the expense of the owner of the batteries.

Recycling of Packaging Materials

Do not dispose of packaging with domestic waste. All packaging materials can be safely recycled according to local regulations.

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