BRAUN

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21318308/VI-25

www.braun.com

PL1

US Helpline:

1-800-BRAUN-11 or

1-800-272-8611

Canada Helpline /

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1-800-387-6657

Warranty US / CAN

For warranty terms and details on money-back satisfaction guarantee please refer to:

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3.2 Warnings 6

3.3 Precautions 7

4 Potential Risks, Side Effects and Skin Reactions 7

5 Benefits 8

6 Intended Use 9

7 How to Use 9

8 Treatment Modes 10

9 Cleaning, Maintenance and Storage 11

10 Travelling with the device 11

11 Troubleshooting 12

12 Labels and Symbols 12

13 Technical Specifications 13

14 Support 13

Our products are designed to meet the highest standards of quality, functionality and design. Using professional light technology, Braun Silk-expert Mini Intense Pulsed Light (IPL) is a hair removal device that can be used comfortably at home.

1 Description

- 1 Treatment window
- 2 SensoAdapt™ skin tone sensor
- 3 Glass filter
- 4 Power indicator light
- 5 Treatment button
- 6 Handle
- 7 Base station with serial number
- 8 Mains cord

2 How Does Silk-expert Mini Work?

Silk-expert Mini is designed to help break the cycle of hair growth. Light energy is transferred through the skin's surface and is absorbed by melanin present in the hair shaft. The absorbed light energy is converted to heat energy (below the surface of the skin), which disables the hair follicle. Treated hairs naturally fall out over the course of a few days to 1-2 weeks.

What to Expect

Immediately post treatment, you should not see any significant side effects from treatment (Refer to section 4 for more information on side effects).

In the first few weeks following the initial treatments, you will still see some hairs growing. These are likely to be hairs that were missed during treatment i.e. they were lying dormant and not in their growth phase (Anagen) during treatment, when IPL is most effective.

Around 6 weeks into the 12 week treatment program you should see a reduction in hair growth. However, many hairs may still not have been treated in their growth phase. It is important to continue the weekly treatments.

After the 12 week program, you should see a significant reduction in hairs within the treated area. Any remaining hairs should be finer and lighter in color. Continued monthly treatments, or treatments as required, should maintain the reduction in unwanted hair.

3 Important Safety Information

3.1 Contraindications

DO NOT use the device if you have very dark skin. The Skin Tone sensor located next to the treatment window of the device will determine the tone of your skin. If the Skin Tone sensor determines your skin is too dark for safe use, operation will be prevented.

DO NOT use the device if your skin has been artificially or naturally tanned recently or is burned from sun overexposure. Your skin may be extra sensitive following sun exposure and particularly susceptible to burning, blistering, discoloring or scarring if you use the device. Wait 1 week until sunburn or tan disappears before using the device. Treating tanned skin with IPL may result in permanent Hyper (darkening) or Hypo (lightening) changes to the treated skin.

DO NOT expose treated areas to the sun. Wait at least 7 days after treatment before exposing treated skin to direct sunlight. Your skin may be extra sensitive following IPL treatment and particularly susceptible to sunburn. Apply sunscreen (SPF 15 or greater) to the treated area or cover the treated area with suitable clothing.

DO NOT use the device directly on the nipples, genitals or around the anus. These areas may have a darker skin color and/or greater hair density. Using the device in these areas may cause discomfort/pain or injure (burn, discolor or scar) your skin.

DO NOT use the device if you have a history of skin cancer or pre-cancerous lesions (e.g. nevi or a large number of moles).

DO NOT use if you are pregnant, lactating, as this device has not been tested with these individuals.

DO NOT use on dark brown or black spots such as birthmarks, moles or warts in the area you wish to treat. Using the device may injure your skin or make existing conditions worse. You may experience side effects such as burns, blisters and skin color changes or scarring.

DO NOT use if you have any skin condition in the treatment area, including psoriasis, vitiligo, eczema, acne, herpes simplex or active infections, wounds or blisters.

DO NOT use if you are under the age of 18 as this device has not been tested on these individuals.

DO NOT use on scars, tattoos or permanent make up in the treatment area.

DO NOT use on beard/facial hair for men.

3.2 Warnings

DO NOT use if you have a known sensitivity to sunlight (photosensitivity) or are taking a medication that causes photosensitivity e.g., certain antibiotics such as tetracycline, NSAIDS such as ibuprofen or retinoids such as Accutane and/or topical retinoids such as Retin A. Always check the instructions that come with your medical products to see if photosensitivity is listed as a side effect. Using the Silk-expert device on photosensitive skin could result in skin damage including blisters and swelling. Seek medical advice prior to application of the appliance as use may harm the skin.

DO NOT use more than once a week on the same area. Damage of skin may occur after prolonged or repeated treatment on one side. It is not recommended to use this device on areas where dermal fillers are present.

Women who are using the device on the face:

DO NOT use if you have had a skin peel treatment on the area of the face you wish to treat. You should wait 30 days after a skin peel before treating.

DO NOT use if you have had a laser skin resurfacing procedure on the area of the face you wish to treat. You should wait 6 months following a laser skin resurfacing procedure before treating.

EYE SAFETY

The device emits flashes of Intense Pulsed Light. Direct exposure is potentially harmful to your eyes. Take care to follow the safety precautions below.

Possible eye injury (potentially leading to loss of vision) or skin injury if instructions are not followed.

Protect the eye from exposure. Do not use the device over the eye lids or close to the eye.

The device only emits a flash if the skin tone sensor (2) detects a valid Skin Tone reading. However:

DO NOT look directly into the treatment window (1) when the device has electrical power and is activated.

DO NOT attempt to release a flash towards the eye.

DO NOT treat areas around the eye (eyebrows or eyelashes).

ONLY press the treatment button when the treatment window is in good contact with the area you wish to treat.

WE RECOMMEND you look away from the device when pressing the treatment button.

ELECTRICAL AND FIRE SAFETY

CAUTION

Always unplug the appliance after use and in case of a power failure to avoid risk or damage to the equipment.

DO NOT subject the equipment to heavy shocks to avoid risk of damage to the lamp.

WARNING: Do not leave the equipment unattended when it is switched on to avoid risk of fire or burns.

DO NOT try to open the device, as this may expose you to dangerous electrical components.

DO NOT use if the device is damaged; for example, cracked housing, cable damage (internal wires visible), cracked or broken glass filter etc. This may result in injury.

DO NOT use if the device becomes uncomfortably hot to touch. This may indicate that the unit is damaged.

DO NOT attempt to use the device on any surface other than skin.

DO NOT operate the device if the mains plug is not easily accessible. You should be able to disconnect at any time.

DO NOT override the safety mechanism inherent to the device.

WARNING: No modification of this equipment is allowed.

Do not use any other component, nor manipulate the device. This may damage the unit and cause it to stop working.

WARNING: Water and electricity are a dangerous combination. To avoid risk of electrocution, - do not use this equipment in wet surroundings (e. g. in a bathroom, near a shower or swimming pool); - do not let water run into the appliance.

WARNING: If the wall socket used to power the equipment has poor connections, the plug of the equipment becomes hot. Make sure you plug the equipment into a properly installed wall socket to avoid the risk of fire and burns.

GENERAL SAFETY

The device should be kept out of the reach of children under 18 years of age. Children should not use, play with, maintain or clean the device. There are a number of potential risks – exposure to the light output causing burns or eye damage, electrocution, strangulation from the cables.

People with reduced physical, sensory or mental capabilities or lack of experience and knowledge may use the device if they have been given supervision or instruction concerning the safe use of the device and understand the hazards involved.

Keep the device out of reach of pets and pests which could damage it.

The glass filter (3) is a critical component that will get hot during use.

DO NOT touch the glass filter during or directly after use.

Regularly inspect the treatment window (1) for damage:

DO NOT use if the glass filter (3) is cracked or missing. Do not use if the built-in reflector is missing or not in correct place.

3.3 Precautions

Hair removal by lasers or intense pulse light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern, and South Asian heritage treated on the face and neck.

DO NOT USE the device on any area where you may want hair to grow back.

The device is not recommended for use on red, grey, white or blonde hair since it is not effective on these hair colors.

4 Potential Risks, Side Effects and Skin Reactions

The use of this device could result in side effects. The table below lists the known skin reactions that may occur following treatment with the device.

Possible Side Effects	How to assess and react
Mild pain / discomfort in the area being treated	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed, and the pain should diminish with continued use.
Warm feeling or tingling sensations during treatment, which typically disappear after a few seconds to a minute and decrease with continued use.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed.

Itchiness in the treatment area	This is quite common for IPL treatments, and should subside after a short period. You can keep on using the device as instructed. Do not scratch the area.
Skin redness during or after treatment which disappears within minutes.	This is expected and is normal for all IPL treatments. You can keep on using the device as instructed once skin redness has disappeared.
Skin redness which does not disappear within minutes and lasts 24-48 hours.	Stop using the device immediately and consult your physician before using it again.
Pain or discomfort that is intense during treatment or persists after a treatment.	Stop using the device immediately and consult your physician before using it again.
In very rare cases	How to assess and react
Swelling and redness that does not disappear within two to three days.	Stop using the device immediately and consult your physician before using it again.
Temporary changes in skin color (lightening or darkening).	If your skin color changes, stop using the device immediately and consult your physician.
Blistering or burning of the skin.	Stop using the device immediately and consult your physician before using it again. Apply a cold pack to the area. Treat with antiseptic or burn cream.

In a clinical trial of 50 subjects, each subject received 12 weekly treatments. The recorded effects were:

- Itchiness: Four subjects experienced mild itchiness after treatment.
- Slight Redness: This was reported by 1 subject.
- Mild Stinging Sensation: Five subjects reported mild discomfort following treatment.
- In-growing Hair: One subject reported experiencing an in-growing hair following treatment.
- Follicle Redness: One subject reported a red spot at the treatment site following treatment.

In all cases, the treated skin returned to normal within 7 days.

5 Benefits

The device is indicated for the permanent reduction of unwanted hair.

A clinical trial was performed to assess both the safety and efficacy of the device. The key features of the clinical trial were:

- All treatments were performed by non-medically trained operators at a single location using devices identical to the one provided in this package;
- All subjects were required to complete questionnaires related to their general health and those who could not be treated, e.g. those who had one or more of the contraindications as listed in Section 3 of this User Manual were excluded. In addition, all subjects were required to provide Informed Consent as per international clinical trial requirements;
- A total of 53 female subjects, aged between 19 and 45 years old, commenced the trial;
- Three subjects left the trial during the treatment period for personal reasons. No subjects left the trial due to problems associated with the use of the device;
- Each subject received 12 weekly treatments to their chosen body location. A typical treatment consisted of shaving the site then applying the device following the same process as outlined in Section 7, Step 3 of this user manual;
- Hair counts were used to measure the change in the amount of hair 6 months after the last treatment and 12 months after the last treatment. The hairs in the treatment sites were counted via high resolution photography and the difference in hair counts was calculated as a percentage change;
- When used as directed, the clinical trial participants showed on average 44% less hair 6 months after their last treatment and 36% less hair 12 months after their last treatment when compared to hair counts taken before treatment (see table below). Actual results did vary from person to person.

Number of subjects at 6 months post treatment	50
Hair reduction at 6 months post treatment	43.9 %
Number of subjects at 12 months post treatment	33
Hair reduction at 12 months post treatment	36.0 %

% of subjects met success (> 30% hair reduction) on all body areas at 12 months post treatment. Subject success is defined as greater than 30 % hair reduction at all treatment sites at 12 months.	66.7 %
---	--------

- The incidence of side effects during the clinical trial was minimal (as outlined in Section 4 above) and the majority of subjects (48 out of 50) would recommend the device to their friends.

6 Intended Use

INDICATION FOR USE

The Braun Silk-expert Mini Hair Removal Device is indicated for the removal of unwanted hair. The Braun Silk-expert Mini is also indicated for the permanent reduction in hair regrowth, defined as the long-term, stable reduction in the number of hairs regrowing when measured at 6, 9 and 12 months after completion of a treatment regime.

BODY AREAS

Women	Men
Legs, arms, underarms, intimate area*, face (below cheek bones), fingers, toes CAUTION DO NOT USE... around the eyes, along the forehead, on scalp, nipples	Shoulders and below (chest, back, underarms, arms, stomach, intimate area*, legs, fingers, toes) CAUTION DO NOT USE... on face, neck, scalp, nipples
*Intimate area: VIO  Bikini area between legs buttocks and skin around anus (peri-anal) DO NOT USE... in mucous membranes such as labia minora, vagina or on the anus	*Intimate area:  Area above penis buttocks and skin around anus (peri-anal) DO NOT USE... on penile shaft, scrotum or anus 

SKIN COLOR

The device is intended for use on light, medium and dark skin tones up to and including skin tone 5 (see Skin tone chart on page 3). The device incorporates a skin tone sensor that measures the color of your skin and will prevent use if your skin tone is too dark.

HAIR COLOR

The device is intended for use on naturally black or brown hair. It may not be as effective on white, grey, blonde or red hair.

7 How To Use

WARNING: Follow each step in the treatment process ensuring you complete all aspects of each step before progressing to the next.

Step 1: Preparing the treatment area

Remove all visible hair on the skin before using the device. You may use your preferred hair removal method. Remove any cosmetics, lotions and cream from the area you wish to treat. Cleanse the treatment

area and pat dry. To keep the treatment window clean and intact and prevent hair burns, make sure there is no more hair remaining on the skin surface.

Step 2: Doing a Test Patch

Before your first treatment on each new body area, we RECOMMEND you test your skin (in that area) for reactions. Perform two flashes next to each other as described below in chapter 8 Treatment Methods, section STAMPING. Wait 15 minutes to examine the test area. If the skin appears normal, begin your first treatment. If a skin reaction develops, wait for another 24 hours. If skin reaction is present after this period, discontinue using the device and consult your physician.

Step 3: Treating with Silk-expert Mini

- a) Use the mains cord (8) to connect the base station to an electrical outlet. Press the treatment button (5) to activate the device. The power indicator (4) will light up once and a fan will start running. This means the device is ready for use.. During operation, the device will heat up slightly.
- b) Press the treatment window (1) firmly against the area to be treated ensuring the skin tone sensor (2) is in contact with the skin. Once the skin tone sensor measures a valid skin tone, the pre-set output power is displayed by the power indicator light color.

If your skin is too dark for treatment with Silk-expert Mini, the skin tone sensor will detect this and prevent the unit from operating. RED light will be shown on the power indicator if your skin is too dark. If a valid skin tone is detected, the power indicator will illuminate pink or white.	
Pink:	high intensity for light skin tones
Light pink:	medium intensity for medium skin tones
White:	low intensity for dark skin tones
Red:	Invalid skin tone / device will not flash

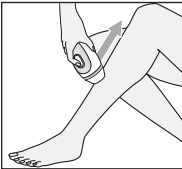
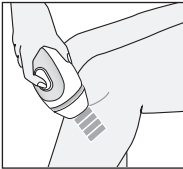
Thanks to the built-in, intelligent SensoAdapt™, which automatically measures the skin color, there is no need for the user to make any manual adjustment for skin tone. The skin tone chart on page 3 is displayed also on the sales box as a guidance to potential purchasers to determine if the device will work on their skin coloration. This does not relate to the Fitzpatrick scale, it is simply a color guide.

- c) Press the treatment button (5). The device will flash and a warm sensation should be felt.
- d) Move the treatment window to the next area to be treated and repeat. The device needs to recharge between flashes. This will take less than a second.
- e) Ensure that the entire area receives treatment, but DO NOT treat an area twice (overlapping flashes). After 10 minutes without use, the device is automatically deactivated. To re-activate, press the treatment button.

If the device is used in high ambient temperatures for extended treatment sessions, short pauses in treatment may be seen to allow cooling. This is normal and operation will resume once the device is cooled.

8 Treatment Modes

The area to be treated has to be hair-free, clean and dry. Before pressing the treatment button, look away from the device to avoid the bright light. In seldom cases, a user may find the flash too bright, however there is no safety risk associated with this. For increased comfort, we recommend treating in a well-lit room to help reduce the perceived brightness.

GLIDING for fast and efficient treatment of large areas (e.g. legs).		STAMPING for treating small, hard-to-reach areas (e.g. face for females, bikini line, underarms, knee, chin, ankles)	
--	---	--	---

1. Place the treatment window (1) firmly on your skin, ensuring the entire front area has skin contact. Power indicator shows pink or white light to confirm that the device is ready to flash.	
2. Press and hold the treatment button (5). The device will begin to flash.	2. Press and instantly release the treatment button (5). The device will emit a flash.
3. Move the device <u>slowly and continuously</u> from the bottom to the top of the treatment area to avoid missing areas. Make sure to always keep full skin contact so that the device flashes continuously.	3. Lift off the device and place it onto the next adjacent spot to be treated. Wait for the power indicator to show white or pink. Then press and release the treatment button.
Once you have completed your first gliding stroke, lift the device off the treatment area. Repeat steps 1 to 3 to treat the next adjacent area in a parallel stroke.	Continue treatment spot by spot until you have completed the entire area you wish to treat.

Make sure to avoid missed areas by placing the treatment window directly next to the previously treated spot without leaving any space between treated areas.

If you feel intense pain, stop treatment immediately and consult your physician before using the device again.

How to treat sensitive areas

Before treating, remember to perform a skin test. Remove any visible hair and cosmetics, lotions or creams. You might want to confirm full hair removal using a mirror.

Face (females only)

When treating the upper lip area, avoid contact with the lips and the nostrils. The central area below the nostrils might be sensitive.

Do not use on permanent / semi-permanent lip liner on the lip margin.

TIP: You can press your lips together to ensure that the light does not activate on the upper lip itself.

Hair growth under the chin can be denser and coarser. This can sometimes lead to slight skin reddening following a treatment. This redness is temporary and should completely disappear after 24 hours.

Intimate area (please see page 9 for suitable treatment areas)

Stand, sit or crouch down to reach all areas to be treated. Treat with the stamping method (spot by spot).

Step 4: Your Treatment Schedule

For optimal results, Silk-expert should be used once a week for 12 weeks to ensure all hair follicles are treated. After the initial treatment phase (up to 12 treatments), we recommend the device is used once every 1 to 2 months, or as required, to maintain results. For the facial area or those with darker skin tones, you may need to extend the initial treatment phase to achieve your best results.

Step 5: Switching Off the Device

When the treatment is finished, unplug the device. Refer to section 9 for details of how best to look after your device.

9 Cleaning, Maintenance and Storage

After treatment, always unplug the device.

After each treatment, inspect the treatment window (1), the skin tone sensor (2), the glass filter (3) and the built-in reflector for damage.

Wipe the treatment window with a dry, lint-free cloth if any marks or black spots are visible.

WARNING: The glass filter (3) can get hot during use. DO NOT touch or clean the filter for at least 5 minutes after use to allow it to cool down.

NEVER use water or other cleaning fluids, as these can damage the device leading to a potential safety hazard.

Store the device in a cool, dry place. Make sure that the treatment window and skin tone sensor is protected from damage.

Regularly check the device (including cords) for visible signs of damage. In case of damage or cracks, stop using the device.

10 Travelling with the Device

The device will work with all commonly encountered mains supplies. Providing you have a suitable adaptor for the mains power cord, the device will operate as normal.

11 Troubleshooting

Problem	Solution
Power indicator (4) does not light up when treatment window is placed on the skin.	Make sure the device is connected to an electrical outlet. Press the treatment button to activate it. Ensure that the entire treatment window area touches the skin. If needed, move the device to improve skin contact.
Flashing sequence slows down.	Device is in heat-protection mode. To get back to normal speed, allow device to cool down.
Red light on power indicator	Invalid skin tone reading. Your skin is either too dark or too light. For safety reasons, the device will not flash.
Red flashing power indicator light	Device error. Device will not operate. Unplug and re-connect. If error persists, contact a Braun Service Center.
Power indicator light flashes alternately red and white.	End of product lifetime. Maximum number of flashes (300,000) has been reached.
The device is broken, cracked or appears defective.	Do not use if the handle, the base station or the mains cord is damaged.

12 Labels and Symbols

The base station is marked with the following symbols:

	Type BF Applied Part
	Class II equipment
	Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures: – Reorient or relocate the receiving antenna. – Increase the separation between the equipment and receiver. – Connect the equipment into an outlet on a circuit different from that to which the receiver is connected. – Consult the dealer or an experienced radio/TV technician for help.
	Follow instructions for use
	The device contains recyclable electronic waste. At the end of its useful life, do not dispose of it in the household waste. To protect the environment, take it to dedicated local collection sites for recycling purposes.
	Conformity mark of the European Union
	Keep Dry

	Optical radiation
	Not suitable for use in a bathtub or shower
	NRTL Certification marking

Other symbols shown on the base station refer to requirements for other regions.

13 Technical specifications

Braun Silk-expert Mini is a filtered broadband Intense Pulsed Light system with the following technical specifications:

Repetition rate:	every 0.63 seconds continuous operation
Max optical output:	4 J/cm ²
Pulse length / profile:	0.5 - 1.0 ms FWHM
Wave length:	510-2000 nm
Power line input:	100 - 240V~, 50-60Hz, 1.7 - 0.91A
Operating temperature:	5°C - 40°C / 41°F - 104°F
Optimum operating temperature:	15°C - 30°C / 59°F - 86°F
Operating humidity:	up to 93% R.H. non-condensing
Operating pressures:	700 hPa to 1060 hPa
Treatment area (spot size):	3 cm ²

When having stored the device outside of the operating temperature range, leave it for 30 minutes at room temperature before using it. The area around the treatment window can reach a maximum temperature of 47.1 °C / 116.78 °F when the device is used for a long period of time.

Transport Conditions:

Temperature range:	- 25°C to 70°C / -13°F to 158°F
Humidity:	up to 93% R.H. non-condensing
Pressures:	500 hPa to 1060 hPa

The device includes 300,000 flashes and is intended for a single user.

Expected service life: 5 years

This device complies with the following EMC standards and latest amendments applicable to Household and Medical Devices. EN 55014-1:2017, EN 55014-2:2015, EN 61000-3-2:2014, EN 61000-3-3:2013, IEC 60601-1-2:2014 (Ed 4.0), EN 62311:2008, EN 62233:2008, FCC CRF 47 Part 15B, 15.107, 15.109, ICES-003 Issue 6.0.

This device also complies with the following safety standards and latest amendments applicable to Household and Medical devices. IEC 60335-1:2010/AMD2:2016, IEC 60335-2-113:2016, IEC 60601-1:2005/AMD1:2012, IEC 62304:2006, IEC 60601-1-11:2015, IEC 60601-1-6:2010/AMD1:2013, IEC 62366-1:2007/AMD1:2014, IEC 60601-2-57:2011, IEC 60601-2-83:2019, IEC 62471:2006, CAN/CSA-C22.2 No. 60601-1:2014, ANSI/AAMI ES 60601-1:2012.

Subject to change without notice

This device is classified as Risk Group Exempt for all hazards

For Spectral Irradiance, please refer to Annex 1 on page 25.

 Manufacturer:
CyDen Ltd., Block A, Bay Studios Business Park, Fabian Way, Swansea SA1 8QB, UK

14 Support

US Helpline: 1-800-BRAUN-11 or 1-800-272-8611

Canada Helpline: 1-800-387-6657

If you are unable to call our Helpline during working hours, we provide online support via our web page (www.braun.com) or Facebook.

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CyDen Ltd, Block A, Bay Studios Business Park, Fabian Way, Swansea SA1 8QB, UK14