

成品尺寸 :85×85mm 展开尺寸 :170×85mm 材质 :80g 白色书纸 ©liudanni 20230315




breco® Neck N5 mini

Designed by Breco. Made in China. Shenzhen Breco Technology Co., Ltd.
TM and © 2023 Breco. All rights reserved.



SIMPLY RELAX

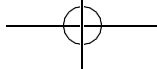

Content


Precautions	24
product configuration	27
Features	28
Product Operation	29
Troubleshooting	37
Product Disposal Requirements	38
Protect	40
Symbol Description	41
FCC Statement	42
Product Specifications	43

EN

PRECAUTIONS

Attention

1. Please read the instructions carefully and do not use this product continuously for a long time to prevent muscle pain caused by excessive massage.
2. Do not use when lying flat or half lying down.
3. Do not use this product during precise operation or driving.
4. Do not use if you have undergone surgery on the massaged part or have skin damage, redness, inflammation, congestion and other symptoms.
5. People with heart disease or medical history, acute disease patients, skin disease patients, malignant tumor patients, bleeding tendency and other blood diseases, osteoporosis, soft tissue diseases, stroke and other physical abnormalities should not use.
6. Do not use it if you have slow movement, physical and mental discomfort or intellectual disability. If you need to use it, please follow the doctor's advice.
7. Pregnant women and menstrual period are not recommended to use the massager. If you need to use it, please follow the doctor's advice.
8. When storing the product, please pay attention to placing it in a place where children are not easy to reach.
9. Children under the age of 8 are prohibited from using it, and children aged 8-12 must use it under the guidance and care of their parents or someone with supervisory ability.
10. Do not scratch or stab this product with sharp objects.

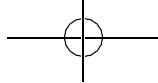

EN

11. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
12. Children should be supervised to ensure that they do not play with the appliance.
13. It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.

Caution

1. This product is not a medical device, and not suitable for diagnosis, treatment and other commercial purposes.
2. Before use, please remove the jewelry on the massaged part or other items that may affect the massage.
3. During use, if there is pain or discomfort in the massage area, please stop using it.
4. When using, if the product is abnormal, please shut down immediately, cut off the power supply and contact the customer service hotline, do not disassemble or repair by yourself.
5. The massage area has a heating surface, people who are not sensitive to heat need to pay attention when using it.
6. When using in other specific environments, please comply with safety regulations.

EN

Product Protection

1. Please use the charging cable configured by our company to charge the product with a standard adapter, and pay attention to the safety of electricity use.
2. Please keep this product away from water, fire and corrosive environment.

Necessary Knowledge

*About frequency of use

1. Use as needed, within 10 minutes at a time. It is recommended to take a proper rest before use after one use.
2. If overused, it may cause physical discomfort and skin injury.
3. Therefore, when you use it for the first time, you need to observe your physical condition and skin changes, and try to shorten the use time as much as possible.

*Cautions Before Use


1. When using this product, please take a personal comfortable posture and completely relax your body.
2. If you feel pain or discomfort on the skin, please stop using it and seek medical attention immediately.

*State After Use

After long-term massage, the massage area may feel red and hot while feeling comfortable, which is normal.

* Cleaning

If you need to clean the appearance, you can use a rag dipped in water or 75% alcohol to wipe the smooth surface, do not use water to rinse.

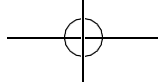

PRODUCT CONFIGURATION



Neck N5 mini



User Manual



Charging Cable

Indications for use:
Breo Neck Massager is intended to
relax neck muscle and relieve minor
neck muscle aches and pains.

27

FEATURES



Bionic Massage




Free Hands




Gentle Heating


28

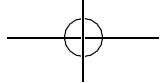



PRODUCT OPERATION

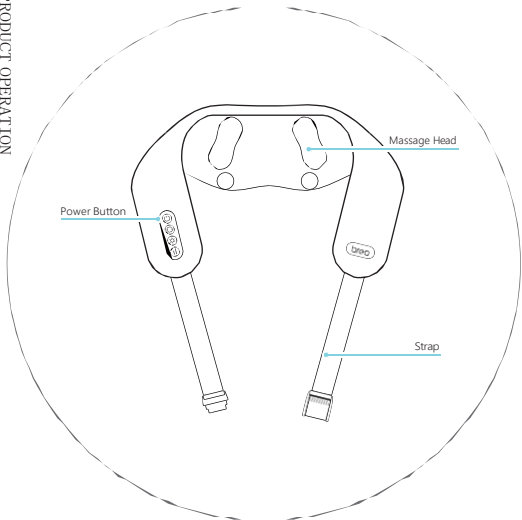

29

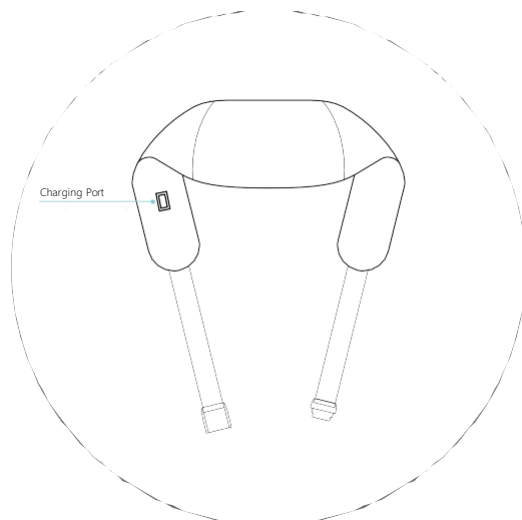

30

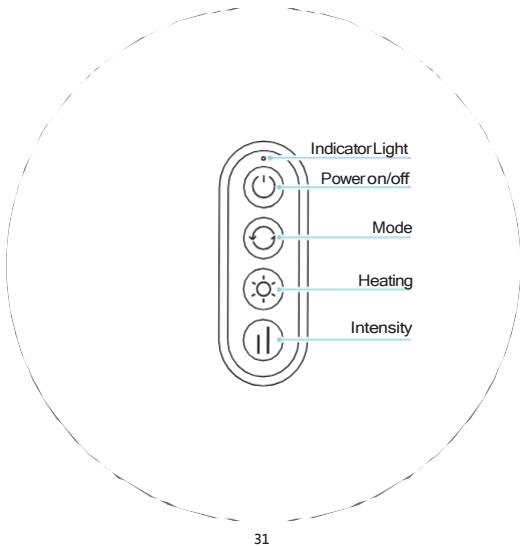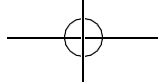

31

EN

Correct Way to Use

Prepare Before Use

Please check whether the product and accessories are complete and free of damage.

Long press the power button for 2 seconds, the product starts to work.

*Mode Selection

The default mode 1 is the default mode when the device is turned on. Short press the mode button to switch between mode 1 and mode 2.

*Heating

The default heating is high-grade when the device is turned on. Short press the heating button to switch gears. The sequence is high heating - low heating - turn off heating, and cycle in this order.

*Intensity



The default intensity is the low intensity when starting up, short press the intensity button to switch between the low intensity and the high intensity.

*Power Off

Automatic shutdown: After 10 minutes of cumulative use, the product stops working.

Manual shutdown: Press and hold the power button for 2 seconds, the product will stop working.

32

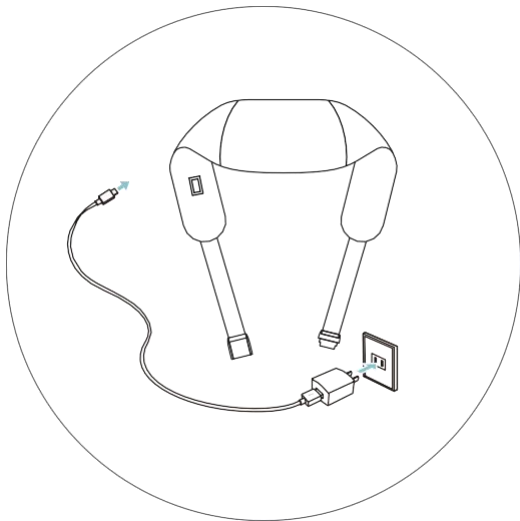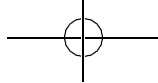


33

EN


Charging

1. This product cannot be used while charging.
 2. Insert the plug of the charging cable into the charging port on the product.
 3. When the product is working, it is connected to the power supply for charging, and the product starts a self-protection program and stops working automatically.
 4. When the battery power is lower than 20%, the indicator light will show a steady red light.
 5. After charging is complete, the indicator light will flash white.
- *It is normal for the device and the power adapter to heat up slightly during charging.
6. After charging, please unplug the power.
 7. This appliance contains batteries that are non-replaceable.

34

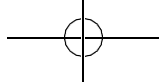

35

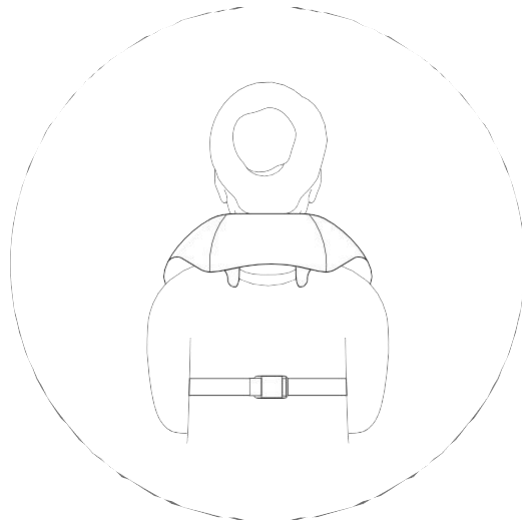

36

EN

TROUBLESHOOTING

Troubleshooting

If there is a failure during operation, please check the following. If the fault cannot be rectified or the product still cannot operate satisfactorily, please cut off the power supply and have it repaired by professional technicians. Do not disassemble the product by yourself!

Problems	Possible Causes	Methods
Fail to power on	Low battery	Please charge it in time.
Intensity is too weak	Low battery	Please charge it in time.
No heating	Heating function is not turned on.	Please turn on the heating function.

If the problem still cannot be solved after troubleshooting, please contact our After-sales Service :

Tel:1-866-546-3587

Mailbox : support@breocare.com

37

EN

PRODUCT DISPOSAL REQUIREMENTS

Notice:

- Before disposing of this product, the battery must be removed as required by regulations;
- When removing the battery, the product must be powered off;
- Batteries should be disposed of safely in accordance with laws and regulations.




Non-professionals do not disassemble by themselves!

Steps to Remove the Battery:

1. Cut off the external power supply and remove the product shell;
2. Unplug the power supply on the battery;
3. Disassemble the battery;
4. Please place waste batteries at a fixed recycling point;
5. Please abide by the local environmental protection laws and disposal guidelines, and care for the environment.

38

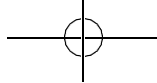

EN

Storage

When the storage period exceeds 3 months, the device should be checked out of the box and recharged before use. Do not place it close to the ground, walls and roofs. The room should be well ventilated to avoid strong sunlight and corrosive gas erosion.

Transportation

The product packaging is equipped with simple shockproof facilities, which are suitable for air, railway, road and ship transportation. Rain, snow, splashing, inversion and collision should be avoided.

Power Check

Check the power supply, if the power supply voltage exceeds the specified range, do not start the machine for use.

39

EN

PROTECT


Maintenance

Take moisture-proof, rodent-proof and anti-pollution measures. Cleaning should be done with the power off. After use, the instrument should be placed in a safe place, and no heavy objects should be placed on the instrument. When a fault occurs, professional technicians should be asked to repair it, and non-professionals are not allowed to disassemble the product by themselves.

Disposal and Recycling

This product should not be thrown away with household waste. When you decide to dispose of this product, please follow local environmental laws and guidelines.

40

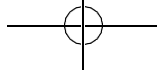

EN

Symbol Description

SYMBOL DESCRIPTION

	Keep Dry		Applied of type BF
	Refer to instruction manual		Serial number
	Manufacturer information: The manufacturer name and address		Disposal in accordance with Directive 2002/96/EC (WEEE)
	Medical device		The device complies with part 15 of the FCC rules
	This device has not been tested for use in an MR environment and should not be used exposed to MR environments while patients are wearing the device. Keep it outside the MRI scanner room.		
	The device complies with compliance with RoHS Directive 2011/65/EU and its amendment directives 2015/863/EU		

41


EN

FCC Statement

FCC STATEMENT

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) this device may not cause harmful interference, and
(2) this device must accept any interference received, including interference that may cause undesired operation.

42

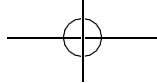

EN

PRODUCT SPECIFICATIONS

Product Name	breo Neck Massager
Model	Neck N5 mini
Rated Voltage	5V 
Rated Power	10W
Power Supply	7.4V  2500mAh rechargeable lithium battery
Net Weight	1200g
Charging Time	≤3.5hours
Working Time	≥60 mins
Product Dimension	L700*W160*H160 (mm)

Shenzhen Breo Technology Co.,Ltd

19/F, Hisense Southern Building, No.1777 Chuangye Road , Nanshan District , Shenzhen , Guangdong,CN,518056

