

CONFIER

Back Support Massage Cushion



User Manual
CF-1313



Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- DO NOT wet, DO NOT use pins, never remove cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- DO NOT use any damaged cables, plugs or loose sockets.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- Misuse or incorrect use excludes any liability for damage.
- Do not use while you sleep.
- Do not use if wet
- Do not use with infants, a helpless person or a person insensitive to heat.
- Overlong application may result in skin burns.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 30 minutes at a time.



PLEASE NOTE !

- DO NOT use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy.
- **The device is not to be used for commercial use or medical purposes.**

Have question ?



Tel: (248) 819-2623
Monday-Friday 9:00AM-4:30PM



Email: supportus@comfier.com

Instruction


Thank you for purchasing Comfier Back Support Massage Cushion
With normal care and proper treatment, it will provide years of reliable service.

Please read all instructions carefully before use this product.

Please retain this manual for future reference.

FEATURES

- 2 pre-programmed massage modes
- 2 levels of Vibration intensity
- 2 heat settings
- Auto 30 minutes timer

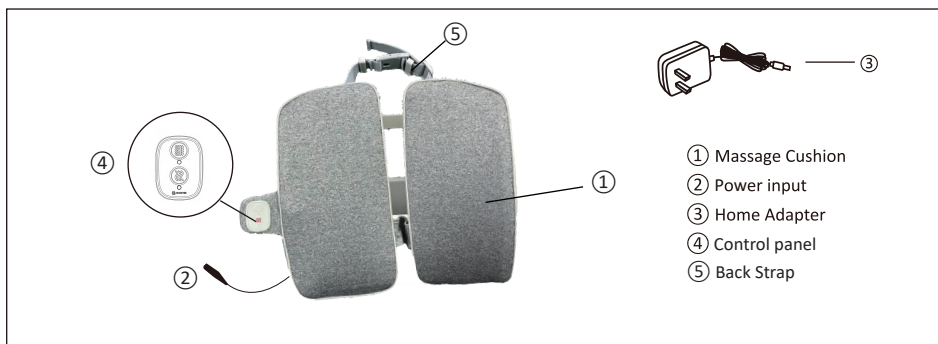

TECHNICAL DATA

Measurements:	16.1x15.4x5.5 inches
Weight:	4.5lbs
Voltage:	Input: AC 100-240V ~50/60Hz Output: 12VDC 1A
Nominal Power:	Max. 12Watt
Automatic runtime:	30 mins

Setup and Operation

1. Place the massage cushion on the chair.
2. Connect the adapter cable to the corresponding cable on the massage cushion.



3. Plug the home adapter into an electric outlet between 100~240V.
 4. Turn on the device by operating controller (refer to controller instructions below.).
 5. When finished unplug the adapter from the outlet.
- The device will automatically switch off after the timer runs out.
 - Do not setup or use the device in a bathroom or similiary wet/damp areas.
 - Ensure the device is not exposed to direct sunlight or placed near sources of heat.

Controller instructions

Vibration Button:

Short press this button to select Vibration Mode 1/ Mode 2/ OFF

The corresponding indicator:

Blue-Mode1


Blue flash-Mode2

Hold the button to cycle between high/low intensity.

The corresponding indicator:

Bright Blue - High intensity

Dark Blue- Low intensity

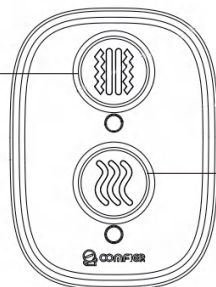

Heat Button:

Short press this button to select High/Low heat level or Heat Off


The heat indicator:

Red -- High level


Dark red-Low level

Installation Steps

Step 1 (Install the Backrest)

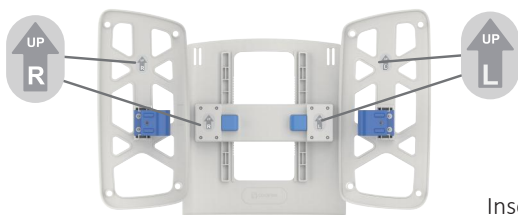

Pay attention to distinguish the left and right, and refer to the sticker label R/ L-correspondence

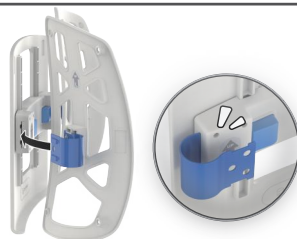


Insert the blue piece into the hole on the frame. Hear a click that means it is firmly inserted (the same operation on the left / right sides). If need, press and hold the blue button to pull out the blue piece.

Step 2 (Install the massage heating pad parts)




The cloth cover connecting strip is on the lower side

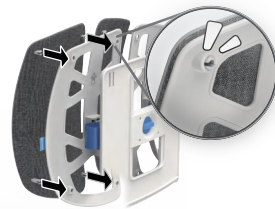


Align the tabs of the cloth cover with the 4 small holes on the frame and press to insert.

Step 3 (Adjust the backrest height)

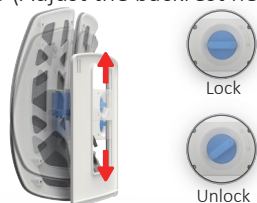

Turn the button to unlocked state and adjust the slider up or down to customize the comfortable backrest height.

Exploded view

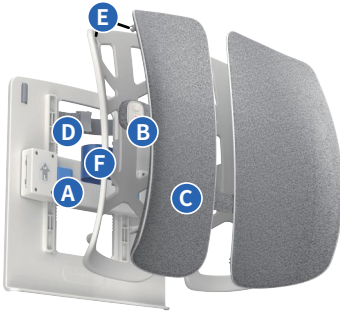

- A** Backrest Removal Button
- B** Control Panel
- C** Massage Heating Pad Part
- D** Elastic Strap
- E** Buckle Connector
- F** Blue Insert

Care and Cleaning Instructions

1. Do not wash in washing machine or immerse the product or the adapter in water.
2. Always switch off before cleaning and remove the adapter from the socket with clean dry hands.
3. The device should be wiped with a slightly moist cloth and dried thoroughly.
Do not wipe the product with alcohol or harsh chemicals such as thinner or bleach.
4. Do not overstretch the elastic straps.
5. Keep away from heat, fire or exposure to direct sunlight.
6. Store the massager in a cool, dry place out of reach of children.
7. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface.

TROUBLESHOOTING GUIDE

Fault	Solutions
Device switches off suddenly	The device will switch off when the timer runs out. Continue using by pressing the 'Power' button. If the device has run for 30 mins., turn back on after resting the unit for 15 mins.
	The device has an overheating protection system be activated, device is switched off for your safety. Please rest device for 15-30mins before restarting.
Device fail to start	Check if device has been used for more than 30 mins. If so, let device rest for 15 mins before continuing. Device may require longer time to cool in warmer environments.
	Check if adapter is plugged into socket.
	Check if adapter output is connected to device power inlet wire.

Warranty

If the product(s) is damaged or defective within a reasonable time, Comfier offers the warranty as below.

Period	Return for refund	Return for replacement	Who paid shipping?
Within 30 Days	Yes	Yes if undamaged	Comfier
31- 90 Days	Yes after being justified	Yes after being justified	Comfier
91 Days - 1 Year	No	Yes after being justified	Buyer

Comfier warrants the product to be free from defects in workmanship and materials, under normal use, for a period of 12 months from the original invoice date, except as noted below. Comfier products can be repaired or replaced within warranty period or returned for refund within 90 days if the products break down in the period of proper use.

The warranty extends only to consumers, NOT to any retailers or resellers. This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. Warranty claims for product defects, unless extended warranty is registered, expire 12 months from the purchase day. Warranty claims CANNOT be processed for items that have exceeded the warranty period. The warranty does NOT cover the replaced, refurbished product. The warranty does have other limitations as below “Warranty Limitations”. To obtain warranty service on Comfier products, follow “Warranty Claim Procedures” below.

Warranty Claim Procedures:

1. Any product the customer believes is defective and seeks for warranty, the customer must provide to Customer Support:
 - a) Valid proof of purchase (see below the definition of valid proof of purchase);
 - b) Pictures and/or description are required to claim defective merchandise;
 - c) Model number and Series Number label attached to each product.
2. If the warranty claim is justified by Comfier Customer Support, customer(s) will obtain from Comfier a Warranty/Return Authorization Number (RMA Number).
3. Comfier reserves the right to specify that items be returned to the designated warehouse for inspection or be inspected by our representative in the field or be retained by customer.
4. Any claim for defective merchandise returns must be packed in original packaging.
5. The item(s) will be repaired or replaced or a refund credit will be issued once the item(s) are delivered. When returning items with RMA# and a prepaid shipping label provided by Comfier, Comfier takes responsibility for any damage or loss incurred in transit.

Valid Proof of Purchase Under this Warranty:

Order numbers from online purchases made through Comfier official store (www.comfier.com) or Comfier's direct-authorized resellers, currently only at Amazon/eBay/Walmart. For product(s) obtained from other Channels, please contact the customer service of the Channel directly for appropriate solution.

Warranty

Warranty Limitations(this warranty does not include):

Normal wear and tear; Any condition resulting from other than ordinary residential wear or any use for which the product was not intended, such as use in rental or contract trade or commercial use; Any condition resulting from incorrect or inadequate maintenance or care; Damage resulting from misuse, abuse, negligence, accidents or shipping damage, exposed to environmental conditions temperatures, water, improper operation,improper installation, improper use of electrical supply; transportation damage; theft,vandalism or;Damage resulting from the attachment or modification of any unauthorized accessories;Fail to show the valid proof of purchase;Lost or stolen products;Items that have expired their warranty period; No quality-related issues (after 30 days of purchase);Free product/ Gift offer by Comfier; unauthorized repairs or modification; Replaced or refurbished product.Dissatisfaction due to buyer's remorse (after 30 days of purchase).

Refund Procedures:

(Refer to <https://www.comfier.com/pages/return-policy> for detail instructions)

Warranty Disclaimer

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. COMFIER SHALL HAVE NO LIABILITY FOR, AND EXPRESSLY DISCLAIMS ANY WARRANTY OR AFFIRMATION OF FACT, EXPRESS OR IMPLIED, OTHER THAN AS SET FORTH IN THIS WARRANTY STATEMENT, INCLUDING, WITHOUT LIMITATION

- (1) THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE;
- (2) ANY WARRANTY OR AFFIRMATION OF FACT RELATED TO MISUSE, IMPROPER SELECTION, RECOMMENDATION, OR MISAPPLICATION OF ANY PRODUCT; AND
- (3) ANY WARRANTY OR AFFIRMATION OF FACT THAT THE CATALOGS, LITERATURE AND WEBSITES IT PROVIDES ACCURATELY ILLUSTRATE AND DESCRIBE PRODUCTS.


Comfier reserves all the right for final interpretation of these Terms of Warranty and Return Policy above and the aright to change, modify, add or remove portions of these terms at any time without prior notification.

After-Sales issues? Warranty Issues?

Please visit:

<http://support.comfier.com>

Or scan the QR code on the right side for tutorial videos,FAQs, manuals, and more information.

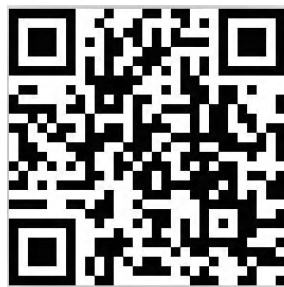

ORDERING REPLACEMENT PARTS

To order additional parts and attachments, visit <http://www.comfier.com>

FCC STATEMENT

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this unit is not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Know more about Comfier



COMFIER TECHNOLOGY CO.,LTD

Address: 573 BELLEVUE RD NEWARK, DE 19713 USA

Tel: (248) 819-2623 (Monday-Friday 9:00AM-4:30PM)

Facebook: www.facebook.com/comfiermassager

Email: supportus@comfier.com

Web: www.comfier.com

*Please provide your order number before contacting customer support