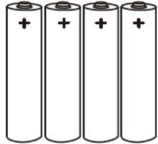
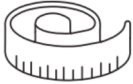



eufy Smart Scale P2 Pro_User Guide (T9149)

What's In the Box

			
Smart Scale P3	AAA Battery (×4)	Soft Measuring Tape (60in / 150cm)	Owner's Manual & Other Documents

Inserting the Batteries

1. Turn the scale over and place it on a flat surface.
2. Open the battery door and insert 4 AAA batteries. Make sure the positive and negative ends are facing the correct polarity directions marked in the battery compartment.
 - The digits will show on the LED display.
3. Press down to close the battery door. Make sure the battery door is fixed in place.

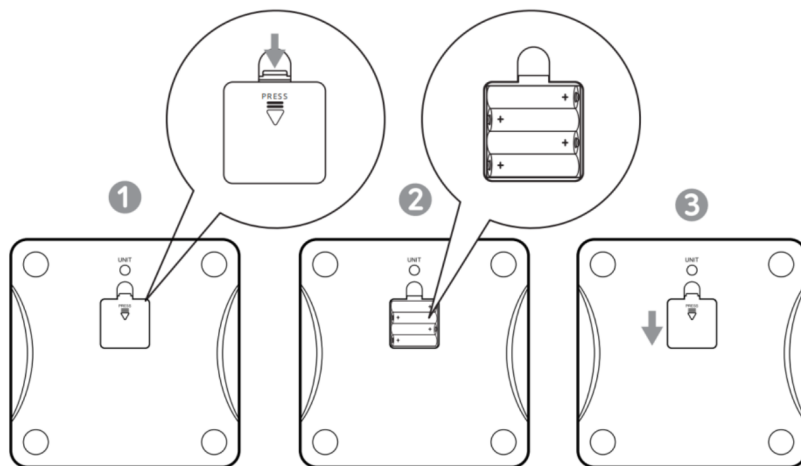


Notes :

- **_Lo** will appear on the LED display when the battery level is low. Replace all batteries at the same time. Do not mix new and used batteries.
- The scale turns off automatically after being idle for 10 seconds.

Selecting the Measuring Unit

1. Press the **UNIT** button at the bottom of the scale to cycle through and select the measuring unit: **lb** (Pound) / **kg** (Kilogram) / **st** (Stone). The measuring unit can be selected and saved within the EufyLife app.

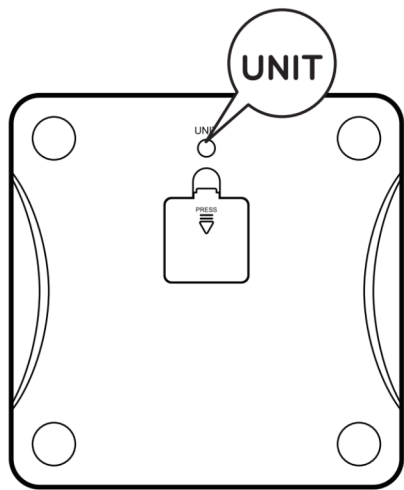


Notes:

- To restore factory settings, press and hold the **UNIT** button at the bottom of the scale for 10 seconds.

Installing the EufyLife App and Pairing

- 1. Download the EufyLife app from the App Store (iOS devices) or Google Play (Android devices).
- 2. Enable Bluetooth on your smartphone or tablet.
- 3. Open the EufyLife app and follow the on-screen instructions to pair the scale with your device.
 - Once paired, the app and scale will automatically connect the next time the app is opened.

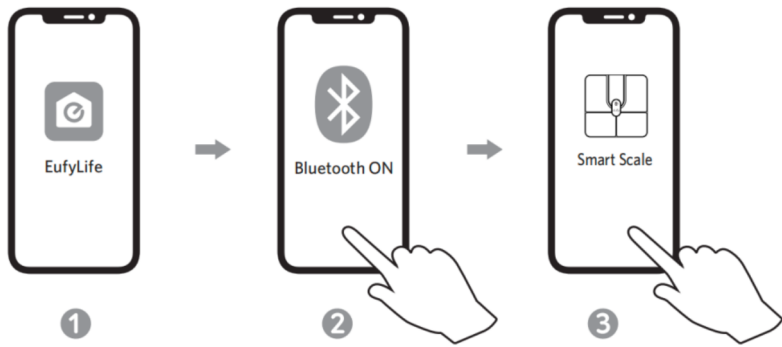

 Bluetooth Status Icon	Status Icon	Indicates
	Flashing	Pairing
	Solid	Connected

 Wi-Fi Status Icon	Status Icon	Indicates
	Off	Not connected
	Flashing slowly	Ready for connecting
	Flashing quickly	Connecting
	Solid	Connected

Notes:

- This scale is compatible with devices running iOS 11.0 and above or Android 5.0 and above.
- To complete the setup, you need to enable Bluetooth on your smartphone or tablet and allow the EufyLife app to use your location.
- Do not pair the scale from your smartphone’s Bluetooth menu. Use the EufyLife app to pair the scale with your smartphone.

Starting a Measurement


Before taking your first measurement, you will need to create and set up a user ID in the EufyLife app.

1. Place the scale on a hard and flat surface.
2. Open the EufyLife app and select " **Add Device** " on the homepage. Follow the in-app instructions to add the smart scale.
3. Step barefoot onto the scale.
 - An analysis of your weight and your body’s composition will begin automatically.
 - After measuring, your weight and other available information will be shown on the LED display. You can find detailed body composition analysis within the EufyLife app.



Using the Measuring Tape

With the measuring tape, you can get accurate measurements according to your requirements.

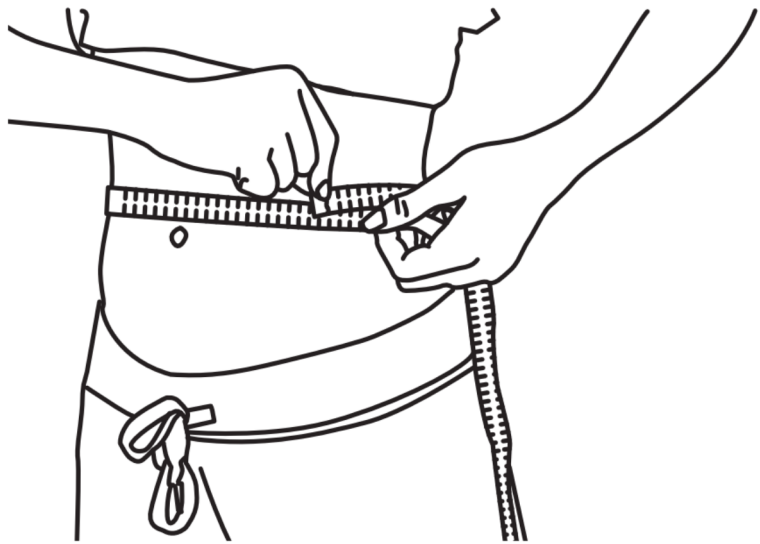

Maintenance

Regularly clean the scale with a dry or slightly dampened soft cloth.



Notes:

- This scale is waterproof, but it is not recommended to wash it with water for a long time or immerse it in water.
- Never use detergents to clean the scale as this may cause it to malfunction or discolor.

Troubleshooting

Error Message

Error Display	Description	Solution
	Overload. The device will power off.	Stop using this scale for measurement.
	Low battery. The device will power off.	Replace all batteries at the same time.

	Wi-Fi setup failure.	• Open the EufyLife app and follow the on-screen instructions to reconnect the scale with your device. • Move the scale to a place near your router before the Wi-Fi setup. • Press and hold the UNIT button for 10 seconds to restore factory settings, then restart the Wi-Fi connection.
---	----------------------	--

When Measuring

Issue	Cause	Solution
Abnormal measurement results: • Too high • Too low • Huge difference between two consecutive measurements	Incorrect posture.	Step onto the scale barefoot and stand still.
	Device is placed on soft carpeting or uneven surface.	Place the device on a hard and flat surface.
LED display does not turn on after stepping onto the scale.	Batteries not installed.	Install 4 new AAA batteries.
	Worn batteries.	Replace all batteries at the same time.
Scale does not proceed to analyse body composition.	Taking a measurement while wearing shoes or socks.	Stand barefoot on the scale.
	Using the weight-only measuring mode in the EufyLife app.	Change the measuring mode in the EufyLife app.

Specifications

Model	T9149
Power	6V (4 × AAA batteries)
LED Display	2.5 × 1.4 in / 64 × 35 mm
Measurement Units	Pound / Kilogram / Stone

Unit Increments	0.1lbs / 0.05kg
Body Fat Unit Increments	0.1 %
Body Water Unit Increments	0.1 %
Muscle Mass Unit Increments	0.1%
Measurement Range	0.2lbs - 397lbs / 0.1kg - 180kg / ost: 0.2lbs - 28st: 5lbs
Activation Weight	3 kg
Working Environment	Temperature: 10°C- 40 °C Humidity: ≤90% RH
Storage Environment	Temperature: -10 °C - 50 °C Humidity: ≤90% RH
Operating Frequency Range	2400 MHz - 2500 MHz
Transmitted Power	Wi-Fi: 18 dBm Bluetooth: 4 dBm
Bluetooth Version	5.0
Wi-Fi Version	2.4 G