User Manual

Product Name:
Core Baby Kit

Recommended Age: 0+ Months

Product Features:
This all-in-one kit covers all baby care essentials, perfect for daily use, keeping your baby healthy, clean, and hygienic.

CONTENTS

Package Contents 01

Cleansing Care 02

Grooming Care 03

Health Care 04


CUSTOMER SERVICE 05

Warranty 06


Warning 07

Package Contents

Cleansing Care

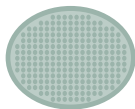Baby Bath Brush

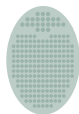Baby Cradle Cap Brush

Finger Toothbrush

Grooming Care

Nail Scissors

Nail Clipper

Nail File

Comb

Health Care

Nasal Aspirator


Gas Relief Tubes

Storage Case

Cleansing Care

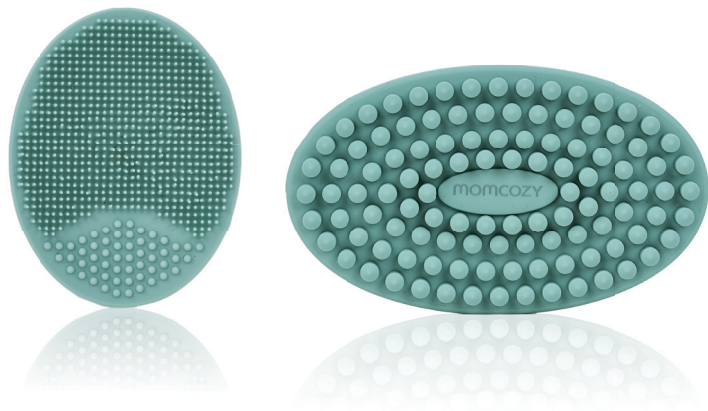Baby Bath Brush

Product Features:

Soft and safe, ideal for cleaning the baby's entire body, removing dirt while protecting delicate skin.

How to Use:

Apply baby body wash on the brush and gently massage the baby's body to clean the skin. Rinse off the body wash, ensuring the baby's skin is clean without residue.

Cradle Cap Brush

Product Features:

Specially designed for the baby's scalp, gently cleans the scalp and helps remove cradle cap, keeping the scalp clean and healthy.

How to Use:

Apply baby shampoo on the brush and gently massage the baby's scalp, being extra gentle around the fontanel. Rinse off the shampoo, ensuring the scalp is thoroughly clean.

Safety

Water Temperature: Ensure the water temperature is around 37° C to avoid overheating or cooling.

Gentle Handling: Always be gentle to avoid irritating the baby's skin or scalp.

Product Cleaning: Clean the brushes with warm water after each use and air dry to prevent bacterial growth.

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Baby Finger Toothbrush

This is a dual-sided finger toothbrush designed with bristles on the front for cleaning your baby's teeth and bumps on the back for tongue cleaning. The thoughtful design allows parents to easily provide full oral care for their baby, ensuring thorough cleaning of both the teeth and tongue.

Product Features:

Front Bristles: Designed for cleaning your baby's teeth and gums. The soft bristles provide gentle oral care, ensuring comfort during use.

Back Bumps: Designed to gently massage and clean the baby's tongue, removing residue and bacteria to maintain oral hygiene.

How to Use:

Wear the Toothbrush: Slide the toothbrush onto your finger, making sure the bristles are covering the fingertip.

Clean Teeth: Gently brush the baby's teeth and gums with the bristled side in soft, circular motions.

Clean Tongue: Use the back bumps to gently massage and clean the baby's tongue.


Clean the Toothbrush: After each use, rinse the finger toothbrush with warm water and mild detergent, and let it air dry.

Safety

Ensure the toothbrush is clean and intact before each use.

Thoroughly clean the toothbrush after each use to prevent bacterial growth.

It is recommended to replace the toothbrush every 3 months to maintain hygiene.

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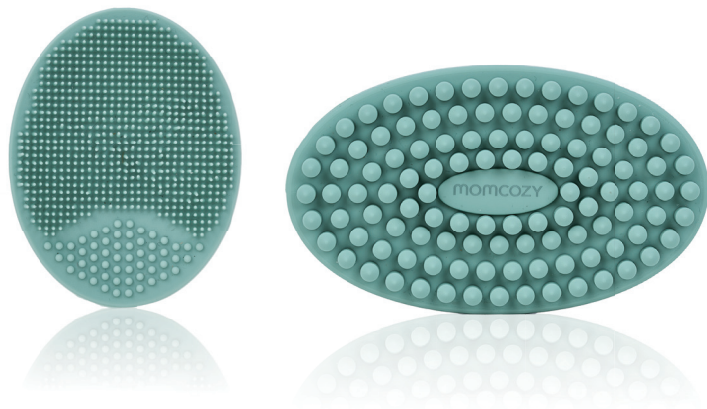Grooming Care

Nail Scissors

Mini blades designed for babies, providing precise nail trimming. Safe and easy to use.

Nail File

Gentle nail file to smooth the edges of the nails, preventing the baby from scratching themselves.

Nail Clipper

Curved design tailored for small baby nails, allowing easy trimming without over-cutting.

How to Use:

- Gently trim the baby's nails using the nail scissors or nail clipper, ensuring smooth cuts without cutting too short
- After trimming, gently use the nail file to smooth the edges of the nails, preventing any rough spots that could scratch.

Safety

Ensure Safety: Be gentle when trimming the baby's nails to avoid cutting the baby's skin.

Clean the Tools: Clean the nail scissors, clippers, and file after each use to prevent bacterial growth.

Avoid Over-Cutting: Don't trim too much of the nail at once to avoid discomfort from overly short nails.

Adult Supervision: Always use these tools under adult supervision to ensure safety

Baby Comb

Curved Design: The comb is designed with a curve to perfectly fit the baby's head, providing a gentle brushing experience that avoids pulling.

Safety

Gentle Handling: Always be gentle when combing, especially around the soft fontanel area.

Safe Storage: Keep the comb out of the baby's reach to avoid accidental use or swallowing.

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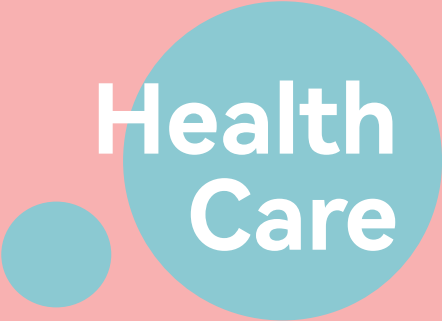
Baby Gas Relief Tubes

Product Features:

Gentle Gas Relief: Designed to help babies release trapped gas from the stomach, alleviating discomfort and fussiness caused by gas buildup.

Safe and Soft Material: Made from medical-grade soft materials to ensure safety and comfort during use.

How to Use:

Clean the gas relief tube before use, rinsing it with warm water.

Apply a small amount of baby-safe lubricant or petroleum jelly to the tip of the tube for easier insertion.

Lay the baby on their back and gently lift their legs as if changing a diaper.

Carefully insert the tip of the gas relief tube into the baby's rectum (no more than 1–2 cm deep).

Hold the tube in place for about 10–15 seconds to allow the gas to pass. Be gentle and avoid forcing the tube.

Gently remove the tube and clean it immediately after use.

Safety

Moderate Use: Use the gas relief tube only when needed. Do not use for more than 10–15 seconds at a time and avoid frequent use within a day.

Hygiene: Clean the gas relief tube thoroughly before and after each use. Regularly disinfect the tube to ensure hygiene.

Use Lubricant: Always apply baby-safe lubricant before use to avoid discomfort or injury during insertion.

Adult Supervision: Always use under adult supervision to ensure the baby's safety.

Emergency: If the baby shows signs of discomfort or if there are any abnormal reactions, stop using the tube immediately and consult a doctor.

Baby Nasal Aspirator

Product Features:

Gentle Nasal Cleaning: The nasal aspirator is designed to gently clean the baby's nasal passages, helping them breathe more easily, especially during colds or congestion.

Soft Silicone Tip: The aspirator comes with a soft silicone nozzle that is safe and non-irritating for the baby's delicate nasal passages.

Easy Squeeze Design: The aspirator features a squeeze bulb that allows parents to control the suction, ensuring gentle and effective removal of nasal secretions.

How to Use:

Before use, clean the aspirator with warm water and mild detergent.

Lay the baby down or keep them in a sitting position, ensuring their head is stable.

Gently insert the soft silicone tip into the baby's nostril, being careful not to insert it too deeply.

Squeeze the bulb to create suction, and release it to gently draw out the mucus.

After each use, clean the aspirator with warm water and mild detergent to maintain hygiene.

Safety

Gentle Handling: Always be gentle when using the aspirator and avoid inserting the tip too deeply to prevent any harm to the baby's nasal passages.

Clean After Each Use: Always clean the aspirator after each use to prevent bacterial growth.

Adult Supervision: Always use the nasal aspirator under adult supervision to ensure safety.

Control Frequency: Use the nasal aspirator as needed and avoid frequent use to prevent irritation to the nasal lining.

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12-Month Limited Warranty
Prolinx Global LTD
Prolinx GmbH
Manufacturer: Hong Kong Lute Technology Co., Limited