momcozy

**ELITE
BABY
KIT**

User Manual

CONTENTS

Package Contents **01**

Cleansing Care **02**

Grooming Care **03**

Health Care **04**


CUSTOMER SERVICE **05**

Warranty **06**


Warning **07**

Package Contents

Cleansing Care

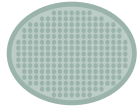Baby Bath Brush

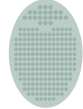Baby Cradle
Cap Brush

Finger
Toothbrush

Comb

Teether

Grooming Care

Electric Nail Trimmer

Electric Nasal Aspirator

Health Care

Heated Tummy Wrap

Thermometer


Medicine
Pacifier

Storage Case

Cleansing
Care

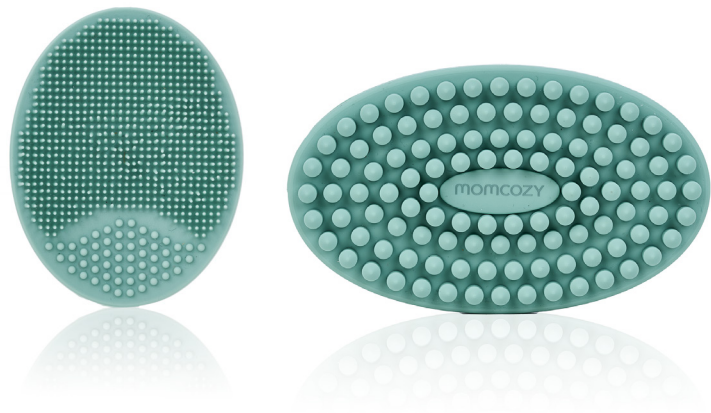Baby Bath Brush

Product Features:

Soft and safe, ideal for cleaning the baby's entire body, removing dirt while protecting delicate skin.

How to Use:

Apply baby body wash on the brush and gently massage the baby's body to clean the skin. Rinse off the body wash, ensuring the baby's skin is clean without residue.

Cradle Cap Brush

Product Features:

Specially designed for the baby's scalp, gently cleans the scalp and helps remove cradle cap, keeping the scalp clean and healthy.

How to Use:

Apply baby shampoo on the brush and gently massage the baby's scalp, being extra gentle around the fontanel. Rinse off the shampoo, ensuring the scalp is thoroughly clean.

Safety

Water Temperature: Ensure the water temperature is around 37° C to avoid overheating or cooling.

Gentle Handling: Always be gentle to avoid irritating the baby's skin or scalp.

Product Cleaning: Clean the brushes with warm water after each use and air dry to prevent bacterial growth.

Adult Supervision: Always use under adult supervision to ensure the baby's safety.

Baby Comb

Curved Design: The comb is designed with a curve to perfectly fit the baby's head, providing a gentle brushing experience that avoids pulling.

Safety

Gentle Handling: Always be gentle when combing, especially around the soft fontanel area.

Safe Storage: Keep the comb out of the baby's reach to avoid accidental use or swallowing.

Baby Finger Toothbrush

This is a dual-sided finger toothbrush designed with bristles on the front for cleaning your baby's teeth and bumps on the back for tongue cleaning. The thoughtful design allows parents to easily provide full oral care for their baby, ensuring thorough cleaning of both the teeth and tongue.

Product Features:

Front Bristles: Designed for cleaning your baby's teeth and gums. The soft bristles provide gentle oral care, ensuring comfort during use.

Back Bumps: Designed to gently massage and clean the baby's tongue, removing residue and bacteria to maintain oral hygiene.

How to Use:

Wear the Toothbrush: Slide the toothbrush onto your finger, making sure the bristles are covering the fingertip.

Clean Teeth: Gently brush the baby's teeth and gums with the bristled side in soft, circular motions.

Clean Tongue: Use the back bumps to gently massage and clean the baby's tongue.

Clean the Toothbrush: After each use, rinse the finger toothbrush with warm water and mild detergent, and let it air dry.

Safety

Ensure the toothbrush is clean and intact before each use.

Thoroughly clean the toothbrush after each use to prevent bacterial growth.

It is recommended to replace the toothbrush every 3 months to maintain hygiene.

Baby Giraffe Teether

Product Features:

Cute Design: Giraffe shape and bright colors attract the baby's attention.

Multiple Chewing Points: 3D chewing points relieve teething discomfort and improve chewing ability.

Safe Material: Made from food-grade rubber, heat-resistant and non-toxic.

UV Sterilization: UV-treated, durable and hygienic.

One-Piece Design: Seamless design, prevents dirt buildup, easy to clean.


Easy to Grip: Anti-slip handle fits baby's hand, easy to hold.

Safety

Clean and sterilize before each use.

Store in a dry place, avoid direct sunlight.

Replace the teether regularly to ensure hygiene.

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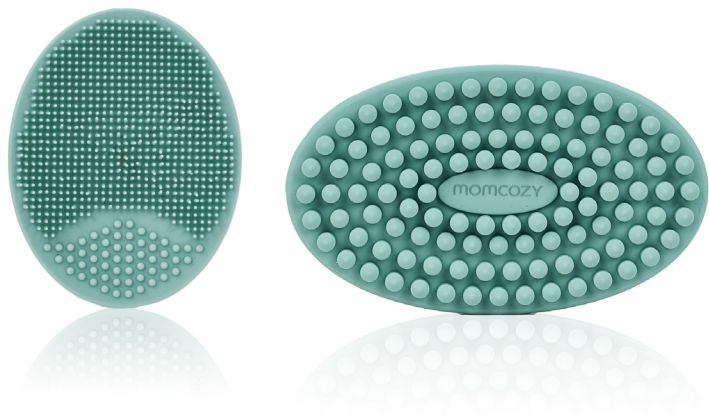Baby Electric Nail File

Product Features:

Rechargeable: Type-C port, no battery replacement needed.

Multiple Heads: 7 grinding heads for different age groups.

Low Noise: Below 35dB, suitable for use while the baby sleeps.

Safe Design: Touch-stop function to protect baby's fingers.

LED Light: 360° illumination, ideal for nighttime use.

How to Use:

Take out a suitable grinding head and install it on the main unit.

Long-press the power button to turn on the main unit, then short-press it to select the appropriate mode.

Place the grinding head close to the nail and adjust the position of the main unit until the nail is trimmed to the desired shape and length.

Operation Instructions:

Product Name: Electric Baby Nail File Kit

Model: NC02

Input voltage: DC 5V $\pm$ 1A

Battery capacity: 600mAh

Charging time: About 1 hour

Working time: About 6 hours

Safety

Clean all parts after each use to ensure hygiene.

Always operate under adult supervision to ensure the baby's comfort.

Do not immerse the product in water to avoid damaging electrical parts.

Only use the provided charging cable, and avoid charging in damp environments.

Baby Electric Nasal Aspirator Instruction

Product Features:

- Suction Levels: Adjustable suction for different nasal needs.
- LED Lights and Music: Provides calming lights and soothing music.
- Safe Material: BPA-free silicone tips, suitable for all ages.
- Removable Storage Cup: Easy to clean, ensuring hygiene.
- Rechargeable: Type-C port with long battery life.
- Portable Design: Includes a storage case for travel.

Operation Instructions:

- Power: Press and hold for 1.5 seconds to turn on/off.
- Suction Level: Short press to switch between levels, long press to pause/resume.
- Light: Short press to toggle the light.
- Music: Short press to turn on music and change tracks, long press for 1 second to turn it off.

How to Use:

- Select the appropriate suction level.
- Gently insert the tip into the baby's nostril to clear mucus.
- Detach and clean the storage cup after use.

Safety

- Clean all parts after each use.
- Always use under adult supervision.
- The device cannot be used while charging.
- Do not submerge the device in water, empty any water residue after cleaning.
- Do not let collected nasal discharge exceed half of the storage cup to avoid overflow.

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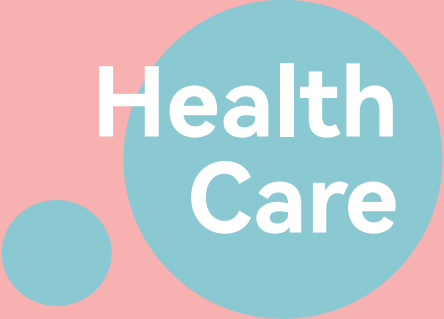
Baby Gas Relief Wrap

Product Features:

The gas relief wrap uses a heated gel pack to apply a warm compress to the baby's abdomen, helping relieve gas discomfort, promoting natural gas release, and providing a comforting warmth.

How to Use:

For Warm Compress: Place the gel pack in the microwave and heat for 15 seconds on medium heat. After heating, take the gel pack out and rub it in your hands to evenly distribute the heat. Insert the appropriately heated gel pack into the designated pocket in the gas relief wrap. Gently wrap the gas relief belt around the baby's waist, ensuring it's comfortably snug but not too tight.

Safety

Temperature Control: Ensure the gel pack is warm, not hot, to avoid burning the baby's delicate skin.
Usage Time: Limit each session to no more than 20 minutes to prevent discomfort from prolonged heat exposure.
Regular Checks: Regularly inspect the gel pack for any damage or leaks to avoid contact with the baby's skin.
Adult Supervision: Always use the product under adult supervision to ensure the baby's safety and comfort.
Storage: After use, store the gel pack in a cool, dry place away from direct sunlight or high heat.

Digital Thermometer

Product Features:

This digital thermometer accurately measures body temperature via oral, axillary, or rectal methods. The LCD display provides clear readings, and the auto-off function conserves battery life.

Measurement:

Oral Measurement: Place the thermometer under the tongue. Approximate measuring time is 1 minute.

Axillary Measurement: Place the thermometer under the armpit. Approximate measuring time is 5 minutes.

Rectal Measurement: Apply water-soluble lubricant to the probe, and gently insert into the rectum. Approximate measuring time is 1 minute.

Reading: After the beep, the final temperature reading will be displayed on the screen.

Turn Off: Press the power button again to turn off, or the device will automatically shut down after about 8 minutes.

How to Use:

Turn On: Press the power button, and the thermometer will beep and begin measuring.

Important Notes

Clean the probe with 75% alcohol before and after each use.

Do not immerse the thermometer fully in water to avoid damage.

Avoid measuring immediately after strenuous exercise, bathing, or sun exposure to ensure accurate readings.

Store the thermometer in a dry, cool place away from direct sunlight.

⚠ WARNING

- **INGESTION HAZARD:** This product contains a button cell or coin battery.
- **DEATH** or serious injury can occur if ingested.
- A swallowed button cell or coin battery can cause
- Internal Chemical Burns in as little as 2 hours.
- **KEEP** new and used batteries **OUT OF REACH** of CHILDREN.
- Seek immediate medical attention if a battery is suspected to be swallowed or inserted inside any part of the body.
- Battery type: LR41 Nominal voltage: 1.5V

Baby Medicine Pacifier

Product Features:

Anti-Choking Design: The dual outlet design ensures a controlled flow rate, reducing the risk of choking during medication administration.

Accurate Dosage Control: The measuring marks on the transparent tube help parents precisely measure the correct amount of medication for safe use.

How to Use:

Fill the feeder with the required amount of medicine, using the measurement marks for precise dosing.

Gently place the nipple in the baby's mouth, allowing the medicine to flow gradually.

Lightly press the silicone bulb at the bottom to ensure a smooth flow of medicine into the baby's mouth.

Safety

Ensure the feeder is clean before use, and wash it immediately after use.

Do not sterilize with hot water or in the microwave. Hand wash with warm water is recommended.

Always administer medicine under adult supervision to ensure the baby's safety.

Regularly inspect the feeder for any damage or leaks before use.

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12-Month Limited Warranty
Manufactured for Hong Kong Lute Technology Co., Limited