NUFACE®

USER MANUAL

ENSURE YOUR BEST RESULTS!

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3) Dual Function Power Button: Turns The Device On/Off

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4) Energizing Hum Indicates The Device Is On (Not Shown)

Energizing Hum Indicates The Device Is On (Not Shown)
4)Zumbido energizante que indica que el dispositivo está encendido (no se muestra)
Betriebsgeräusch (Brummen, nicht dargestellt)
震动声表示仪器已开启 (未显示)
기기에서 느껴지는 진동은 기기가 켜져 있음을 나타냅니다 (표시되지 않음)

5) Protective Cap

Capuchon protecteur
Cubierta de protección
Schutzkappe
보호 캡
보호 캡

6) NuFACE FIX® Line Smoothing Serum 0.5 fl oz / 15 mL

NuFACE FIX® Sérüm Lissant 0.5 oz / 15 mL
Serum NuFACE Line Smoothing FIX 0.5 fl oz / 15 mL
NuFACE Faltenglättendes, Fixierendes Serum 0.5 oz / 15 mL
NuFACE FIX® 淡纹精华 0.5 oz / 15 mL
NuFACE FIX® 세럼 0.25 oz / 7.4 mL

7) Micro USB Cable (Not Shown)

Câble micro USB
Cable micro USB
Micro-USB-Kabel
Micro-USB 连接线 (未显示)
마이크로 USB 케이블

INTENDED USE

The NuFACE FIX® Device is intended for facial skin stimulation and is indicated for over-the-counter cosmetic use.

IMPORTANT INFORMATION

- Use your FIX Device only as described within this User Manual. The risks and dangers of using the FIX Device in any way that does not follow this User Manual might be unknown.
- Your NuFACE FIX® Device has a unique serial number located on the back of the device. Record your serial number and register your device online at www.MuNuFACE.com. A registration link can be found at the bottom of the webpage.

CONTRAINDICATIONS FOR USE

- Do not use this device if you have a history of epilepsy or seizures.
- Do not use if you have had recent surgery or other surgical procedures.
- Do not use if you have an electronic implanted device such as a defibrillator, neurostimulator, pacemaker, or ECG Monitor. If you use a body-worn medical device, such as an insulin pump - consult your healthcare provider.
- Do not use if you have any active cancer or over any suspicious or cancerous lesions.
- Do not use this device if you are pregnant. The safety of electrical stimulation during pregnancy has not been established.
- Do not use on infants or children under the age of 18.

*Consult your healthcare provider if any of the contraindications above apply to you.

⚠ WARNINGS

- If you are under the care of a physician, consult your physician before using this device.
- Apply stimulation only to normal, intact, clean, healthy skin. Do not use over red, swollen, infected, or inflamed areas or skin eruptions, e.g., phlebitis, thrombophlebitis, broken capillaries, varicose veins, etc. or near cancerous or open wounds.
- Do not operate FIX Device in the presence of electronic monitoring equipment (e.g., heart monitors, heart tracing or in close proximity (e.g. 1 m) to shortwave, microwave or electromagnetic therapy equipment as it may cause the medical equipment to malfunction. Such equipment may produce instability in FIX Device output and cause burns at the stimulation site.

Consult your healthcare provider:

- if you experience an adverse reaction from using the device and discontinue use.
- before using the device if you have painful areas.

Do not apply stimulation/use the device:

- across your chest as the introduction of electrical current into the chest may cause rhythmic disturbances to your heart, which could be lethal.
- in the eye area (circular muscle within the orbital rim), on the temples, or crossing over the heart.
- over areas of skin that lack normal sensation or painful areas.
- while bathing or showering.
- when there is a tendency to bleed following acute injury or breaks in bones.

⚠ WARNINGS (CONTINUED)

- while driving, operating machinery, or during any activity in which uncontrollable muscle twitching may put the user or the public at unnecessary risk of injury.
- at extreme temperatures.
- if it's not functioning properly due to damage caused by dropping or accidental submersion in water.

Discontinue use:

- if you suspect the device to be malfunctioning, or not operating properly please contact NuFACE.
- Never drop or insert any object into any openings on the device.
- Keep all electrical appliances (including the FIX Device and any NuFACE device) away from water (including baths, showers, and sinks).
- If the device begins to overheat, turn the device OFF immediately.
- Use of accessories, transducers and cables other than those specified, except for transducers and cables sold by the manufacturer of the FIX Device as replacement parts for internal components, may result in increased emissions or decreased immunity of the system.
- Keep this device out of the reach of children and pets.
- Use of controls or adjustments, or the performance of procedures other than those specified herein, may result in injury.
- Do not disassemble the device as this may cause damage, malfunction, electrical shock, fire, or injury. There are no user-serviceable parts inside. **Note: Disassembly voids the warranty.**

PRECAUTIONS

The NuFACE FIX® Device is designed for cosmetic use only and for individuals in good health. If this is not the case for you, please consult your healthcare provider before use.

- The long-term effects of transcutaneous electrical stimulators are unknown.
- If you have suspected or diagnosed heart disease, you should follow precautions recommended by your healthcare provider.
- Some users may experience skin irritation or hypersensitivity due to the electrical stimulation or Microcurrent Activator. Irritation can usually be reduced by using an alternate NuFACE Microcurrent Activator.
- Skin irritation and burns beneath the electrodes have been reported with the use of other powered muscle stimulators. If you experience skin irritation or burning, discontinue use immediately.
- If using the FIX Device causes persistent pain or discomfort, discontinue use immediately.
- A slight tingling sensation may occur while using the FIX Device. If you feel the current, apply more NuFACE Microcurrent Activator. Lowering the intensity may reduce or eliminate this sensation.
- There may be a perception of flashing lights during the use of your FIX Device, resulting from the stimulation of your optic nerve. If you notice this condition constantly without the use of your device, consult your healthcare provider.
- The battery of the device is not removable by the user. If the device is to be disposed of do not put into a fire. Dispose of the device per local regulations.
- Inspect the device and do not use the device if you see any signs of internal battery leaks.

PRECAUTIONS (CONTINUED)

- The FIX Device is intended to be stored and used at room temperature. Do not store FIX Device in the sun or on a hot surface. High temperatures may be damaging. The device should be operated, transported and stored in the manner specified in this manual. Failure to use and maintain the device in accordance with the instructions outlined in this manual will void the warranty.
- Only use the Micro USB Cable provided with your FIX Device.

The NuFACE FIX® Device is intended for cosmetic use on the face as described. Any consequence resulting from misuse, application to other parts of the body or face, connection to improper voltage sources, dirty spheres, expired or non-NuFACE conductive alternative, or any other improper application is neither the responsibility of NuFACE nor its affiliates.

BEFORE USING YOUR DEVICE

1) Remove your FIX Device and Micro USB cable from the box.

2) Connect your FIX Device to a USB outlet to begin charging, safely away from any source of water, on a table or countertop. A single short pulse will occur when the device is properly connected. The white LED light will blink, showing that your FIX Device is charging. The light will stay illuminated when device is fully charged.

Note: The FIX Device cannot be used while plugged in and will not turn ON while charging.

3) **Important: Charge your FIX Device for 4 hours before the first use.** The white LED light flashes to indicate that the beauty device is charging. When charging is complete, the white light will light up continuously.

4) **Important: Test your skin to make sure you are not sensitive to the NuFACE Microcurrent Activator.** To do so, rub a dime-sized amount of Microcurrent Activator on your forearm and leave for 15-30 minutes before washing it off. If an allergic reaction (rash, redness, etc) develops, immediately contact NuFACE for recommendations on alternative NuFACE formulas. You should not initiate treatment with the FIX Device until your reaction to the Microcurrent Activator has been resolved.

5) When your FIX Device is ON it will emit an energizing hum and the white LEDs will illuminate.

Note: There is only one intensity level to the device.

6) **Important: Always ensure that the two (2) spheres remain in direct contact with the skin.** Avoid pulling or stretching the skin during treatment.

7) The FIX Device energizing hum will emit 3 rapid pulses during treatment when the battery needs recharging.

Note: The FIX Device will NOT contract the muscles. The current is soft, soothing, and mild when used correctly.

DEVICE FEEDBACK — INDICATES

- | | |
|----------------------|--|
| 3 Rapid Pulses ————— | Low Battery. Need to charge device. |
| 2 Long Pulses ————— | Device is charging. |
| 1 Long Pulse ————— | Device is still attached to charger. |
| | Disconnect in order to turn device ON. |

NUFACE MICROCURRENT ACTIVATOR


The essential first step to your microcurrent treatment to ensure optimal results. **Important:** Always apply Microcurrent Activator when using the FIX Device. Refer to packaging for use by date.

TREATMENT AREAS

The FIX Device is designed to smooth, and tighten the appearance of fine lines & wrinkles around the eyes, lips and forehead. Follow instructions below for best results, and refer to your Quick Start for more detailed images.

TREATMENT TECHNIQUES

Ensuring that both spheres are in contact with the skin, apply light to medium pressure while using the FIX Device feathering technique. This quick erasing motion delivers bursts of NuFACE® microcurrent energy to help firm, smooth and tighten skin.

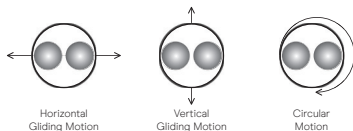HOW TO USE YOUR FIX DEVICE

Get ready to Firm, Smooth and Tighten.

The instructions below are designed to provide the most effective results when followed correctly.

1) Cleanse your face with an oil-free cleanser. Men should shave before treatment — facial hair can interfere with the conductivity of the spheres. **Note: Any cream, lotion, self tanner, and/or oil-based products on the skin can block the microcurrent.**

2) **Important: Apply a mask-like layer of the NuFACE Microcurrent Activator section-by-section as you treat.**

Note: Always apply a NuFACE Microcurrent Activator when using your FIX Device. Reapply if skin becomes dry during treatment to avoid any uncomfortable sensations.

3) Unplug the Micro USB Cable from your FIX Device and turn it ON by pressing the ON/OFF button; the white LEDs will illuminate and the energizing hum will begin when the device is ON.

4) Begin treatment by using the recommended technique shown. Refer to your Quick Start Guide for additional illustrations.

5) You can treat one or more treatment areas for the full 3-minute cycle. After treatment is complete, tap remaining NuFACE Microcurrent Activator into skin.

6) Turn OFF the FIX Device after treatment. See the "CLEANING" section for proper cleaning instructions.

Your FIX Device will automatically turn OFF after 3 minutes.

TREATMENT FREQUENCY

The FIX Device may be used up to 2 times a day per treatment area, 5 days a week.

Can be used as the finishing step to your NuFACE Trinity or Mini treatment.

RESULTS

Results may vary based on skin type and skin condition. Best results are achieved with regular use.

CARING FOR YOUR FIX DEVICE

Cleaning: Slightly dampen a soft, lint-free cloth with water or rubbing alcohol to clean your device after every treatment. **Do not immerse device in water at any time.**

Charging: Before first use, charge your FIX Device for 4 hours. After being fully charged, the device shall operate for up to 40 full 3-minute treatment cycles.

Battery Care: For ideal battery life, charge your FIX Device when the low battery indicator begins (3 rapid pulses during treatment). Note: The rechargeable battery of the FIX Device is not removable by the user. Use caution prior to using the device if it has been stored or not used for an extended period of time. Inspect the device and do not use it if there are any signs of internal battery leaks.

Important: If the device is to be disposed of, do not put into a fire. Dispose of the device and its batteries per local regulations.

Troubleshooting: If your FIX Device is not turning ON, plug it in and recharge device for 4 hours. If your device is still not working, if you need assistance, if you need to report unexpected events, or if you have any questions; please call NuFACE Client Services at 888-360-9171.

WARRANTY

NuFACE® warrants that all NuFACE® FIX devices are free from manufacture defects in materials and workmanship for a period of one (1) year from the date of your original purchase. Our obligation under this Warranty is limited to replacement of the product which will not extend the period of your Warranty. This Warranty does not cover damage covered by abuse, improper care, application of harmful chemicals, product that may have been altered in any way, loss of the device and product purchased from an unauthorized seller. NO OTHER WARRANTY EXPRESSED OR IMPLIED, WHETHER OF FITNESS FOR A PARTICULAR PURPOSE OR OF MERCHANTABILITY OR OF ANY OTHER KIND, SHALL EXIST IN RESPECT TO SUCH PRODUCT.

For Devices Purchased in the US: For more information about your warranty or to initiate a warranty claim and receive a return authorization number, contact our Client Services Team via our website at www.MyNuFACE.com/warranty or by calling (888)360-9171. Must show receipt from an authorized seller as proof of purchase.

For Devices Purchased Outside of the US, International Returns and International Warranties: For all product purchased outside of the United States, warranty and returns are per your place of purchase—spa, retailer or distributor. Please consult your receipt or place of purchase for their return and warranty policy.

Conforms to EN IEC/UL: 60601-1, 60601-1-6, 60601-1-11, 60601-2-10, 60601-1-2, 60601-1-4, 60335-1, 60335-2

SYMBOLS

One or more of the following markings may appear on the FIX Device, User Manual or packaging.

European Conformity

European Representative

Part Number

Refer to Instruction Manual

Manufacturer

Comply with European Restriction of Hazardous Substances Directive 2002/95/EC

Product must be disposed of in accordance with the WEEE directive

Type BF applied part, which is equipment that provides a particular degree of protection against electric shock. Particularly regarding allowable leakage current and of the protective earth connection.

Class II double insulated electrical device

Ingress protection rating: Protected from touch by fingers and objects greater than 12 millimeters. Protected from water spray less than 15 degrees from vertical.

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All Rights Reserved. / Tous droits réservés. / Todos los derechos reservados. / Alle Rechte vorbehalten. / 版权所有 / 모든 권리 보유.
Medical Device Safety Service GmbH