NŪFACE TRINITY®

EFFECTIVE LIP AND EYE ATTACHMENT

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NuFACE TRINITY® EFFECTIVE LIP AND EYE ATTACHMENT

Accessoire efficace pour lèvres et yeux Trinity
Trinity Effektiver Lippen- und Augenaufsatz
Trinity 眼唇提拉替换头 / Trinity 트리니티 엘르



Note: Refer to NuFACE Trinity® User Manual for more device information.

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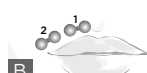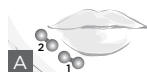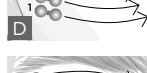
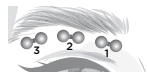
LIFT - Glide/hold both spheres sequentially over each treatment area position as shown for 1 beep (5 seconds). Complete each sequence 3 times, then repeat on the other half of face. Do not place spheres directly on or into the nose. Do not place spheres directly into tear duct or on eyelid.

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HELLO FROM NuFACE®!

Congratulations on receiving your Trinity Effective Lip and Eye Attachment and achieving a new level of anti-aging precision.

We're excited for you to experience the accuracy of its dual-wand design that delivers targeted microcurrent around the eyes and lips. Simply swap out your Trinity Facial Trainer Attachment and let the Trinity Effective Lip and Eye Attachment do the rest.

As always, don't forget to take before and after photos and tag us at @MyNuFACE. We'd love to follow your skincare journey!

Have an uplifting day,
Your NuFACE Family

INTENDED USE

The Trinity Effective Lip and Eye Attachment is intended for facial stimulation and is indicated for over-the-counter cosmetic use.

IMPORTANT INFORMATION

- Use the Trinity Effective Lip and Eye Attachment only as described within this User Manual. The risks and dangers of using it in any way that does not follow this User Manual might be unknown.
- Your Effective Lip and Eye Attachment has a unique serial number located on the underside of the attachment. Record your serial number and register your device online at www.MyNuFACE.com. A registration link can be found at the bottom of the webpage.
- In addition to guidance within this manual, also refer to your NuFACE Trinity® Device manual.

CONTRAINDICATIONS FOR USE

- Do not use this device if you have a history of epilepsy or seizures.
- Do not use if you have had recent surgery or other surgical procedures.
- Do not use if you have an electronic implanted device such as a defibrillator, neurostimulator, pacemaker, or ECG Monitor. If you use a body-worn medical device, such as an insulin pump – consult your healthcare provider.
- Do not use if you have any active cancer or over any suspicious or cancerous lesion.
- Do not use this device if you are pregnant. The safety of electrical stimulation during pregnancy has not been established.
- Do not use on infants or children under the age of 18.

*Consult your healthcare provider if any of the contraindications above apply to you.

⚠ WARNINGS

- If you are under the care of a physician, consult your physician before using this device.
- Discontinue using the device and consult your healthcare if you experience an adverse reaction from using the device.
- Apply stimulation only to normal, intact, clean, healthy skin. Do not use over red, swollen, infected, or inflamed areas or skin eruptions, e.g., phlebitis, thrombophlebitis, broken capillaries, varicose veins, etc.
- Consult your healthcare provider before use if you have any underlying skin conditions.
- Do not use the Effective Lip and Eye Attachment on the groin, or crossing over the heart.
- Do not apply stimulation across your chest as the introduction of electrical current into the chest may cause rhythmic disturbances to your heart, which could be lethal.
- Do not operate the Effective Lip and Eye Attachment in the presence of electronic monitoring equipment (e.g., heart monitors, heart tracing) or in close proximity (e.g. 1 m) to shortwave, microwave or electromagnetic therapy equipment as it may cause the medical equipment to malfunction. Such equipment may produce instability in the device output and cause burns at the stimulation site.
- Do not use this device when there is a tendency to bleed following acute injury or breaks in bones.
- Do not apply stimulation over areas of skin that lack normal sensation.
- Device should not be used at extreme temperatures.
- Do not apply stimulation over painful areas. If you have painful areas, you should consult with your healthcare provider before using this device.
- Keep all electrical appliances (including the Trinity Effective Lip and Eye Attachment) away from water (including baths, showers, and sinks).
- Do not apply stimulation while bathing or showering.

⚠ WARNINGS (CONTINUED)

- Do not apply stimulation while driving, operating machinery, or during any activity in which uncontrollable muscle twitching may put the user or the public at unnecessary risk of injury.
- If the device begins to overheat, turn the device OFF immediately.
- Keep this device out of the reach of children and pets.
- Never operate the device if it is not functioning normally, or not operating properly due to damage caused by dropping or accidental submersion in water.
- Discontinue using the device immediately if you suspect the device to be malfunctioning, or not operating properly. Return the device to NuFACE for examination.
- Never drop or insert any object into any openings on the device.
- Use of accessories, transducers and cables other than those specified, except for transducers and cables sold by the manufacturer of the Effective Lip and Eye Attachment as replacement parts for internal components, may result in increased emissions or decreased immunity of the system.
- Use of controls or adjustments, or the performance of procedures other than those specified herein, may result in injury.
- Do not take the device apart as this may cause damage, inappropriate function, electrical shock, fire or injury. The device parts cannot be replaced or repaired by users. **Note: Disassembly voids the warranty.**

PRECAUTIONS

The Effective Lip and Eye Attachment is designed for cosmetic use only and for individuals in good health. If this is not the case for you, please consult your healthcare provider before use.

- The long-term effects of transcutaneous electrical stimulators are unknown.
- If you have suspected or diagnosed heart disease, you should follow precautions recommended by your healthcare provider.
- Some patients may experience skin irritation or hypersensitivity due to the electrical stimulation or conductive gel. Irritation can usually be reduced by using an alternate NuFACE conductive medium.
- Skin irritation and burns beneath the electrodes have been reported with the use of other powered muscle stimulators. If you experience skin irritation or burning, discontinue use immediately.
- If using the Effective Lip and Eye Attachment causes persistent pain or discomfort, discontinue use immediately.
- A slight tingling sensation may occur while using the Effective Lip and Eye Attachment. Lowering the intensity may reduce or eliminate this sensation.
- There may be a perception of flashing lights during the use of your Effective Lip and Eye Attachment, resulting from the stimulation of the eye area. If you notice this condition constantly without the use of your device, consult your healthcare provider.
- Some users may experience a metallic taste in their mouth during the use of your device, resulting from stimulation and varying sensitivity around the mouth. If you notice this condition constantly without the use of your device, consult your healthcare provider.
- The Trinity Effective Lip and Eye Attachment is intended to be stored and used at room temperature. Do not store the attachment in the sun or on a hot surface. High temperatures may be damaging. The device should be operated, transported and stored in the manner specified in this manual. Failure to use and maintain the device in accordance with the instructions outlined in this manual will void the warranty.
- Do not store the device in the sun or on a hot surface. High temperatures may be damaging. The device should be operated, transported and stored in the manner specified in this manual. Failure to use and maintain the device in accordance with the instructions outlined in this manual will void the warranty.

PRECAUTIONS (CONTINUED)

The NuFACE Trinity Effective Lip and Eye Attachment is intended for cosmetic use on the face as described. Any consequence resulting from misuse, application to other parts of the body or face, connection to improper voltage sources, dirty spheres, expired or non-NuFACE conductive alternative, or any other improper application is neither the responsibility of NuFACE nor its affiliates.

GETTING STARTED

Follow these important steps before using your new attachment.

- 1) Follow steps 1-5 of the Trinity Manual "Getting Started" section.
- 2) Always firmly press the attachment release button on the back of the device to attach or detach the treatment attachment. Refer to #4 in the "Getting To Know Your Device" section. Guide the attachment in place and then release the button to secure the attachment to your Trinity device. The device will not turn on if the attachment is not securely in place.
- 3) Remember to take a before photo of your cleansed face in a well-lit room. After your first treatment, take an after photo and repeat the process in the same well-lit room every eight weeks to track your results.

DEVICE SOUND — INDICATES

Ascending Beeps	Powered ON
Descending Beeps	Powered OFF
Single Beep	Intensity Adjusted
Dual-Tone Beep	End of Glide Interval
Low Tone Beep	Battery Needs To Be Charged

TREATMENT AREAS

The Trinity Effective Lip and Eye Attachment is designed to stimulate smaller surface areas of the face, such as the lip and eye areas. The instructions below are designed to provide users with the most effective results when followed correctly.

TREATMENT TECHNIQUES

Always ensure that you are correctly holding the device in either a vertical or horizontal position, as required per the specified treatment area. Reference the Treatment Key in the Quick Start on how to properly position the spheres during treatment.

HOW TO USE YOUR EFFECTIVE LIP AND EYE ATTACHMENT

Get ready to reduce the look of expression lines around the eyes and lips.

- 1) Cleanse your face with an oil-free cleanser. Men should shave before treatment – facial hair can interfere with the conductivity of the spheres. **Note: Any cream, lotion, self-tanner, and/or oil-based products on the skin can block microcurrent.**
- 2) Apply a NuFACE Gel Primer (sold separately) on your clean, dry skin section-by-section as you treat. **Note: Always apply a NuFACE Gel Primer when using your device and attachment.** Reapply if skin becomes dry during treatment to avoid any uncomfortable sensations.

HOW TO USE YOUR EFFECTIVE LIP AND EYE ATTACHMENT (CONTINUED)

3) Remove your device from the charging cradle and turn it on by pressing the ON/OFF button. An ascending audible beep is emitted and the lights will turn on notifying you that your device is ON.

4) Adjust the intensity level to your comfort by pressing the intensity adjustment buttons on the front of your device, either up (+) or down (-). As you increase or decrease intensity, the number of illuminated blue lights will change. At a maximum intensity, all five lights will be illuminated. At a minimum level, only one light will be illuminated.

5) Place the Effective Lip and Eye Attachment in direct contact with the treatment area. Always keep both spheres on your face at the same time and glide or hold with light to medium pressure.
Note: You may not feel any sensation during treatment, depending on the intensity level and amount of Gel Primer used. This is normal and the device is still working.

6) The device emits a dual-tone beep every five (5) seconds to indicate when to move to the next treatment area. This beep will change to two (2) low tones when the battery is in need of recharging.
Note: The Trinity Device will automatically shut off after 20 minutes of use.

7) After treatment is complete, massage remaining NuFACE Gel Primer into skin.

TREATMENT FREQUENCY

- Use on targeted areas of concern (3 minutes for the lips, 5 minutes for the eyes).
- For the first 60 days, use 5 days per week.
- After 60 days, use 2-3 times per week for maintenance.
- Allow 24 hours between treatments.

RESULTS

Results may vary based on skin type and skin condition. Optimal results are achieved with regular use.

CARING FOR YOUR EFFECTIVE LIP AND EYE ATTACHMENT

CLEANING

Slightly dampen a soft, lint-free cloth with water or rubbing alcohol to clean your device and treatment attachment. Do not immerse device or treatment attachment in water at any time.
Do not immerse device or treatment attachment in water at any time.

TROUBLESHOOTING

If your Trinity Device is not turning on, please make sure all of the following are true:

- The Effective Lip and Eye Attachment is completely secured to the device and fits snugly. Make sure to firmly press the attachment release button on the back of the device when attaching and detaching the Effective Lip and Eye Attachment.
 - The Trinity Device battery is fully charged. If not, please charge your device for a full 12 hours.
- If your device is still not working, please call NuFACE client services at 888-360-9171. We're here to help!**

WARRANTY

For Devices Purchased in the US: NuFACE warrants that the NuFACE® devices and attachments are free of Manufacture defects in materials and workmanship for (1) year from the date of purchase for the NuFACE Trinity® Device, (2) years for the Trinity PRO Device. Abnormal wear and tear, submersions and drops are not covered by this warranty. This warranty is void if the product has been subject to accident, misuse, abuse, improper maintenance or repair, unauthorized modification, or any other use not found in the printed instructions. **Must show receipt as Proof of Purchase.**

For more information on your warranty, visit
www.MyNuFACE.com/pages/warranty-returns

Procedure for Warranty Repairs or Replacement: If warranty service is necessary you can download a Warranty Form on our website at: www.MyNuFACE.com/warranty or call Client Services at 888-360-9171 to receive further instructions. Please be sure to return all parts (device, cradle and charger) and include a completed Warranty Form.

NuFACE - Returns / Repairs
1325 Sycamore Ave, Suite A
Vista, CA 92081
USA
1-888-360-9171

For Devices Purchased Outside of the US, International Returns and International Warranties:

All product returns outside of the United States are per your place of purchase - spa, retailer or distributor. Please consult your receipt or place of purchase for their return and warranty policy.

Conforms to EN IEC/UL : 60601-1, 60601-1-6, 60601-1-11, 60601-2-10, 60601-1-2, 60601-1-4, 60335-1, 60335-2

SYMBOLS

One or more of the following markings may appear on the Trinity Device, User Manual or packaging.

	European Conformity		Comply with European Restriction of Hazardous Substances Directive 2002/95/EC
	European Representative		Product must be disposed of in accordance with the WEEE directive
	Part Number		FCC-approved equipment authorization mark
	Refer to Instruction Manual		Class II double insulated electrical device
	Manufacturer		
	Type BF applied part, which is equipment that provides a particular degree of protection against electric shock. Particularly regarding allowable leakage current and of the protective earth connection.		

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All Rights Reserved. / Tous droits réservés. / Alle Rechte vorbehalten. / 版权所有 / 모든 권리 보유.
NuFACE® is a registered trademark of the Carol Cole Company.
Medical Device Safety Service GmbH