

GLEE



English

IMPORTANT SAFEGUARDS

Your GLEEM toothbrush has been carefully designed to offer you a unique brushing experience that is both safe and effective. When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

Read all instructions before using.

DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

- Close supervision is necessary when this product is used by, on, or near children or persons with disabilities.
- Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer. Only use charger provided with your appliance.
- Keep cord away from heated surfaces.
- Never drop the appliance or insert any object into any opening of the appliance.
- Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.
- This appliance contains batteries that are non-replaceable.
- Do not disassemble the product except when disposing of the battery. When taking out the battery for disposal of the unit, use caution so as not to short the positive (+) and negative (-) terminals.
- Cleaning and maintenance shall not be performed by children.
- Children shall not play with the appliance.


IMPORTANT

- Periodically check the entire product/cord for damage. Avoid sharp bending, jamming, and squeezing of cord. A damaged or non-functioning unit should no longer be used. If product/cord is damaged, contact GLEEM Customer Service.
- Do not touch the power plug with wet hands. This can cause an electric shock.
- When unplugging, always hold the power plug instead of the cord.
- Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact GLEEM Customer Service.
- The appliance has no user serviceable parts. Do not open and disassemble device except when disposing of the battery. Follow handling and disposal instructions of the battery/device described in this manual.
- Do not use on children under age 3.
- This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
- It is recommended to use individual brush handles for each family member.
- This appliance contains magnets and generates magnetic fields. If you use a medical device or have an implanted medical device and you have questions about the safe usage of this toothbrush, consult the medical device manufacturer or a physician.
- For tooth brushing only. Do not chew or bite on brush. Do not use as a toy or without adult supervision.
- To avoid brush head breakage which may generate a small-parts choking hazard or damage teeth:
 - Prior to each use, ensure the brush head fits well. Discontinue use of the toothbrush if brush head no longer fits properly. Never use without a brush head.
 - If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
 - Replace brush head every 3 months or sooner if brush head becomes worn.
 - Clean brush head properly after each use (see CLEANING). Proper cleaning ensures the safe use and functional life of the toothbrush.


CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.





SAVE THESE INSTRUCTIONS

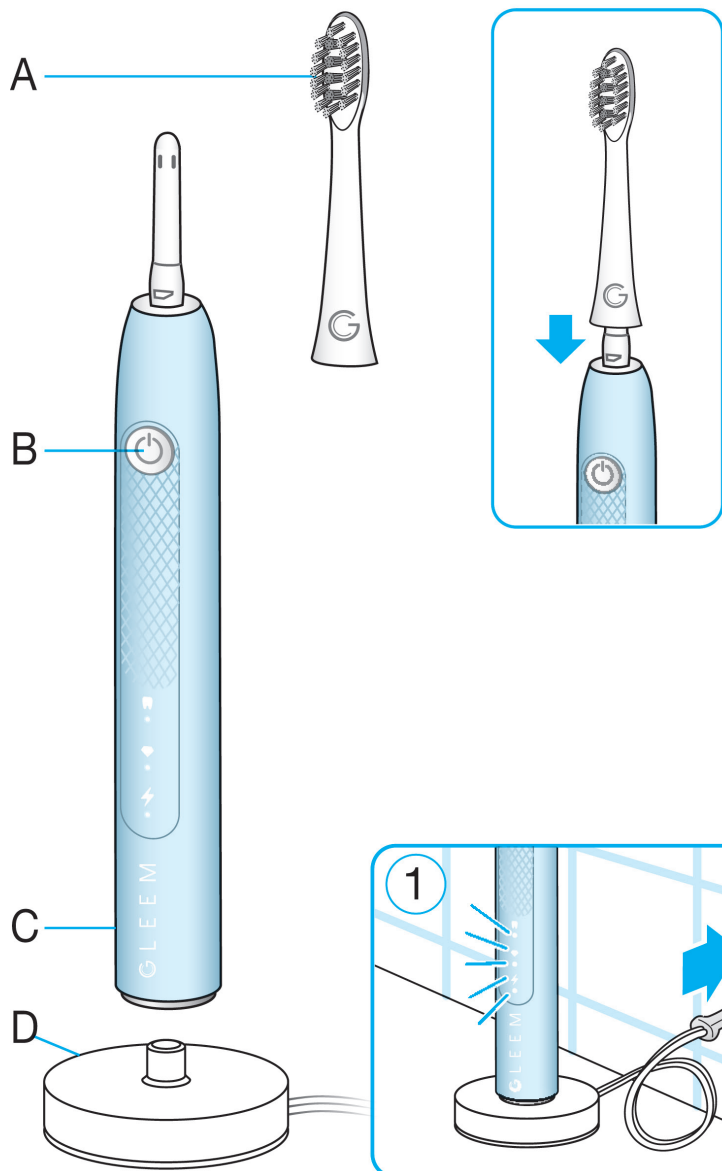

2

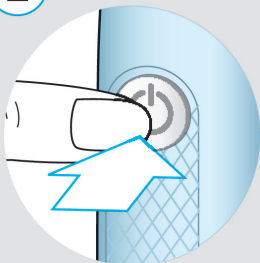

I / on

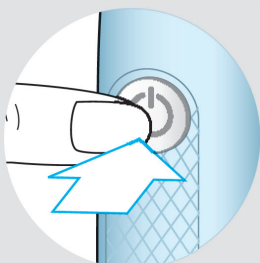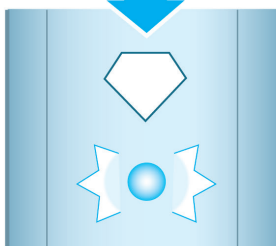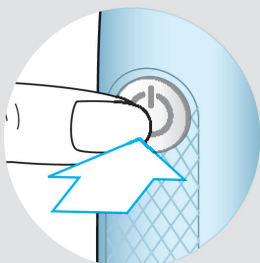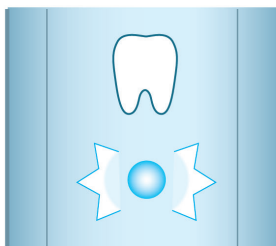

O / off

DESCRIPTION

- A Brush head
- B On/off button
- C Handle
- D Charging unit

Accessories:

E Travel Case

Note: Content may vary based on the model purchased.

Specifications

For voltage specifications please refer to the bottom of the charging unit. Noise level: ≤ 65 dB (A)

CHARGING AND OPERATING

- Plug the charging unit (D) into an electrical outlet and place the toothbrush handle (C) on the charging unit. A full charge takes about 22 hours and enables up to 14 days of regular brushing (twice a day for 2 minutes).
- Recharge the toothbrush when brush battery is running low. Overcharging is impossible.
- To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

USING THE TOOTHBRUSH

- Wet brush head and apply any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching the appliance on (B).
- When brushing your teeth with one of the GLEEM brush heads, guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface.
- Apply light pressure and start brushing in back and forth motions, just like you would do with a manual toothbrush.
- Start brushing the outsides, then the insides and finally the chewing surfaces. Brush all four quadrants of your mouth equally. Do not press too hard, simply let the brush do all the work. Brush at least 2 minutes to thoroughly remove plaque. You may also consult your dentist or dental hygienist about the right technique for you.
- Turn off the handle by pressing the on/off switch.

MODES

Your GLEEM Rechargeable Toothbrush has 2 modes: Classic Clean and GLEEM Clean. Classic Clean mode provides sonic bristle vibrations for daily use. GLEEM Clean mode provides additional sonic bristle vibrations for a brilliant smile.

To turn on Classic Clean mode, press the on/off button once; the LED light below the tooth icon will turn on. To change to GLEEM Clean mode, press the on/off button again (while brush is on Classic Clean mode); the LED light below the diamond icon will turn on. To turn off completely, press the on/off button again (while brush is on GLEEM Clean mode); the LED light will turn off. An incorporated timer memorizes the elapsed brushing time, which allows you to pause or resume during your brushing cycle. To pause, press and hold the on/off button during the brushing cycle until the motor stops; the brush will stop vibrating. To resume, press the on/off button again (within 2 minutes) to continue where you left off during the brushing cycle. GLEEM timer: A short stuttering sound at 30-second intervals reminds you to brush equally all four quadrants of your mouth. A long stuttering sound indicates the end of the professionally recommended 2-minute brushing time.

INDIVIDUAL BRUSH HEADS

GLEEM brush heads can be used with any GLEEM brush handle, including GLEEM Rechargeable Toothbrushes and GLEEM Power Toothbrushes (battery-powered). GLEEM brush heads are designed for everyday cleaning and are made with soft, Nylon bristles for a gentle clean.

REPLACEMENT BRUSH HEADS

GLEEM does not recommend the use of non-GLEEM replacement brush heads.

- GLEEM has no control over the quality of non-GLEEM replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-GLEEM replacement brush heads.
- GLEEM cannot ensure a good fit of non-GLEEM replacement brush heads.
- GLEEM cannot predict the long-term effect of non-GLEEM replacement brush heads on handle wear.

All GLEEM replacement brush heads carry the GLEEM logo and meet GLEEM's high quality standards. GLEEM does not sell replacement brush heads or handle parts under any other brand name.

CLEANING


After brushing, rinse your brush head for several seconds under running water with the handle switched on. Switch off the handle and remove the brush head. Rinse both parts separately under running water then wipe them dry before reassembling the toothbrush. From time to time also clean the charging unit with a damp cloth; never place the charging unit in water.

ENVIRONMENTAL NOTICE

Product contains batteries and/or recyclable electric waste. For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations. When taking out the battery for disposal, use protective clothing (gloves, safety glasses, etc.) to prevent skin contact with leaking batteries.

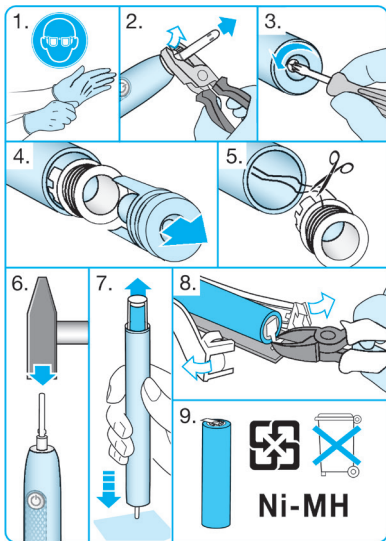

Caution! Opening the handle will destroy the appliance. Subject to change without notice.

SATISFACTION GUARANTEED OR YOUR MONEY BACK:

Try it with no risk for 60 days! If you are not 100% satisfied, call 1-855-210-8802 to receive a full refund. Return the product in its original packaging along with the original cash receipt postmarked within 60 days of purchase to receive a pre-paid card in the amount of purchase. Terms of prepaid card apply. Your GLEEM toothbrush must be returned in its original packaging and be postmarked no later than 60 days after the date of purchase. Please allow 6-8 weeks for handling. In the U.S., the shipment cost to return the product is the responsibility of the consumer. Sorry: No refunds will be made for any GLEEM product damaged by accident, neglect, unreasonable use or lost in the post. Terms and conditions apply. For further details please contact GLEEM Customer Service.

FOR USA ONLY

CUSTOMER SERVICE

For questions or to speak with a GLEEM Customer Service Representative, call 1-855-210-8802.

