WATER FLOSSER ADVANCED

powered by **BRAUN**

91449346

WATER FLOSSER ADVANCED

Internet:

www.oralb.com

www.oralb.ca











www.braun.com

www.service.braun.com

Braun GmbH
Frankfurter Strasse 145
61476 Kronberg,
Germany/Allemagne

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91449346

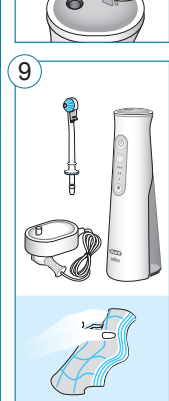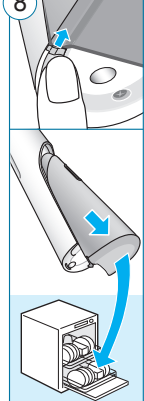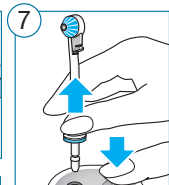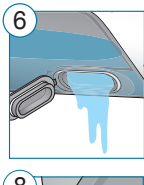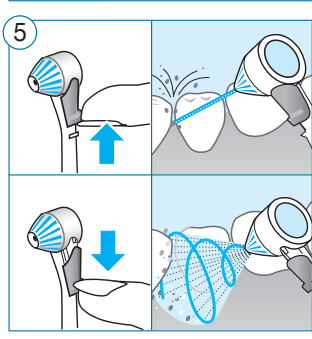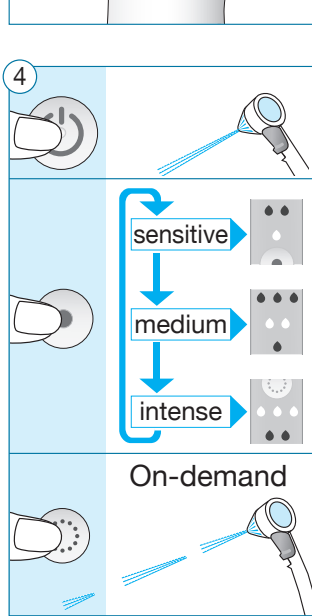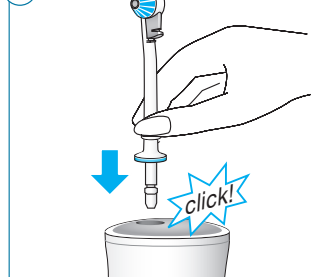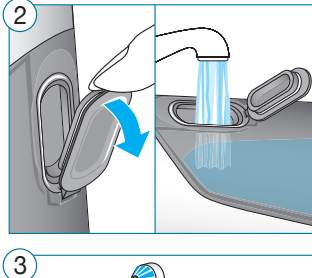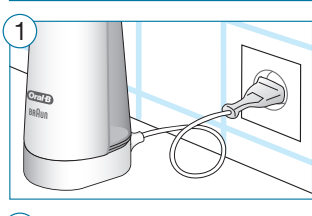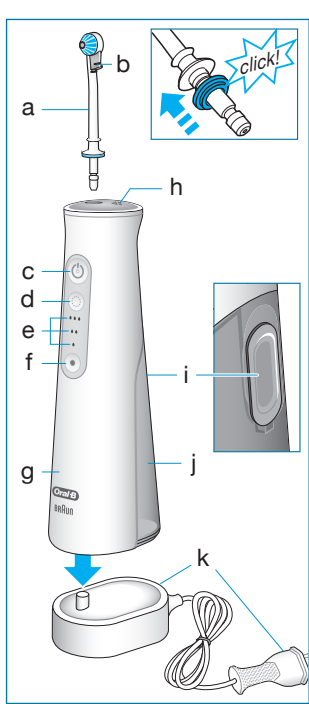IMPORTANT SAFEGUARDS

Your Oral-B® irrigator has been carefully designed to offer you and your family a unique flossing experience that is both safe and effective.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

Read all instructions before using.

DANGER

In order to reduce the risk of electrocution:

1. Do not use while bathing.
2. Do not place or store the charger where it can fall or be pulled into a tub or sink.
3. Do not place the charger in water or other liquid.
4. Do not reach for a charger that has fallen into water. Unplug immediately.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by, on, or near children or persons with disabilities.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.
6. Do not use on self when drowsy or on another who is drowsy or sleeping.
7. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.

IMPORTANT

1. Periodically check the entire product/cord for damage. If the product/cord is damaged, take it to an Oral-B Service Centre. A damaged or non-functioning unit should no longer be used. Do not modify or repair the product. This may cause fire, electric shock or injury.
2. The oral irrigator is not intended for use by children under age 8. Children aged from 8 years and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge can use the product, if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
3. Children shall not play with the appliance.
4. Cleaning and user maintenance shall not be performed by children.
5. Oral irrigators shall be used only with liquids defined in the manufacturer's instructions.
6. Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in the appliance. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
7. Only use the charger provided with your appliance. Do not charge any other device with this charger, as this may cause damage.
8. This appliance contains batteries that are non-replaceable. The appliance has no user serviceable parts. Do not open and disassemble device except when disposing of the battery. Follow handling and disposal instructions of the battery/device described in this manual.
9. When unplugging, always pull the power plug instead of the cord. Do not touch the power plug with wet hands. This can cause an electric shock.
10. Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and

may be fatal if directed into these areas.

11. Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
12. This appliance is a personal care device and is not intended for use on multiple persons.
13. Small parts may come off, keep out of reach of children.

To avoid irrigator nozzle breakage which may generate a small parts choking hazard or damage teeth:

- If the product is dropped, the irrigator nozzle should be replaced before the next use even if no damage is visible. Clean the appliance properly after each use (see "Cleaning recommendations" paragraph).
- Proper cleaning ensures the safe use and functional life of the appliance.

CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

SAVE THESE INSTRUCTIONS

Description

- a Nozzle
- b Switch
- c On/off button (continuous jet)
- d On-demand button
- e Intensity level
- f Intensity selection button
- g Handle
- h Release button
- i Water reservoir cap
- j Water reservoir (150 ml)
- k Charger

Note: Content may vary based on the model purchased.

Specifications:

For voltage specifications refer to the bottom of the charging unit.

Charger Type 3726

Handle Type 3720

Charging and Operating

Plug the charger (k) into an electrical outlet and place the handle (g) on the charger (picture 1). While charging the on/off button (c) light is fading in and out. Once it is fully charged the light turns off. A full charge typically takes 14 hours.

If the battery is running low the on/off button is flashing in red.

First Use

Attach the nozzle (a) to the handle until it locks into place indicated by a click sound (picture 3). Fill the water reservoir

(j) with some fresh and clean water by lifting the reservoir cap (i) on the back of the handle (picture 2).

Hold the handle over a washbasin, press the on/off button (c) and allow the entire water volume to run through the nozzle. When the device has not been used for more than one week repeat this step to ensure flushing.

Jet functions:

On applicable nozzles, you can choose between different jet stream types by sliding the jet switch (b) down or up (before doing so, turn off the water stream with the on/off button (c) on the handle (g)):

Rotating stream type for general cleaning and gum massage (Jet switch in lower position): The rotation of the water stream with its micro-bubbles are designed for easy penetration and gentle gum massage (picture 5).

Focused stream type (Jet switch in upper post on): The focused water stream with its micro-bubbles are designed for targeted cleaning (picture 5).

Press the on/off button (c) for a continuous jet stream.

Press the on-demand button (d) to trigger a single jet stream. Upon releasing the on-demand button the single jet stream stops. Press the intensity selection button (f) successively to adjust the intensity level (sensitive - regular/medium - intense) (depending on model). You may pre-select the intensity level before you turn on the handle. The intensity level selected last is memorized (picture 4). Note: In the beginning, we recommend using the sensitive level for your teeth and gums to get accustomed to the usage.

Using the irrigator

Fill the water reservoir (j) with lukewarm water by lifting the

reservoir cap (i) on the back of the handle (g).

Lean over the sink and direct the nozzle in your mouth. Close your lips slightly and turn on the handle by pressing the on/off button (c).

Guide the jet along the gumline and pause briefly between teeth.

Alternatively, you may also point the nozzle tip at a desired area and use the on-demand button (d) for a single jet stream.

Keep your mouth slightly open to allow the water to run into the sink.

Refill water reservoir as many times as needed to complete your cleaning routine. It may take multiple refills to clean your entire mouth.

After use

Make sure the water reservoir is completely empty. Press the on/off button (c) and leave it turned on until no more water is running through the nozzle, to prevent any deposit build-up or bacteria from collecting in the standing water. Turn off the handle by pressing the on/off button again once the water is empty. A longer operation of the handle without water may damage the motor.

To drain any water out of the handle hold it over the washbasin and tilt it from side to side (picture 6).

Release the nozzle (a) by pressing the release button (h) and pulling the nozzle at the same time (picture 7).

Lift up the cap (i) and let the water reservoir air-dry. Always wipe the handle and nozzle dry (picture 8).

Note: Make sure you always empty the water reservoir before storing or travelling with the irrigator.

In general we recommend to replace the nozzle every 3 months.

Cleaning recommendations

Clean and dry handle after each use and leave the water reservoir cap (i) open.

The water reservoir (j) is dishwasher-safe and can be detached via the snap-fit on the bottom of the handle.

Move the lever and slide the reservoir down (picture 8). Regularly clean the charger and handle with a damp cloth. Store the appliance in a cool and dry place.

Subject to change without notice.

Environmental Notice


Product contains batteries and/or recyclable electric waste.

For environment protection do not dispose in household waste, but for recycling take to electric waste collection points provided in your country.

Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

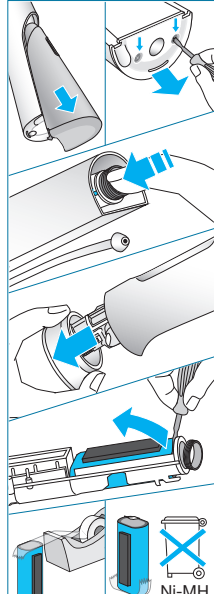Caution! Opening the handle will destroy the appliance and invalidate the warranty.

Satisfaction Guaranteed or Your Money Back:

Try it with no risk for 60 days! If you are not 100% satisfied, call 1-877-593-4199 to receive a full refund. Return the product in its original packaging along with the original cash receipt postmarked within 60 days of purchase to receive a pre-paid card in the amount of purchase. Terms of prepaid card apply.

Your Oral-B irrigator must be returned in its original packaging and be postmarked no later than 60 days

after the date of purchase. Please allow 6-8 weeks for handling. In the U.S., the shipment cost to return the product is the responsibility of the consumer. Sorry: No refunds will be made for any Oral-B product damaged by accident, neglect, unreasonable use or lost in the post.

Terms and conditions apply. For further details please contact our customer service.

FOR USA ONLY

Limited 2-Year Warranty

Oral-B warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail. If the appliance exhibits such a defect, Oral-B will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B. This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B replacement parts and accessories:

Contact an authorized Oral-B Service Center.

Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To order original Oral-B replacement parts and accessories (including replacement nozzles) by phone:

Please call 1-800-566-7252.

To speak to an Oral-B Consumer Service Representative: Please call 1-800-566-7252.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

Oral-B reserves the right, without any prior notice, to make any technical modifications or changes to system implementation of device

features, interface, and menu structure changes that are deemed necessary to ensure that Oral-B systems function reliably.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages. The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B replacement parts & accessories:

Contact an authorized Oral-B Service Center.
Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To speak to an Oral-B Service Representative:

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