PROFESSIONAL

powered by **BRAUN**

PRECISION 5000/5500

TRIACTION 5000/5500

+ SmartGuide

English	6
Español	14
Français	21

Internet:

www.oralb.com
www.oralb.ca
www.braun.com
www.service.braun.com

Braun GmbH
Frankfurter Strasse 145
61476 Kronberg
Germany/Allemagne/Alemania

Charger type 3757

Handle type 3762

Display type 3742

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







Tipo de mango 3762







Tipo de pantalla 3742

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Type de manche 3762

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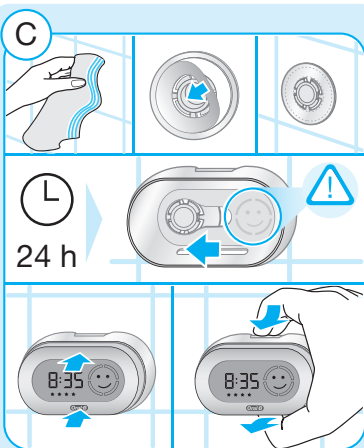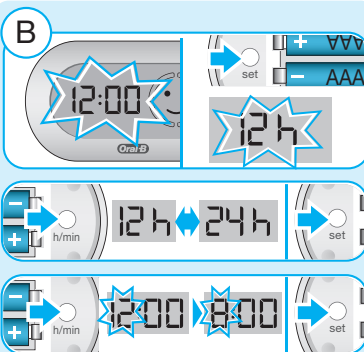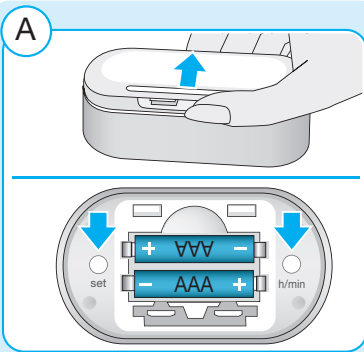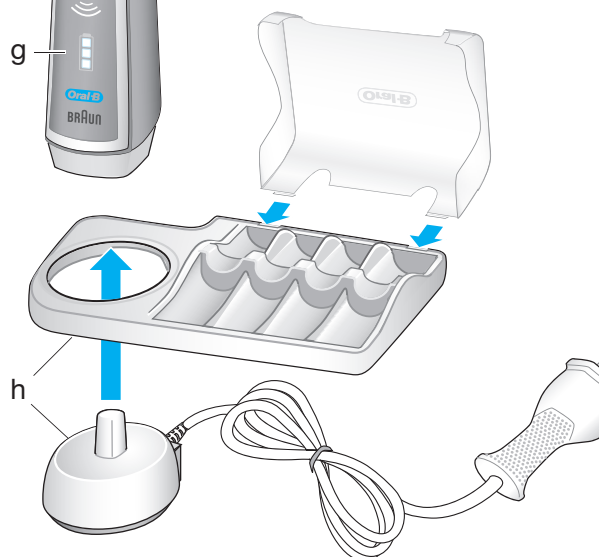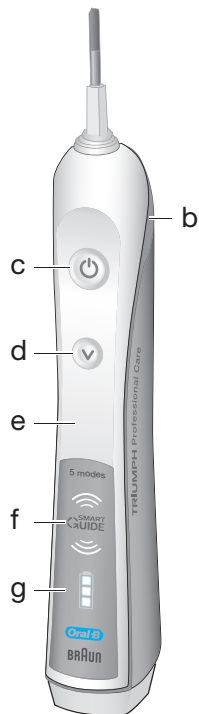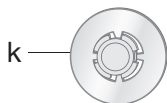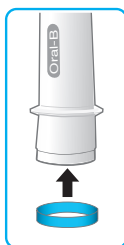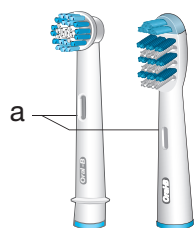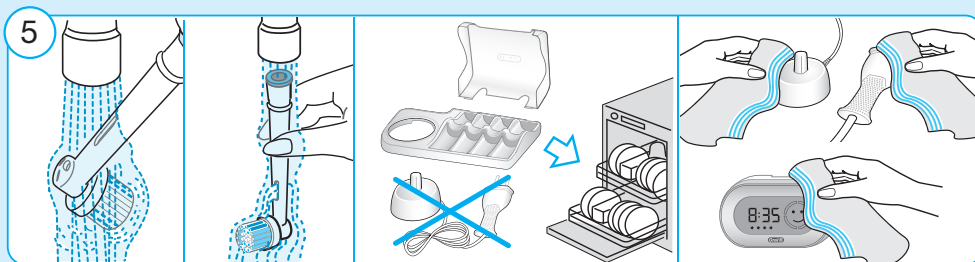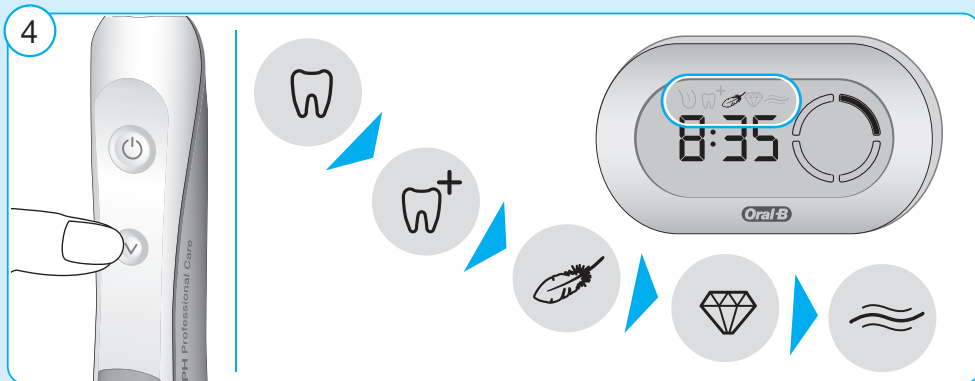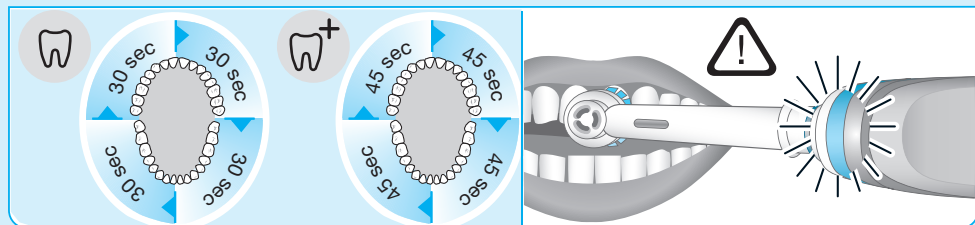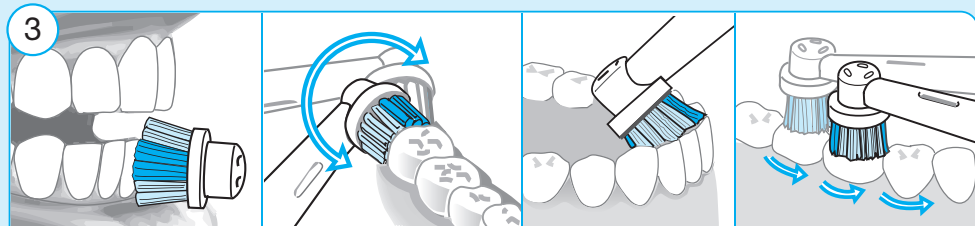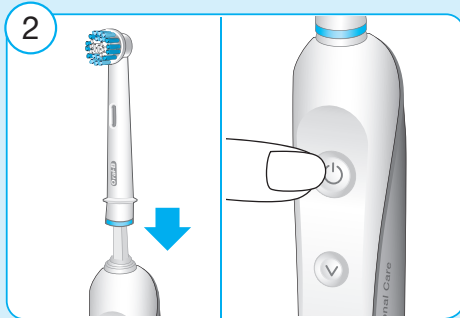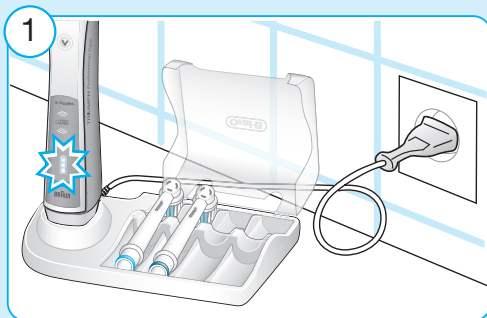IMPORTANT SAFEGUARDS

Oral-B® Professional™ has been carefully designed to offer you and your family a unique brushing experience that is both safe and effective.

Read all instructions before using.

DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by, on, or near children or invalids.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.

IMPORTANT

1. Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Center. A damaged or non-functioning unit should no longer be used.
2. If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
3. Do not touch the power plug with wet hands. This can cause an electric shock.
4. When unplugging, always hold the power plug instead of the cord.
5. Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service Center. Do not disassemble the product except when disposing of the battery.
6. When taking out the battery for disposal, use caution not to short the positive (+) and negative (-) terminals.
7. This product is not intended for use by children under age three.

CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

SAVE THESE INSTRUCTIONS

English

Important Information

- Your Oral-B Triaction/Precision 5000/5500 uses an interactive wireless display. To avoid electromagnetic interference and/or compatibility conflicts, deactivate the radio transmission before using it in restricted surroundings, such as aircrafts or specially marked areas in hospitals. Deactivate radio transmission by simultaneously pressing the on/off and mode button for 3 seconds until the radio transmission display (f) goes off. Follow the same procedure to reactivate transmission again.
- Persons with pacemakers should always keep the toothbrush more than six inches (15 centimeters) from the pacemaker while turned on. Anytime you suspect interference is taking place, deactivate the radio transmission of your toothbrush.

Description

- a Brush head
- b Pressure sensor light
- c On/off button
- d Brushing mode button
- e Handle
- f Radio transmission display
- g Charge indicator
- h Charging unit (basic charger and brush head compartment with protective cover)
- i SmartGuide
- j SmartGuide holder
- k Wall mount
- l Adhesive ring

Connecting and charging

Your toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charger into an electrical outlet. For voltage specifications refer to the bottom of the charging unit.
- The flashing green lights indicate that the toothbrush is charging (picture 1, page 5); once it is fully charged, the lights turn off. A full charge can take up to 24 hours and enables up to ten days of regular brushing (twice a day, 2 minutes).
Note: After deep discharge, the green lights might not flash immediately; it can take up to 10–15 minutes.
- If the rechargeable battery is low the red charge indicator light is flashing for some seconds when turning your toothbrush on/off.
If the rechargeable battery is empty the motor will stop; it will need a 40-minute charge before you can use it for a 2-minute brushing.

- For everyday use, the toothbrush handle can be stored on the charging unit to maintain it at full power. Overcharging is impossible. To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and fully discharge the handle with regular use at least every 6 months.

SmartGuide adjustments

The SmartGuide (i) comes with a demonstration mode. Deactivate this mode by removing the compartment cover in the back (picture A, page 4) and push the «set» or «h/min» button.

A flashing «12:00» indicates that you can now adjust the clock time.

Note: If the SmartGuide display is blank before first use (no active demonstration mode / no batteries inserted) the batteries will need to be replaced/inserted before continuing.

Batteries are supplied.

Clock Time

Press and hold the «set» button until «12h» is displayed. By briefly pressing the «h/min» button, you can toggle between «12h» or «24h» daytime format. Push «set» to confirm your choice. Then the hour digit starts flashing. Press the «h/min» button until the right hour appears and confirm with «set». Proceed same way to set both minute digits (picture B, page 4).

SmartGuide Icons

Displayed Icon	Meaning
	«Daily Clean» mode
	«Deep Clean» mode
	«Sensitive» mode
	«Whitening» mode
	«Massage» mode
	Too much pressure is applied. Reduce the brushing force.
	«Quadrant circle» Move to the next quadrant of your mouth when the relevant segment flashes.
	You reached the dentist recommended 2-minute brushing time (3-minute brushing time within the «Deep Clean» mode).
	After use you will receive a star with each completed quadrant.
	SmartGuide Batteries are low. Change batteries of the SmartGuide.

Using your toothbrush

Brushing technique

Wet brush head and apply any kind of toothpaste. To avoid splashing, guide the brush head to your teeth before switching on the appliance (picture 2, page 5). When brushing your teeth with one of the rotating brush heads, guide the brush head slowly from tooth to tooth, spending a few seconds on each tooth surface (picture 3, page 5). When using the Deep Sweep brush head place the toothbrush bristles against the teeth at a slight angle towards the gumline. Apply light pressure and start brushing in back and forth motions, just like you would do with a manual toothbrush. With any brush head start brushing the outsides, then the insides and finally the chewing surfaces.

Brush all four quadrants of your mouth equally. You may also consult your dentist or dental hygienist about the right technique for you.

In the first days of using any electric toothbrush, your gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist. If you have sensitive teeth and gums, Oral-B recommends using the «Sensitive» mode (optional in combination with the Oral-B Sensitive Gum Care brush head).

Timer

You can choose between «Professional» or «2-Minute» timer and «Count up» or «Count down» timer.

The «2-Minute» timer signals with four beeps and a Smiley on the SmartGuide that the dentist recommended 2-minute brushing time has been reached.

In addition the «Professional» timer signals with short beep sounds at 30-seconds intervals to move to the next quadrant of your mouth (45 seconds intervals within «Deep Clean» mode) (picture 3, page 5). A single beep indicates completion of the first quadrant, two beeps the second one, three beeps the third one and four beeps sound when you reach the recommended 2-minute brushing time (3-minute brushing time within «Deep Clean» mode).

The timer memorizes the elapsed brushing time, even when the handle is briefly switched off during brushing. The timer resets if pausing longer than 30 seconds or briefly pressing the mode button (d) during pause.

The «Count up» timer displays your actual brushing time. The «Count down» timer indicates the remaining brushing time.

Personalizing the timer

Your toothbrush comes with the «Professional» and «Count up» timer activated. To change this, press and hold the brushing mode button for 3 seconds until the SmartGuide shows «:30» and «timer». Briefly press the on/off button to switch to the «2-Minute» timer. Push the mode button to confirm. The SmartGuide now shows «up» and «timer». Press the on/off button to switch to the «Count down» timer and confirm with the mode button. Leave the setup menu by pressing the mode button.

Pressure sensor

If you put excessive pressure on your teeth, the red pressure sensor light (b) will light up (picture 3, page 5). Additionally, when brushing in «Daily Clean»

or «Deep Clean» mode, the brushing mode will change to the «Sensitive» mode. Periodically check the operation of the pressure sensor by pressing moderately on the brush head during use.

Brushing modes

Your toothbrush offers different brushing modes:

- «Daily Clean» – Standard mode for daily mouth cleaning
- «Deep Clean» – Mouth cleaning with extended brushing time of 45 seconds per quadrant (3 minute total brushing time)
- «Sensitive» – Gentle, yet thorough cleaning for sensitive areas
- «Whitening» – Polishing for occasional or everyday use
- «Massage» – Gentle massage of gums


Your toothbrush automatically starts in the «Daily Clean» mode. To switch to the other modes, successively press the mode button (picture 4, page 5). To return from any mode back to the «Daily Clean» mode, press and hold the mode button.

The brushing mode is memorized even when the handle is briefly switched off during brushing. When pausing longer than 30 seconds or briefly pressing the mode button (d) during pause the brushing mode resets to the «Daily Clean» mode.

Brush heads


Oral-B also offers you a variety of different brush heads that fit your Oral-B Professional toothbrush handle and can be used for tooth-by-tooth precision cleaning.

Oral-B Precision Clean brush head

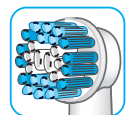
provides effective and gentle cleaning action for both teeth and gums. The brush head cups each tooth with longer and denser bristles and reach deep between teeth at the same time.

Oral-B FlossAction brush head

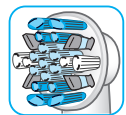features micropulse bristles that allow for an outstanding plaque removal in the interdental areas.

Oral-B ProWhite brush head

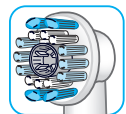features a special polishing cup to naturally whiten your teeth. Please note that children under 12 years old should not use the Oral-B ProWhite brush head.

Oral-B Sensitive Gum Care brush head

provides thorough cleaning for sensitive teeth and gums. It has an extra soft bristle texture that is gentle on teeth and gums.

Oral-B Ortho brush head

with special tufting pattern thoroughly cleans and gently removes plaque from teeth with braces.


Oral-B Power Tip brush head

The Power Tip is designed for cleaning between teeth as well as around bridges, crowns and implants. Position the Power Tip at a 45 degree angle to the gumline. Slowly guide the Power Tip along the gumline, pausing between teeth for excellent results.

Our dynamic Deep Sweep brush head can be used with all common brushing techniques.

Oral-B Deep Sweep brush head

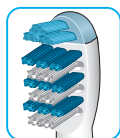features triple zone cleaning action for outstanding plaque removal, even in between teeth.

Oral-B brush heads feature light blue INDICATOR® bristles to help you monitor brush head replacement need. With thorough brushing, twice a day for two minutes, the blue color will fade halfway down approximately in a 3-month period. If the bristles splay before the color recedes, you may exert too much pressure on teeth and gums.

We do not recommend using the Oral-B FlossAction or Oral-B ProWhite brush head with braces. You may use the Oral-B Ortho brush head, specifically designed to clean around braces and wires.

Synchronizing your toothbrush

To avoid interference of display messages when using more than one Oral-B Triaction/Precision 5000/5500 toothbrush with SmartGuide in the same household, you need to assign the handles to their respective SmartGuide.

You can synchronize a maximum of 2 handles with one SmartGuide. To do so, push the «h/min» button

in the battery compartment for at least 3 seconds. The SmartGuide shows «L---». Turn on your handle to start the synchronization process. It is finished when «L-1-» is displayed. Turn off the handle. To assign a second handle to the same SmartGuide turn on the other handle until «L-2-» is displayed. After having synchronized the second handle, «Full» will be displayed. Exit the synchronization mode by pressing the «h/min» or «set» button, otherwise it will end automatically after 30 seconds.

Cleaning recommendations

After brushing, rinse your brush head under running water with the handle switched on. Switch off and remove the brush head. Clean handle and brush head separately; wipe them dry before reassembling. Disassemble the charging unit before cleaning. Brush head compartment and protective cover are dishwasher-safe. Both basic charger and SmartGuide should be cleaned with a damp cloth only (picture 5, page 5).

Wall mount

Use the adhesive ring (l) with the wall mount (k) if you prefer to mount the SmartGuide on a wall, mirror or other surface. Make sure the chosen surface is clean and dry before adhering. Stick the adhesive ring on the even side of the wall mount and press it firmly to the wall. Wait 24h before sliding the SmartGuide holder (j) into the fixed wall mount. Place the SmartGuide into its holder (picture C, page 4). The adhesive ring will not work on dirt-repellent surfaces.

Subject to change without notice.

Environmental Notice


This product contains rechargeable batteries. In the interest of protecting the environment, please do not dispose of the product in the household waste at the end of its useful life. Disposal can take place at an Oral-B Service Center or at appropriate collection points provided in your country.

Removing batteries

Handle battery

At the end of the product's useful life, open the handle as shown, remove the rechargeable battery and dispose of it according to local environmental regulations.

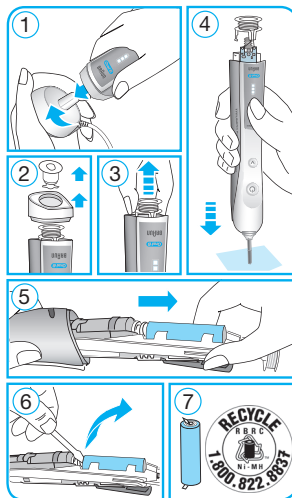Caution! Opening the handle will destroy the appliance and invalidate the warranty.

SmartGuide batteries

Batteries may leak if empty or not used for a long time. Avoid skin contact whilst handling leaking batteries.

Keep away from children.

Do not swallow. If swallowed, seek medical advice at once.

Remove fully discharged batteries immediately.

Replace all batteries at the same time.

Do not mix used and new batteries, different brands or types. Do not disassemble, recharge or dispose of in fire.

Satisfaction Guaranteed or Your Money Back!

Try your Oral-B Professional brush for 60 days, starting from the day of purchase. If you are not 100% satisfied, receive a full refund. No questions asked. For a full refund, please return your product along with your name, address, contact phone number and the original cash register receipt indicating purchase price and date to the address below.

U.S. only:

60-Day Challenge /
Returns Processing Center Suite 2 /
4608 Appliance Drive /
Belcamp, MD 21017

Canada only:

Attention: 60-Day Challenge
Braun MBG
PO Box 12060
Saint John, NB
E2L 5E7

Your Oral-B Professional toothbrush must be returned in its original packaging and be postmarked no later than 60 days after the date of purchase. Please allow 6–8 weeks for handling. The shipment cost to return the product is the responsibility of the consumer.

Sorry: No refunds will be made for any Oral-B Professional product damaged by accident, neglect, unreasonable use or lost in the post.

Terms and conditions apply. For further details please contact our customer service.

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

FOR USA ONLY

Limited 2-Year Warranty

Oral-B warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase at retail.

If the appliance exhibits such a defect, Oral-B will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B. This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B replacement parts and accessories:

Contact an authorized Oral-B Service Center.

Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To order original Oral-B replacement parts and accessories (including replacement brush heads) by phone: Please call 1-800-566-7252.

To speak to an Oral-B Consumer Service Representative: Please call 1-800-566-7252.

Declaration of conformity

This appliance complies with part 15 of the FCC rules. Operation is subject to the following two conditions: (1) this appliance may not cause harmful interference, and (2) this appliance must accept any interference received, including interference that may cause undesired operation.

Display unit: Tested to comply with FCC standards. FOR HOME OR OFFICE USE.

To change or modify any part of this appliance could void the authority to operate the equipment.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages. The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice. Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B replacement parts & accessories:

Contact an authorized Oral-B Service Center.

Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To speak to an Oral-B Service Representative:

Please call 1-800-566-7252.

Declaration of conformity

Operation is subject to the following two conditions: (1) this appliance may not cause harmful interference, and (2) this appliance must accept any interference received, including interference that may cause undesired operation.

This Class B digital apparatus complies with Canadian ICES-003.

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