





Kids' Power Toothbrush

Charger Type 3757, Handle Type 4733

powered by **BRAUN**

97544926

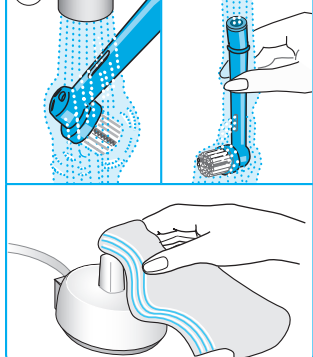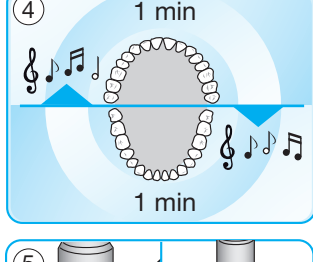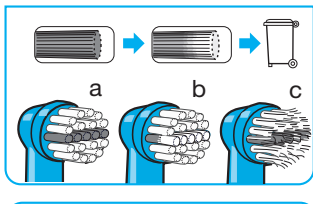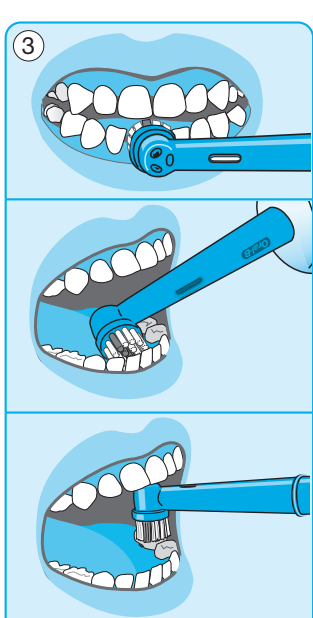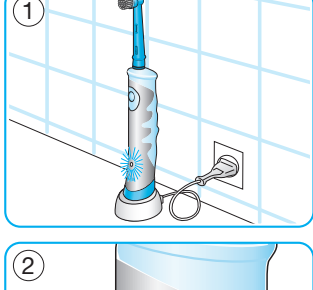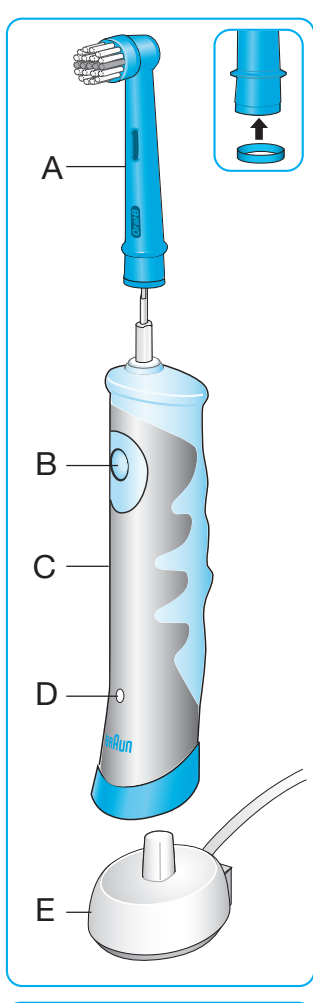IMPORTANT SAFEGUARDS

Dear parents,

The Braun Oral-B Kids' PowerToothbrush has been developed with leading dental professionals to effectively remove plaque from children's teeth. It will make it easier for your child to brush properly and thoroughly. We hope your child will enjoy his/her new Braun Oral-B appliance.

This appliance is suitable for children 3 years or older. Please read the use instructions carefully and completely before showing your child how to use the appliance. We recommend that parents supervise their children when they start using the appliance in order to help them become familiar with it.

Read all instructions before using.

DANGER

In order to reduce the risk of electrocution:

1. Do not place or store the charger where it can fall or be pulled into a tub or sink.
2. Do not place the charger in water or other liquid.
3. Do not reach for a charger that has fallen into water. Unplug immediately.
4. Do not use while bathing.

WARNING

To reduce the risk of burns, electrocution, fire or injury:

1. Close supervision is necessary when this product is used by, on, or near children or handicapped.
2. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
3. Keep cord away from heated surfaces.
4. Never drop the appliance or insert any object into any opening of the appliance.
5. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged or if it has been dropped into water. If any of these occur, please contact a service center and return the product for examination and repair.

IMPORTANT

1. Periodically check the cord for damage. If cord is damaged, take the charging unit to an Oral-B Service Center. A damaged or non-functioning unit should no longer be used.
2. If the toothbrush handle is dropped, the brush head should be replaced before the next use even if no damage is visible.
3. Do not touch the power plug with wet hands. This can cause an electric shock.
4. When unplugging, always hold the power plug instead of the cord.
5. Do not modify or repair the product. This may cause fire, electric shock or injury. Consult your dealer for repairs or contact an Oral-B Service Center. Do not disassemble the product except when disposing of the battery.
6. When taking out the battery for disposal, use caution not to short the positive (+) and negative (–) terminals.
7. This product is not intended for use by children under age three.
8. This toothbrush is a personal care device and is not intended for use on multiple patients in a dental practice or institution.
9. For tooth brushing only. Do not use as a toy or without adult supervision.
10. Do not chew or bite on brush.
11. Replace brush (head) every 3 months or sooner if brush head becomes worn.

CAUTION

If you are undergoing treatment for any oral care condition, consult your dental professional prior to use.

SAVE THESE INSTRUCTIONS

Description

- A Brush head
- B On/off switch
- C Squish Grip handle
- D Charge indicator
- E Charging unit

Specifications

For voltage specifications please refer to the bottom of the charging unit.

Connecting and charging

This toothbrush has a waterproof handle, is electrically safe and designed for use in the bathroom.

- Plug the charging unit (E) into an electrical outlet and place the toothbrush handle (C) on the charging unit.
- The charge indicator light (D) flashes while the handle is being charged, once it is fully charged the light turns off (1). A full charge can take up to 22 hours and enables up to seven days of regular brushing (twice a day, 2 minutes).
Note: If the charging light does not flash immediately, continue charging. It will appear in 10–15 minutes.
- For everyday use, the handle can be stored on the plugged-in charging unit to maintain the appliance at full power. Overcharging is impossible.

Battery maintenance

To maintain the maximum capacity of the rechargeable battery, unplug the charging unit and discharge the handle by regular use at least every 6 months.

Brushing with a music timer

Use the following steps to show and explain to your child how to use the appliance properly:

- When the handle is taken from the plugged-in charging unit, a short melody will sound.
- Put a small amount of toothpaste onto the brush head (we recommend kids toothpaste).
- First guide the brush head to the teeth, only then turn on the handle by pressing the on/off switch (B) (2).
- Start with the lower teeth: Slowly guide the brush head from tooth to tooth;
 - first the outsides,
 - then the insides,
 - finally the chewing surfaces (3).Do not scrub or press too hard. The toothbrush automatically provides the proper cleaning action.
- After one minute of operation, the brush movement will be interrupted by a short melody. This tells your child to move on to the upper teeth. Clean them – again tooth by tooth – until the next melody plays (after another minute of operation). This way you can make sure that your child brushes his/her teeth for at least 2 minutes as recommended by dental professionals (4).
- Turn off the handle by pressing the on/off switch (B) while the brush head is still in the mouth.

We recommend to replace the brush head every three months. The brush head is provided with a row of blue INDICATOR® bristles (3a). With proper brushing twice a day for two minutes, the color will disappear half-way down these bristles in a 3-month period (3b). If bristles begin to splay before the color starts to recede, this indicates that your child exerts too much pressure on teeth/gums (3c).

Replacement brush heads: Stages Power (EB10)

Available at your retailer or at Oral-B Braun Service Centers.

In the first few days of using any electric toothbrush, gums may bleed slightly. In general, bleeding should stop after a few days. Should it persist after 2 weeks, please consult your dentist or dental hygienist.

Cleaning


After use, rinse the brush head for several seconds under running water with the handle switched on. Then switch off the handle and take off the brush head. Rinse both parts separately under running water and wipe them dry before putting the handle back onto the charging unit. From time to time, also clean the charging unit with a damp cloth (5).

Environmental Notice

PRODUCT CONTAINS NICKEL-METAL HYDRIDE RECHARGEABLE BATTERIES. MUST BE RECYCLED OR DISPOSED OF PROPERLY.

Battery removal at the end of the product's useful life

Open the handle as shown, remove the battery and dispose of it according to local environmental regulations.

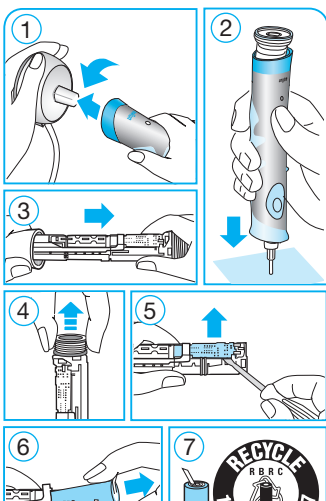Caution! Opening the handle will destroy the appliance and invalidate the warranty.

Subject to change without notice.

FOR USA ONLY

Limited 2-Year Warranty

Oral-B warrants this appliance to be free of defects in material and workmanship for a period of two years from the date of original purchase.

If the appliance exhibits such a defect, Oral-B will, at its option, repair or replace it without cost for parts and labor. The product must be carried in or shipped, prepaid and insured (recommended), to an authorized Oral-B Service Center. Proof of date of original purchase is required.

The cost to ship the product to an authorized Service Center is the responsibility of the consumer, the cost to return the product to the consumer is the responsibility of Oral-B.

This warranty does not cover finishes or normal wear, nor does it cover damage resulting from accident, misuse, dirt, water, tampering, unreasonable use, service performed or attempted by unauthorized service agencies, failure to provide reasonable and necessary maintenance, including recharging, or units that have been modified or used for commercial purposes.

ALL IMPLIED WARRANTIES, INCLUDING ANY IMPLIED WARRANTY OF MERCHANTABILITY OR FITNESS FOR ANY PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO TWO YEARS FROM DATE OF ORIGINAL PURCHASE. IN NO EVENT WILL ORAL-B BE RESPONSIBLE FOR CONSEQUENTIAL DAMAGES RESULTING FROM THE USE OF THIS PRODUCT.

Some states do not allow the exclusion or limitation of incidental or consequential damages, so the above limitation may not apply to you. This warranty gives you specific legal rights, and you may have other legal rights which vary from state to state.

How to obtain service in the USA

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

A. Carry the product in to the authorized Oral-B Service Center of your choice or,

B. Ship the product to the authorized Oral-B Service Center of your choice.

Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B replacement parts and accessories:

Contact an authorized Oral-B Service Center.

Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To order original Oral-B replacement parts and accessories (including replacement brush heads) by phone:

Please call 1-800-566-7252.

To speak to an Oral-B Consumer Service Representative:

Please call 1-800-566-7252.

FOR CANADA ONLY

Limited 2-Year Warranty

In the event this appliance fails to function within the specified warranty period because of defects in material or workmanship, and the consumer returns the unit to an authorized service center, Braun Canada will, at its option either repair or replace the unit without additional charge to the consumer. This warranty does not cover any product which has been damaged by dropping, tampering, servicing performed or attempted by unauthorized service agencies, misuse or abuse. This warranty excludes units which have been modified or used for commercial purposes.

To the extent allowed by law, Braun Canada shall not be responsible for loss of the product, loss of time, inconvenience, commercial loss, special or consequential damages. The provisions of this warranty are in addition to and not a modification of or subtraction from the statutory warranties and other rights and remedies contained in applicable provincial legislation.

How to obtain service in Canada

Should this product require service, please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To obtain service:

- A. Carry the product in to the authorized Oral-B Service Center of your choice or,
- B. Ship the product to the authorized Oral-B Service Center of your choice.
Pack the product well. Ship the product prepaid and insured (recommended). Include a copy of your proof of purchase to verify warranty coverage, if applicable.

To obtain original Oral-B replacement parts & accessories:

Contact an authorized Oral-B Service Center.

Please call 1-800-566-7252 to be referred to the authorized Oral-B Service Center closest to you.

To speak to an Oral-B Service Representative:

Please call 1-800-566-7252.

FOR USA AND CANADA

Warranty replacement brush heads

Oral-B warranty will be voided if electric rechargeable handle damage is found to be attributed to the use of non-Oral-B replacement brush heads.

Oral-B does not recommend the use of non Oral-B replacement brush heads.

- Oral-B has no control over the quality of non-Oral-B replacement brush heads. Therefore, we cannot ensure the cleaning performance of non-Oral-B replacement brush heads, as communicated with the electric rechargeable handle at time of initial purchase.
- Oral-B cannot ensure a good fit of non-Oral-B replacement brush heads.
- Oral-B cannot predict the long-term effect of non-Oral-B replacement brush heads on handle wear.

All Oral-B replacement brush heads carry the Oral-B logo and meet Oral-B's high quality standards. Oral-B does not sell replacement brush heads or handle parts under any other brand name.

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Braun GmbH