


PMD Beauty Clean Acne

PMD Clean Acne is a smart facial cleansing device using SonicGlow™ Technology to ensure a deep, daily cleanse. Clean Acne removes skin's impurities by breaking down dirt and oil from within the pores while firming and smoothing the skin. Furthermore, Clean Acne utilizes blue light therapy to kill acne-causing bacteria leaving you with clear, glowing skin. Suitable for all skin types, the gentle silicone brush is hygienic, antibacterial, waterproof and - best of all - never needs replacing.



Benefits

- ✓ Removes makeup and impurities from the skin
- ✓ Treats mild to moderate acne
- ✓ Minimizes the appearance of blackheads
- ✓ Reduces oil congestion
- ✓ Lifts, firms and tones the skin

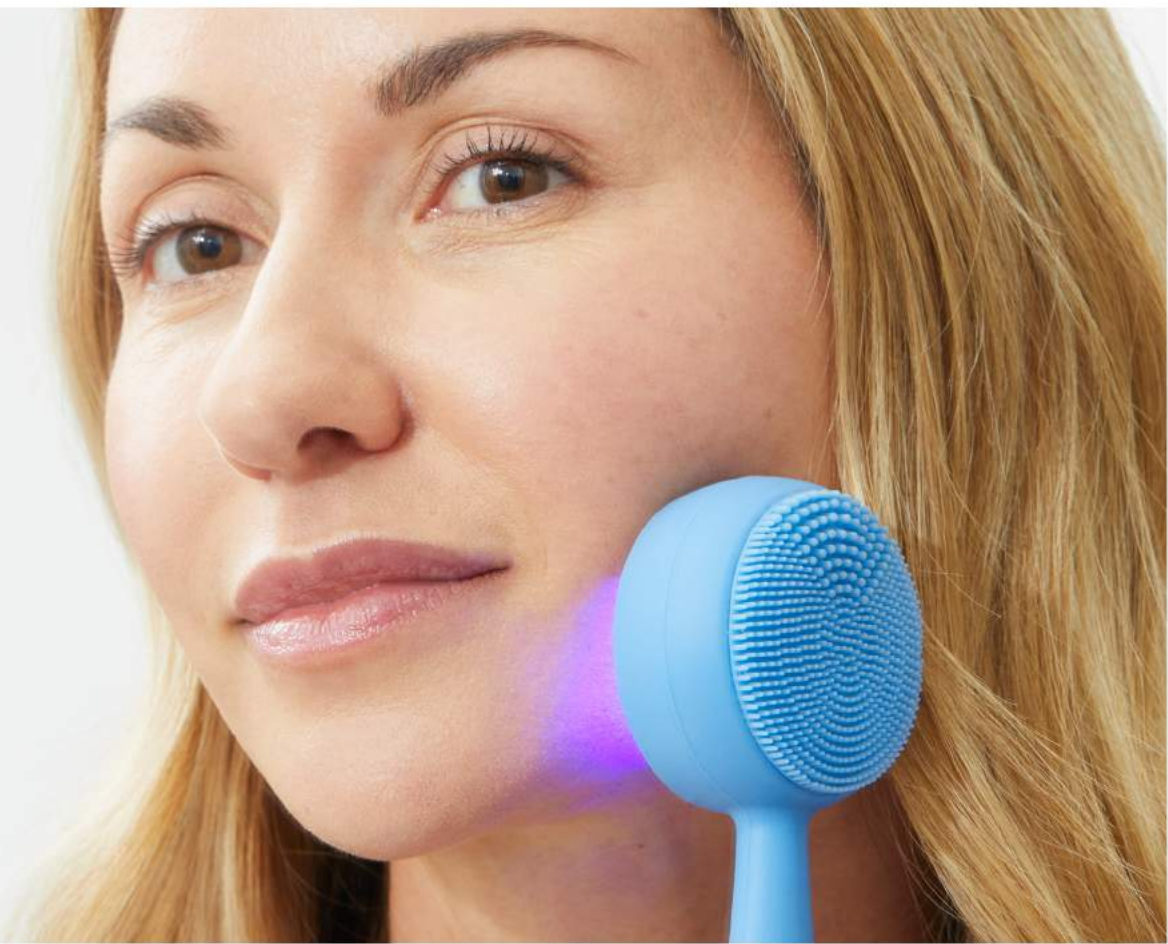

Technical Information

SonicGlow™ technology:
over 7,000 vibrations
per minute

Rechargeable,
USB charging
cord included

Waterproof,
can be used
in the shower

2 modes for cleansing

3-minute Blue Light
Therapy Spot Treatment:
powerful medical grade
LEDs used specifically for
the treatment of mild to
moderate acne

Includes
goggles

2-year
Warranty

Includes a
storage case

How to Use



1: Rinse face with warm water and gently rub in cleanser.



2: Rinse PMD Clean Acne under water and turn it on.



3: Pair with mode 1 or 2 for best cleansing results.



4: Move brush head in circular motions, cleansing the entire face. It is safe to spend more time in areas that may be of personal concern.



5: Rinse face and PMD Clean Acne of any excess cleanser and pat dry.



6: Turn on blue light therapy and apply to areas with mild to moderate acne. The PMD Clean Acne will vibrate to signify the beginning of the 3 minute spot treatment. Repeat as needed.

Use only as directed. Not following instructions may lead to rash, redness, itching or even injury. If this occurs, discontinue use and consult a medical professional. Wear goggles and do not look directly into the LED while using. Product is intended for external use only. Avoid sensitive areas, including any area with broken, irritated, or itching skin. Keep product out of reach of children. PMD is not responsible for malfunction or injury if the product is not used as intended. Not intended for clinical or commercial use. If you feel any discomfort or skin redness after use that persists for more than two days, please stop using this device. If you are sensitive to light or have other conditions after use, please contact customer service. The effects of the PMD Clean Acne are not always immediately visible and results may vary.