
RENPHO

USER MANUAL

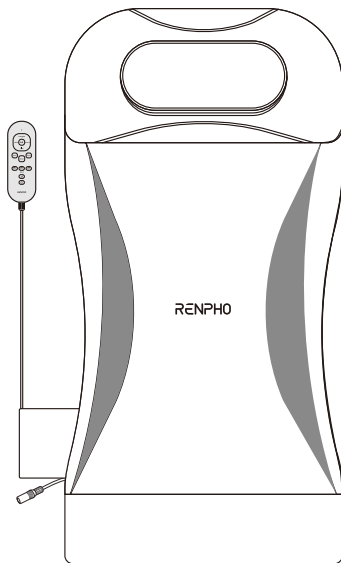MASSAGE CUSHION

Model: R-E002

CONTENTS

• ENGLISH	01-11
-----------	-------

• DEUTSCH	12-22
-----------	-------

• FRANÇAIS	23-33
------------	-------

• ITALIANO	34-44
------------	-------

• ESPAÑOL	45-55
-----------	-------

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

1. Consult a doctor before using the device, if you fall into any of the following user groups:
 - Pregnant or lying-in women
 - Children
 - Elderly (over 65 years old)
 - Persons using a pacemaker or other implantable medical device
 - Persons with heart disease, malignant tumors, or other serious health conditions
 - Persons receiving medical treatment or feel discomfort during use
 - Persons suffering from diabetes, or thrombosis, varicose veins

- Persons with high or abnormal blood pressure
 - Persons with skin infections, sensitivity, redness, wound, and/or burns on the skin
2. DO NOT use this device for more than 30 minutes at a time. Extensive use could lead to the product's excessive heating and shorter life. Should this occur, discontinue use and allow the unit to cool before operating.
 3. Never use this device while driving, sleeping, or ingesting substances that reduce reaction time (e.g. painkiller medications, alcohol).
 4. Never use this device when your skin is wet after a shower, sweating after or during exercise, slippery after applying lotion or oil, etc.
 5. NEVER use this device on an uneven or unstable floor or places that is smaller than the product size.
 6. This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
 7. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
 8. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
 9. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 10. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
 11. Do not carry this appliance by supply cord or use cord as a handle.
 12. Keep the cord away from heated surfaces.

13. Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
14. Never drop or insert any object into any opening.
15. Do not use outdoors.
16. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
17. To disconnect, turn all controls to the off position, then remove plug from outlet.
18. The appliance is only to be used with the power supply unit provided with the appliance.
19. The appliance must only be supplied at safety extra low voltage corresponding to the marking on the appliance.


SAVE THESE INSTRUCTIONS

WHAT'S IN THE BOX


Massage Cushion x1

Power Adapter x1

ABOUT THE PRODUCT

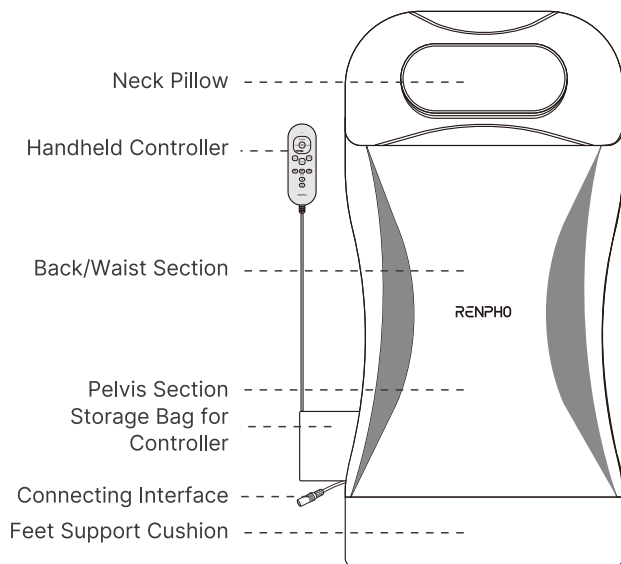
Handheld Controller

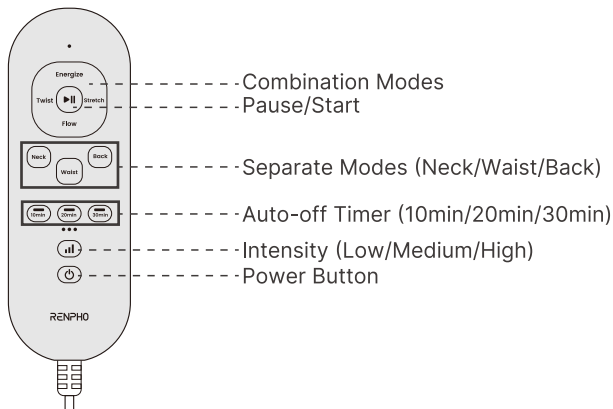
USING THE PRODUCT

1. Preparation


- Unzip the cushion and lay it on a flat and solid floor with the logo facing upwards.

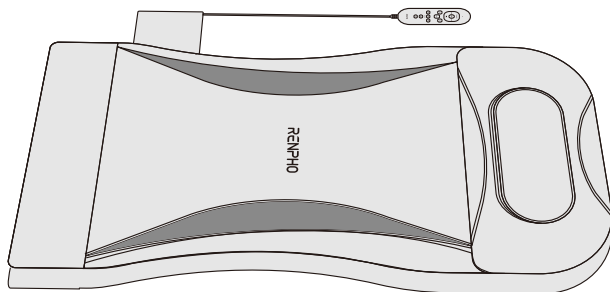- b. Plug the adapter into the appliance's Connecting Interface and into an outlet.

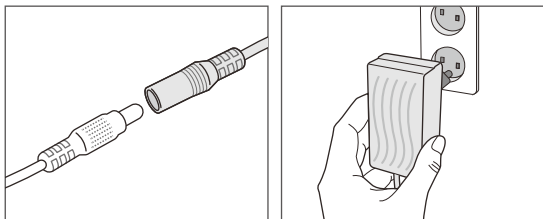- c. Lay down with your shoulders 5cm away from the neck pillow.

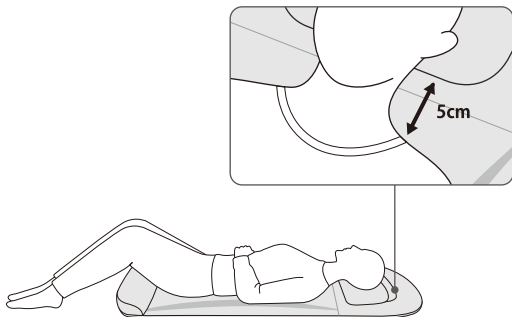2. Operation


Handheld Controller	Operation
	Press to turn the appliance on.
Energize	A dynamic combination of shoulder, back, and waist massage, followed by targeted neck massage. (Default Mode)
Flow	A wave-like motion from hips to neck, transitioning into a combined shoulder, back, and waist massage.
Twist	Twist your whole body diagonally from side to side.
Stretch	Gradually stretches your waist, back, and shoulders, then eases into a full-body stretch with neck relief.
Neck	Neck and shoulders combined stretch.
Waist	Back and waist combined stretch.
Back	Waist and hip combined stretch.
	Press to Pause/Start the process. The pause time will not be counted into the timer. The device will automatically turn off after 10 minutes of pause.
10min 20min 30min	Press one of three buttons to choose timer among 10 min (default) / 20 min / 30 min.
	Press to cycle through 3 intensity level: Low / Medium (default) / High.

Note: When a function/mode is selected, corresponding indicator will light up.

CLEANING & MAINTENANCE

1. Storage

- a. Roll up gently from both ends to the middle.
- b. Zip up slowly.

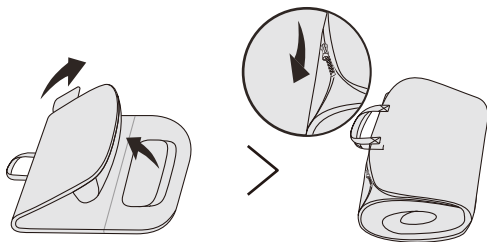Note:

- Always unplug the device after use.
- Place the device in a safe, cool, and dry place when not in use.
- DO NOT store under direct sunlight or in a wet environment.

2. Cleaning

- Always unplug the device prior to cleaning.
- After use, wipe the cushion with a soft, damp cloth using non-abrasive cleaners. If necessary, clean with a sterilizer.
- DO NOT submerge the device in water or any other liquid. If liquid is spilled on the device, wipe it lightly with dry cloth.
- Keep away from all solvents and harsh detergents.

3. Maintenance

- DO NOT attempt to repair this device. There are no user-serviceable parts.
- Abide by local discard and recycle regulations when discarding the device.

TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
Device doesn't turn on.	Disconnected power or loose adapter connection causing poor contact.	Reconnect the adapter and re-plug to the power source.
	Power is off.	Press the Power button to turn on the device.
Device shuts off automatically.	Auto-off timer is activated.	Press the Power button to restart the device if needed.
	The device is overloaded, activating overload protection.	Restart the device after alleviating the load.
	The motor auto shuts-off due to long time use, activating overheating protection.	Use the device after the motor cools off (approximately 30 minutes).

TECHNICAL SPECIFICATIONS

Product Name: Massage Cushion

Model: R-E002

Input: 12V  2A

Rated Power: 24W

Auto-off Timer: 10/20/30min

Dimensions: 47.24 × 23.62in

Product Weight: 7.94lbs

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 Tel: **+1(844) 417 0149 (US&CA)**
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)
+44(0)7434-666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM(GMT)

 Email: support@renpho.com (US&CA)
support-au@renpho.com (AU)
support-eu@renpho.com(UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2014/35/EU, Directive 2011/65/EU, The Electrical Equipment (Safety) Regulations 2016, The Electromagnetic Compatibility Regulations 2016, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:
<https://renpho.eu/pages/compliance-data>

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FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1: This device may not cause interference.
- 2: This device must accept any interference, including interference that may cause undesired operation of the device.

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