


ARTEMIS

LED Light Therapy Mask

USER MANUAL

Model: E100B

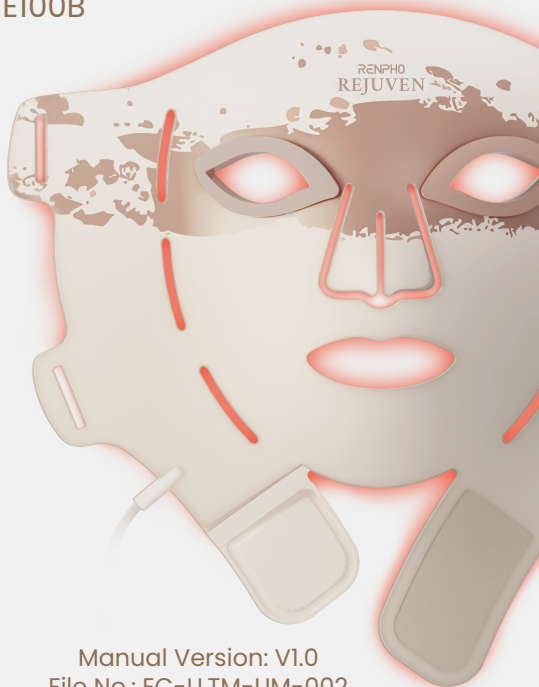

Manual Version: V1.0
File No.: EC-LLTM-UM-002

RENPHO
REJUVEN

CONTENTS

1. RECOMMENDATION OF USE	1
1.1 Warnings.....	1
1.2 Contraindications (Do not use if any of the following occur).....	3
1.3 Adverse Reactions	4
1.4 Precautions.....	4
2. INTENDED USE & INDICATIONS FOR USE	5
3. PRODUCT DESCRIPTION	6
3.1 Specification	7
3.2 Product List	9
3.3 Description of LCD Screen of Controller	10
3.4 Description of Controller Buttons	11
4. USING INSTRUCTIONS	12
4.1 Check for These Situations Before the Use	12
4.2 Light Sensitivity Test	13
4.3 First Use and Preparation	13
4.4 Operation Instructions	14
5. CLEANING, DISPOSAL AND TROUBLE SHOOTING	18
5.1 Cleaning	18
5.2 Storage	19
5.3 Trouble Shooting	19
5.4 Disposal of Electrical and Electronic Equipment Waste (E-waste)	21
6. SERVICE	22
7. ESSENTIAL PERFORMANCE	23
8. SAFETY AND EMC	23
9. LABEL SYMBOLS	29
10. MANUFACTURER INFORMATION	31

1. RECOMMENDATION OF USE

Please read the complete instruction manual before using the LED light therapy mask. It will give you a better understanding of how the product works to enjoy the best treatment results.

1.1 Warnings

- Do not leave the device unattended when turned on to avoid the risk of fire or burns.
- This device is not designed for use by persons with limited physical, sensory, or mental abilities, or lack of experience and knowledge, unless they are trained, supervised, or used together to avoid fire hazards or burns.
- This device is not intended for use by children. Children should be supervised to ensure that they do not play with the device to avoid the risk of fire and burns.
- Do NOT wash electrical parts with water or other liquids.
- Water and electricity are a dangerous combination, and in order to avoid the risk of electric shock:
 - Do not use the device in a humid environment (e.g., in a bathroom or shower or pool attachment)
 - Do not allow water to enter the machine
- If the main cord is damaged, it may be replaced by Shenzhen Eyco Technology Co., Ltd, a service centre authorized by Shenzhen

Eyco Technology Co., Ltd or similarly qualified persons in order to avoid risk of electrocution.

- Do not look directly at the light-emitting diodes when using this device.
- Do not use the device if it appears to be operating incorrectly, and contact customer support.
- Do not use the device in combination with muscle heating creams or other medications which may result in damage to skin tissue.
- Do not expose any parts of the device to chemical solvents or harsh cleaning fluids.
- The device contains silicone and plastics that come into contact with the skin during use. Do not use the device if you have a known allergic reaction to contact with silicone or plastics.
- Avoid fluids on the light panel and controller. Discontinue use if there is water in the device.
- Opening or altering the device may damage the operating circuits and will result in voiding the warranty.
- Please keep the unit out of reach of children, as it may cause an entanglement or strangulation risk.
- Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.
- Use of controls or adjustments or performance of procedures other than those specified herein may result in HAZARDOUS radiation exposure.

- The controller of this device contains lithium batteries, and personnel without adequate training should not replace the batteries without authorization to avoid related risks.
- This device contains a Lithium-ion Battery. Before disposal of this device, please review the local laws and requirements related to Lithium-ion Battery disposal at your area. The preferred method of disposal is recycling.

1.2 Contraindications (Do not use if any of the following occur)

- Do not use the device if you are under the age of 22.
- Do not use this device if you have any photosensitive disorders or genetic conditions of the eye.
- Do not use the device if you are pregnant, breastfeeding, or preparing to get pregnant.
- Do not use the device on infants or children.
- Do not use the device if you have a medical history of seizures that were triggered by light.
- Do not use the device if you are taking photosensitive drugs.
- Do not use the device over malignant and nonmalignant tumors.
- Do not use the device over open or bleeding wounds.
- Do not use the device if you have or may have heart disease.
- Do not use the device over area with phlebitis,

thrombophlebitis and varicose veins.

1.3 Adverse Reactions

- This product has no side effects, but if you are sensitive to light, please stop using immediately.
- If you experience any discomfort, skin discoloration, or any other concerns develop, stop using immediately and consult your physician.
- The material of this mask is silicone, and the strap is neoprene. If you are allergic to any materials found in this device, please discontinue use immediately.

1.4 Precautions

- Discontinue use if the device overheats or becomes uncomfortably hot.
- Use only the supplied USB cable. The use of an incorrect power cord could damage the device.
- Contact your doctor if you have any medical questions regarding the use of this device.
- Clean the device regularly to avoid cross-infection.
- It will void your warranty if you do not operate the device per the manual's instructions.
- Store the device according to the protocol worked out in this manual.
- The normal function of the product may be disturbed by solid electromagnetic

interference. If so, reset the product to resume normal operation by following the instruction manual. Please use the product in another location if the function cannot resume.

- The device should only be used with manufacturer-approved accessories.
- Opening or altering the device may damage the operation circuits and will result in voiding the warranty.
- If your device requires repair or service, please contact customer support.
- Use a specification (input: 100-240Va.c., 50Hz/60Hz, output: DC 5V, 1A) adapter to charge the device before it can be charged normally.
- Do not subject the controller and package to shocks, such as dropping it onto the floor.

2. INTENDED USE & INDICATIONS FOR USE

- Mode 1 (Red + Infrared light): treatment of full face wrinkles.
- Mode 2 (Blue light): Treatment of mild to moderate inflammatory acne.
- Mode 3 (Mixed light): Treatment of mild to moderate inflammatory acne.

3. INTENDED USE & INDICATIONS FOR USE

LED phototherapy masks are over-the-counter light-emitting diode (LED) devices that emit energy and are used in dermatology to treat acne and wrinkles. The device emits red (630nm) and infrared (850nm) light to treat wrinkles and blue (415nm) light to treat mild to moderate acne. The device has to be fitted with eye cups and strapped to the treatment area.

1. The mode of power supply powers the device and controls its operation through the lithium battery inside the controller. The specifications for the adapter that can be used are as follows: input AC 100–240V, 50/60Hz; output DC 5V, 1A, or 2A to charge the controller.

2. The corresponding equipment of E100B, is soft case, and the material is silicone.

3. The controller has three buttons to adjust the brightness level, switch modes, set time, and display. The display shows the battery level, treatment mode, time, and brightness information of the current controller. The controller is powered by an internal rechargeable lithium battery.


3.1 Specification

Specification	
Product name	LED light therapy mask
Model	E100B
Power source	DC 5V/1A. Rechargeable Li-Ion batteries
Weight	Mask: about 210.6g Controller: about 102g
Dimensions	Mask: 296*214*4.5mm Controller: about 118*55*28.2mm
Location for Use	Face
Environment of Use	OTC
Wavelengths	Red: 630nm Blue: 415nm NIR: 850nm
Energy Source	LEDs
Mode	3
Irradiance	R+NIR: 6~8 mW/cm ² Blue light: 2.0~4.0 mW/cm ² Mixed light: 9.0~12.0 mW/cm ²
Irradiance corresponding to each brightness level	R+NIR: Level 1: 6~6.4 mw/cm ² Level 2: 6.4~6.8 mw/cm ² Level 3: 6.8~7.2 mw/cm ²

	Level 4: 7.2–7.6 mw/cm² Level 5: 7.6–8.0 mw/cm² Blue : Level 1: 2.0–2.4 mw/cm² Level 2: 2.4–2.8 mw/cm² Level 3: 2.8–3.2 mw/cm² Level 4: 3.2–3.6 mw/cm² Level 5: 3.6–4.0 mw/cm² R+B+NIR: Level 1: 9.0–9.6 mw/cm² Level 2: 9.6–10.2 mw/cm² Level 3: 10.2–10.8 mw/cm² Level 4: 10.8–11.4 mw/cm² Level 5: 11.4–12.0mw/cm²
Recommended treatment time	10 minutes each time. 3–4 treatment a week, reduce to 1–2 treatment a week once the results shown.
Number of LEDs	108
Service life	About 2 years
Safety and EMC	IEC 60601-1, IEC 60601-1-2, IEC 60601-1-11, IEC 60601-2-57, IEC 62471, IEC 62133-2
Biocompatibility	ISO 10993-5, ISO 10993-10, ISO 10993-23
Operating Environment	Temperature: 0–35℃ Humidity: 30–95%RH Atmospheric pressure. 700hPa–1060hPa

Storage Environment	Temperature: -20-50℃ Humidity: 10-80%RH Atmospheric pressure: 700hPa-1060hPa
---------------------	---

3.2 Product List

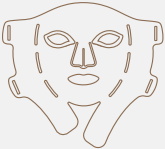Component	Picture	Quantity
E100B Mask		1pc (According to the model you buy)
Controller		1pc
USB cable		1pc
Strap		2pc

Component	Picture	Quantity
Eye cups		2pc
User manual	/	1pc


3.3 Description of LCD Screen of Controller.

LCD screen to display treatment duration, treatment mode, brightness level and power status. Please refer to table below for details:



- Battery symbol: indicates battery level
- LEVEL: indicates brightness level
- Color: Adjusts the mode. R NIR is displayed as mode 1, B is displayed as mode 2, and R NIR B is displayed as mode 3
- Time: displays the time left for the current session

3.4 Description of Controller Buttons

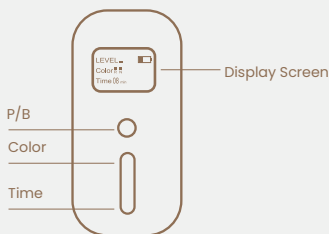

Explanation

- **P/B:** long press the button to turn on/off the device, default red + infrared mode, default 10 minutes; Short press the button to adjust brightness level and there are 5 levels. The five levels range from 1 to 5, with level 1 being the dimmest and level 5 being the brightest.
- **Color:** press the button to select treatment mode from R+NIR, B, R+NIR+B.
- **Time:** press the button to switch between 10, 20, 30-minute sessions.
- **Display Screen:** display battery level.

4. USING INSTRUCTIONS

4.1 Check for These Situations Before Use

Please refer to the table below before use, which will provide more effective instructions for your choice.

Level	Clinical features	Pictures
Mild Acne	One or several pimples and redness in small areas	
Moderate Acne	More pimples and more redness over a wider area	
Severe Acne	Widespread and redness with many pimples and pustules. May have hard bumps under the skin	
Note: The subject device is not suitable for severe acne. Please do not use the device if you have severe acne.		

4.2 Light Sensitivity Test

When first using the device, the brightness should be set to level 1 to test for photosensitivity. Carefully follow these instructions:

- 1) IR emitted from this device. Do not stare at the light source, since it may harm your eyes. Caution is advised.
- 2) If your skin feels hot, uncomfortable, or itchy during the treatment stop immediately.
- 3) Wait until the treatment is fully completed. It will take 10 minutes to complete a full treatment.
- 4) The mask will automatically turn off after the treatment is completed. Remove the mask and wait for at least 10 hours.

Note: Check your skin regularly over the ten-hour period. If your skin appears red and/or itchy, the mask is unsuitable for you. If you experience photophobia, please contact your medical advisor or doctor.

4.3 First Use and Preparation

- 1) Unpack the package and check for all the parts described in section 3.2.
- 2) Before use it is important to fully charge the controller.
 - Connect the adapter to the USB inlet of the controller.
 - Plug the adapter into the power outlet (Please note, Use certified adapters, adapter

input: 100-240Va.c., 50Hz/60Hz, output: DC 5V, 1A) to charge the controller.

- During charging, the LCD screen display the battery capacity.
- It takes about 3~4 hours to fully charge the device.

Note:

1. The recommended AC adapter should be certified by IEC/EN 60601-1 or IEC/EN 62368-1.
2. There are two USB ports on the controller which are only used to connect the Mask and power it. The interfaces are for power only and do not transmit data.
3. The device cannot be used while the controller is being charged.

4.4 Operation Instructions

4.4.1 Clean your skin

Before starting the treatment, thoroughly cleanse your skin, remove any cosmetics that contain a chemical or acidic element that may affect light output, and lead to a decreased or in-effective treatment.



4.4.2 Connect controller and mask


Insert the cable of mask into USB ports of controller. Both USB ports can be plugged in without any difference.

4.4.3 Place mask on your face

- Make sure that the inside of the mask (Transparency is more visible to the side of the treatment lamp bead) is facing you, and place the eye cups into the mask according to below picture.




- Slide the strap through the slot.
- Attach the mask to your face, adjusting the position of the straps according to the size of your head, and secure it to the back of your head.

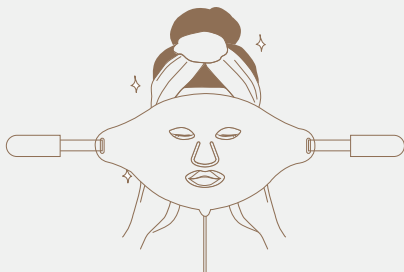


4.4.4 Turn on the device

Long press P/B button on controller to turn on the device

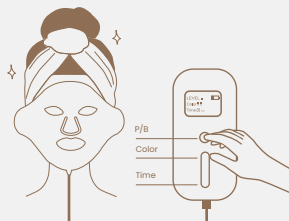

4.4.5 Select the treatment mode

Press Color button to select mode from R+NIR, B, R+NIR+B, default mode is R+NIR. LCD displays the corresponding light mode.

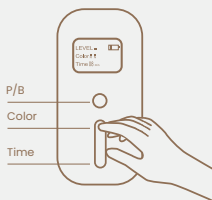

4.4.6 Adjust brightness level


Short press P/B button on controller to adjust brightness level. LCD screen displays corresponding brightness level. There are 5 levels. Increasing from level 1 to level 5, the brightness in the current mode will be brighter,

with level 1 being dark and level 5 being the brightest.

Note: In addition, for initial use, we recommend starting treatment at the lowest level of 1 and gradually increasing to the brightest level of 5.

4.4.7 Set treatment time

Press Time button to set treatment time. Default time is 10 mins.

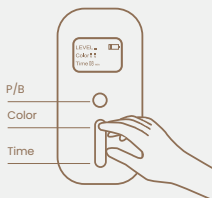

Note: Our recommended time is 10 minutes each time. 3-4 treatment a week, reduce to 1-2 treatment a week once the results shown.

4.4.8 Turn off the device

Long press P/B button to turn off the device or it will automatically shut down when the treatment time is over.

4.4.9 Take off the mask

Remove the mask from your face.

4.4.10 Clean the mask (Please refer to Section 5.1 for details)

5. CLEANING, DISPOSAL AND TROUBLE SHOOTING

5.1 Cleaning

The cleaning methods of the device are instructed as follows:

- 1) The device must be shut down before cleaning.
- 2) Wipe the mask (the inner shell contacting your skin) with soft and clean wipes, and ensure that no liquid enters the interior of the device.
- 3) Allow the device to air dry naturally.
- 4) Use a clean wipe moistened with alcohol (70%) to clean the mask, ensuring that no liquid enters the interior of the device.
- 5) Visually check the cleanliness of the mask and make sure there are no dust, impurities, or any other pollutant; otherwise, please repeat the above step again.
- 6) Check the device to see if any visible contaminants remain. If there are visible contaminants on the device, repeat the above steps.
- 7) Store the mask in your cloth bag until you use it next time.

Note: Do not use any other cleaning products as these may damage the plastics.

- 8) Clean the device regularly to avoid cross-infection.
- 9) After removing the device from storage at minimum or maximum temperature/humidity, please wait at least 1h before use.

Note:

- a) DO NOT immerse the product in liquid.
- b) DO NOT dismantle the housing of the unit by yourself.
- c) DO NOT drop the device or apply external pressure onto the unit since this may damage or weaken the LEDs and other components.
- d) Once the cord or device is damaged or does not operate normally, contact the distributor. This device only is serviced, repaired, and disassembled by Shenzhen Eyco Technology Co., Ltd. All units should be returned to the manufacturer for repair, including Warranty repair and Out-Of-Warranty repair, please do not disassemble or repair the returned units by yourself.

5.2 Storage

Store the mask in the supplied packaging.
Please store the device in Temperature:
-20-50°C Humidity: 10-80%, atmospheric
pressure: 700hPa-1060hPa.

5.3 Trouble Shooting

In this event that the device fails to perform as intended, the following notes will help to identify potential problems with the device and its setup.

Problem	Cause	Solution
The controller can't turn on the mask.	Connection between controller and main unit is loose; The controller has insufficient power.	Reconnect controller to main unit by using USB cable; Charge the controller; If it still does not work, please contact the manufacturer/ distributor.
The controller cannot be charged.	Adapter is not proper.	Ensure that you are using a 5V adapter. If not, the controller will not be able to charge. If you use 5V adapter for charging and still cannot charge successfully, please contact the manufacturer.
	Charging cable is not properly inserted into the power outlet.	Insert the charging cable properly into the power outlet.
	Charging cable is not properly inserted into the controller.	Insert the charging cable properly into the controller.
	Power outlet is not working.	Try a different power outlet.
During use, the mask does not shut off after the end of treatment time.	The device may be damaged.	1. Manually press the P/B button or turn off the device. 2. If the device does not shut off after the previous step, please contact Customer Service.

During use, the mask shut off before the end of treatment time.	The device is not charging; The device may be damaged or there's something wrong with the programming.	1. Turn off the device and restart the device. 2. If the device is still automatically turned off before the end of the treatment time after the previous step, please record the device situation and contact customer service.
The straps on both sides of the mask were unsteadily fixed.	This may impact the effect of treatment or prevent the normal use of the equipment.	1. Check whether there are stains on both sides of the straps that may affect wearing. 2. Contact customer service for assistance.

5.4 Disposal of Electrical and Electronic Equipment Waste (E-waste)

As each city has its regulation of the disposal of e-waste, please learn about recycling or disposal programs in your local area by searching “E-waste recycling or disposal” online. You can also contact the relative department in your city directly to consult the policy.

6. SERVICE

- 1) Customers are provided a one-year warranty service from the date of purchase, but damage caused by time and use is not within the scope of the service guarantee.
- 2) The product warranty period is calculated from the date of purchase by the customer. Please keep the certificate and stamp the official seal of the dealer after receiving the product to ensure your rights.
- 3) After the warranty period expires, the consumer must bear the cost of component replacement and repair.
- 4) The warranty scope of this product does not include failures caused by any factors other than the defects of the product itself.
 - Natural disasters, other unpredictable factors, damage caused by any force majeure factors.
 - Product damage caused by improper use, product damage caused by improper storage, consumption of the product.
 - The power adapter and USB data cable are normal losses and generally not covered by the warranty.
 - The damage caused by disassembling, modifying, and repairing this product by personnel not authorized by the company.
 - The company does not bear the logistics costs of product returns.
- 5) After-sales service or any questions, please consult the dealer or manufacturer.

TEL: +1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

E-MAIL: support@renpho.com

Note: Thank you for choosing our company's products, this product has undergone strict quality

management and has been shipped from the factory. The previous inspection is qualified, and if the product fails under normal use and is not damaged by humans or animals, please refer to the items in the warranty card for repair services.

7. ESSENTIAL PERFORMANCE

The essential performance of the system may be lost or degraded because of electromagnetic disturbances. For expected degradations and instructions of the basic safety and essential performance maintenance in the case of electromagnetic disturbances, see the table below:

Essential Performance	Degradation Caused by Electromagnetic Disturbances	Essential Performance Maintenance
Product Function	No degradation	Not Applicable
Note: The main functions such as the power on or off function, the button function, wavelength energy output function, the function of automatically stopping treatment when the treatment time ends and treatment modes of the device can keep working normally.		

8. SAFETY AND EMC

Before using the device, keep an appropriate distance from radio frequency (RF) sources whenever possible. The sources include but not limit to:

- Radio and TV stations
- Portable and mobile RF communication devices (cell phones, two-way radios, base station, etc.)
- High-frequency surgical units, such as diathermy, electrocautery, argon beam coagulators, etc.
- X-ray, CT, or MRI devices
- Security systems (e.g., electromagnetic anti-theft systems (EAS), and metal detectors)

These devices are also possible sources of interference as they may emit higher levels of electromagnetic radiation.

WARNING:

- Use of this equipment adjacent to or stacked with other equipment should be avoided because it could result in improper operation. If such use is necessary, this equipment and the other equipment should be observed to verify that they are operating normally.
- Portable RF communications equipment (including peripherals such as antenna cables and external antennas) should be used no closer than 30 cm (12 inches) to any part of the LED light therapy mask device or system, including cables specified by the manufacturer. Degradation of the performance of this equipment could result.
- Use of accessories and cables other than those specified or provided by the Shenzhen Eyco Technology Co., Ltd manufacturer of this equipment

could result in increased electromagnetic emissions or decreased electromagnetic immunity of this equipment and result in improper operation.

● Magnetic and electric fields may interfere with the performance of the device or system. Make sure that all peripheral components operated in the vicinity of the device comply with the relevant EMC requirements. X-ray equipment, MRI devices, radio systems (cellular phones) and so forth, are possible sources of interference because they may emit higher levels of electromagnetic radiation. Verify the performance of the system before use.

Guidance and manufacturer's declaration – electromagnetic emissions	
This equipment is intended for use in the electromagnetic environments specified below, and the purchasers or users shall ensure that it is used in these electromagnetic environments.	
Emissions	Compliance
RF emissions(Radiated) CISPR 11	Group 1 Class B
RF emissions (Conducted) CISPR 11	Group 1 Class B
Harmonic emissions IEC 61000-3-2	Class A
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Applicable

Guidance and manufacturer's declaration – electromagnetic immunity

The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Environments
Electrostatic discharge (ESD) IEC 61000-4-2	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	±8 kV contact ±2 kV, ±4 kV, ±8 kV, ±15 kV air	HOME HEALTHCARE environment
Electrical fast transient/ burst (EFT) IEC 61000-4-4	±2 kV for power port	±2 kV for power supply lines	HOME HEALTHCARE environment
Surge IEC 61000-4-5	± 0.5kV, ± 1 kV line to line	± 0.5kV, ± 1 kV line(s) to lines	HOME HEALTHCARE environment
Voltage dips, short interruptions and Voltage variations on power supply input lines IEC 61000-4-11	0% UT, 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% UT, 1 cycle at 0° 70% UT, 25 / 30 cycle at 0° 0% UT, 250 / 300 cycle	0% UT, 0.5 cycle at 0°, 45°, 90°, 135°, 180°, 225°, 270° and 315° 0% UT, 1 cycle at 0° 70% UT, 25 / 30 cycle at 0° 0% UT, 250 / 300 cycle	HOME HEALTHCARE environment

Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	30 A/m	30 A/m	HOME HEALTHCARE environment
Conducted RF IEC 61000-4-6	3 Vrms 0.15 MHz to 80MHz 6 Vrms in ISM bands between 0.15 MHz and 80MHz	3 Vrms 0.15 MHz to 80MHz 6 Vrms in ISM bands between 0.15 MHz and 80MHz	HOME HEALTHCARE environment
Radiated RF IEC 61000-4-3	10 V/m 80 MHz to 2.7GHz 80% AM at 1KHz	10 V/m 80 MHz to 2.7GHz 80% AM at 1KHz	HOME HEALTHCARE environment

NOTE

1. UT is the a.c. mains voltage prior to application of the test level; 25 / 30 and 250 / 300 cycle means 25 / 250 for 50 Hz system and 30 / 300 for 60Hz system.

2. The ISM (industrial, scientific and medical) bands between 0,15 MHz and 80 MHz are 6,765 MHz to 6,795 MHz; 13,553 MHz to 13,567 MHz; 26,957 MHz to 27,283 MHz; and 40,66 MHz to 40,70 MHz. The amateur radio bands between 0,15 MHz and 80 MHz are 1,8 MHz to 2,0 MHz, 3,5 MHz to 4,0 MHz, 5,3 MHz to 5,4 MHz, 7 MHz to 7,3 MHz, 10,1 MHz to 10,15 MHz, 14 MHz to 14,2 MHz, 18,07 MHz to 18,17 MHz, 21,0 MHz to 21,4 MHz, 24,89 MHz to 24,99 MHz, 28,0 MHz to 29,7 MHz and 50,0 MHz to 54,0 MHz.

3. Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radio, AM and FM radio broadcast, and TV broadcast cannot, theoretically, be predicted with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, consider conducting an electromagnetic site survey. If the measured field

strength in the location the system is used exceeds the applicable RF compliance level listed in this table, observe the system to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the system.

4. At 80 MHz and 800 MHz, the higher frequency range applies.

5. These guidelines may not apply in all situations. Electromagnetic propagation is affected by the reflection from structures, objects, and people.

6. Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.

7. This test is applicable only to ME EQUIPMENT and ME SYSTEM intended for use in the HOME HEALTHCARE ENVIRONMENT.

8. The carrier shall be modulated using a 50% duty cycle square wave signal.

9. r.m.s., before modulation is applied.

Test specifications for ENCLOSURE PORT IMMUNITY to RF wireless communications equipment

Test frequency (MHz)	Band ^{a)} (MHz)	Service ^{a)}	Modulation ^{b)}	IMMUNITY TEST LEVEL (V/m)
385	380-390	TETRA 400	Pulse modulation ^{b)} 18Hz	27
450	430-470	GMRS 460, FRS 460	FM ^{c)} ±5 kHz deviation 1 kHz sine	28
710	704-787	LTE Band 13, 17	Pulse modulation ^{b)} 217Hz	9
745				
780				

810	800-960	GSM 800/900, TETRA 800, iDEN 820, CDMA 850, LTE Band 5	Pulse modulation ^{b)} 18Hz	28
870				
930				
1720	1700-1900	GSM 1800; CDMA 1900; GSM 1900; DECT; LTE Band 1, 3, 4, 25; UMTS	Pulse modulation ^{b)} 217Hz	28
1845				
1970				
2450	2400-2570	Bluetooth, WLAN, 802.11 b/g/n, RFID 2450, LTE Band 7	Pulse modulation ^{b)} 217Hz	28
5240	5100-5800	WLAN 802.11a/n	Pulse modulation ^{b)} 217Hz	9
5500				
5785				

9. LABEL SYMBOLS

1		Manufacturer
2		The symbol indicates that the product should not be discarded as unsorted waste but must be sent to separate collection facilities for recovery and recycling.

3		CONSULT INSTRUCTIONS FOR USE
4		The date of manufacture
5		Type BF Applied Part
6		Batch code
7		Symbol for "SERIAL NUMBER". This symbol shall be accompanied by the manufacturer's serial number
8		Unique device identifier
9		Recycle
10		Temperature limitation
11		Humidity limitation
12		Atmospheric pressure limitation
13		Keep dry

14		Protects against fixed foreign objects with a diameter of more than 12.5 mm and is waterproof against drips with an inclination of 15°.
15		Medical device
16		MR Unsafe
17		Warning
18	RISK GROUP EXEMPT This product does not emit UV light and does not cause any serious harm.	Lamp classification group in IEC 62471:2006

10. MANUFACTURER INFORMATION



Manufacturer: Shenzhen Eyco Technology Co., Ltd
 2 /F, Building 1, Jiashitai Science Park, 164# Fengtang
 Ave Fuhai Street, Baoan District, Shenzhen 518103,
 China Guangdong, CN 518103

RENPHO
REJUVEN