
RENPHO

USER MANUAL

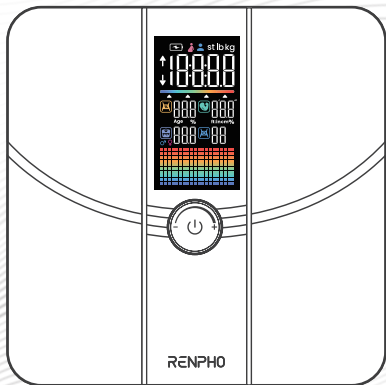ELIS 3X

Model: R-MSB04

CONTENTS

• ENGLISH	01-12
-----------	-------

• DEUTSCH	13-25
-----------	-------

• FRANÇAIS	26-37
------------	-------

• ITALIANO	38-49
------------	-------

• ESPAÑOL	50-61
-----------	-------

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. The scale is for household use only and not for medical or commercial purposes. It should not be used in diagnosing a disease, disorder or abnormal physical state.
4. Children under 8 or individuals who are pregnant, have osteoporosis, a pacemaker, artificial joints, portable ECG or other metal implants should use the scale for weighing purposes only and make sure to wear shoes to prevent any electric current when stepping onto the scale.
5. All body composition data provided is for reference and should not be used as a substitute for advice from a licensed healthcare professional.
6. Do not overload the scale. (Maximum Capacity: 400 lbs / 180 kg)
7. Inspect the scale for any damage before each use. Do not use a damaged scale to avoid inaccuracy or injury.
8. Only use the provided charging cable and charger (DC 5V 0.5A, not included) to charge the scale.
9. To ensure accuracy, use the scale on a hard, flat and dry surface.
10. To avoid interference or danger, do not use the scale in an environment near equipment with strong electromagnetic fields.
11. To prevent injury, do not stand on the edges of the scale, jump on or stomp it.
12. The scale is a delicate and precise measuring device. Handle with care to avoid breakage.
13. The scale is not waterproof. Keep it dry and never submerge it in water. To clean the surface, use a damp cloth or glass cleaner. Do not use soap or chemical cleaning agents.
14. Do not expose the scale to direct sunlight or use it in high temperatures or extremely humid areas. Store the scale in a cool, dry place.
15. Do not attempt to disassemble the scale for any reason, as this may cause damage or void the warranty.

About the Battery

1. Do not expose the appliance or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated appliance or battery being short- circuited by metal objects.
3. This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
10. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100 °C (212 °F).

SAVE THESE INSTRUCTIONS

WHAT'S IN THE BOX

Elis 3X x1

USB-C Charging Cable x1

TECHNICAL SPECIFICATIONS

Product Name: Body Fat Scale

Model: R-MSB04

Maximum Capacity: 400 lbs / 180 kg

Division: 0.2 lb / 0.05 kg


Dimension: 11.02 x 11.02 x 0.98in / 280 x 280 x 25mm

Battery Capacity: 500mAh 1.85wh

Input: 5V $\overline{\text{---}}$ 0.5A

Rated Power: 2.5W

ABOUT YOUR PRODUCT

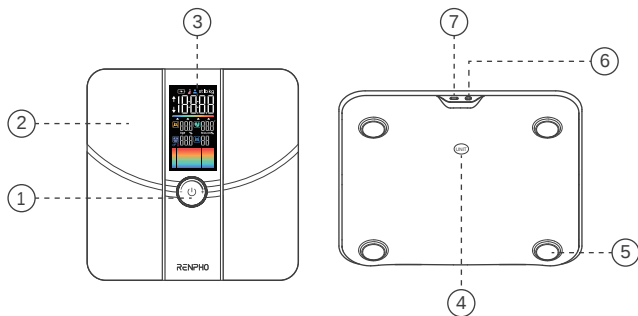1. Dial (For inputting personal information)
2. ITO Coated Platform
3. VA Display
4. Unit Conversion Key
5. Foot Pads
6. Reset Button
7. USB-C Charging Port

Icons

Icons	Description
	Thin / Standard / Overweight / Severely Overweight
	Body Fat Percentage
	Muscle Percentage
	BMI
	Visceral Fat Index
	Male / Female
	Weight Trend Each section represents the most recent weight measurement taken over individual days.
	Weight Increase / Weight Decrease
	Battery Indicator ⚡ Charging 🔋 Fully Charged
	Pregnancy Mode
	User Recognition Indicates existing user profile history matches the current weight reading.


Icons	Description
st bkg / ft:incm	Units Tap and hold Unit Conversion Key to change the unit.
Err	Overloaded Maximum Capacity: 400 lbs / 180 kg.
Lo	Low Battery Charge the scale.
CAL	Calibrating Calibrate the scale before each use.

CHARGING YOUR PRODUCT

Charging Time: approx. 3 hours


DC 5V 0.5A Charger (Not Included)

USB-C Charging Port

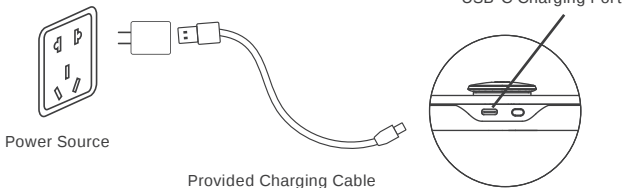
INITIALIZING/RESETTING YOUR PRODUCT

1. Before first use, press the Reset Button or connect the scale to a power source using the provided charging cable and charger (not included) to initialize the scale.
2. To reset the scale, press the Reset Button.

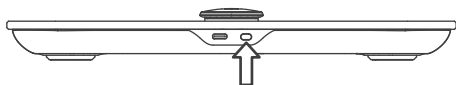CALIBRATING YOUR PRODUCT

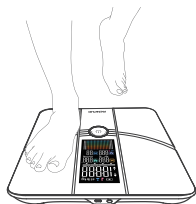CAL → 0.0

or

CAL → 0.00

Calibrate the scale before each use to ensure accuracy.

1. Place the scale on a flat, hard surface. Avoid any carpeting.
2. Step on the scale with a single foot to activate it and then step off.
3. The display will show “CAL” and then “0.0” or “0.00” to indicate calibrating is complete.

USING YOUR PRODUCT

1. Make sure the scale is placed on a flat, hard surface.
2. All features in the Icon section can be accessed with or without using the RENPHO Health App.

Method 1: Use with dial

A. Press the dial to turn on the scale.

Note: The display will automatically turn off after 10 seconds of inactivity. Avoid extended pauses during setup to prevent restarting the process.

B. Create a Profile

Step 1: Turn the dial to select “User 1-8” or “GUEST”. Press the dial to enter the profile settings. Please note that any changes made will not be saved in “GUEST” upon auto-off.

Step 2: Turn the dial to select gender. Then press the dial to confirm.

Step 3: Turn the dial to select age. Then press the dial to confirm.

Step 4: Turn the dial to select height. Then press the dial to confirm. When selecting height, press and hold the unit conversion key at the back to switch between cm/ft.in.


Step 1

Step 2

Step 3

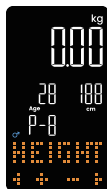Step 4

- C. After inputting all information and calibrating the scale, step on the scale with bare and dry feet, keep still and wait for the final measurement to be shown on the scale's display.

Select user

Turn the dial to cycle through users and press to confirm. Select “YES” to enter weighing, select “no” to revise your information.

Note:

- User Identification: After profiles are set, the scale automatically identifies the user for each weight reading. If weights are similar, select the user before weighing.
- Guest Mode: Data in "GUEST" mode will not be saved after auto-off.

Method 2: Use with App

A. Download the RENPHO Health App

Scan the QR code below or search for “RENPHO Health” from the Apple Store or Google Play to download the app.

System requirement: iOS 12.0 or later, Android 7.0 or later

Note: Due to continuous updates and improvements, the “RENPHO Health” App may appear slightly different.

B. Log in / Sign up



- Log in directly with an existing RENPHO App account.
- If you do not have any RENPHO App account, tap "Sign up" using your email address and fill in the personal information; this step is required for calculating your health measurements accurately.

Note: You may also log in using Guest Mode if you do not wish to register. However, your measurement data will not be saved if you log it in on another device.

C. Pair the Scale

Before pairing, enable the wireless setting on your smartphone.

- Open the RENPHO Health App, place the scale on a hard and flat surface, and step onto it with bare and dry feet.
 - Tap "+" in the upper-right corner of the "Home" page and select your body fat scale from the device list.
 - Follow the in-app instructions to complete pairing.
- D. After calibrating the scale, step on the scale with bare and dry feet, keep still and wait for the final measurement to be shown on the scale's display.
- E. After measuring, the scale will transmit the data to the app. Open the app to sync and view your detailed measurements, track your progress, and access additional features and insights.

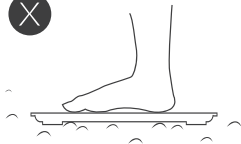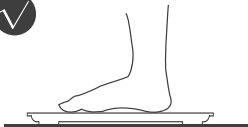TROUBLESHOOTING GUIDE

1. Display Issues

Display	Measurement	Description / Cause
	Weight, BMI, body fat percentage, muscle percentage, visceral fat index	• User identification was detected by the scale. • This measurement was taken with bare feet.
	Weight, BMI	• User identification was detected by the scale. • This measurement was taken with socks or shoes.
	Weight only	• User identification was not detected by the scale. • This measurement was taken with bare feet.
		• User identification was not detected by the scale. • This measurement was taken with socks or shoes.

2. Other Issues

Problem	Solution
The scale cannot detect body fat and / or other data.	a. Ensure your feet are dry and bare. b. Stand and keep still on the scale for about 10s to finish the weighing process. If still not available, try applying a little moisturizer to your feet and measure again or ask another person to repeat the 2 steps above, please be careful not to slip on the weighing platform.
I registered an account on the RENPHO website but cannot use it to log in on the RENPHO Health App.	The RENPHO website and RENPHO Health app are two different systems. Create an account on the RENPHO Health App.
I can get my data of BMI, body fat, muscle mass, body water, bone mass and metabolic age from the scale's display, but my family or friends failed.	Your family or friends must complete their profile information on scale via dial or on your app using the "add a member" function or log in to the RENPHO Health App with their own account.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not void the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns.

RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

✉ **Email:** **support-eu@renpho.com (UK&EU)**
support@renpho.com (US&CA)
support-au@renpho.com (AU)

☎ **TEL:** **+44(0)7434 666088 (UK ONLY)**
Monday-Friday 9:00AM-5:00PM (GMT)
+1(844) 417 0149 (US&CA)
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)
+49 40 23969285 (DE ONLY)
Monday-Friday 9:00AM-6:00PM (CET)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

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