
RENPHO

User Manual

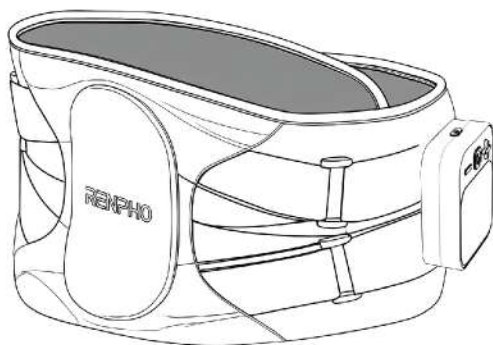Model: R-HAP65

Emberace Dorsi

CONTENTS

- ENGLISH 01-11

- DEUTSCH 12-22

- FRANÇAIS 23-33

- ITALIANO 34-44

- ESPAÑOL 45-55

IMPORTANT:

Read this manual carefully before first use and keep it in a convenient place for future reference.

	DANGER
	Hot Surfaces.

DANGER		
• Do not sit on the pad.	• Do not use while sleeping. • Never put the pad under the body.	• Do not fold the pad. • Avoid sharp folds, as it may damage the product.

ICONS

This guide contains important information about using, operating, and storing the product. Please read it through and note that non-observance of these instructions can result in serious injury or damage to the product.

WARNING

These notes must be observed to prevent any damage to the product.

CAUTION

These notes must be observed to prevent any injury to the user.

NOTE

These notes provide useful additional information on the installation or operation.

CAUTION

To reduce the risk of burns, electric shock and fire, this pad must be used in accordance with the following instructions.

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING THIS DEVICE.

1. This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the device.
3. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
4. The device must not be located immediately below a socket-outlet.
5. When storing the device, allow it to cool down before folding.
6. Do not crease the device by placing items on top of it during storage.
7. Examine the device frequently for signs of wear or damage.
8. The device is for home use only and is not intended for any medical, diagnostic, or therapeutic purposes. Its design and function are solely focused on daily comfort and relaxation, rather than any medical use.
9. Do not use the device if wet.
10. This device must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
11. The device is not to be used by young children over the age of three unless the controls have been pre-set by a parent or guardian, or the child has been adequately instructed on how to operate the controls safely.
12. Do not use if folded or creased.
13. Do not use this device with an infant, an immobile person, or a person with poor blood circulation.
14. Do not tuck in the wired area of this device.
15. Do not use dry-cleaning fluid on this device, cleaning solvents may have a deteriorating effect on the insulation of the heating element.
16. Do not route electrical cords on the device. Do not allow the cords to be pinched.
17. When in use, keep the control away from drafty areas that may become damp or wet such as an open window.
18. Do not poke the device with pins - They may damage the electric wiring.

19. Be sure that this device is used on an alternating-current supply circuit of the proper voltage.
20. Turn off or unplug when not in use.
21. Excessive heating may result from folding device when in actual use.
22. If improper operation of this device is observed, discontinue its use immediately and consult RENPHO Customer Service regarding its repair.
23. Do not tuck the device or supply cord in a sofa or other furniture. All portions of the device should be visible at all times.
24. Wrapping of the cord around the controller may damage the cord. Loop loosely when storing.

DANGER - to reduce the risk of burns, electric shock and fire, this device must be used in accordance with the following instructions.

1. Do not destroy.
2. Do not use while sleeping.
3. Do not use on an infant.
4. This device is not to be used on or by an invalid. Sleeping or unconscious person, or a person with poor blood circulation unless carefully attended.
5. Do not use on areas of insensitive skin.
6. Burns can occur regardless of control setting, check skin under the device frequently.
7. Do not use in an oxygen atmosphere (applies to hospital devices).
8. Do not use pins or other metallic means to fasten this device in place.
9. Do not sit on or crush the device - avoid sharp folds.
10. Never pull this device with the supply cord.
11. Do not use the cord as a handle.
12. Carefully examine the inner cover before each use. Discard the device if the inner covering shows any sign of deterioration.

About the Battery

1. Do not expose the device or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated device or battery being short-circuited by metal objects.
3. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product

- manual for proper charging instructions.
8. Do not leave a battery on prolonged charge when not in use.
 9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
 10. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

SAVE THESE INSTRUCTIONS

ABOUT THE PRODUCT

1. WHAT'S IN THE BOX

Emberace Dorsi × 1

Extension Strap × 1

USB-C Charging Cable × 1

2. TECHNICAL SPECIFICATIONS

Product Name: Heating Pad

Model: R-HAP65

Input: 5V $\equiv$ 2A

Rated Power: 10 W

Battery Voltage: DC 3.7V

Battery Capacity: 4000mAh 14.8Wh

Auto-off Time: 20 min



Extension Strap Length: 13.8 in / 35 cm

Product Size: 43.3 × 8.5 in / 110 × 21.5 cm

Temperature Range: 104°F - 140°F / 40°C - 60°C

Applicable Waist Circumference: 25.6 - 55 in / 65 - 140 cm

3. Product Overview

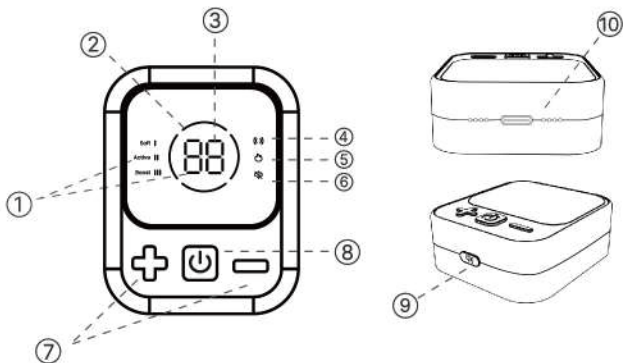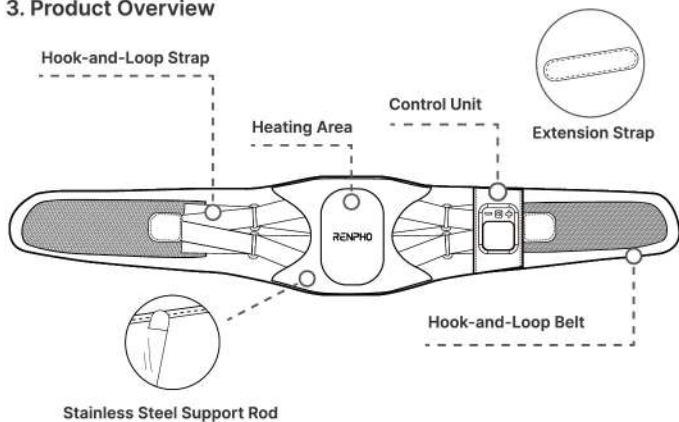1. Vibration Level Indicator
2. Battery Indicator
3. Heat Level Indicator
4. Vibration On Indicator
5. Heat On Indicator

6. Audio Off Indicator
7. Heat Button
8. Power / Vibration Button
9. Audio Button
10. USB-C Charging Port


USING THE DEVICE

Before Use



- * Wear thinner clothing and place the device on the clothing to use it instead of using it directly on your skin.
- * Keep the control unit side facing outward during use.

1. Wear the Device

Step 1: Position the heating area over your spine at the lower back.

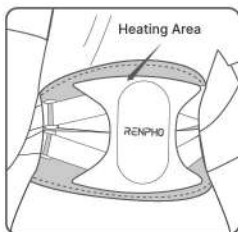
Step 2: Firmly and evenly pull both sides of the hook-and-loop belt and the hook-and-loop strap around your waist, overlapping them securely across your abdomen.

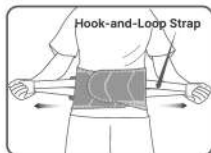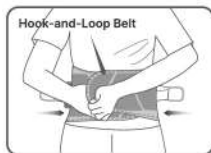2. Use the Extension Strap (For Larger Waist Circumference)

Connect the extension strap to the hook-and-loop belt, then follow the instructions in "1. Wear the Device".

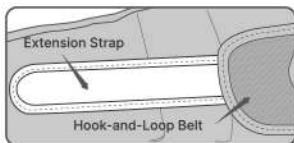BUTTON INTRODUCTION

Button	Operation
	Press and hold to turn On / Off the device. Press to increase / decrease the vibration level: L1 (Soft) (default), L2 (Active), L3 (Boost) and L0 (Off).
	Press to increase the heat level: H1, H2, H3 (default), H4, H5 and H0 (Off).
	Press to decrease the heat level: H4, H3 (default), H2, H1, and H0 (Off).
	Press and hold to turn the audio On (default) / Off.

BATTERY INDICATOR

Battery Level	Indicator
70% - 100%	Solid Green
40% - 70%	Solid Yellow
5% - 40%	Solid Red
<5%	Flashing Red

CLEANING & MAINTENANCE

1. Use a wet towel or brush to clean the surface of the device.
2. To store the device, loosely roll it toward the straps, loop the straps over the rolled ends, and then store the device in a cool, dry, dark place after each use. Avoid sharply folding the device.

BUTTON INTRODUCTION

Problem	Cause	Solution
The device stops working suddenly.	Batteries are weak or depleted.	Charge the device and restart it.
	20-min auto-off timer is activated.	Restart the device if needed.
No heating or insufficient heating	No temperature / low heat is selected.	Adjust the heat level.
	The device needs time to warm up.	Allow up to 1 minute to heat up.
	The device is overheating.	Built-in overheat protection. Allow the device to cool down and restart it.
Fail to turn on the device.	Batteries are weak or depleted.	Charge the device and restart it.
	Power Button was pressed for only a short duration.	Press and hold the Power Button to turn on the device.
Slow charging / no charging.	Use other types of charging cables.	Use the provided charging cable only.
Feels too hot.	Higher heat level is selected.	Adjust to a lower heat level.
	The device was used directly on the skin.	Wear thinner clothing and place the device over the clothing.
The device does not conform comfortably to your abdominal curve.	The support rods are applying excessive pressure to your body.	Remove the support rods to relieve pressure.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty.

If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

-  **TEL: +44(0)7434 666088 (UK ONLY)**
Monday-Friday 9:00AM-5:00PM (GMT)
+1(844) 417 0149 (US&CA)
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)
+49 40 23969285 (DE ONLY)
Monday-Friday 9:00AM-6:00PM (CET)
-  **Email: support-eu@renpho.com (UK&EU)**
support@renpho.com (US&CA)
support-au@renpho.com (AU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2011/65/EU, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:
<https://renpho.eu/pages/compliance-data>

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