
RENPHO

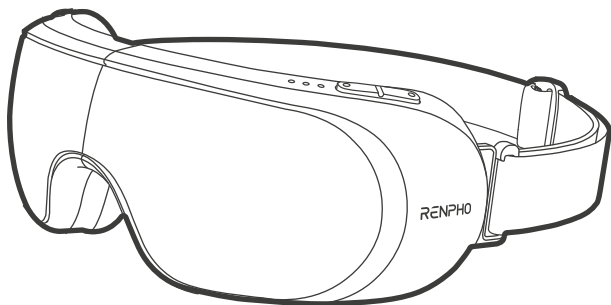Eyeris Thermacool

Eye Massager

User Manual

Model: R-WME06

Contents

English..... 01~15

Deutsch..... 16~30

Français..... 31~45

Italiano..... 46~59

Español..... 60~73

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.

1. This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.
4. Cleaning and maintenance should not be conducted by children.
5. Please consult your doctor before use if you fit into any of the following user groups:
 - Persons who had undergone eye surgery or eye conditions such as glaucoma, cataract and retinal detachment
 - The elderly
 - Pregnant or postpartum women
 - Persons who are allergic to leather material
 - Persons with skin infections, sensitivity, redness and/or burns on the skin
 - Persons with any swelling or wound around the eye area
 - Persons using a pacemaker or other implantable medical devices
 - Persons with a heart condition, malignant tumors, or other serious health conditions
 - Persons undergoing medical treatment
 - Persons with high or abnormal blood pressure
6. If you experience any discomfort during use, stop using the device immediately.
7. Remove contact lenses, glasses, false eyelashes, hairpins, and any other accessories before use.
8. Do not use the device while driving, walking, or running.
9. Avoid contact with sharp objects.
10. High temperatures may cause burns to occur regardless of the

control setting. Do not use the device on infants, persons who are with disabilities, sleeping or unconscious. Check frequently where the skin is in contact with the heated area to reduce the risk of blistering.

DANGER - To reduce the risk of electric shock:

1. Never use pins or other metallic fasteners with this device.
2. Do not use the device with wet hands or submerge the device in water or any liquid.

CAUTION - To reduce the risk of burns:

1. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
2. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

1. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
2. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
3. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this device if it is not working properly, has been damaged, or dropped into water.
5. Never drop or insert any object into any opening.
6. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

About the battery:

1. Do not use while bathing or in a shower.
2. Do not place or store appliance where it can fall or be pulled into a tub or sink.
3. Do not place in or drop into water or other liquid.
4. Do not expose the appliance or battery to excessive temperatures.

5. Be aware of the risk of terminals of the battery-operated appliance or battery being short-circuited by metal objects.
6. This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.
7. Do not dismantle, open or shred battery.
8. Do not subject cells or batteries to mechanical shock.
9. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
10. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
11. Do not leave a battery on prolonged charge when not in use.
12. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
13. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

SAVE THESE INSTRUCTIONS



WHAT'S IN THE BOX

Eyeris Thermacool x 1

Storage Bag x 1


USB-C Charging Cable x 1

ABOUT YOUR EYE MASSAGER

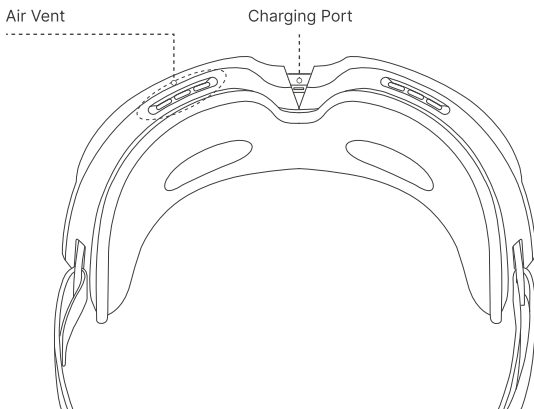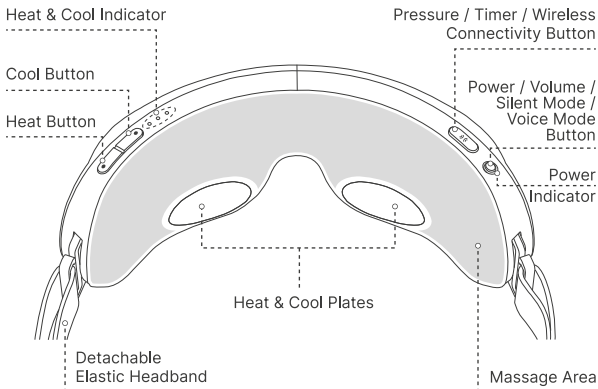CHARGING YOUR EYE MASSAGER



Fully charge the device before initial use.

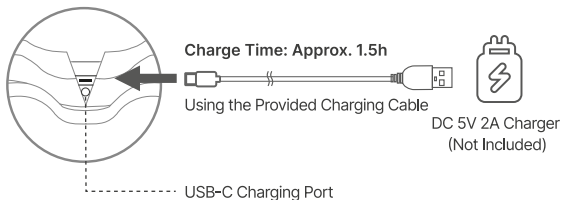Note: The device cannot be turned on while charging.

WEARING YOUR EYE MASSAGER

1. Before wearing the device, remove all makeup and any eye accessories (e.g. eyeglasses, false eyelashes, etc.).
2. Wear the device and ensure it fits over your eyes comfortably. The heat & cool plates are meant to be applied on the eye bag areas. Adjust the elastic headband with the buckle at the back to find the right fit for your head.

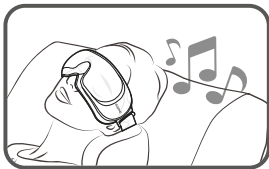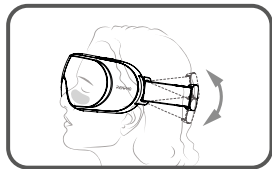Note:

- On first use, it is recommended for you to start with soft pressure before switching to a stronger pressure once you are comfortable.
- People with low heat sensitivity or delicate skin are recommended to start on the low heat setting.

- If you feel any pain during your session, it is recommended for you to switch back to a softer pressure or loosen the headband.
- It is not recommended that the device be used for more than 20 minutes at a time.
- Your vision may become blurry after using the device. This is an expected occurrence, and the effect is only temporary. In the case of blurred vision, it is recommended that you close your eyes for a few minutes until your vision goes back to normal.
- After use, it is recommended that you wipe any sweat, eye cream, or cosmetic residue left on the on the protein leather surface and heat & cool plates. Keeping the device clean will prolong the lifespan of the device.

USING YOUR EYE MASSAGER

Tip: To use voice control, turn on the device first, say “Hi Eyeris” and proceed with your clear voice command within 8 seconds once you hear “Welcome”. If you fail to give a clear voice command within 8 seconds, repeat the above steps, excluding turning on the device.

Turn on the eye massager
and say “Hi Eyeris”

The eye massager will
respond with “Welcome”

Give a clear voice command to
the eye massager within 8s, such
as “Soft Pressure”

The eye massager will respond
with the corresponding voice
broadcast before execution

Button	Manual Control	Voice Control	Voice Broadcast
/	/	Hi Eyeris <i>(Turn on the device and wake it up when you want to use the voice-control function)</i>	Welcome
	a. Press and hold to turn on / off the device	/	Welcome to the RENPHO Eye Massager
		Turn Off	Goodbye
	b. Press to cycle the music volume: level 0 (default) / 1 / 2 / 3 / 4 / 5	Increase Volume Decrease Volume	Maximum Volume <i>(Level 5)</i> Mute <i>(Level 0)</i> <i>(No level 0 voice prompts during button-controlled operation)</i>
	c. Double press to switch between voice mode and silent mode	Sound On Sound Off	Voice Mode Silent Mode
	a. Press to adjust the air pressure level: Soft pressure: default / Medium pressure / Strong pressure / No pressure	Soft Pressure Medium Pressure Strong Pressure No Pressure	Soft Pressure Medium Pressure Strong Pressure No Pressure
	b. Press and hold to select the auto-off timer: 10 min / 15 min (default) / 20 min	Change Time	Ten Minutes Fifteen Minutes Twenty Minutes

Button	Manual Control	Voice Control	Voice Broadcast
	c. Double press to turn off wireless connectivity <i>(Only applicable when wireless connectivity is connected)</i>	/	Bluetooth Disconnected
	a. Press to enable the heat function and adjust the heat level: Low heat: 97°F (36°C) / Medium heat: 100°F (38°C), default / High heat: 104°F (40°C)	Low Heat Medium Heat High Heat	Low Heat Medium Heat High Heat
	b. Press and hold to switch back to the previous music track	/	Previous Track
	a. Press to enable the cool function adjust the cooling level: Low cool: 66°F (19°C) / Medium cool: 63°F (17°C), default / High cool: 59°F (15°C)	Low Cool Medium Cool High Cool	Low Cool Medium Cool High Cool
	b. Press and hold to switch to the next music track	/	Next Track
	Press and hold simultaneously to disable the heat & cool function <i>(Press heat/cool button to enable the heat/cool function)</i>	No Temperature	No Temperature

Power Indicator

Status	Indicator
Normal operation.	Solid white.
Wireless connectivity is connected.	Solid blue.
Low battery.	Solid red.
Charging.	Breathing red.
Fully charged.	Solid white.
Charging at a voltage higher than 5V.	Continuous quick red flashing for 1 minute followed by turning off, which indicates an incompatible charger was used.

Note:

- The heating and cooling functions cannot be used simultaneously.
- Memory function: Recalls last settings (music volume excluded) when turned on.
- Use a DC 5V 2A charger (not included). The battery indicator turns breathing red during charging and solid white when fully charged.

CONNECTING TO YOUR SMARTPHONE

1. Turn on the device and enable wireless connectivity on your smartphone.
2. Search for "R-WME06" via your wireless connectivity of your smartphone and connect to your smartphone to play your music.
3. Once the device has been connected successfully, the power indicator will turn solid blue, followed by a voice prompt.

Note:

If the initial pairing is successful, the device will automatically connect to your smartphone in future sessions, as long as the wireless connectivity has been enabled.

CARE & MAINTENANCE

1. Turn off the device and allow it to cool down before cleaning or storing it away.
2. Use a soft cotton cloth that has been slightly dampened with mild detergent to wipe the surface of the device. Do not use cloth that has been submerged in water, bleach, or other corrosive cleaners.
3. Ensure the device is clean and dry before storing it in a cool, dry place. Keep it away from direct sunlight, moisture, heat, and corrosive gases.
4. Charge the device periodically to prolong its lifespan. If the device has not been used for over 3 months, charge the device before use.
5. Do not disassemble or repair the device by yourself, as this may damage the device and void the manufacturer's warranty.

TECHNICAL SPECIFICATIONS

Product Name: Eye Massager

Model: R-WME06

Input: 5V $\overline{=}$ 2A

Battery Capacity: 1700mAh 6.46Wh

Battery Voltage: DC 3.8V

TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
Cannot turn on the device.	The battery is weak or has been depleted.	Charge the device and restart it.
	The device is charging.	Pull the charging cable / plug out and restart the device.
	The Power Button was pressed for only a short duration.	Press and hold the Power Button to turn on the device.
Insufficient heat / cooling or no heat / cooling.	No temperature / low heat / low cool was selected.	Press the Heat / Cool Button to adjust the level.
	The device takes time to heat up gradually to avoid discomfort caused by rapid heat.	Wait several minutes after turning on the device.
	The device has overheated, and the built-in overheat protection has turned the device off.	Allow the device to cool down before restarting it.
The compressions are causing pain.	The elastic band has been strapped too tightly.	Adjust the headband with the buckle and find the right fit for your head before continuing to use the device.
	Strong Pressure has been selected.	Switch to Soft Pressure.
	The device hasn't been properly aligned with the eyes and temples.	Refer to the "Wearing Your Eye Massager" section.
Cannot connect to the device via wireless connectivity.	The wireless connectivity setting has not been enabled on your mobile device.	Enable the wireless connectivity setting on your mobile device.
	The device is currently connected to another mobile device.	Double press the Pressure Button to disconnect, then connect to the new device.

Problem	Possible Cause	Solution
Charging failed or charging is very slow.	The charger used was incompatible (voltage too high).	Use a DC 5V 2A charger (not included).
	Another type of charging cable was used.	Use the provided charging cable only.
No music or the music sounds faint.	The music streamed via wireless connectivity is either at volume level 0 or a low volume level.	Adjust the volume level on your mobile device.
	The device is set to either volume level 0 or a low volume.	Press the Power Button to adjust the volume level.
No voice prompts broadcasted.	Silent Mode is currently selected.	Double press the Power Button to switch to voice mode.
No pressure is being applied.	No Pressure is currently selected.	Press the Pressure Button to switch the pressure level.
The device stopped working suddenly.	The battery is weak or has been depleted.	Charge the device.
	The auto-off timer was activated and turned off the device.	Restart the device if needed.
Air remains inside the device's chamber, even when powered off.	The device is not expelling the residual air.	Turn the device off and wait 10 seconds for the device to deflate before removing it from the head. Restart after the device has fully deflated.
The internal fan continues running after switching from cool to heat function or turned off.	Cool the device to ensure safe operation and component longevity.	No action is needed. The fan will stop automatically in 1 minute.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one year limited manufacturer warranty from the date of delivery. For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not void the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

-  **Tel:** **+1(844) 417 0149 (US&CA)**
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)
+44(0)7434-666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM(GMT)
-  **Email:** **support@renpho.com (US&CA)**
support-au@renpho.com (AU)
support-eu@renpho.com(UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dis-pose of any product parts as they may be required for inspection/repair.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1: This device may not cause interference.
- 2: This device must accept any interference, including interference that may cause undesired operation of the device.

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1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

3. Personal data

All end-users of WEEE are reminded that you are responsible for deleting personal data from the WEEE to be disposed of.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/53/EU, Directive 2011/65/EU, Radio Equipment Regulations 2017, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:
<https://renpho.eu/pages/compliance-data>

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