RENPHO

USER MANUAL

MASSAGE CUSHION

Model:RF-BM066

CONTENTS

• ENGLISH	01-13
-----------	-------

• DEUTSCH	14-25
-----------	-------

• FRANÇAIS	26-38
------------	-------

• ITALIANO	39-50
------------	-------

• ESPAÑOL	51-62
-----------	-------

Thank you for your purchase of RENPHO Massage Cushion. The massager offers a powerful massage. While enjoying muscle and body relief, you should use it with care and safety in mind. To ensure safety and avoid any property damage, please read this manual carefully and use as indicated. Non-compliance may cause serious injuries or death and/or severe property damage.

IMPORTANT INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. **WARNING** – To reduce the risk of burns, fire, electric shock, or injury to persons.
 - a) This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
 - b) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
 - c) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
 - d) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 - e) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.

- f) Do not carry this appliance by supply cord or use cord as a handle.
- g) Keep the cord away from heated surfaces.
- h) Never operate the appliance with the air openings blocked.
Keep the air openings free of lint, hair, and the like.
- i) Never drop or insert any object into any opening.
- j) Do not use outdoors.
- k) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- l) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 4. The appliance is only to be used with the power supply unit provided with the appliance.
- 5. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
- 6. **WARNING:** Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- 7. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.

WARNING

- 1. This appliance should not be used for any medical purposes. It is designed for household use only.
- 2. **WARNING** -To reduce the risk of fire, electric shock, or injury to persons, please follow the instructions below:
Consult a doctor before using the device, if you fall into any of the following user groups:
 - Pregnant or lying-in women
 - Children
 - Elderly (over 65 years old)

- Persons using a pacemaker or other implantable medical device
 - Persons with heart disease, malignant tumors, or other serious health conditions
 - Persons receiving medical treatment or feel discomfort during massages
 - Persons suffering from diabetes, or thrombosis, varicose veins
 - Persons with high or abnormal blood pressure
 - Persons with skin infections, sensitivity, redness, wound, and/or burns on the skin
3. DO NOT use this device for more than 30 minutes at a time.
 4. DO NOT sit or stand on the massage (back) part of the device.
The device should only be used when attached to a chair with the back in the vertical position.
 5. When the device is working, do not touch the moving massage mechanism with any body part except your back.
 6. To avoid pinching, DO NOT lean on the back massage mechanism in the cushion when adjusting your body position.
 7. Caution is needed when using it on upholstered furniture.
 8. The appliance must not be located immediately below a socket-outlet.

SAVE THESE INSTRUCTIONS

About the RENPHO Massage Cushion

1. What's in the box

RENPHO Massage Cushion x 1

Adapter x 1

User Manual x 1

2. Product Introduction

3. Using the RENPHO Massage Cushion


3.1 Important tips before use

- It is not recommended to use this Massage Cushion in the driver's seat of a vehicle.
- It is not recommended to lay Massage Cushion completely horizontal when in operation.

3.2 Preparation

- Place the massager on a recliner, sofa seat or office chair, fasten the straps and secure the massager to the chair.
- **When the adapter is plugged in, the power button on the controller will flash red. The back massage nodes will move and return to the default position and once complete, the red light will go out and enter standby mode.** You can press the power button now to turn the massager on or you can just press the power button and use all functions after powered on without waiting for the reset process to be complete.

3.3 Operation Instruction

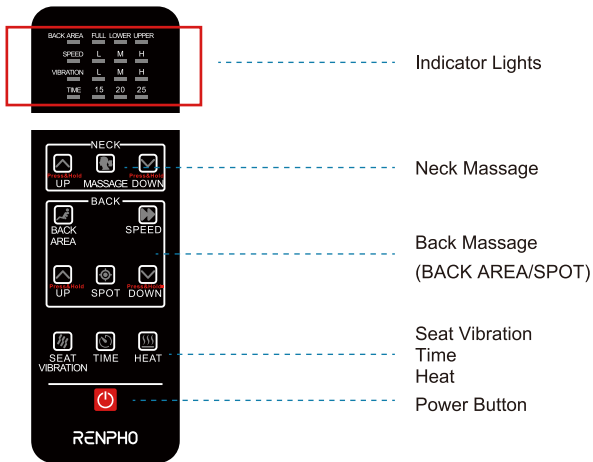1. Press the “Power” Button to turn on the massager; the indicator light will turn red indicating the massager is on. The Neck Kneading, Full Back Massage functions and the Heating function of neck and back are on by default. Seat Vibration function is off. The default timer is 15 minutes. The indicators on the controller will remain in the “on” status accordingly.
2. Activate the functions according to your preferences.

Power

“”: Press to turn the massager on (indicator light will turn on). Press again to turn off (flashing light indicates that the massager nodes are returning to its reset configuration; the lights will turn off once reset is complete). During the reset process, there is no need to wait until the reset is complete to turn the massager on. You can press this button to turn it on immediately.

Neck Massage

“MASSAGE” button: Press this button to turn on the neck massage nodes (indicator light will turn on). Press this button again to turn the neck massage nodes off (indicator light will turn off).

“Up/Down” button: Press and hold this button to have the neck massage nodes move up/down (indicator light will turn on). The message nodes will move upward/downward to the highest/lowest point, and then stop (indicator light will flash). Release the button to have the neck massage nodes stop moving and start massaging at the specified point (indicator light will turn off).

Regardless of whether the neck massage is on or off, the position of the neck massage nodes can be adjusted.

Back Area Massage

“BACK AREA” button: Press this button to turn on/off the back massage function and switch among FULL/LOWER/UPPER back massage mode. (Corresponding indicator light will turn on).

The back massage nodes will stop once the reset is complete. During the reset process, the back massage could be powered on by pressing “BACK AREA” or “SPOT” button.

Back Spot Massage

“SPOT” button: Press this button to activate spot massage mode (the message nodes will move up and down on the specified spot; indicator light will turn on). Press it again to turn off the back spot massage mode (indicator light will turn off).

The back massage nodes will stop once resetting is complete. During the reset process, the back massage could be powered on by pressing this button.

“UP” button: When the spot massage mode is on, press and hold this button to have the back massage nodes move up.

“DOWN” button: When the spot massage mode is on, press and hold this button to have the back massage nodes move down.

Note: When “UP/DOWN” button is released, the massage nodes will massage the specified area. When the massage nodes reach the lowest/highest point, this button will flash to indicate that the lowest/highest point has been reached. They will work only when the spot massage mode is on.

“Speed” button: Press this button to select and cycle through L/M/H massage speed (corresponding indicator light will turn on). When powering off the back massage function, the speed indicator light will turn off. This button will work only when the Back Massage or Spot Massage is on.

Seat Vibration

“Seat Vibration” button: Press this button to turn on/off the seat vibration function and switch among L/M/H. The default vibration intensity is L.

Timer Function

“TIME” button: Press this button to select and cycle through 15 (Default) /20/25 minutes of massage (corresponding indicator light will turn on).

Heating Function

“HEAT” button: Press to activate the heating function of the neck massage and the back massage (indicator light will turn on); press again to turn it off (indicator light will turn off).

Note: Heating function will work only when the Neck or Back Massage is on.

Product Care and Maintenance

- Place the massage cushion in a safe, cool, and dry place when not in use.
- Always unplug the massager cushion prior to cleaning.
- After use, wipe the massager cushion with a soft, damp cloth using non-abrasive cleaners. If necessary, clean with a sterilizer.
- DO NOT submerge the massager in water or any other liquid.
- Keep away from all solvents and harsh detergents.
- DO NOT attempt to repair this massager. There are no user-serviceable parts.

Troubleshooting

PROBLEM	CAUSE	SOLUTION
Massager doesn't turn on.	Disconnected power or loose adapter connection causing poor contact.	Reconnect the adapter and re-plug to the power source.
	Power is off.	Press the Power button to turn on the massager.
Massager shuts off automatically.	Massager shuts-off when time is up.	Press the Power button to turn on the massager.
	The product is overloaded, activating overload protection.	Restart the massager after alleviating the load.
	The motor auto shuts-off due to the device being on too long, activating overheating protection.	Use the massager after the motor cools off (approximately 30 minutes).
The height is not ideal.	Neck massage nodes have not been adjusted to a proper height.	Adjust the height of the neck massage nodes following according to the "Operation Guide".

PROBLEM	CAUSE	SOLUTION
Pain is felt during the use.	Direct skin contact with the product directly.	Avoid direct contact with the skin. Do not lean on the massager too closely.
Heating function does not work.	1. Heating function hasn't been activated. 2. The heat function is manufactured according to International safety standards, which limits how fast the product can heat up. It needs a few minutes to warm up after activating heating function.	1. Press the Heating button to activate this function (indicator light will turn on). 2. Please wait 2-3 minutes after activating this function.
Unsatisfied experience of using this product in the car.	This product is not designed for car seat.	Don't use it in the car.

Specifications

Product Name: RENPHO Massage Cushion

Input: 12V  5A

Rated Power: 60W

Dimension (Folded): 30.7in (L) x 16.9in (W) x 7.5in (H)
78.0cm (L) x 43.0cm (W) x 19.0cm (H)

Net Weight: 11.24 lbs (5.1 Kg)

Auto-off Timer: 15/20/25 Minutes

Symbol	Definition
	Direct Current

Warranty Policy

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

Customer Service

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

 **Tel: +1(844) 417 0149 (US&CA)**

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+44(0)7434-666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM(GMT)

 **Email:** support@renpho.com (US&CA)

support-au@renpho.com (AU)

support-eu@renpho.com(UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1: This device may not cause interference.
- 2: This device must accept any interference, including interference that may cause undesired operation of the device.

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1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2014/35/EU, Directive 2011/65/EU, The Electrical Equipment (Safety) Regulations 2016, The Electromagnetic Compatibility Regulations 2016, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:

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