
RENPHO

User Manual

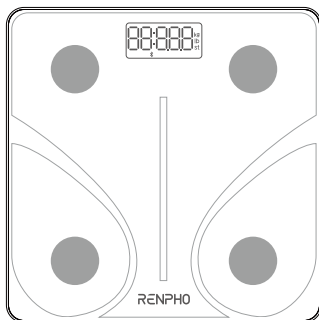Elis 1L

Model: R-A019

Contents

English..... 01-11

Deutsch..... 12-22

Français..... 23-33

Italiano..... 34-44

Español..... 45-55

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. The scale is for household use only and not for medical or commercial purposes.
4. Children under 8 or individuals who are pregnant, have osteoporosis, a pacemaker, artificial joints, portable ECG or other metal implants should use the scale for weighing purposes only and make sure to wear shoes to prevent any electric current when stepping onto the scale.
5. All body composition data provided is for reference and should not be used as a substitute for advice from a licensed healthcare professional.
6. Do not overload the scale. (Maximum Capacity: 500 lbs / 220 kg)
7. Inspect the scale for any damage before each use. Do not use a damaged scale to avoid inaccuracy or injury.
8. In case of any malfunction, check the battery power first. If necessary, replace the batteries.
9. To ensure accuracy, use the scale on a hard, flat and dry surface.
10. To avoid interference or danger, do not use the scale in an environment near equipment with strong electromagnetic fields.
11. To prevent injury, do not stand on the edges of the scale, jump on or stomp it.
12. The scale is a delicate and precise measuring device. Handle with care to avoid breakage.
13. The scale is not waterproof. Keep it dry and never submerge it in water. To clean the surface, use a damp cloth or glass cleaner. Do not use soap or chemical cleaning agents.
14. Do not expose the scale to direct sunlight or use it in high temperatures or extremely humid areas. Store the scale in a cool, dry place.
15. Do not attempt to disassemble the scale for any reason, as this may cause damage or void the warranty.

About the battery

1. Do not expose the appliance or battery to excessive temperatures.
2. Replace all batteries of a set at the same time.
3. Different types of batteries or new and used batteries are not to be mixed.
4. Exhausted batteries are to be removed from the appliance and safely disposed of.
5. If the appliance is to be stored unused for a long period, the batteries should be removed.
6. Do not use modified or damaged batteries.
7. Be aware of the risk of terminals of the battery-operated appliance or battery being short-circuited by metal objects.
8. This appliance contains non-rechargeable batteries, these batteries are not to be recharged.
9. Always insert batteries correctly according to the polarity (+ and -) marked on the battery and scale.
10. Do not heat, short-circuit, force discharge, encapsulate and/or modify the batteries. Do not deform or dismantle the batteries.
11. Keep batteries out of the reach of children.
12. Do not charge the batteries.
13. Do not weld or solder directly to batteries.
14. Do not dispose of batteries in fire.
15. Do not allow children to replace batteries without adult supervision.
16. Store unused batteries in their original packaging away from metal objects. If already unpacked, do not mix or jumble batteries.

SAVE THESE INSTRUCTIONS

ABOUT THE PRODUCT

1. Technical Specifications

Product Name: Body Composition Scale

Model: R-A019


Capacity: 0.4-500lbs / 0.2-220kg

Division: 0.2lb / 0.05kg

Dimension: 12 x 12 x 1inch / 30 x 30 x 2.5cm

Batteries: 3 x 1.5V AAA (included)


2. Product Overview

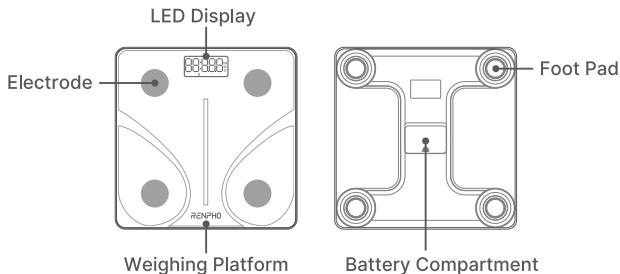BEFORE USE

1. Install the Batteries

Install the provided batteries into the scale before use. Always insert batteries correctly according to the polarity (+ and -) marked on the battery and scale.

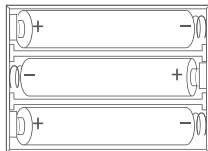Note:

- Always select the correct size and type of battery most suitable for the intended use. Do not mix new and used batteries or batteries of different types or brands. Do not charge the batteries.
- Clean the battery contacts and also those of the scale prior to battery installation.
- If "Lo" appears on the display, please replace the batteries and make sure to replace all batteries in a set at the same time.

2. Pair the Scale with the App

1) Download the App

Scan the QR code below or search for “Renpho Health” from the Apple App Store or Google Play to download the Renpho Health App.

Renpho Health

Note: Due to continuous updates and improvements, the "Renpho Health" App may appear slightly different.

2) Log in/Sign up

- Log in directly with an existing RENPHO App account.
- If you don't have any RENPHO App account, tap **“Sign up”**.
Register a new account for Renpho Health App using your email address and complete your profile information.

3) Add Device

Before adding the device, enable Bluetooth on your smartphone.

For iOS (iOS 12.0 or above), Bluetooth permission is required.
For Android system (Android 7.0 or above), location is also required.

-Method 1


- Step onto the scale with bare and dry feet on the scale's conductive areas to activate the scale.
- Open the App, and a window will pop up and tap **“Connect”**.

-Method 2:

- Tap "+" on the top right corner of the "Device" page to pair manually.
- Tap "Body Fat Scale" and follow the in-the-app instruction.
- Tap "+".

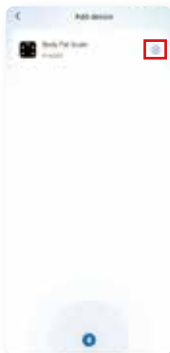
3. Calibrate the Scale

Calibrate the scale before each use to ensure accuracy.

- Place the scale on a flat, solid and even surface. Avoid any carpeting or uneven flooring.
- Step on the scale to activate it and then step off.
- The scale will show "CAL" and then "0.0" or "0.00" to indicate the calibrating is complete.



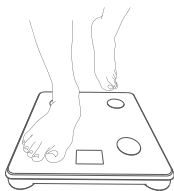
CAL → 00
or
CAL → 000

USING THE PRODUCT

Measure

- Ensure the scale is placed on a flat, solid and even surface.
- After calibrating the scale, step on the scale's conductive areas with bare and dry feet.
- Stand and keep still until the measurement flashes twice and then the display shows "0000" to indicate the final measurement. If the scale is connected to the Renpho Health App, you can get additional body composition indexes (e.g. BMI, body fat, etc.) on the App once the measuring is complete.

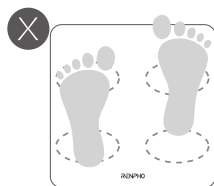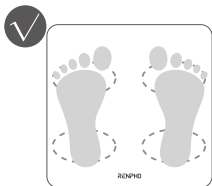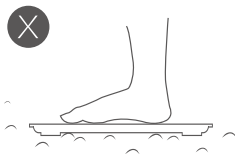Tips:

- It is recommended to use the scale at the same time and under the same condition each time to ensure accuracy.
- To avoid falling off the scale, never stand on the edges of the scale.

FREQUENTLY ASKED QUESTIONS

1. Why isn't the scale detecting body fat and/or other data?

- a. Ensure your feet are dry and bare.
- b. Stand and keep still for 5-10s on the scale's conductive areas to get the measurement, facing the LED display.
 - If still not, try applying some moisturizer to your feet and measure again. Alternatively, ask another person to repeat the 2 steps above.

2. I have registered an account on the RENPHO website, why can't I log into the Renpho Health App?

Because the RENPHO website and Renpho Health App are two different systems, you must create an account for the Renpho Health App. If the Bluetooth connection fails, please try the following steps:

1) Reset the scale.

Remove one battery from the scale, wait for 5s and reinstall it.

2) Reset the Bluetooth connection.

Ensure Bluetooth and GPS (Android only) are enabled on your phone.

- a. Open the "Renpho Health" App.
- b. Step one foot on the scale to activate it.
- c. After a few seconds of loading, a window will pop up on the App and tap "Connect".

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty.

If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

TEL: +44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

+1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+49 40 23969285 (DE ONLY)

Monday-Friday 9:00AM-6:00PM (CET)

Email: support-eu@renpho.com (UK&EU)

support@renpho.com (US&CA)

support-au@renpho.com (AU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

1: This device may not cause interference.

2: This device must accept any interference, including interference that may cause undesired operation of the device.

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1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

3. Personal data

All end-users of WEEE are reminded that you are responsible for deleting personal data from the WEEE to be disposed of.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/53/EU, Directive 2011/65/EU, Radio Equipment Regulations 2017, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:

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