RENPHO

USER MANUAL

Shiatsu Foot & Calf Massager

Model: R-D001R

CONTENTS

• ENGLISH	01 ~ 13
-----------	---------

• DEUTSCH	14 ~ 27
-----------	---------

• FRANÇAIS	28 ~ 41
------------	---------

• ITALIANO	42 ~ 54
------------	---------

• ESPAÑOL	55 ~ 67
-----------	---------

Thank you for choosing the RENPHO Foot and Calf Massager. Please use this product with the utmost care and safety in mind. To ensure safety and avoid any accidental damage or injury, please read this manual carefully and use the product only as indicated. Non-compliance may cause serious injuries, death and/or severe property damage.

IMPORTANT INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. **WARNING** – To reduce the risk of burns, fire, electric shock, or injury to persons.
 - a) This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
 - b) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
 - c) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
 - d) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 - e) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
 - f) Do not carry this appliance by supply cord or use cord as a handle.
 - g) Keep the cord away from heated surfaces.

- h) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
 - i) Never drop or insert any object into any opening.
 - j) Do not use outdoors.
 - k) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
 - l) To disconnect, turn all controls to the off position, then remove plug from outlet.
4. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.
 5. WARNING: Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
 6. WARNING: Do not stand on or in appliance. Use only while seated.
 7. DANGER-To reduce the risk of electric shock: The appliance should never be left unattended when plugged in. Unplug from the outlet when not in use, and before cleaning, putting on or taking off parts.
 8. WARNING: Do not carry this appliance by supply cord or use cord as a handle.
 9. WARNING: Keep the cord away from heated surfaces.
 10. DANGER: Do not reach for an appliance that has fallen into water. Unplug immediately.
 11. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
 12. The switch or controlled unit must not be allowed to get wet during washing and that during drying the cord must be positioned to ensure that water does not flow into the switch or control.
 13. SERVICING OF DOUBLE-INSULATED PRODUCTS. In a double-insulated appliance, two systems of insulation are provided instead of grounding. No grounding means is provided on a double-insulated product, nor should a means for grounding be added to the appliance. Servicing a double-insulated product

requires extreme care and knowledge of the system, and should be done only by qualified service personnel. Replacement parts for a double-insulated appliance must be identical to the parts they replace.

14. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.

ABOUT THE BATTERY

1. Remove and immediately recycle or dispose of used batteries according to local regulations and keep away from children.
2. Do NOT dispose of batteries in household trash or incinerate.
3. Even used batteries may cause severe injury or death.
4. Call a local poison control center for treatment information.
5. Clean the battery contacts and also those of the equipment prior to battery installation.
6. Non-rechargeable batteries are not to be recharged.
7. Do not force discharge, recharge, disassemble, heat above (manufacturer's specified temperature rating) or incinerate. Doing so may result in injury due to venting, leakage or explosion resulting in chemical burns.
8. Ensure the batteries are installed correctly according to polarity (+ and -).
9. Do not mix old and new batteries, different brands or types of batteries, such as alkaline, carbon-zinc, or rechargeable batteries.
10. Remove and immediately recycle or dispose of batteries from equipment not used for an extended period of time according to local regulations.
11. Always completely secure the battery compartment. If the battery compartment does not close securely, stop using the product, remove the batteries, and keep them away from children.
12. Do not expose the appliance or battery to excessive temperatures.

13. Always select the correct size and type of battery most suitable for the intended use. Information provided with the equipment to assist correct battery selection should be retained for reference.
14. In case of ingestion of a cell or battery, seek medical assistance promptly. Swallowing lithium coin cells or batteries can cause chemical burns, perforation of soft tissue, and in severe cases can cause death. They must be removed immediately if swallowed.
15. Do not allow children to replace batteries without adult supervision.
16. Do not force discharge batteries. Do not deform batteries. Do not dispose of batteries in fire.
17. Exhausted batteries should be immediately removed from equipment and properly disposed of. When discharged batteries are kept in the equipment for a long time, electrolyte leakage can occur causing damage to the equipment and/or personal injury.
18. A lithium battery with a damaged container should not be exposed to water.
19. Store unused batteries in their original packaging away from metal objects. If already unpacked, do not mix or jumble batteries.
20. Do not use modified or damaged batteries.
21. Be aware of the risk of terminals of the battery-operated appliance or battery being short-circuited by metal objects.

WARNING

1. Do Not Crush – Avoid sharp folds.
2. Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.

SAVE THESE INSTRUCTIONS

ABOUT THE FOOT AND CALF MASSAGER

1. Specifications

Model: R-D001R

Input: AC 100-240V 50/60Hz

Rated Power: 60W

Auto-off Timer: 15min/20min/30min

Product Size: L584*W395*H246mm/23*15.6*9.7 inches

Product Weight: 7.5kg/16.5lb

2. Package Included

1 x Foot and Calf Massager

1 x Remote Control

3. Overview

Control Panel

Detachable foot cover with zipper

Magnetic area designed to place the remote control

Remote Control

4. Button Introduction

Button Name	Function
Power	Power ON/OFF
Air Mode	Full feet(Default), Forefeet, Rear feet, OFF
Air Intensity	1 (Default) / 2 / 3 / 4
Knead Intensity	L(Default) / M / H / OFF
Time (Auto-off)	15min (Default) / 20min / 30min
Heat	Instep+Bottom (Default), Instep, Bottom, OFF

Tips: Use your thumb to press buttons to get better touch sensitivity.

USING THE PRODUCT

Massage Postures

When massaging calves, it is recommended to turn on only the Air pressure or Kneading function, and start from the lowest intensity.

Massage Feet

Massage Calves

Massage Arms

We recommend you limit massaging to 15-30 minutes. Longer massage time may cause mild pain.


Operation

1. Connect the plug to the power source. The Power indicator flashes White accompanied by a beep sound and the massager enters standby mode.

2. Press the  Button to turn on the massager. Check the default Mode and Intensity in **Button Introduction** chart.
3. Press the  Button to switch the auto-off time among 15, 20, or 30 minutes. The timer will re-time each time you press the Time Button.
4. Press  to turn ON/OFF the Mode function or select Full foot, Fore foot or Rear foot area.
5. Press  Button to adjust the Air Intensity among 1/2/3/4.
6. Press the  Button to adjust the Knead Intensity among L/M/H/OFF.
7. Press the  Button to turn ON/OFF the heating function or select Full (Bottom+Instep), Instep or Bottom heating area.

Note:


- A: Heating, Kneading and Air pressure function can work independently.
 - B: Each press is accompanied by a beep sound and corresponding indicator light (When off, the LED goes out).
 - C: If none of the Heating, Kneading and Air pressure function is working, the Time and Power button indicator lights are on, and the massager will enter the standby state if there is no further operation within 10 seconds.
8. The foot massager fits up to US Men's size 14 (EU Men's size 48.5).
 9. Removable mesh: Unzip to remove and clean

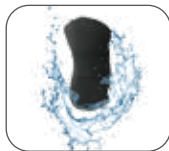REPLACE THE BATTERY

1. For AU Only


- Use a screwdriver to unscrew the fixing screw on the battery compartment.
 - Push the locking tab on the bottom edge of the remote control and hold it.
 - Pull the battery door out and install a CR2025 battery.
 - Push the battery door back in and tighten the fixing screw.
- Tip: Keep the "+" symbol on the battery facing upwards when placing the battery in the battery door.

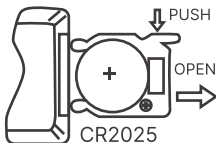2. For Not AU

- Push the locking tab on the bottom edge of the remote control and hold it.
- Pull the battery door out and install a CR2025 battery.
- Push the battery door back in and tighten the fixing screw.

Tip: Keep the "+" symbol on the battery facing upwards when placing the battery in the battery door.

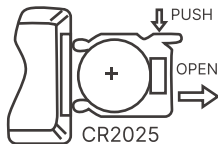PRODUCT CARE AND MAINTENANCE

1. The use of non-rated voltage is prohibited.
2. It is prohibited for non-professional personnel to repair the massager or replace parts.
3. Do not use or store the appliance in wet environments or high temperatures.
4. Unplug the power cord from the outlet after using.
5. Ensure that the operation of this product does not exceed more than 30 minutes at any given time.
6. If you feel unwell or in any unusual discomfort, please stop using the product immediately.
7. If a high temperature or foul odor is experienced while using the product, stop using immediately and unplug the device.
8. Any Individual who may be pregnant or recently gave birth, has a heart condition or a pacemaker, suffers from diabetes or other sensory impairment, phlebitis and/or thrombosis, joint dysfunction, inflammation or injuries, weak bones, osteoporosis, a high fever, an abnormal or curved spine, is at an increased risk of developing blood clots, has pins/screws / artificial joints or other medical devices implanted in his / her body, is being treated for a medical condition or physical injury or is under a doctor's care should consult with a medical professional before using a massaging device designed for home use.

9. Do not use sharp objects to handle the product in order to avoid damage.
10. Please keep the product clean. You can use soft cotton cloth and neutral detergent to clean it.
11. The product is a temperature-protection device; it will automatically shut off when the temperature exceeds the safe limit.

TROUBLESHOOTING GUIDE

Common problem	Possible Cause	Solution
The product does not work.	1. Disconnected power or loose plug causing poor contact. 2. The power switch is off or time is up. 3. The motor auto shuts off due to the device being on too long activating overheating protection.	1. Re-plug to the power source. 2. Turn the power switch on again. 3. Use the massager after the motor cools off naturally.
The massage intensity is too high or too low.	1. Your feet is positioned at the wrong angle. 2. Didn't choose intensity level which is right for you.	1. Adjust the contact angle to get moderate intensity. 2. The massager has 4 intensity levels, choose the intensity you prefer. We recommend low intensity to begin with.

Common problem	Possible Cause	Solution
No heating or insufficient heating.	1. The heating function isn't turned on. 2. The working heat area might be bottom, not full/Instep. 3. The product structure and working principle limits how fast the bottom can heat up. Therefore, the heating-up may seem a bit slow, wait a few minutes longer. 4. The ambient temperature or the internal temperature of the product is too high, the heating function will be turned off because of the over-heated protection.	1. Press the "Heat" button. The white indicator light should turn on. 2. Choose the Full Heat area to get more sufficient heating. 3. It takes about a few minutes to fully heat up. Please be patient as this is to avoid the temperature rising too fast. 4. After the internal temperature drops, the heating function will be auto-turned on.
The touch panel is not working properly.	1. Fingers are wet. 2. Insufficient touch area.	1. Dry your fingers before touching the button. 2. Use your thumb to press buttons to get better touch sensitivity.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

TEL: +1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

Email: support@renpho.com (US&CA)

support-au@renpho.com (AU)

support-eu@renpho.com (UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a retail store or other collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2014/35/EU, Directive 2011/65/EU, The Electrical Equipment (Safety) Regulations 2016, The Electromagnetic Compatibility Regulations 2016, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address: <https://renpho.eu/pages/compliance-data>

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FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions: (1) This device may not cause interference. (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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