RENPHO

USER MANUAL

EYE SPA PODS

Model: R-G009

CONTENTS

• ENGLISH	01-09
-----------	-------

• DEUTSCH	10-17
-----------	-------

• FRANÇAIS	18-26
------------	-------

• ITALIANO	27-35
------------	-------

• ESPAÑOL	36-44
-----------	-------

Thank you for purchasing the RENPHO Eye Spa Pods! Please read this instruction manual carefully and use the product only as indicated. RENPHO accepts no responsibility for any cause of injuries due to misuse or negligence to follow the instructions outlined. Please keep this manual for future reference.

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.

1. This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.
4. Cleaning and maintenance should not be conducted by children.
5. Please consult your doctor before use if you fit into any of the following user groups:
 - Persons who had undergone eye surgery or eye conditions such as glaucoma, cataract and retinal detachment
 - The elderly
 - Pregnant or postpartum women
 - Persons who are allergic to leather material
 - Persons with skin infections, sensitivity, redness and/or burns on the skin
 - Persons with any swelling or wound around the eye area
 - Persons using a pacemaker or other implantable medical devices
 - Persons with a heart condition, malignant tumors, or other serious health conditions
 - Persons undergoing medical treatment
 - Persons with high or abnormal blood pressure
6. If you experience any discomfort during use, stop using the device immediately.
7. Remove contact lenses, glasses, false eyelashes, hairpins,

and any other accessories before use.

8. Do not use the device while driving, walking, or running.
9. Avoid contact with sharp objects.
10. High temperatures may cause burns to occur regardless of the control setting. Do not use the device on infants, persons who are with disabilities, sleeping or unconscious. Check frequently where the skin is in contact with the heated area to reduce the risk of blistering.

DANGER - To reduce the risk of electric shock:

1. Never use pins or other metallic fasteners with this device.
2. Do not use the device with wet hands or submerge the device in water or any liquid.

CAUTION - To reduce the risk of burns:

1. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
2. The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

1. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
2. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
3. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this device if it is not working properly, has been damaged, or dropped into water.
5. Never drop or insert any object into any opening.
6. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

About the battery:

1. Do not use while bathing or in a shower.
2. Do not place or store appliance where it can fall or be pulled into a tub or sink.
3. Do not place in or drop into water or other liquid.

4. Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
5. Do not expose the appliance or battery to excessive temperatures.
6. Be aware of the risk of terminals of the battery-operated appliance or battery being short-circuited by metal objects.
7. This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.
8. Do not dismantle, open or shred battery.
9. Do not subject cells or batteries to mechanical shock.
10. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
11. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
12. Do not leave a battery on prolonged charge when not in use.
13. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
14. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

SAVE THESE INSTRUCTIONS

ABOUT THE PRODUCT

1. What's in the Box

1 x Eye Spa Pods

1 x USB-C Charging Cable

2. Specifications

Product Name: Eye Spa Pods


Model: R-G009

Input: 5V $\overline{=}$ 1A

Rated Power: 5W

2 x Battery Capacity: 600 mAh
2.22Wh


3. Overview

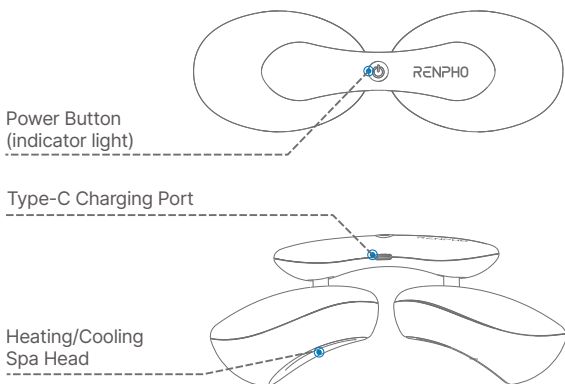USING THE PRODUCT

Remove any eye accessories (e.g., eyeglasses, contact lenses, false eyelashes, etc.) before using.

- Slide and rotate the spa heads to fit your eyes or put them a little lower under your eyes area.
- Lie down and turn on the device. Keep your eyes closed while using. You can also choose to hold the device and place it over your eyes for use.

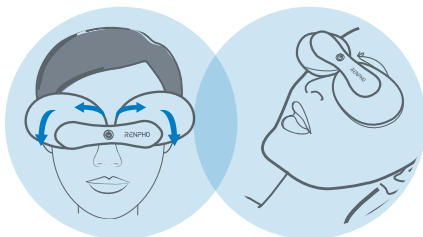Tip:

- It is not recommended to use the device for more than 20 minutes at a time.
- Your vision may become blurry for a while, as expected after using the device. It is advised to close your eyes for a few minutes until your vision goes back to normal.

Power button:

Press and hold for 2s to turn on/off the device.
Short press to cycle through 3 spa modes
(default: cooling mode).

- Each press is accompanied by a beep sound.

Spa Mode	Indicator Light	Temperature	Timer	Support use upon a full charge
Cooling	Blue	20°C / 68°F	3mins	about 6 sessions
Heating	Red	39°C / 102°F	10mins	about 6 sessions
Cooling & Heating Cycle	Purple	20-39°C / 68-102°F	8mins	about 4 sessions

Note:

- As the device needs to dissipate the heat to support the cooling function, it is not recommended to operate the cooling mode in succession. It is recommended to operate this mode every 30mins.

CHARGING

- a. Charge the device using the provided charging cable and a DC 5V 1A mobile charger (not included).
- b. It takes about 2.5h to fully charge the device.

Status	Indicator Light
Low battery / Charging	Flashing red
Fully charged	Solid red
Overcharging (at the voltage of 6-30V)	Flashes red quickly for 1 minute

Note:

- Do not turn on the device while charging.

CLEANING AND MAINTENANCE

- Turn off the device before cleaning. Do not rinse the device. Use a soft cotton cloth dampened with mild detergent to wipe it clean.
- Keep away from all corrosive cleaners and harsh detergents. Keep away from humidity, heat and/or corrosive gas.
- If the device has not been used for over 3 months, please recharge it before use.
- Do not attempt to disassemble the device by yourself.

TROUBLESHOOTING GUIDE

Problem	Cause	Solution
The device is not charging	Using improper charger	Please use the DC 5V 1A charger
	Loose plug causing poor contact	Check the connection and re-insert into the device
	The charger or charging cable is damaged	Check and replace the charger or charging cable if damaged
Power-on failure	Didn't turn on the device as instructed	Press and hold the power button for 2s to turn on the device
	Built-in over-heat protection is being activated after use	It is recommended to rest 30mins before resuming use

WARRANTY POLICY

Your RENPHO product purchase is covered by a one year limited manufacturer warranty from the date of delivery. For warranty terms and conditions, please visit:

<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not void the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

Tel: +1(844) 417 0149 (US&CA)
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)
+44(0)7434-666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM(GMT)

Email: support@renpho.com (US&CA)
support-au@renpho.com (AU)
support-eu@renpho.com(UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dis-pose of any product parts as they may be required for inspection/repair.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation. Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1: This device may not cause interference.
- 2: This device must accept any interference, including interference that may cause undesired operation of the device.

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1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2011/65/EU, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:

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