RENPHO

Heating Pad

Instant Soothing Warmth

User Manual

CONTENTS

ENGLISH	01~07
---------	-------

- - - - -

DEUTSCH	08~15
---------	-------

- - - - -

FRANÇAIS	16~23
----------	-------

- - - - -

ITALIANO	24~31
----------	-------

- - - - -

ESPAÑOL	32~39
---------	-------

IMPORTANT:

Read this manual carefully before first use and keep it in a convenient place for future reference.

	 DANGER
	Hot Surfaces.

 DANGER		
 • Do not sit on the pad.	 • Do not use while sleeping. • Never put the pad under the body.	 • Do not fold the pad. • Avoid sharp folds, as it may damage the product.

ICONS

This guide contains important information about using, operating, and storing the product. Please read it through and note that non-observance of these instructions can result in serious injury or damage to the product.

WARNING

These notes must be observed to prevent any damage to the product.

CAUTION

These notes must be observed to prevent any injury to the user.

NOTE

These notes provide useful additional information on the installation or operation.

CAUTION

To reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions.


IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING THIS PRODUCT.

1. This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the product.
3. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
4. The product must not be located immediately below a socket-outlet.
5. When storing the product, allow it to cool down before folding.
6. Do not crease the product by placing items on top of it during storage.
7. Examine the product frequently for signs of wear or damage.
8. The product is for home use only and is not intended for any medical, diagnostic, or therapeutic purposes. Its design and function are solely focused on daily comfort and relaxation, rather than any medical use.
9. Do not use the product if wet.
10. This product must not be used by persons insensitive to heat and other very vulnerable persons who are unable to react to overheating.
11. The product is not to be used by young children over the age of three unless the controls have been pre-set by a parent or guardian, or the child has been adequately instructed on how to operate the controls safely.
12. Do not use if folded or creased.
13. Do not use this product with an infant, an immobile person, or a person with poor blood circulation.
14. Do not tuck in the wired area of this product.
15. Do not use dry-cleaning fluid on this pad, cleaning solvents may have a deteriorating effect on the insulation of the heating element.
16. Do not route electrical cords on the pad. Do not allow the cords to be pinched.
17. When in use, keep the control away from drafty areas that may become damp or wet such as an open window.
18. Do not poke the product with pins - They may damage the electric wiring.
19. Be sure that this product is used on an alternating-current supply circuit of the proper voltage.
20. Turn off or unplug when not in use.
21. Excessive heating may result from folding pad when in actual use.
22. If improper operation of this product is observed, discontinue its use immediately and consult RENPHO Customer Service regarding its repair.
23. Do not tuck the product or supply cord in a sofa or other furniture. All portions of the product should be visible at all times.
24. Wrapping of the cord around the controller may damage the cord. Loop loosely when storing.

DANGER - to reduce the risk of burns, electric shock and fire, this product must be used in accordance with the following instructions.

1. Do not destroy.
2. Do not use while sleeping.
3. Do not use on an infant.
4. This product is not to be used on or by an invalid. Sleeping or unconscious person, or a person with poor blood circulation unless carefully attended.
5. Do not use on areas of insensitive skin.
6. Burns can occur regardless of control setting, check skin under the product frequently.
7. Do not use in an oxygen atmosphere (applies to hospital pads).
8. Do not use pins or other metallic means to fasten this product in place.
9. Do not sit on or crush the product - avoid sharp folds.
10. Never pull this product with the supply cord.
11. Do not use the cord as a handle.
12. Carefully examine the inner cover before each use. Discard the product if the inner covering shows any sign of deterioration.

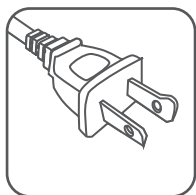
This product has a polarized plug (one blade is wider than the other) as a safety feature. This plug can only be inserted into a polarized outlet in one way. If the plug does not fit fully into the outlet, reverse the plug. If it still does not fit, contact a qualified electrician to install a proper outlet. Do not attempt to modify the plug in any way.

SAVE THESE INSTRUCTIONS

WHAT'S IN THE BOX

- 1 x Heating Pad
- 1 x Controller with Power Cord

ABOUT THE PRODUCT


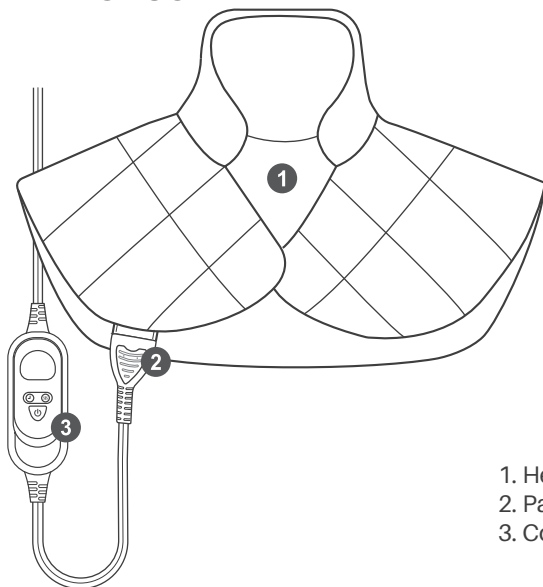
1. Heating Pad
2. Pad Connector
3. Controller with Power Cord

USING THE PRODUCT

1. Before Use

Reminder: Check each part before use and ensure there are no damages.
Push the controller plug into the pad connector and ensure the inner prong is locked in place. (To remove the controller, lift the locking tab with your finger and pull the controller plug out directly.)

Plug the cord into an AC 120V (US&CA) / AC 230V (EU&UK&AU) outlet.

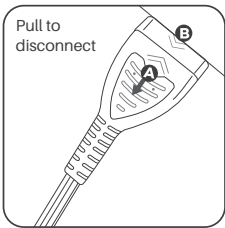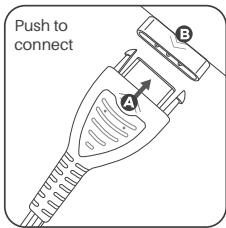A. Controller plug
B. Pad connector

Note:
If there is no power, make sure the electrical outlet is functional and all connections are secure. Do not turn on the pad until it has been fully assembled.

2. Operation

Buttons	Functions
	Press to set a timer among 2(Default) / 1.5 / 1 / 0.5 hours. The timer will reset each time it's switched.
	Press to switch the heating levels among 1 / 2 / 3 / 4 / 5 / 6.
	Press to turn ON/OFF the product.

Note:
1. Always unplug the product from the electrical outlet whennot in use!
2. It's recommended to start at the lowest heat level and gradually increase as desired.

CLEANING & MAINTENANCE

1. Before cleaning

Disconnect the power cord from the electrical outlet BEFORE detaching the controller from the pad.

Caution: DO NOT dry clean. DO NOT bleach. DO NOT use wringer. DO NOT immerse control or cord in water or other liquids. DO NOT use mothballs or sprays of any kind. DO NOT dry with a hair dryer or heater.

2. Cleaning

Machine wash using a neutral detergent or mild soap in cold water on a gentle/delicate cycle or hand wash. Do not wash with hot water as the pad may shrink.

Tip: It is recommended to wash once every 1-2 months.

3. Air drying

Drape the pad over a clothesline or hang over a shower rod. Do not use clothespins. Do not tumble dry.

Caution: The pad must be COMPLETELY dry before reconnecting the controller plug or power plug.

4. Storage

Fold the pad neatly and gently but not tightly without any sharp bends in the heating elements, and store it in a cool, dry place where there is no any heavy object placed upon it.

TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
The product does not work.	The power plug is not connected properly to the electrical outlet.	Re-plug into the outlet.
	The pad and the controller are not connected properly.	Follow the “Before Use” instructions.
The Controller Displays F1.	PTC disconnection error.	Unplug the pad connector and re-start the product.

The Controller Displays F2.	NTC short circuit error.	Contact customer service.
The Controller Displays F3.	NTC disconnection error.	Contact customer service.
The Controller Displays F4 / F5.	Optocoupler controllable silicon short circuit error.	Contact customer service.
The Controller Displays F6.	PCBA temperature detection system short circuit or open circuit protection.	Contact customer service.
The Controller Displays F7.	Mains power abnormality protection.	When the voltage returns to the normal range of 180V - 275V, the product can be restarted normally.

Tip:

Please do not attempt to disassemble or repair the product on your own to avoid any potential danger. If you are unable to resolve the above fault symptoms, please contact our after-sales customer service.

SPECIFICATIONS

Product Name: Heating Pad

Model: R-HAP07

Input: AC 220-240V~ 50Hz (EU&UK&AU)

AC 120V 60Hz (US&CA)

Rated Power: 90W (UK&EU&AU)

130W (US&CA)

Power Cord Length: 8.9 ft / 2.7 m

Dimension: 19 × 24 in / 48.3 × 61 cm

Temperature Range: 104 °F - 140°F (40 °C - 60 °C)

WARRANTY (1-YEAR PERIOD)

The manufacturer's warranty is void under the following circumstances but not limited to:

- Improper or inadequate maintenance or modification.
- Accident, misuse, abuse, contamination, or other external causes.
- The use of unspecified adapter and accessories.
- Loss or damage in transit.
- Damages that occur as a result of one's failure to follow the instructions.
- This warranty does not apply to expendable or consumable parts and does not extend to any product from which the serial number has been removed.

CUSTOMER SERVICE

We have an award-winning warranty, exchange, and customer service team that guarantees hassle-free solutions to any issue you might have within the 24 business hours.

TEL: +1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

Email: support@renpho.com (US&CA)

support-au@renpho.com (AU)

support-eu@renpho.com (UK&EU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2014/35/EU, Directive 2011/65/EU, The Electrical Equipment (Safety) Regulations 2016, The Electromagnetic Compatibility Regulations 2016, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:

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FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory compliance

This device complies with Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions: (1) This device may not cause interference, and (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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