
# RENPHO
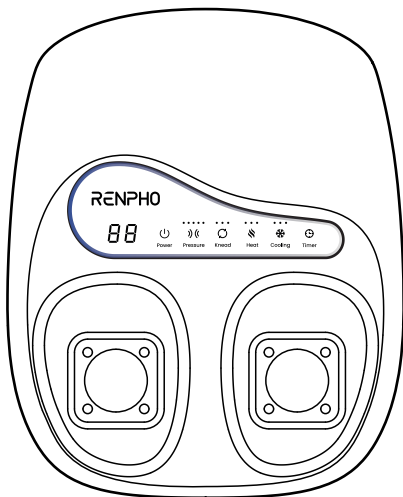
## Solheel Thermacool

Foot Massager

## User Manual

Model: R-WMF08

# Contents

---

-----  
· ENGLISH 01 ~ 13

-----  
· DEUTSCH 14 ~ 25

-----  
· FRANÇAIS 26 ~ 37

-----  
· ITALIANO 38 ~ 50

-----  
· ESPAÑOL 51 ~ 63  
-----

# IMPORTANT SAFETY INSTRUCTIONS

## READ ALL INSTRUCTIONS BEFORE USING.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. **WARNING** – To reduce the risk of burns, fire, electric shock, or injury to persons.
  - a) This appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
  - b) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
  - c) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
  - d) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
  - e) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
  - f) Do not carry this appliance by supply cord or use cord as a handle.
  - g) Keep the cord away from heated surfaces.
  - h) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
  - i) Never drop or insert any object into any opening.
  - j) Do not use outdoors.
  - k) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.

- l) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 4. DANGER: Do not use while bathing or in a shower.
- 5. DANGER: Do not place or store appliance where it can fall or be pulled into a tub or sink.
- 6. DANGER: Do not place in or drop into water or other liquid.
- 7. WARNING: Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
- 8. The appliance is only to be used with the power supply unit provided with the appliance.
- 9. The appliance has a heated/cooling surface. Persons insensitive to heat/cooling must be careful when using the appliance.
- 10. WARNING: Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
- 11. WARNING: Do not stand on or in appliance. Use only while seated.
- 12. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.

## About the Battery

1. Remove and immediately recycle or dispose of used batteries according to local regulations and keep away from children.
2. Do NOT dispose of batteries in household trash or incinerate.
3. Even used batteries may cause severe injury or death.
4. Call a local poison control center for treatment information.
5. Clean the battery contacts and also those of the equipment prior to battery installation.
6. Non-rechargeable batteries are not to be recharged.
7. Do not force discharge, recharge, disassemble, heat above (manufacturer's specified temperature rating) or incinerate. Doing so may result in injury due to venting, leakage or explosion resulting in chemical burns.
8. Ensure the batteries are installed correctly according to polarity (+ and -).
9. Do not mix old and new batteries, different brands or types of batteries, such as alkaline, carbon-zinc, or rechargeable batteries.
10. Remove and immediately recycle or dispose of batteries from equipment not used for an extended period of time according to local regulations.
11. Always completely secure the battery compartment. If the battery compartment does not close securely, stop using the product, remove the batteries, and keep them away from children.
12. Do not expose the appliance or battery to excessive temperatures.
13. Always select the correct size and type of battery most suitable for the intended use. Information provided with the equipment to assist correct battery selection should be retained for reference.
14. In case of ingestion of a cell or battery, seek medical assistance promptly. Swallowing lithium coin cells or batteries can cause chemical burns, perforation of soft tissue, and in severe cases can cause death. They must be removed immediately if swallowed.

15. Do not allow children to replace batteries without adult supervision.
16. Do not force discharge batteries. Do not deform batteries. Do not dispose of batteries in fire.
17. Exhausted batteries should be immediately removed from equipment and properly disposed of. When discharged batteries are kept in the equipment for a long time, electrolyte leakage can occur causing damage to the equipment and/or personal injury.
18. A lithium battery with a damaged container should not be exposed to water.
19. Store unused batteries in their original packaging away from metal objects. If already unpacked, do not mix or jumble batteries.
20. Do not use modified or damaged batteries.
21. Be aware of the risk of terminals of the battery-operated appliance or battery being short-circuited by metal objects.

# **SAVE THESE INSTRUCTIONS**

# ABOUT THE DEVICE

## 1. What's in the Box

Solheel Thermacool × 1

Remote Control × 1

Power Adapter × 1

## 2. Technical Specifications

Product Name: Foot Massager

Model: R-WMF08

Input: 12V  $\equiv$  3A

Rated Power: 36W

Rated Voltage: AC 100-240V, 50/60Hz


Auto-off Timer: 15 / 20 / 30 Min

Product Weight: 6.6 lbs / 3 kg

Product Size: 16.4 x 13.2 x 6.9 in

/ 41.7 x 33.5 x 17.6 cm

## 3. Overview

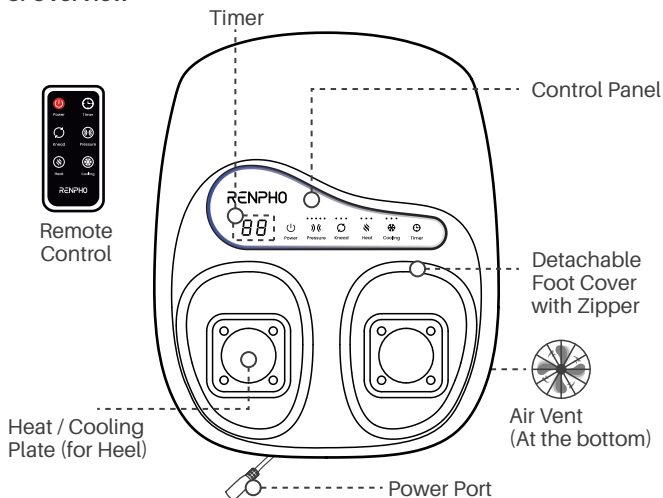
## USING THE DEVICE

- \* Do not use the device for more than 30mins each time.
  - \* Do not stand up while massaging.
  - \* To prevent blocking the air vent, use only on a flat, hard floor or on low-pile carpets and rugs.
1. Connect the device to the power adapter and plug it into a power outlet. Once the device is plugged in, the timer "— —" will flash, indicating that the device is in standby mode.
  2. The device must be used with the power adapter (12V 3A) provided.
  3. Use the Control Panel or Remote Control (CR2025 battery preinstalled) to operate the device, ensuring the remote is pointed directly at the panel. Remove the isolation strip from the remote control before the first use.

Remote Control

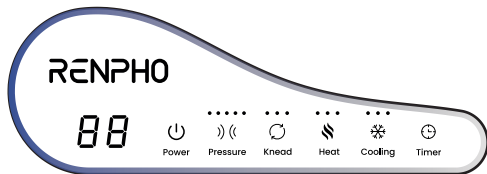Control Panel

## Button Instruction




| Button                                                                           | Operation                                                                                                                                                                                                    |
|----------------------------------------------------------------------------------|--------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
|  | Press to turn ON /OFF the device.                                                                                                                                                                            |
|  | Press to cycle through 15 (default) / 20 / 30 min.<br>* The remaining time is shown by the count down display on the left.  |
|  | Press to cycle through 3 kneading levels.                                                                                                                                                                    |
|  | Press to cycle through 5 pressure levels.                                                                                                                                                                    |
|  | Press to cycle through 3 heat levels<br>* Heating areas: Insteps and heels.                                                                                                                                  |
|  | Press to cycle through 3 cooling levels.<br>* Cooling area: Heels                                                                                                                                            |

### Note:

1. The heating and cooling functions cannot be used simultaneously.
2. Unplug the device after each use.
3. Knead, pressure, heat, and cooling functions can be used separately or in combination. If these functions are all turned off, the device will automatically enter standby mode after 30 seconds of inactivity.
4. Memory Function: The device automatically recalls the last selected settings when you press the Power Button again. However, the memory function will not work if the power is cut off.

## MAINTENANCE AND CLEANING

1. It is prohibited for non-professional personnel to repair the device or replace its parts.
2. Place and store the device in a cool and dry environment. Do not expose the device to direct sunlight. Keep away from humidity, heat, and/or corrosive gas.
3. Do not place any heavy objects on the device.
4. Keep the device clean. Unplug the device before cleaning.
5. Use a soft, damp cloth to wipe the main body for daily cleaning. To wipe out the stains caused by oil or cosmetics, use a soft, damp cloth with neutral detergent or special detergent for leather. Never immerse the device in water or other liquids.
6. After extended use, unplug the device and check the air vent. Then, clear any debris, dust and piles of carpets and rugs by using tweezers or by blowing it out with compressed air (e.g., from a hair dryer set to a cool, low setting).
7. The detachable foot covers are machine washable. To wash them, unzip and remove the covers. For your safety, do not operate the device without the covers attached. To attach, align the silicone with the protruding dots of the device.
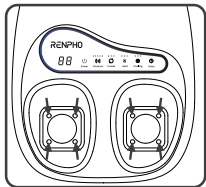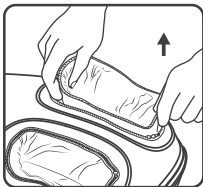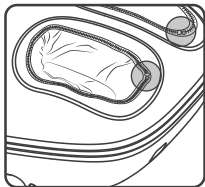
# TROUBLESHOOTING GUIDE

| Problems                                          | Causes                                                                                          | Solutions                                                                           |
|---------------------------------------------------|-------------------------------------------------------------------------------------------------|-------------------------------------------------------------------------------------|
| The device does not work.                         | The power plug is not connected properly.                                                       | Reconnect the plug to the outlet and check the power indicator.                     |
|                                                   | Not using the supplied adapter, or the interface is wrong, or the power supply is insufficient. | Use the adapter provided with the device.                                           |
|                                                   | Applying excessive pressure could cause nodes to stall.                                         | Do not tread on nodes.                                                              |
| The device shuts off suddenly.                    | The plug is loose, causing poor contact.                                                        | Re-plug the adapter and restart the device.                                         |
|                                                   | The timer runs out.                                                                             | Restart the device for another session.                                             |
| Insufficient heat / cooling or no heat / cooling. | The heat / cooling function is disabled during massage.                                         | Press Heat / Cooling Button enables the corresponding function.                     |
|                                                   | The device takes time to heat up gradually to avoid discomfort caused by rapid heat.            | Wait several minutes after turning on the device.                                   |
|                                                   | Thick socks are worn during massage, reducing sensation.                                        | Wear thinner socks to improve foot sensation and the overall effect of the massage. |

| Problems                                          | Causes                                                                 | Solutions                                                                                                                                                                                                              |
|---------------------------------------------------|------------------------------------------------------------------------|------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------------|
| Insufficient heat / cooling or no heat / cooling. | The air vent is blocked by debris, dust, or piles of carpets and rugs. | Only use on a flat, hard floor, or on low-pile carpets and rugs. Or unplug the device, then clear any blockage using tweezers or blow it out with compressed air (e.g., from a hair dryer set to a cool, low setting). |

## WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:

**<https://renpho.com/pages/warranty-terms-and-conditions>**

Note: Product registration is not required for the warranty.

If you choose not to register your product, it will not diminish the product warranty.

## CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

TEL: +1(844) 417 0149 (US&CA)

Monday-Friday 9:00AM-4:30PM

1800 MY RENPHO (1800 69 7367) (AU ONLY)

Monday-Friday 9:00AM-6:00PM (AEST)

+44(0)7434 666088 (UK ONLY)

Monday-Friday 9:00AM-5:00PM (GMT)

Email: [support@renpho.com](mailto:support@renpho.com) (US&CA)

[support-au@renpho.com](mailto:support-au@renpho.com) (AU)

[support-eu@renpho.com](mailto:support-eu@renpho.com) (UK&EU)

\*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

### 1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

### 2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a retail store or other collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2014/35/EU, Directive 2011/65/EU, The Electrical Equipment (Safety) Regulations 2016, The Electromagnetic Compatibility Regulations 2016, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:

<https://renpho.eu/pages/compliance-data>

## FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

## ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- (1) This device may not cause interference.
- (2) This device must accept any interference, including interference that may cause undesired operation of the device.

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