
RENPHO

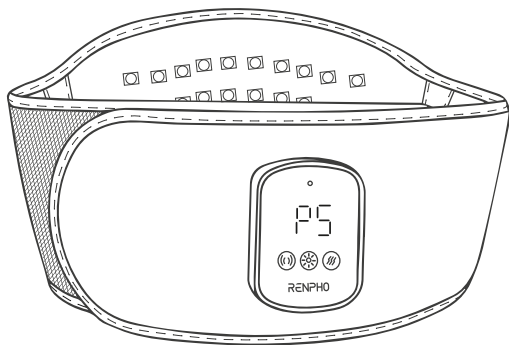**Embrace
Dorsi Red**

User manual

Model: R-WTR01

CONTENTS

● ENGLISH	01~11
-----------	-------

● DEUTSCH	12~22
-----------	-------

● FRANÇAIS	23~33
------------	-------

● ITALIANO	34~44
------------	-------

● ESPAÑOL	45~55
-----------	-------

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING.

1. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.
2. This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
3. Children should be supervised to ensure that they do not play with the device.
4. **WARNING** - To reduce the risk of burns, fire, electric shock, or injury to persons.
 - a) This device should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
 - b) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
 - c) Close supervision is necessary when this device is used by, on, or near children, invalids, or disabled persons.
 - d) Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
 - e) Never operate this device if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
 - f) Do not carry this device by supply cord or use cord as a handle.
 - g) Keep the cord away from heated surfaces.
 - h) Never operate the device with the air openings blocked. Keep the air openings free of lint, hair, and the like.

- i) Never drop or insert any object into any opening.
 - j) Do not use outdoors.
 - k) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
 - l) To disconnect, turn all controls to the off position, then remove plug from outlet.
5. **DANGER** : Do not use while bathing or in a shower.
 6. **DANGER** : Do not place or store device where it can fall or be pulled into a tub or sink.
 7. **DANGER** : Do not place in or drop into water or other liquid.
 8. **WARNING**: Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
 9. The device has a heated surface. Persons insensitive to heat must be careful when using the device.
 10. **WARNING**: Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.

About the Battery:

1. Do not expose the device or battery to excessive temperatures.
2. Be aware of the risk of terminals of the battery-operated device or battery being short-circuited by metal objects.
3. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.
4. Do not dismantle, open or shred battery.
5. Do not subject cells or batteries to mechanical shock.
6. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
7. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging

instructions.

8. Do not leave a battery on prolonged charge when not in use.
9. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
10. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

SAVE THESE INSTRUCTIONS

ABOUT THE DEVICE

1. What's in the Box

Emberace Dorsi Red x 1

Extension Strap x 1

USB-C Charging Cable x 1

2. Technical Specifications

Product Name: Red Light Waist Band

Model: R-WTR01

Input: 5V $\equiv$ 2A

Battery Voltage: DC 3.7V

Battery Capacity: 5000mAh 18.5Wh

Charging Time: 4 h


Auto-off Timer: 30 Min

Extension Strap Length: 13.8 × 3.2 in / 35 × 8 cm



Product Size: 46.1 × 8.3 in / 117 × 21 cm

Applicable Waist Circumference: 27.6 - 55 in / 70 - 140 cm

3. Overview

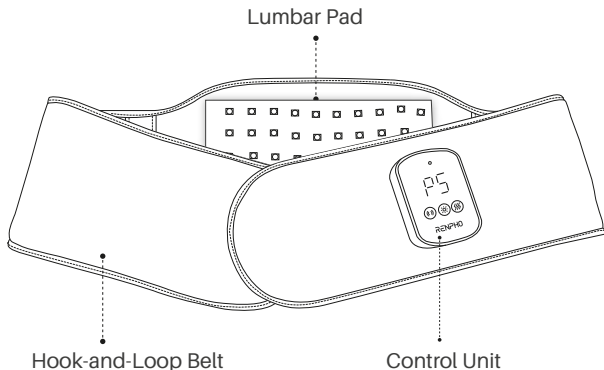Control Unit

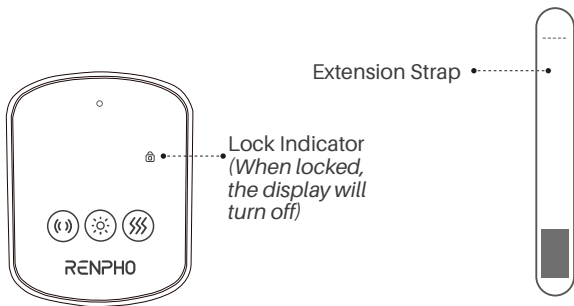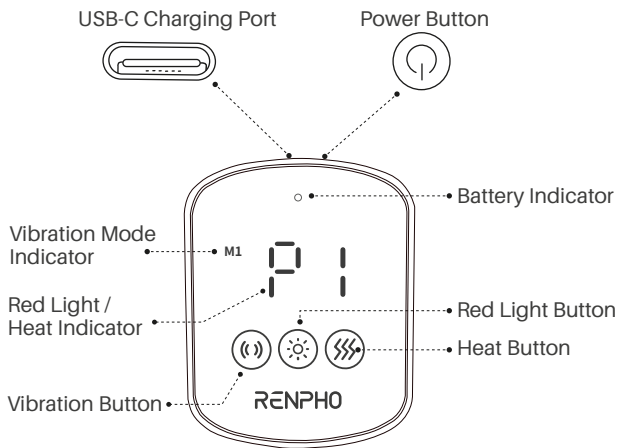
USING THE DEVICE

Before Use


- * Place the device on your skin to use.
- * Keep the control unit side facing outward during use.

1. Wear the Device

Step 1: Position the application area at your lumbar region.

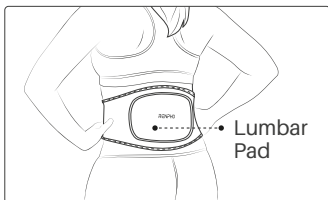
Step 2: Firmly and evenly pull both sides of the hook-and-loop belt and overlap them securely across your abdomen.

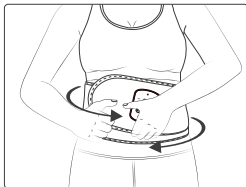2. Use the Extension Strap (For Larger Waist Circumference)

Connect the extension strap to the hook-and-loop belt, then follow the instructions in "1. Wear the Device".

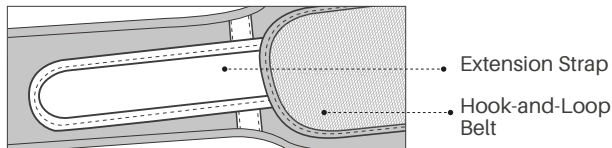Note:

- * The device will shut off automatically after 30 minutes of use.
- * **Do not use the device during sleep.**
- * Do not expose eyes directly into the red light.
- * It is not recommended to use the device for more than 30 minutes at a time.
- * Avoid the risk of skin damage by beginning with the lowest settings of the red light and heat when both are enabled simultaneously.

Button Introduction

Button	Operation
	Press and hold to turn On / Off the device. Press to lock / unlock (default) the device.
	Press to select the vibration mode: M1, M2, M3, and Off.
	Press to select the red light intensity: P1, P2, P3, P4, P5, and P0 (Off).
	Press to select the heat level: L1, L2, L3, and Off.

Battery Indicator

Battery Level/Status	Indicator
70% - 100% / Fully Charged	Solid Green
40% - 70%	Solid Yellow
5% - 40%	Solid Red
<5% / Charging	Flashing Red

CLEANING & MAINTENANCE

1. Do not submerge the device in water or wash it in a machine.
2. Use a wet towel or brush to clean the surface of the device.
3. To store the device, loosely roll it toward the straps, loop the straps over the rolled ends, and then store the device in a cool, dry, dark place after each use. Avoid sharply folding the device.
4. Avoid placing items on the device to prevent damage.

TROUBLESHOOTING GUIDE

Problem	Cause	Solution
The device stops working suddenly.	Batteries are weak or depleted.	Charge the device and restart it.
	30-min auto-off timer is activated.	Restart the device if needed.
No heating or insufficient heating.	No temperature / low heat is selected.	Adjust the heat level.
	The device needs time to warm up.	Allow up to 5 minutes to heat up.
	The device is overheating.	Built-in overheat protection. Allow the device to cool down and restart it.
Fail to turn on the device.	Batteries are weak or depleted.	Charge the device and restart it.
	Power Button was pressed for only a short duration.	Press and hold the Power Button to turn on the device.
Slow charging / no charging.	Use other types of charging cables.	Use the provided charging cable only.
The red light fails to light up or is blinking.	The light beads are malfunctioning.	Contact our Service Center.
	Batteries are weak or depleted.	Charge the device and restart it.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:
<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty.
If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

- ☎ TEL: +44(0)7434 666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM (GMT)
+1(844) 417 0149 (US&CA)
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)
+49 40 23969285 (DE ONLY)
Monday-Friday 9:00AM-6:00PM (CET)
- ✉ Email: support-eu@renpho.com (UK&EU)
support@renpho.com (US&CA)
support-au@renpho.com (AU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1: This device may not cause interference.
- 2: This device must accept any interference, including interference that may cause undesired operation of the device.

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1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2011/65/EU, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:
<https://renpho.eu/pages/compliance-data>

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