
RENPHO

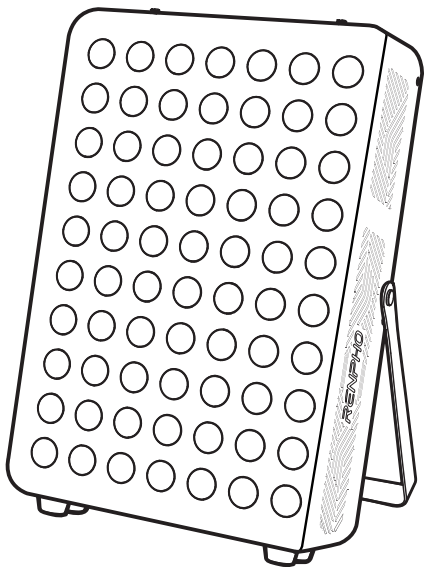**Red Light
Panel**

User manual

Model: QB70-M

CONTENTS

● ENGLISH	01~09
-----------	-------

● DEUTSCH	10~17
-----------	-------

● FRANÇAIS	18~25
------------	-------

● ITALIANO	26~33
------------	-------

● ESPAÑOL	34~41
-----------	-------

IMPORTANT SAFETY INSTRUCTIONS

READ ALL INSTRUCTIONS BEFORE USING.

1. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the appliance.
3. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.
4. If the product continues operating after the preset timer ends, stop using it immediately and check whether the controller is functioning properly. Improper use may cause the product to overheat, smoke, deform, or even catch fire.
5. Do not use the device for more than 30 minutes per session.
6. Do not look directly at the light source without wearing the protective goggles.
7. Individuals with photosensitivity are advised to conduct a skin tolerance test prior to use.

WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:

1. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
2. Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
3. Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.

4. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water.
5. Do not carry this appliance by supply cord or use cord as a handle.
6. Keep the cord away from heated surfaces.
7. Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
8. Never drop or insert any object into any opening.
9. Do not use outdoors.
10. Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
11. To disconnect, turn all controls to the off position, then remove plug from outlet.

DANGER - To reduce the risk of electric shock:

1. The appliance should never be left unattended when plugged in. Unplug from the outlet when not in use, and before cleaning, putting on or taking off parts.
2. WARNING: Do not carry this appliance by supply cord or use cord as a handle.
3. WARNING: Keep the cord away from heated surfaces.
4. DANGER: Do not reach for an appliance that has fallen into water. Unplug immediately.
5. If the supply cord is damaged, it must be replaced by the manufacturer, its service agent or similarly qualified persons in order to avoid a hazard.
6. The switch or controlled unit must not be allowed to get wet during washing and that during drying the cord must be positioned to ensure that water does not flow into the switch or control.
7. SERVICING OF DOUBLE-INSULATED PRODUCTS. In a double-insulated appliance, two systems of insulation are provided instead of grounding. No grounding means is


provided on a double-insulated product, nor should a means for grounding be added to the appliance. Servicing a double-insulated product requires extreme care and knowledge of the system, and should be done only by qualified service personnel. Replacement parts for a double-insulated appliance must be identical to the parts they replace.

SAVE THESE INSTRUCTIONS

WHAT'S IN THE BOX

Red Light Panel x1
Goggles x1
Door Hook x1
Power Cord x1
Extended Rope x1
Sling x1

ABOUT YOUR PRODUCT

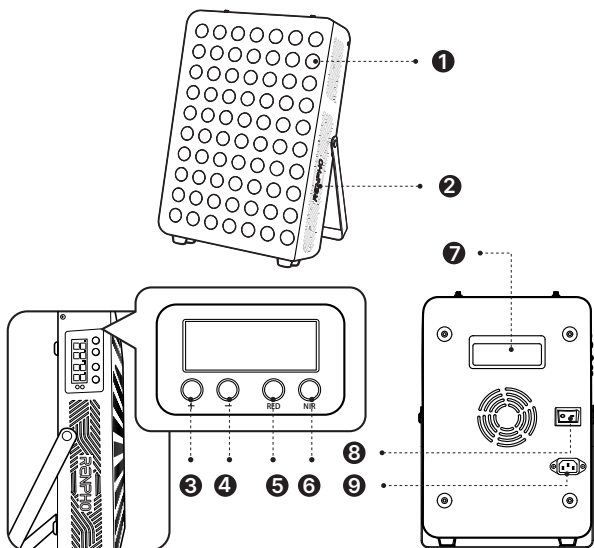- 1. LED Lights
- 2. Stand
- 3. Timer Up
- 4. Timer Down
- 5. RED Light

- 6. NIR Light
- 7. Recessed Handle
- 8. Power Switch
- 9. Power Port

USING YOUR PRODUCT

1. Before use

Place the appliance on a firm, even surface, free from obstructions. Then, plug into an appropriate electrical outlet. Insert the end of the power cord into the power port of the appliance. Turn on the power switch.

2. Operation

	Press to increase the auto-off timer in 5-min increments (5 mins to 30 mins).
	Press to decrease the auto-off timer in 5-min increments (30 mins to 5 mins).
	Press to enable / disable the red light (660 nm). This is for a visibly radiant complexion.
	Press to enable / disable the infrared light (850 nm). This is for deep relaxation and comfort.

Note:

- The infrared light is invisible to the naked eye.
- When there is no light available, the display will show "OFF" after 8s, then shut off.

TECHNICAL SPECIFICATIONS

Product Name: LED Red Light Panel

Model: QB70-M

Wavelength: 660 nm / 850 nm

Input: AC 100-240V 50 / 60Hz

Size: 13 x 9.1 x 2.4 in

Product Weight: 5.95 lbs

Working Temperature: -5°C-40°C

TROUBLESHOOTING GUIDE

Problem	Possible Cause	Solution
Power-on failure.	The power switch hasn't been turned on.	Turn on the power switch behind the appliance.
The appliance shuts off suddenly.	The power plug hasn't been connected.	Ensure the pins of the power plug are snugly fitted in the socket.
	Auto-off timer is activated.	Restart the device if needed.
The appliance doesn't light on when press "NIR".	The infrared light is invisible to the naked eye.	As long as the timer is on, there's no need to do anything.

WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery.

For warranty terms and conditions, please visit:
<https://renpho.com/pages/warranty-terms-and-conditions>

Note: Product registration is not required for the warranty.
If you choose not to register your product, it will not diminish the product warranty.

CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

☎ TEL: +44(0)7434 666088 (UK ONLY)
Monday-Friday 9:00AM-5:00PM (GMT)
+1(844) 417 0149 (US&CA)
Monday-Friday 9:00AM-4:30PM
1800 MY RENPHO (1800 69 7367) (AU ONLY)
Monday-Friday 9:00AM-6:00PM (AEST)
+49 40 23969285 (DE ONLY)
Monday-Friday 9:00AM-6:00PM (CET)

✉ Email: support-eu@renpho.com (UK&EU)
support@renpho.com (US&CA)
support-au@renpho.com (AU)

*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

FCC Regulatory Compliance

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

ISED Regulatory Compliance

This device contains licence-exempt transmitter(s)/receiver(s) that comply with Innovation, Science and Economic Development Canada's licence-exempt RSS(s). Operation is subject to the following two conditions:

- 1: This device may not cause interference.
- 2: This device must accept any interference, including interference that may cause undesired operation of the device.

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1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/30/EU, Directive 2014/35/EU, Directive 2011/65/EU, The Electrical Equipment (Safety) Regulations 2016, The Electromagnetic Compatibility Regulations 2016, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address: <https://renpho.eu/pages/compliance-data>

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