
# RENPHO

## USER MANUAL
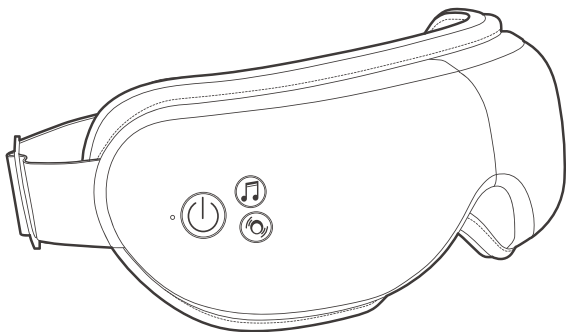
# EYERIS SMART

Model: RF-EM001S

# Contents

---

English..... 01-14

Deutsch..... 15-29

Français..... 30-44

Italiano..... 45-59

Español..... 60-73

Thank you for purchasing the RENPHO Eye Massager!  
Please read this instruction manual carefully and use the product only as indicated. RENPHO accepts no responsibility for any cause of injuries due to misuse or negligence to follow the instructions outlined. Please keep this manual for future reference.

## **IMPORTANT SAFETY INSTRUCTIONS**

### **READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.**

1. This device is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the device by a person responsible for their safety.
2. Children should be supervised to ensure that they do not play with the device.
3. The device is intended for household use only and not intended for medical use. It should not be used in treatment or alleviation of disease or abnormal physical state.
4. Cleaning and maintenance should not be conducted by children.
5. Please consult your doctor before use if you fit into any of the following user groups:
  - Persons who had undergone eye surgery or eye conditions such as glaucoma, cataract and retinal detachment
  - The elderly
  - Pregnant or postpartum women
  - Persons who are allergic to leather material
  - Persons with skin infections, sensitivity, redness and/or burns on the skin
  - Persons with any swelling or wound around the eye area
  - Persons using a pacemaker or other implantable medical devices
  - Persons with a heart condition, malignant tumors, or other serious health conditions

- Persons undergoing medical treatment
  - Persons with high or abnormal blood pressure
6. If you experience any discomfort during use, stop using the device immediately.
  7. Remove contact lenses, glasses, false eyelashes, hairpins, and any other accessories before use.
  8. Do not use the device while driving, walking, or running.
  9. Avoid contact with sharp objects.
  10. High temperatures may cause burns to occur regardless of the control setting. Do not use the device on infants, persons who are with disabilities, sleeping or unconscious. Check frequently where the skin is in contact with the heated area to reduce the risk of blistering.

**DANGER - To reduce the risk of electric shock:**

1. Never use pins or other metallic fasteners with this device.
2. Do not use the device with wet hands or submerge the device in water or any liquid.

**CAUTION - To reduce the risk of burns:**

1. Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
2. The device has a heated surface. Persons insensitive to heat must be careful when using the device.

**WARNING - To reduce the risk of burns, fire, electric shock, or injury to persons:**

1. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
2. Close supervision is necessary when this device is used by, on, or near children, invalids, or disabled persons.
3. Use this device only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this device if it is not working properly, has been damaged, or dropped into water.
5. Never drop or insert any object into any opening.
6. Do not operate where aerosol (spray) products are being

used or where oxygen is being administered.

### **About the battery:**

1. Do not use while bathing or in a shower.
2. Do not place or store device where it can fall or be pulled into a tub or sink.
3. Do not place in or drop into water or other liquid.
4. Never operate on a soft surface such as a bed or couch where the air openings may be blocked.
5. Do not expose the device or battery to excessive temperatures.
6. Be aware of the risk of terminals of the battery-operated device or battery being short-circuited by metal objects.
7. This device contains batteries that are non-replaceable. When the battery is at end of life, the device shall be properly disposed of.
8. Do not dismantle, open or shred battery.
9. Do not subject cells or batteries to mechanical shock.
10. In the event of a cell leaking, do not allow the liquid to come in contact with the skin or eyes. If contact has been made, wash the affected area with copious amounts of water and seek medical advice.
11. Secondary cells and batteries need to be charged before use. Always use the correct charger and refer to the manufacturer's instructions or product manual for proper charging instructions.
12. Do not leave a battery on prolonged charge when not in use.
13. After extended periods of storage, it may be necessary to charge and discharge the cells or batteries several times to obtain maximum performance.
14. CAUTION - The battery used in this device may present a risk of fire or chemical burn if mistreated. Do not recharge, disassemble, heat above 100°C (212°F).

# **SAVE THESE INSTRUCTIONS**

## ABOUT THE PRODUCT

### 1. What's in the Box

1 × Eyeris Smart

1 × USB-C Charging Cable

### 2. Technical Specifications

Product Name: Smart Eye Massager

Model: RF-EM001S

Input: 5V  $\equiv$  1A

Rated Power: 5W

Battery Capacity: 1500mAh 5.55Wh

Battery Voltage: DC 3.7V

Charging Time: about 3.5 hours

Working Time: 90-240mins

## USING THE PRODUCT


### 1. Charging

- a. Fully charge the device before the first use. Charge the device using the provided USB-C charging cable and a DC 5V 1A charger (not included).
- b. Low power warning: The indicator flashes red slowly, accompanied by a voice prompt "Low power, please charge".
- c. It takes about 3.5 hours to fully charge the device. The indicator flashes red slowly while charging. Once charging is complete, the indicator will turn solid red.

**Note:** Do not use the device while charging.


## 2. Wearing

- a. Remove any eye accessories (e.g. eyeglasses, contact lenses, false eyelashes, etc.) that may hinder you from using the device.
- b. Wear the device and ensure it fits your eyes properly. Adjust the elastic band to fit your head.

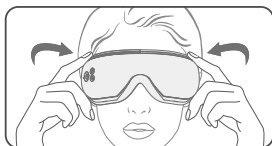**Note:** Your vision may become blurry for a while, as expected after using the device. It is recommended to close your eyes for a few minutes until your vision goes back to normal.

## 3. Operation
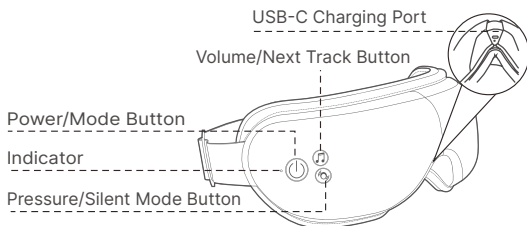
## Indicator

| Status                                                   | Indicator                                                                                      |
|----------------------------------------------------------|------------------------------------------------------------------------------------------------|
| Device is on                                             | Turns solid white                                                                              |
| Successful Bluetooth music connection/app connection     | Flashes blue slowly                                                                            |
| Successful Bluetooth music connection and app connection | Turns solid blue                                                                               |
| Low battery/charging                                     | Flashes red slowly                                                                             |
| Fully charged                                            | Turns solid red                                                                                |
| Charging at a voltage higher than 5V                     | Keeps flashing red quickly for 1min and then goes out, indicating an improper charger is used. |

**a. Press and hold to turn on/off the device.**

Once the device is turned on, the indicator will be solid white and the default setting is Air Pressure Level: soft, Heat: on, Volume Level: 3, Voice Mode: on, Bluetooth: on.

**b. Press to cycle through 5 modes:**

M1: Air Pressure + Heat

M2: Air Pressure + Heat + Vibration

M3: Air Pressure

M4: Heat Only

M5: Vibration

a. Press to cycle through soft pressure (Default)/strong pressure/no pressure.

b. Press and hold to switch between silent mode and voice mode.

Once silent mode is activated, the voice broadcast and built-in music will be disabled. The music from your phone or RENPHO Health App will be paused but can be resumed via your phone or app.

a. Press to cycle the music volume through Level 3 (Default)/Level 4/Level 5 (Max.)/Level 0/ Level 1/Level 2.

b. Press and hold to enter the next track when the device is playing music from your phone or RENPHO Health App.

**Note:** The volume of the voice broadcast cannot be adjusted.

## 4. Bluetooth Music Connection

a. Turn on the device and enable Bluetooth on your smartphone.

b. Search for “**RF-EM001S Music**” and connect it to play your own music playlist through the device.

**Note:** Once paired successfully, the device will be automatically connected to Bluetooth every time it is turned on.

## 5. Auto-off Timer

The device will automatically shut off after 15 minutes of continuous use unless you select another timer via the app.

## 6. Memory Function

- Except in DIY mode, the device will recall the last selected mode and pressure level once turned on again.
- Except in DIY mode, the device will recall the last selected pressure level and volume while switching modes.

## APP INSTRUCTIONS

### 1. Download the App

Scan the QR code below or search for “**RENPHO Health**” from the Apple App Store or Google Play to download the RENPHO Health App.

RENPHO Health

**Note:** Due to continuous updates and improvements, the "RENPHO Health" App may appear slightly different.

## 2. Log in/Sign up

- a. Log in directly with an existing RENPHO App account.
- b. If you don't have any RENPHO App account, tap **“Sign up”**. Register a new account for RENPHO Health App using your email address and complete your profile information.

## 3. Connect to the App

Before connecting to the RENPHO Health App, enable Bluetooth on your smartphone (Android System also requires location access) and turn on the eye massager.


### -Method 1

Open the app, and a window will pop up and tap **“Connect”**.

## -Method 2

- Open the app, and tap “+” on the top right corner of the “**Device**” page to pair manually.
- Select “**Message**” and tap “**Eyeris Smart**”.
- A window will pop up and tap “+”.
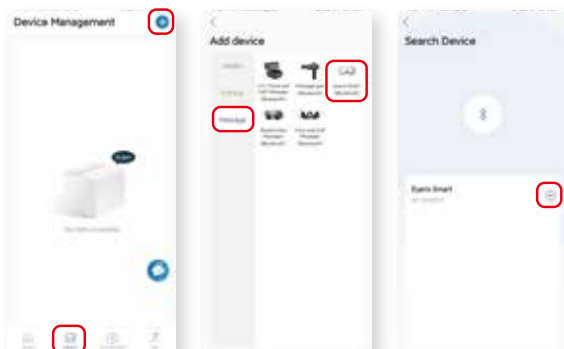
### Tips:

- When the eye massager is connected to the app, you will hear a voice prompt “app connected”.
- Once the device has started operating the mode set via the app, please do not touch any button on the device while wearing it to avoid changing the settings. Alternatively, you can lock the buttons via the app before wearing the device.

Want to get more out of Eyeris Smart? Scan the QR code or enter the URL to unlock more app features!

<https://renpho.com/pages/help-and-manual>

Discover how to lock the device buttons before use, preset the standard mode, set up a customized massage session and enjoy meditation MV or music. Take your eye relaxation to the next level with these awesome features!

## CLEANING & MAINTENANCE

1. Turn off the device or unplug the charging cable and allow the device to cool down before cleaning or storing.
2. Wipe the surface with a soft cotton cloth slightly dampened with mild detergent. Do not use abrasives, solvents, chemical cleaners, or harsh detergents.
3. Do not expose the device to direct sunlight or high temperatures. If the device has not been used for over 3 months, please charge it from time to time to prolong its lifespan.
4. It is recommended to wipe off any sweat, eye cream, or cosmetic residue left on the protein leather surface after use to prolong the lifespan of the device.
5. Do not disassemble the device by yourself. Otherwise, it may cause damage and void the manufacturer's warranty.

## TROUBLESHOOTING GUIDE

| Problem                                     | Possible Cause                                               | Solution                                                                                  |
|---------------------------------------------|--------------------------------------------------------------|-------------------------------------------------------------------------------------------|
| Power-on failure.                           | Out of power.                                                | Charge the device.                                                                        |
|                                             | Short press the power button.                                | Long press the power button to turn on the device.                                        |
| Low air pressure.                           | Low battery.                                                 | Charge the device.                                                                        |
| No air pressure.                            | No air pressure is selected.                                 | Switch to soft/strong pressure.                                                           |
|                                             | The mode without air pressure is selected.                   | Select the mode with air pressure function.                                               |
| The device shuts off suddenly.              | Out of power.                                                | Charge the device.                                                                        |
|                                             | The built-in auto-off timer is activated.                    | Restart the device for another massage session if needed.                                 |
| The device cannot be charged.               | Use other types of charging cables.                          | Use the provided Type-C charging cable only.                                              |
|                                             | Use a charger with a voltage higher than 5V.                 | Use a DC 5V 1A charger.                                                                   |
| Fail to connect the device to the app.      | Your phone's Bluetooth is not enabled.                       | Enable Bluetooth on your smartphone and pair the device following the "APP INSTRUCTIONS". |
|                                             | The eye massager has been paired with another mobile device. | Unbind the eye massager on another mobile device.                                         |
| Fail to play Bluetooth music on the device. | Bluetooth music connection is disabled.                      | Enable Bluetooth on your smartphone and connect to "RF-EM001S Music".                     |
|                                             | Silent Mode is activated.                                    | Long press the Pressure button to switch to voice mode.                                   |

| Problem                                                                                       | Possible Cause | Solution         |
|-----------------------------------------------------------------------------------------------|----------------|------------------|
| The app is not working.                                                                       | App error.     | Restart the app. |
| <b>If you encounter any other problems, please feel free to contact our customer service.</b> |                |                  |

## WARRANTY POLICY

Your RENPHO product purchase is covered by a one-year limited manufacturer warranty from the date of delivery. For warranty terms and conditions, please visit:

**<https://renpho.com/pages/warranty-terms-and-conditions>**

Note: Product registration is not required for the warranty. If you choose not to register your product, it will not diminish the product warranty.

## CUSTOMER SERVICE

Please feel free to contact us if you have any questions or concerns. RENPHO Customer Service Team guarantees a quick response and hassle-free solutions to any issue you may have within business hours.

✉ **Email:** **support-eu@renpho.com (UK&EU)**  
**support@renpho.com (US&CA)**  
**support-au@renpho.com (AU)**

☎ **TEL:** **+44(0)7434 666088 (UK ONLY)**  
Monday-Friday 9:00AM-5:00PM (GMT)  
**+1(844) 417 0149 (US&CA)**  
Monday-Friday 9:00AM-4:30PM  
**1800 MY RENPHO (1800 69 7367) (AU ONLY)**  
Monday-Friday 9:00AM-6:00PM (AEST)

\*For defective products or the return of items, please contact us with your order number within the specified warranty period. DO NOT dispose of any product parts as they may be required for inspection/repair.

### 1. Meaning of the crossed-out wheeled bin symbol

The crossed-out wheeled bin symbol means that you are legally obliged to dispose of these appliances separately from unsorted municipal waste, i.e. old appliances must be disposed of separately from household waste. Batteries or rechargeable batteries that are not permanently installed must be removed beforehand and disposed of separately.

### 2. Options for returning old appliances

Owners of WEEE can return or collect it within the framework of the possibilities set up and available by public waste management authorities to ensure proper disposal.

### 3. Personal data

All end-users of WEEE are reminded that you are responsible for deleting personal data from the WEEE to be disposed of.

1. The symbol of the crossed-out garbage can on batteries or accumulators indicates that they must not be disposed of in household waste at the end of their life.
2. You are legally obliged to return old batteries and accumulators after use. You can do this free of charge at a collection point in your area. Addresses of suitable collection points can be obtained from your city or local government.
3. Batteries may contain substances that are harmful to the environment and human health. The separate collection and recycling of old batteries and accumulators is intended to avoid negative effects on the environment and human health.
4. Please avoid the generation of waste from old batteries as far as possible, e.g. Please avoid littering public spaces by not carelessly leaving batteries or electrical and electronic equipment containing batteries.

Hereby, Shenzhen Ruiyi Business Technology Co., Ltd. declares that this equipment complies with Directive 2014/53/EU, Directive 2011/65/EU, Radio Equipment Regulations 2017, Restriction of the Use of Certain Hazardous Substances in Electrical and Electronic Equipment Regulations 2012, and other relevant provisions. The full text of the Declaration of Conformity can be obtained at the following internet address:

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