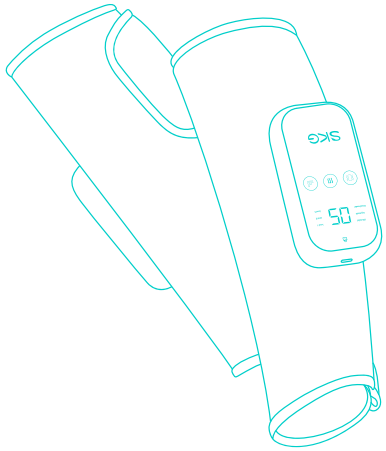




The diagram illustrates the control panel of the SKG K5 Pro smart heating pad. It features a central LED digital display showing '50'. Above the display is a lock screen indicator. Below the display are three circular buttons: Message mode (with a message icon), Heating mode (with a wavy line icon), and Power level (with a bar chart icon). The SKG logo is prominently displayed below the buttons. At the top of the device is a Type-C charging port, and at the bottom is a power button/lock screen button.

Labels and their corresponding features:

- Type-C charging port
- Lock screen indicator
- LED digital display one
- Message mode
- Heating mode
- Power level
- Power button/lock screen button

1. Before using the massager, please check whether it can be started. If not, please charge it with the charging cable supplied with the appliance.
2. Separate the Velcro tapes on the massager, so that the product is fully unfolded.
3. Bundle the massager on your calf. Make sure that the power button of the main unit is facing upwards and adjust the tightness appropriately.
4. After pressing and holding the power button to start up this product, the buzzer will beep twice; the product will be powered on by default; the unlock screen indicator will be on simultaneously in green. Tap the power button to unlock/lock the screen, and the lock screen indicator will be on in red.
5. When this product is started, press the mode button to switch between modes: lift massage mode, leg comfort mode and kneading relaxation mode. This product will be in the lift massage mode by default once it is started and at the level 1 power; press the power key to change the power from level 1 to level 3 in cycle. When the mode/power is switched, the buzzer will beep once.
6. Once the product is switched on, each time you press the heating mode button, the temperature of the product will be switched once in the order of 40°C, 45°C, 50°C (temperature shown on the display). The buzzer beeps once when the heating temperature is changed.
7. After turning on this product, the run time set is 30 minutes by default. After the time set is up, the buzzer will beep twice; all loads will be off, and this product will be powered off.
8. The product will not work when it is being charged. When it is connected with the USB power supply, the lock screen indicator will be green and blinking, and the green indicator will always be on when it is fully charged.

	Operating environment	Storage and transportation
Temperature range	0°C~35°C	-20°C~50°C
Humidity range	15%~90%RH	0%RH~90%RH
Atmospheric pressure	70KPA~106KPA	70 KPA ~110KPA

