



SKG

Eye Massager
E7-EN

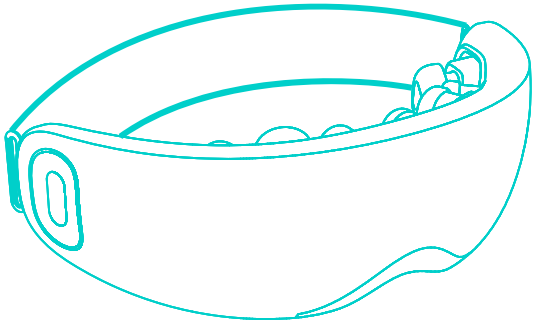

Manual

Product diagram

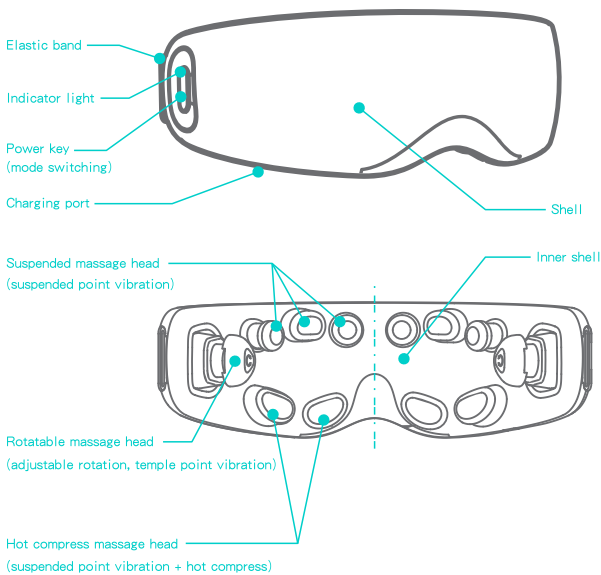

Operation guidance

* Before using this product, please adjust the elastic band to the appropriate length, and then wear it to the eye area. Properly adjust the position and tightness of the massage head according to the wearing situation, and then start the host to experience the massage effect.

Operation steps:

1) Long press the "power key" of the host for about 2 seconds, and the host starts up (with voice prompt, the prompt light is on). The first boot defaults to the vibrant mode, and the mode can be switched by short pressing the "power key". In the power on state, press and hold the "power key" for about 2 seconds to turn off the product (with voice prompt, the indicator light turns off)

2) When you press the "power key" once, the mode will be switched to one gear, and the four modes will be switched in a cycle (with voice prompts: "vibrant mode", "soothing mode", "sleep mode" and "eye exercises mode"; background music will be played when you are in vibrant mode, soothing mode and eye exercises mode). Under different modes, the massage techniques are different, and each mode has the function of hot compress.

3) The host is timed for 15 minutes ($\pm$ 30 seconds). When the host is powered on, it will start timing. It will automatically shut down after 15 minutes ($\pm$ 30 seconds), and there will be a voice prompt with the shutdown (if it is automatically shut down in sleep mode, the shutdown prompt tone will not be played)

4) When the mobile phone is turned on, the Bluetooth searches 'S1511HB'. For the first connection, you need to manually click the connection. After pairing and saving, each time you turn on the Bluetooth, it will connect automatically. After the connection is successful, the voice prompts "Bluetooth connected". After the mobile phone disconnects Bluetooth, the host voice prompts "Bluetooth disconnected". When the mobile phone connects the device with Bluetooth and plays mobile music, the music brought by the device stops and plays mobile music.

5) When the product is under low battery, it will be accompanied by a low power prompt, and the indicator light will change from blue light to orange light; The product will shut down automatically in case of power shortage. When the machine is shut down, the massage mode before the shutdown will be recorded. After the next startup, the machine will return to the massage mode before the last shutdown. When you start up for the first time, enter the default mode (vibrant mode).

Basic product parameters

Device Name	E7-EN	Model	S1511HB
Rated Input Voltage	5V	Rated Input Power	5W
Battery Capacity	1500mAh/5.55Wh		
Note: This product is rechargeable and is suitable for charging with type-C interface. Users need to charge with a standard charger with an output of 5V , 1A.			

Recommendations:

- Please use it in normal environment: temperature: 0 $^{\circ}$ C ~ 40 $^{\circ}$ C; Relative humidity: < 80% RH, no condensation; Atmospheric pressure: 50kPa ~ 106kpa.
- Please select normal transportation and storage conditions: ambient temperature range: - 20 $^{\circ}$ C ~ 55 $^{\circ}$ C; Relative humidity range: $\leq$ 95% RH, no condensation; Atmospheric pressure: 50kPa ~ 106kpa.
- It is recommended that you use it no more than twice a day, no more than half an hour each time, according to your personal situation.

Daily maintenance

- When the massage device surface is dirty, it can be wiped with a clean soft cloth dipped in a little water. It cannot be wiped with volatile liquids such as alcohol. The product cannot be disassembled and rinsed in water.
- When not in use for a long time, please cover the product to avoid dust.
- In the process of cleaning the machine, prevent water or other liquids from entering the shell to avoid short circuit during use.

SKG Health Technologies Co.,Ltd.
Website: <https://global.skg.com/>
Address: 23A Floor,Building 3,Zhongke R & D Park,No. 009,Gaoxin South 1st Road,High-tech Zone Community,Yuehai street,Nanshan District,Shenzhen,Guangdong,China
*If you have any questions, please feel free to send an email to support@skg.com

Notes for Product Scrap

Please recycle the rechargeable battery (lithium ion battery)!

 Danger	The battery in this product is a dedicated rechargeable battery, which is not intended for other products. Please do not charge the removed rechargeable battery. - Do not throw it in the fire, or heat it up to high temperature - Do not pierce it with nails, impact it by gravity or decompose it for modification - Do not connect the positive and negative poles of the battery directly with metal - Do not store it together with necklaces, hairpins, etc. - Do not charge, use or put it in the places in high temperature such as fireside or in the blazing sun - The situations above may cause heating, fire and breakage.
 Warning	Keep the removed rechargeable battery out of the reach of infants and young children in case that it is swallowed by mistake - When the device is mistakenly swallowed, please consult a doctor immediately - When the battery fluid leaks, do not touch the fluid by hand, and it can be processed in the following way - Do not rub your eyes with your hands when battery electrolyte enters your eyes in case of blindness. Rinse your eyes with clean water immediately and seek medical advice. When your body or clothing is contaminated by the battery electrolyte, it may cause skin irritation or injury. Please rinse with water. If you still feel discomfort, seek medical advice immediately. - Only the scrapped products can be decomposed. Otherwise, the accident of fire, electric shock, injuries may be resulted in. When you need to repair the battery (not to replace the battery), please consult the dealer where you purchased this product

The lithium-ion battery in this product is a valuable resource that can be recycled and reused, which is not allowed to be replaced. Please hand it over to a professional agency that can recycle the rechargeable batteries when this product is scrapped. Before scrapping this product, remove the battery from this appliance. Before removing the battery, the appliance must be powered off. The battery should be disposed of safely, which can be taken out in the steps below:

- With professional tools, open the upper and lower cases.
- Remove the circuit board and the motor.
- Take out the battery, cut the battery wire with professional pliers, and then take out the battery.

Matters needing attention

Attention

- Before using this product, please carefully read and abide by the operation instructions in the operation manual. The company is not responsible for the consequences caused by incorrect use or non-compliance with the instructions of this product.
- As we constantly improve our products, the technical specifications are subject to change without notice.
- Do not use this product after facial makeup, which will damage your makeup and pollute the contact surface of the product. At this time, please wipe the contaminated surface with clean water.
- Do not use the eye massager while driving.
- Do not use the eye massage instrument in places with excessive high temperature and humidity.
- Do not use in patients with stroke, hypertension and hypotension.
- Do not use when there are symptoms of eye diseases such as glaucoma, cataract, retinal detachment and red eye disease.
- Please avoid using eye massage after drinking.
- It is forbidden for patients who have performed eye surgery and have heart disease.
- Before use, please make sure there are no contact lenses in the eyes, and remove the hairpin, hair ornament, etc.
- This product is a personal care massage product and cannot replace drugs and medical devices.
- Children, pregnant women and women in menstrual period are forbidden to use eye massage instrument.
- There is a lithium battery inside. Please do not put it in a high temperature environment.
- For sensitive skin, facial allergic dermatitis and plastic materials on the face, use this product under the guidance of a doctor.
- This product has heating function. Please pay attention when using it for those who are not sensitive to heat.
- If you feel unwell during use, please stop using it immediately.
- Do not get massage in other parts of the body except the eyes. It is forbidden to disassemble this product by yourself.
- The product has no waterproof function. Please do not wash it or use it in high humidity environment. If the lining is stained, you can wipe it gently with a wet towel



§ 15.19 Labelling requirements.

This device complies with part 15 of the FCC Rules. Operation is subject to the condition that this device does not cause harmful interference.

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) This device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

§ 15.21 Information to user.

Any changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

§ 15.105 Information to the user.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.



This device was tested for typical body - worn operations. To comply with RF exposure requirements, a minimum separation distance of 0cm must be maintained between the user's body and the handset, including the antenna. Third - party belt - clips, holsters, and similar accessories used by this device should not contain any metallic components. Body - worn accessories that do not meet these requirements may not comply with RF exposure requirements and should be avoided. Use only the supplied or an approved antenna. This device in compliance with the essential requirements and other relevant provisions of Directive 1999/5/EC. All essential radio test suites have been carried out.

1. CAUTION - RISK OF EXPLOSION IF BATTERY IS REPLACED BY AN INCORRECT TYPE. DISPOSE OF USED BATTERIES ACCORDING TO THE INSTRUCTIONS.
2. The device complies with RF specifications when the device used at 0mm from your body please visit <https://www.skg.com/userManual> to find our CE Declaration of Conformity(DoC) info. and Importer contact info.



Old electrical appliances must not be disposed of together with the residual waste, but have to be disposed of separately. The disposal at the communal collecting point via private persons is for free. The owner of old appliances is responsible to bring the appliances to these collecting points or to similar collection points. With this little personal effort, you contribute to recycle valuable raw materials and the treatment of toxic substances.

If you experience discomfort or skin irritation when wearing your smart watch, then we recommend you try cleaning your device. Sometimes residue or foreign materials build up around your device and may aggravate your skin. It's also possible that you are not wearing the watch correctly. We recommend making sure to clean and adjust your watch regularly for a more comfortable fit.

Caution:

- If you experience skin irritation when wearing your watch, please refrain from wearing, and wait two to three days to see if your symptoms ease. If symptoms persist or worsen, please consult a physician.
- If you have eczema, allergies or asthma, you may be more likely to experience a skin irritation or allergy from a wearable device.

IC Caution

This device complies with Industry Canada licence-exempt RSS standard(s). Operation is subject to the following two conditions: (1) This device may not cause interference, and (2) This device must accept any interference, including interference that may cause undesired operation of the device. Le présent appareil est conforme aux CNR d'Industrie Canada applicables aux appareils radio exempts de licence. L'exploitation est autorisée aux deux conditions suivantes:

- (1) L'appareil ne doit pas produire de brouillage, et
- (2) l'utilisateur de l'appareil doit accepter tout brouillage radioélectrique subi, même si le brouillage est susceptible d'en compromettre le fonctionnement.

- The appliance it must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision.

90*100mm 128g C2S, 印刷颜色3252C+Black 风琴4折