
SKG

Neck & Shoulder Massager H5 Mini

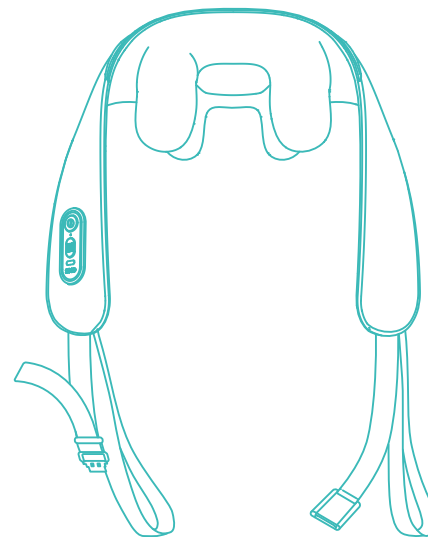SKG Health Technologies Co., Ltd.

Address: 23A Floor, Building 3, Zhongke R&D Park, No. 009,
Gaoxin South 1st Road, High-tech Zone Community, Yuehai Street,
Nanshan District, Shenzhen, Guangdong, China.

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IMPORTANT SAFEGUARDS

EN

WHEN USING ELECTRICAL PRODUCTS, ESPECIALLY WHEN CHILDREN ARE PRESENT, BASIC SAFETY PRECAUTIONS SHOULD ALWAYS BE FOLLOWED, INCLUDING THE FOLLOWING: READ ALL INSTRUCTIONS BEFORE USING

DANGER – TO REDUCE THE RISK OF ELECTRIC SHOCK:

Always unplug the device from the electrical outlet immediately after using and before cleaning.

- DO NOT reach for a device that has fallen into water. Unplug it immediately.
- DO NOT use the device while bathing or in a shower.
- DO NOT place or store the device where it can fall or be pulled into sink or tub.
- DO NOT place or drop the device in water or other liquid.
- DO NOT plug or unplug the charging cable with wet hands.

WARNING – TO REDUCE THE RISK OF BURNS, ELECTRIC SHOCK, FIRE OR INJURY TO PERSONS:

- The device should NEVER be left unattended when plugged in. Unplug from outlet when not in use and before putting on or taking off parts or attachments.
- Close supervision is necessary when this device is used by, on, or near children, invalids or disabled persons.
- Children should not use this device.
- Use this device only for its intended use as described in this manual. DO NOT use attachments not recommended by SKG; specifically any attachments not provided with the device.
- NEVER operate this device if it has a damaged power button, if it is not working properly, if it has been dropped or damaged, or dropped into water.
- NEVER drop or insert any object into any opening of the device.
- DO NOT operate the device where aerosol (spray) products are being used or where oxygen is being administered.
- DO NOT operate the device under a blanket or pillow. Never cover the device when it is in operation, as excessive heating can occur, causing fire, electric shock or injury to persons.
- To disconnect, turn all controls to the "off" position.
- NEVER operate the device with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- NEVER operate the device on a soft surface such as a bed or couch where the air openings may be blocked.
- Use heated surfaces carefully. DO NOT use over insensitive skin areas or in the presence of poor circulation. The unattended use of the heated device by children or incapacitated persons may be dangerous.
- DO NOT use this massager in close proximity to loose clothing, jewelry, or long hair.

SAVE THESE INSTRUCTIONS

CAUTION—PLEASE READ ALL INSTRUCTIONS CAREFULLY BEFORE USE.

- This product is not intended for medical use. It is intended only to provide a luxurious massage.
- Consult your doctor prior to using this product if you are pregnant, you have a pacemaker or have any concerns regarding your health.
- This product is not recommended for use by diabetics.
- NEVER leave the device unattended, especially when children are present.
- DO NOT use this product for more than 15 minutes at a time.
- NEVER use this product directly on swollen or inflamed areas, or on skin eruptions.
- DO NOT use this product as a substitute for medical attention.
- NEVER use this product while lying flat or in bed.
- NEVER use this product with a necklace on.
- This product should NEVER be used by any individual suffering from any physical ailment that would limit the user's capacity to operate the controls.
- This product should not be used by children or invalids without adult supervision.
- This product is not intended for use in automobiles.
- This product is intended for household use only.

I-Instructions for Use

The Neck Massager is designed for temporary relief of pain associated with sore and aching muscles in the back of the neck due to strain from exercise or normal household and work activities.

II-Basic Product Parameters

Model: H5 Mini

Product Name: Neck & Shoulder Massager

Rated Input Voltage: 5V ===

Battery Capacity: 2000mAh/14.8Wh

Rated Input Power: 10W

Operating environment:

Temperature range: 0°C-40°C

Relative Humidity: <80% RH, non-condensing

Atmospheric Pressure: 50Kpa-106Kpa

Storage and transportation environment:


Temperature range: -20°C to +55°C

Relative Humidity range: <95% RH, non-condensing

Atmospheric Pressure: 50Kpa-106Kpa

III-Device Description

1. General Description

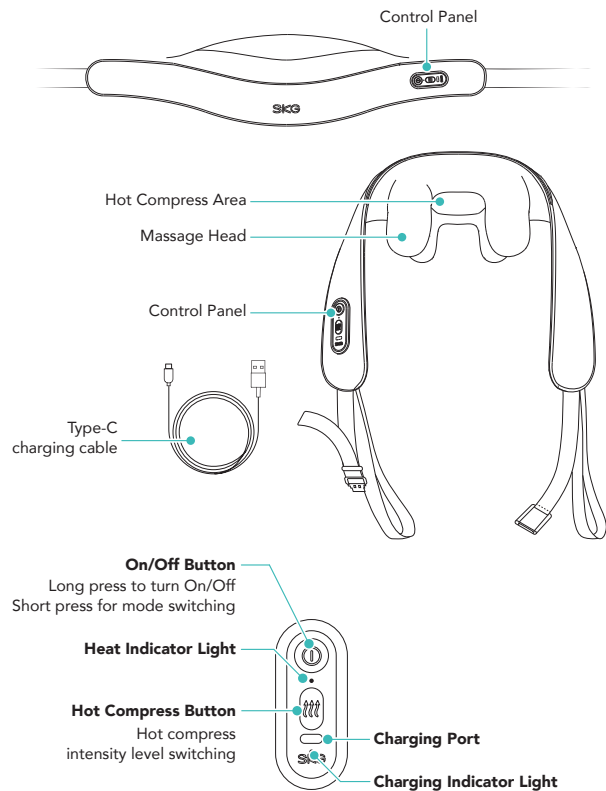2.Power On/Off

Press and hold the power button to turn on or off.

3.Massage Mode Control

This massager offers three speed settings:

- Mode 1: High speed (25 times/min)
- Mode 2: Medium speed (31 times/min)
- Mode 3: Low speed with vibrations

How to use:

- Press and hold to power on at high speed.
 - Press once for medium speed.
 - Press twice for low speed with vibrations.
 - Press three times to return to high speed.
- Repeat this sequence to cycle through the modes.

4.Heat Mode Control

This massager offers two heat settings:

- Low heat (orange light)
- High heat (red light)

The indicator light will be off when the heat mode is deselected.

To operate:

- Begin at low heat level.
- Press again for high heat.
- Press a third time to turn off heat.

Repeat the cycle as desired.

5.Power/Charge Indicator


Light Color	Light Color	Light Color
Blue light	Steady on	Battery level >25%
Orange light	Steady on	Battery level <24%
Orange light	Breathing	Battery charging
Blue light	Breathing	Full battery

6.Caution

All servicing of this massager must be performed by authorized SKG service personnel only.

Note: The device will automatically shut off after 10 minutes for your safety.


IV-Method of Wearing

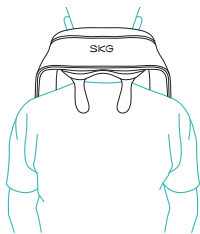
Place the massager around your neck


Pull the port to the back (avoid discomfort caused by overtightening)

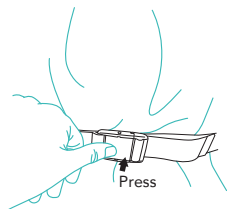After use, press the buckle to release the shoulder strap

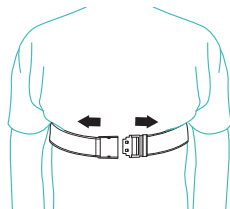Release the back clips

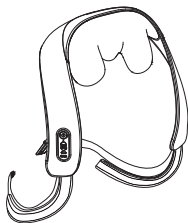Carefully remove the product from the shoulder to prevent falling

V-Maintenance

Storage

Place the massager in its box or in a safe, dry, cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the fabric surface. Keep the product out of reach of infants, children, and anyone who is fully or partially unconscious.

Cleaning

Turn off and allow it to cool down before cleaning. Spot clean only with a soft, slightly damp cloth. Never allow water or any other liquids to come into contact with the device. DO NOT immerse in any liquid to clean. NEVER use abrasive cleaners, brushes, gasoline, kerosene, glass / furniture polish or paint thinner for cleaning.

VI-Troubleshooting

Refer to the table for solutions to any problems. If the problem persists, please return the product to the dealer for repair.

Problem	Possible Cause	Solution
Unable to start	Low battery	Charge the product
	It can't be turned on	Long-press the power button to turn it on
	The product is still connected to the power socket	Remove the charging cable from the device
Unable to charge or the charging light is not on	The charging cable is not properly inserted into the power socket	Replug the charging cable
	Wrong charging cable used	Use a Type-C connector cable

⚠ WARNING

No modification of this device is allowed.
Do not modify this device without authorization of the manufacturer.
Any modification requires proper inspection and testing to ensure the device remains safe for use.

VII-Notes for Product Scrap

DANGER

Please support the recycling of rechargeable battery (lithium-ion battery).

The battery included in this product is a dedicated rechargeable battery, not intended in other products. Please do not charge the removed rechargeable battery.

- DO NOT throw it into a fire, or heat it up to high temperature.
- DO NOT pierce it with nails, impact it by gravity or disassemble it for modification.
- DO NOT connect the positive and negative poles of the battery directly with metal.
- DO NOT store it together with necklaces, hairpins, etc.
- DO NOT charge, use or place the battery in areas with high temperatures, such as near a fireplace or under direct sunlight.

• The above situations may cause heating, fire and breakage.

WARNING

Keep the removed rechargeable batteries out of reach of children to prevent accidental ingestion.

If a battery is swallowed, consult a doctor immediately.

In the event of a battery fluid leak, avoid touching the fluid and handle according to safety procedures.

- If battery electrolyte enters the eyes, rinse with clean water and seek medical advice without rubbing your eyes.
- If battery electrolyte contacts the skin, rinse the area with water and get medical help if needed.

Only disassemble scrapped products to prevent accidents such as fires or electric shocks. For battery repair, contact the dealer where you purchased the product. This does not include battery replacement.

The lithium-ion battery in this product is recyclable and should not be replaced; surrender it to a professional recycling agency when the product is discarded.

Always turn off the device before removing the battery for disposal.

To safely dispose of the battery, follow these simplified steps:

Use specialized tools to detach the device's casing.

Remove the internal circuit board and motor.

Extract the battery, sever its wires with pliers, and remove it completely.

VIII-Certification Statement

Note: This product has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when used in a residential environment. This product generates, uses, and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications.

However, there is no guarantee that interference will not occur in a specific installation. If this product does cause harmful interference to radio or television reception, which can be determined by turning the product off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the distance between the device and the receiver.
- Connect the device to an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This device in compliance with the essential requirements and other relevant provisions of Directive (EU)2023/988. CAUTION: THERE IS A RISK OF EXPLOSION IF THE BATTERY IS REPLACED WITH AN INCORRECT

TYPE. DISPOSE OF USED BATTERIES ACCORDING TO THE INSTRUCTIONS.

EU Declaration of Conformity:

SKG Health Technologies Co., Ltd. hereby declare that this device complies with the essential requirements and other relevant provisions of Directive (EU)2023/988. A copy of the EU Declaration of Conformity is available online at <https://www.skg.com/pag-es/skg-declaration-of-conformity>.

Old electrical appliances must not be disposed of with the residual waste; they must be disposed of separately. Disposal at public collection points is free for private individuals. The owner of old devices is responsible for bringing them to these collecting points or similar facilities. By making this small personal effort, you contribute to the recycling of valuable raw materials and the proper treatment of toxic substances.

IX-Register & Support

Register your SKG product at support@skg.com to maximize support and visit <https://www.skg.com/> for expert tips, tutorials, and FAQs to optimize your device use.

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