


Use & care manual

FOOT MASSAGER

SL-594

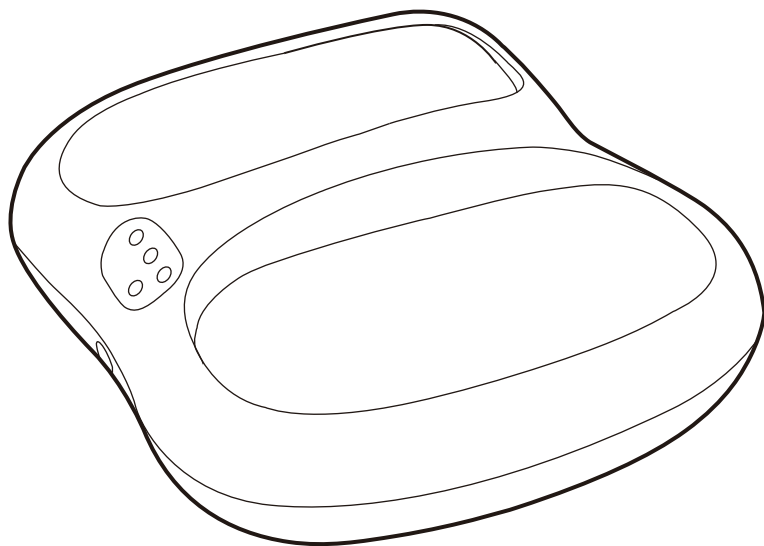

 **SNAILAX[®]**

Important Safety Instruction



Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these instructions for further reference!

- The massage device conforms to recognized technical principles and the latest safety regulations.
- Do not wet, do not use pins, never remove cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Never use any power sources other than the home adapter provided as original equipment with this device.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- Never operate under blankets where the air opening may be blocked.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use this device while driving.
- Do not use while you sleep.
- Do not use with infants, a helpless person or a person insensitive to heat.
- Overlong application can result in skin burns.
- Do not use this product for medical treatment.
- For indoor use only
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 15 minutes at a time.

Important Safety Instruction

- We equip the massager with a protection against over-heating for safe use.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy.

Failure to observe the above instructions may constitute misuse of the product and can cause serious injury or burns.

SAVE THESE INSTRUCTIONS



The device is not to be used for commercial use or medical purposes. If you have health concerns, consult your doctor before using this product.



Only use indoors!



Do not use the device in damp rooms
(e.g. when bathing or showering)





• Do not wash!

• Do not dry clean!

Product Parts

Configuration Figure of Functions

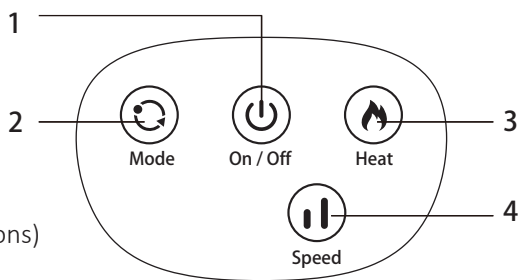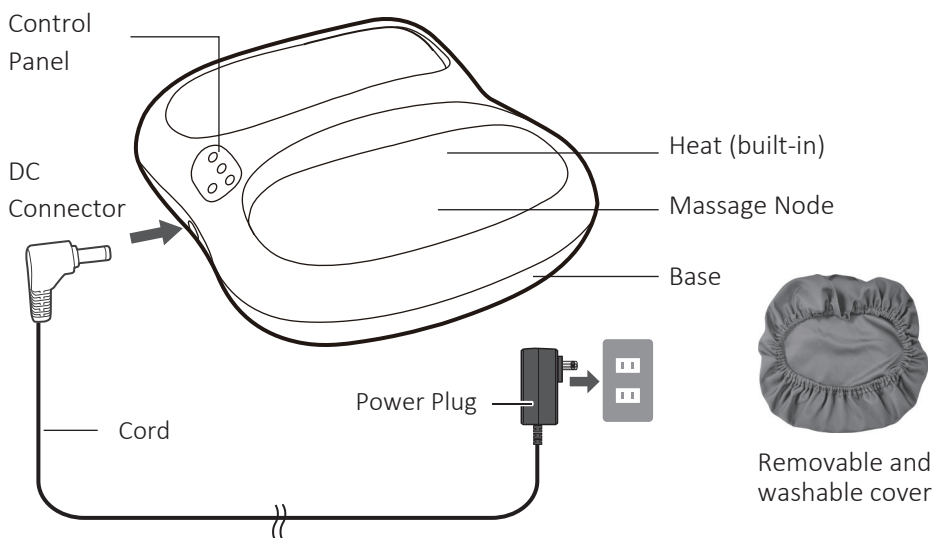

1. On/Off Button: ON/OFF
2. Mode Button : (2 modes for options)
3. Heat Button: (Low / High)
4. Speed Button: (Low / High)



ATTENTION

Please use the slightly lower level massage function when starting up the appliance. Stop using it if you feel anything abnormal.

Operations for Use

After the unit is powered on, the whole machine enters the standby state.
The detailed functions of each control button are as follows



Power Button

Press to turn on the machine, the 'rolling & scraping' vibration' and 'heating' functions are turned on. The lights on the control buttons light up. The initial state of the boot runs in 'Mode I' and 'low speed', 'high heating'. 15 minutes automatic timing, after running for 15 minutes, the massager will automatically shut down;
Press this button again to shut down.



Mode Button

Press this button to switch the massage modes.

Mode I: The rolling, scraping and vibration functions are running simultaneously, the Mode button lights up green.

Mode II: Only vibration massage, the Mode button lights up red.



Heat Button

Press this button to cycle low/high/off settings of heating function.

Low Level: The heat button lights up green

High Level: The heat button lights up red



Speed Button

Press this button to cycle low/high speed or off settings of Mode I or Mode II massage function.

Low speed: The speed button lights up green

High speed: The speed button lights up red

Model	SL-594
Name	Foot Massager
Power input	DC 12V  4A
Power consumption	48 W
Default Time	15 mins

Cleaning & Maintenance Instructions

- Before cleaning, make sure the power supply is turned off and the plug has been pulled out from the socket.
- It is forbidden to dip the whole product in water or other liquid when cleaning.
- Use a soft cloth to wipe the surface of the product and control panel.
- Use a piece of cloth that has been dipped with cleaning agents to wipe the main unit of the foot Massager, and then use a piece of dry and soft cloth to wipe it.
- Do not pull the power plug by damp hands.
- Do not use water washing product, do not spray the water or other liquids to the interior of the product.
- Do not use corrosive cleaning agents like benzene or thinner to clean the Foot Massager, otherwise, the color of upholstery will fade or the Foot Massager will crack or be eroded.
- Do not iron the upholstery.
- When product is wet, let the product dry naturally and thoroughly.

Maintenance

- Make sure the machine is turned off and the plug is pulled out from the socket.
- Sweep off the dust. Keep the foot massager clean and away from damp places.
- Store the massager in a cool, dry place out of reach of children.
- Avoid contact with sharp edges or pointed objects which might cut or puncture the cover of massager.

TROUBLESHOOTING GUIDE

Fault	Solutions
There is a clattering noise in the movement	This is perfectly normal caused by the working of product. It is normal status.
The massage device switches off suddenly	Device may have been in use for 15 minutes. Auto-shut off function activated. Continue using by pressing the 'power' button after resting the unit for 15 minutes.
The massage device fails to start	Check if device has been used for more than 15 mins. If so, let device rest for 15 mins before continuing. Device may require longer time to cool in warmer environments.
	Check if adapter is plugged into socket.
	Check if adapter output is connected to device power inlet.

Warranty

If you have any issue about the product, please feel free to contact us by sending email to support@snailax.com. We will strive to provide the best service possible within 24 hours.

30 days unconditionally Return

Snailax product can be returned to receive a full refund for any reason within 30 days. Please contact our customers service (support@snailax.com), our staff will contact you within 24 hours.

90 days refund/replace

Snailax product can be returned / replaced within 90 days if the product break down in the period of proper use.

12 months Warranty

If the product break down within 12 months in the period of proper use, customers can enjoy the relevant warranty to get them replaced.

Attention!

No warranty will be given to any force majeure and man-made causes to a defective product, such as improper care, personal tear down and willful damage, etc.

Extend Warranty for Free

- 1) Enter the following URL or scan the QR code below to find Snailax facebook page and like it, enter "Warranty" to messenger to extend your warranty from 1 year to 3 years.



<https://www.facebook.com/snailax>

- OR 2) Send message "Warranty" and email us support@snailax.com to extend your warranty from 1 year to 3 years.

Have question ?

Tel: 734-709-6982

Monday-Friday 9:00AM-4:30PM

Email: support@snailax.com

FCC STATEMENT

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this unit is not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.



2750 Carpenter Road,
Suite 1 B, Ann Arbor, 48108
Tel: 734-709-6982
support@snailax.com
www.snailax.com