

Use & care manual

NECK & SHOULDER SHIATSU MASSAGER WITH HEAT SL-632U



Snailax **SNAILAX**

Safety Instructions



Please read the following instructions carefully before putting your massage device into use to ensure trouble-free operation and optimum efficiency. Please retain these operating instructions for further use!

- The massage device conforms to recognized technical principles and the latest safety regulations.
- Do not wet, do not use pins, never remove cover.
- This item is NOT A TOY. Close supervision is necessary when this appliance is used by, on, or near children or disabled persons.
- This appliance should never be left unattended when plugged in.
- Never use any power sources other than a standard 110V wall outlet or the 12V Car adaptor provided as original equipment with this unit.
- Any possible repairs may only be carried out by authorized expert staff. Improper use and unauthorized repairs are not permitted for safety reasons and lead to loss of warranty.
- Never touch the power plug with wet hands.
- Please avoid contact of the device with water, high temperatures and direct sunlight.
- Do not use any damaged cables, plugs or loose sockets.
- Never operate under blankets where the air opening may be blocked.
- If plugs or cords are damaged, they must be replaced by the manufacturer, a service representative or by qualified personnel.
- In case of malfunction, disconnect immediately from the mains.
- Do not use if you have skin disorders, open wounds, or swollen or inflamed areas.
- Misuse or incorrect use excludes any liability for damage.
- Do not use this device while driving.
- Do not use while you sleep.
- In order to avoid excessively stimulating the muscles and nerves, the recommended massage time should not be exceeded 15 minutes at a time.
- Each massage – even a hand massage – must be refrained from during pregnancy or if one or more of the following complaints are present in the massage area: recent injuries, thrombotic diseases, all kinds of inflammations and swellings, and cancer. It is recommended to consult a doctor before having a massage for the treatment of ailments and illnesses.
- If you are dependent on electrical aids e.g. pacemakers, please be sure to consult your doctor for medical advice before taking a massage.
- The packaging material supplied should not be used as a toy.
- The device is not to be used for commercial use or medical purposes
- Do not use with infants, a helpless person or a person insensitive to heat.
- Overlong application can result in skin burns.

Have question ?



Tel: 734-709-6982
Monday-Friday 9:00AM-4:30PM




Email: support@snailax.com

Instruction

Congratulations on your purchase of the Snailax Massager. It combines lighted heat therapy together with 8 massage nodes which rotate in unison to provide you with a deep-penetrating warmth and rejuvenating massage at home or away. Your body will feel more relaxed and rejuvenated.

FEATURES

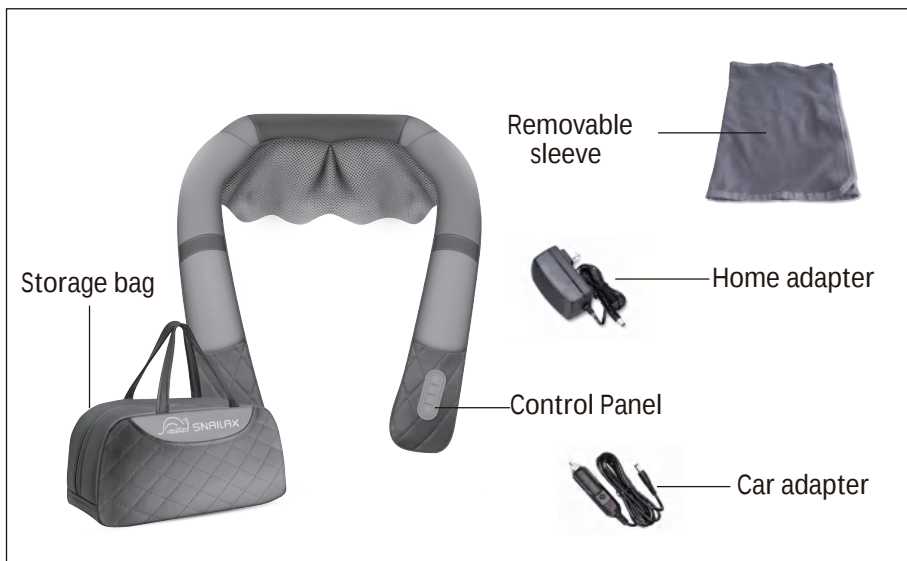

CONTENTS

- Neck & Shoulder Shiatsu massager with heat
- Home adapter
- 12V car adapter
- Removable sleeve
- Storage bag

TECHNICAL DATA

Measurements:	19.81 x 15 inches
Weight:	3.5lbs
Voltage:	Input: AC 100-240V ~50/60Hz Output: 12VDC 1500mA
Nominal Power:	max. 18Watt
Automatic runtime:	15 minutes

Setup and Operation

1. Connect the adapter cable to the corresponding port on the neck massager.
2. Plug the home adapter into an outlet or the car adapter into your car's power receptacle.
3. Turn on the massager and select your desired settings through below controller panel
4. Additional dust-proof sleeve is provided to cover the massage heads for softening purpose when needed.

Power Button:

Press once to activate the massager and shiatsu functions at medium intensity, with automatic reverse rotation and heating feature enabled. Press again to deactivate the massager.

Intensity Button:

Defaults to medium intensity, indicated by fast blinking. Press again to cycle through intensity levels: high (steady light), low (slow blinking), and back to medium (fast blinking).




Direction Button:

By default, the device operates in automatic bidirectional shiatsu mode, switching directions every minute (steady light). Press to manually switch between clockwise and counterclockwise directions (blinking light). Press a third time to return to automatic reversing mode, cycling through this sequence.

Heat Button

When the massager is turned on. Press to turn on or off the heat.


Please use in the manner shown below

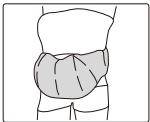


Abdomen




Neck&Shoulder

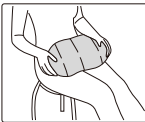


Thigh

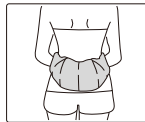

Waist

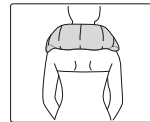

Neck

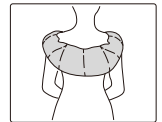

Back



- Remove the adapter from the outlet when not in use.
- Do not setup or use the massager in a bathroom or similiary wet/damp areas.

- For your own safety we equip the massager with a protection against over-heating. Heat is created in the motors when placed under great stress and excessive pressure. Before this could create any risk and after a runtime of 15 Minutes The massager will switch off automatically and will only allow switching on again after an appropriate downtime.
- Likewise, you should allow your body some rest periods. To avoid straining your muscles too much, we recommend to not exceed a continuous massage of 15 Minutes.

Care and Cleaning Instructions

CARE AND MAINTENANCE

Always unplug the unit before cleaning it. Never touch the power plug with wet hands. Before cleaning, allow unit to cool.

Wipe the unit with a soft, dry cloth. Do not use cloths containing any kind of chemical or alcohol and solvent liquids.

Never submerge any part of the unit in liquid.

If the unit is particularly dirty, use the following instructions:

1. Dampen a cloth in water or a 3%–5% mild detergent solution.
2. Wipe the dirty areas with the wet cloth.
3. Clean the cloth in water and wring it out, then wipe away any remaining detergent off the surface of the unit.
4. Wipe the areas with a soft, dry cloth to dry.
5. Wait for the unit to become fully dry before use.

Place massager in a safe, dry, and cool place. Avoid contact with sharp edges or pointed objects which might cut or puncture the surface.

Removable Sleeve

- Sleeve wraps around unit, covering massage heads, and fastens with loop strip.
- Wash in cold water with mild detergent. Allow to air dry.

TROUBLESHOOTING

Problem	Solution
Massage heads rotating at reduced speed.	The load on the massage heads is excessive. Reduce pressure and try again.
Device switches off suddenly.	Device may have been in use for 15 minutes. Auto-shutoff function activated. Continue using by pressing the 'On' button after resting the unit for 15 minutes.
	Over-heat protection function activated. Device is switched off for your safety. Please rest device for 15-30 minutes before restarting.
Device fails to start.	Check if device has been used for more than 15 minutes. If so, let device rest for 15 minutes before continuing. Device may require longer time to cool in warmer environments.
	Check if adaptor is plugged into socket.
	Check if adaptor output is connected to device power inlet wire.
Jerky motion of massage heads/Vibrations from device.	Excessive pressure may have been applied to the device. Reduce pressure and try again.

WARRANTY

If the product(s) is damaged or defective within a reasonable time, Snailax offers the warranty as below.

Period	Return for refund	Return for replacement	Who paid shipping?
Within 30 Days	Yes	Yes if undamaged	Snailax
31- 90 Days	Yes after being justified	Yes after being justified	Snailax
91 Days - 1 Year	No	Yes after being justified	Buyer
1 Year - 3 Years	No	Only to buyers who have received extended warranty	Buyer

Snailax warrants the product to be free from defects in workmanship and materials, under normal use, for a period of 12 months from the original invoice date, except as noted below. Snailax products can be repaired or replaced within warranty period or returned for refund within 90 days if the products break down in the period of proper use.

The warranty extends only to consumers, NOT to any retailers or resellers. This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. Warranty claims for product defects, unless extended warranty is registered, expire 12 months from the purchase day. Warranty claims CANNOT be processed for items that have exceeded the warranty period. The warranty does NOT cover the replaced, refurbished product. The warranty does have other limitations as below "Warranty Limitations". To obtain warranty service on Snailax products, follow "Warranty Claim Procedures" below.

Warranty Claim Procedures:

1. Any product the customer believes is defective and seeks for warranty, the customer must provide to Customer Support:
 - a) Valid proof of purchase (see below the definition of valid proof of purchase);
 - b) Pictures and/or description are required to claim defective merchandise;
 - c) Model number and Series Number label attached to each product
2. If the warranty claim is justified by Snailax Customer Support, customer(s) will obtain from Snailax a Warranty/Return Authorization Number (RMA Number)
3. Snailax reserves the right to specify that items be returned to the designated warehouse for inspection or be inspected by our representative in the field or be retained by customer.
4. Any claim for defective merchandise returns must be packed in original packaging
5. The item(s) will be repaired or replaced or a refund credit will be issued once the item(s) are delivered. When returning items with RMA# and a prepaid shipping label provided by Snailax, Snailax takes responsibility for any damage or loss incurred in transit.

Valid proof of purchase under this Warranty:

Order numbers from online purchases made through Snailax official store (www.snailax.com) or Snailax's direct-authorized resellers, currently only at Amazon/eBay/Walmart. For product(s) obtained from other Channels, please contact the customer service of the Channel directly for appropriate solution.

WARRANTY

Warranty Limitations(this warranty does not include):

Normal wear and tear; Any condition resulting from other than ordinary residential wear or any use for which the product was not intended, such as use in rental or contract trade or commercial use; Any condition resulting from incorrect or inadequate maintenance or care; Damage resulting from misuse, abuse, negligence, accidents or shipping damage, exposed to environmental conditions temperatures, water, improper operation,improper installation, improper use of electrical supply; transportation damage; theft,vandalism or;Damage resulting from the attachment or modification of any unauthorized accessories;Fail to show the valid proof of purchase;Lost or stolen products;Items that have expired their warranty period; No quality-related issues (after 30 days of purchase);Free product/ Gift offer by Snailax; unauthorized repairs or modification; Replaced or refurbished product.Dissatisfaction due to buyer's remorse (after 30 days of purchase).

Refund Procedures:

(Refer to <https://www.snailax.com/pages/refund-policy> for detail instructions)

Warranty Disclaimer

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. SNAILAX SHALL HAVE NO LIABILITY FOR, AND EXPRESSLY DISCLAIMS ANY WARRANTY OR AFFIRMATION OF FACT, EXPRESS OR IMPLIED, OTHER THAN AS SET FORTH IN THIS WARRANTY STATEMENT, INCLUDING, WITHOUT LIMITATION

- (1) THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE;
- (2) ANY WARRANTY OR AFFIRMATION OF FACT RELATED TO MISUSE, IMPROPER SELECTION, RECOMMENDATION, OR MISAPPLICATION OF ANY PRODUCT; AND
- (3) ANY WARRANTY OR AFFIRMATION OF FACT THAT THE CATALOGS, LITERATURE AND WEBSITES IT PROVIDES ACCURATELY ILLUSTRATE AND DESCRIBE PRODUCTS.

Snailax reserves all the right for final interpretation of these Terms of Warranty and Return Policy above and the aright to change, modify, add or remove portions of these terms at any time without prior notification.

Extend Warranty for Free

1) Extend your warranty from 1 year to 3 years by visiting our Snailax online service center at <http://support.snailax.com> or by scanning the QR code below:



2) Alternatively, like our Snailax Facebook page at <https://facebook.com/snailax> or scan the QR code below, to message us with your purchase proof



Your application for a free warranty extension must be completed within 60 days of purchase.

ORDERING REPLACEMENT PARTS

To order additional parts and attachments ,visit <http://www.snailax.com>

FCC STATEMENT

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio / TV technician for help.

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

WARNING: Changes or modifications to this unit is not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Know more about Snailax



SNAILAX CORPORATION



3914 Bestech Dr, Suite A, Ypsilanti, MI 48197, USA

Tel: 734-709-6982 (Monday-Friday 9:00AM-4:30PM)

Email: support@snailax.com

Web: www.snailax.com

*Please provide your order number before contacting customer support