soocas

Contents

Soocas NEOS II Ultra Electric Toothbrush User Manual(EN)	01
Soocas NEOS II Ultra Elektrische Zahnbürste Bedienungsanleitung(DE)	13
Manuel de l'utilisateur de la brosse à dents électrique NEOS II Ultra de Soocas(FR)	25
Soocas NEOS II Ultra Cepillo de Dientes Eléctrico Manual de Usuario(ES-ES)	37
Soocas NEOS II Ultra Cepillo Dental Eléctrico Manual de Usuario(ES-LA)	49
Manuale d'uso dello spazzolino elettrico Soocas NEOS II Ultra(IT)	61
Instrukcja obsługi szczoteczka elektryczna Soocas NEOS II Ultra(PL)	73
Soocas NEOS II Ultra Elektrische Tandborstel Gebruikershandleiding(NL)	85
97	Soocas NEOS II Ultra(AR) قىيؤا بىره كىلا نازىسألا قاشرف مءءءسم لىلد
Soocas NEOS IIウルトラ電動歯ブラシユーザーマニュアル(JP)	109
Soocas NEOS II 울트라 전동칫솔 사용 설명서(KR)	121
Soocas NEOS II Ultra 電動牙刷使用說明書(ZH-TW)	133
Soocas NEOS II Ultra 電動牙刷使用說明書(ZH-HK)	145

Thank you for choosing the Soocas NEOS II Ultra Electric Toothbrush.

Before using your Soocas NEOS II Ultra Electric Toothbrush for the first time, please read this user manual carefully to ensure proper operation, achieve optimal cleaning results, and keep the manual for future reference.

Product Name: Soocas NEOS II Ultra Electric Toothbrush

Model: T04 Ultra

Package Contents:

- Toothbrush Handle
- Travel Case
- MultiClean Ultra Clean Brush Head x1
- MultiClean Ultra Soft Brush Head x1
- Brush Head Covers
- Magnetic Charger
- Quick Start Guide
- Cleaning Brush for Reservoir

Support E-mail: support@soocas.com

Soocas Official Website: soocas.com

Manufacturer: Shenzhen Soocas Technology Co., Ltd.

Address: Room 1101, 1102, 1103, 1104, Building 2, Chongwen Area, Nanshan iPark, No.3370 Liuxian Avenue, Fuguang Community, Taoyuan Subdistrict, Nanshan District, Shenzhen, China

Made in China

Important Safety Instructions

It is normal to experience mild discomfort or gum bleeding the first time you use this product. This is because your gums may not yet be fully accustomed to the high-frequency vibration and sweeping cleaning motion. For beginners, we recommend starting with the Gentle Mode. This condition may persist for 1–2 weeks after initial use. If bleeding continues for more than two weeks without noticeable improvement, stop using this product immediately and consult a dentist.


WARNINGS


- If you have undergone oral or gum surgery within the past two months, or if you have dental braces, crowns, implants, or similar conditions, consult your dentist before using this toothbrush.
- This product is not intended for use by persons (including children) with reduced physical, sensory, or mental capabilities, or lack of experience and knowledge, unless they are supervised or instructed by a responsible person regarding the safe operation of the product.
- Supervise children to ensure they do not play with this product.
- Children must not perform cleaning or maintenance without supervision.
- Do not immerse the magnetic charger in water or any other liquid.
- Use a 5V power adapter for charging.
- Ensure the charging port and power cord are completely dry before charging, as moisture may cause electric shock or short circuit.
- Do not place metal objects on the charger surface while it is powered on.
- Place the product on a stable surface during charging.
- Charge in environments between 0°C and 45°C. The product cannot be charged in excessively low or high temperatures.
- This product contains a built-in rechargeable battery. Keep it away from fire and avoid charging in high-temperature conditions.
- Do not strike or drop this product. If the product is damaged due to collision or falling, stop using it immediately.
- If the toothbrush (brush head, handle, or charger) is damaged in any way, stop using it immediately.
- Do not clean the electric brush head, handle, or charger in a dishwasher.
- The Soocas NEOS II Ultra Electric Toothbrush complies with electromagnetic safety standards. If you have a cardiac pacemaker or other implanted device, consult your physician or the device manufacturer before use.
- If you have any medical condition, consult your physician before using the Soocas NEOS II Ultra Electric Toothbrush.
- This appliance is designed solely for cleaning teeth, gums, and tongue. Do not use it for any other purpose. If you experience discomfort or pain, stop using it immediately and contact your physician.
- Do not use brush heads with bristles that are frayed or bent. It is recommended to replace brush heads every 3 months or sooner.
- If your toothpaste contains peroxide, baking soda, or bicarbonate (common in whitening toothpastes), thoroughly clean the brush head with soap and water after each use to help prevent possible plastic degradation.

Notes


- The product and accessory images shown in this manual are for reference only. The actual product may differ due to updates or upgrades. Please refer to the actual product.
- This product is intended for household use only and must not be used in commercial settings.

Your NEOS II Ultra Electric Toothbrush

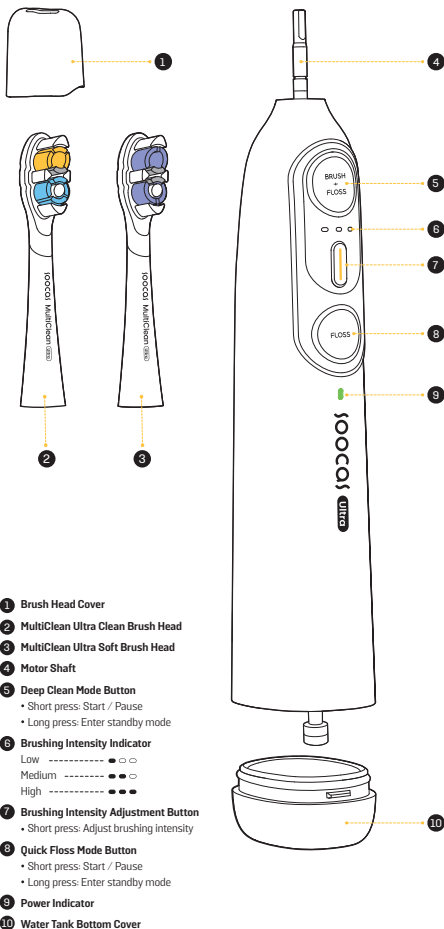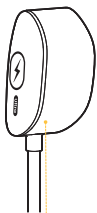
11

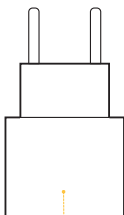12

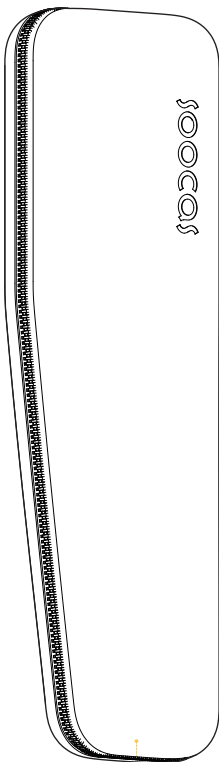13

11 Magnetic Wireless Charging Base

12 Adapter (specific models only)

Note: The box contents may vary depending on the model purchased.

13 Travel Case

Soocas NEOS II Ultra Electric Toothbrush

Name	Soocas NEOS II Ultra Electric Toothbrush
Model	T04 Ultra
Rated Power	12W
Rated Voltage	3.6V
Rated Input	5V ---1A
Charging Time	≈8 hours
Power Consumption	≈ 0.3W
Time To Enter the Condition Off Mode	≈ 20min
IP Rating	IPX8
Battery	Rechargeable
	3.6V --- , 3000mAh, 10.8Wh
Operating Time (full to empty)	30 days+
Battery Type	Lithium-ion
Number of Battery	1
Weight of Battery	53g Max
Capacity of Battery	3000mAh

Magnetic Wireless Charger

Name	Soocas Magnetic Wireless Charger
Model	T04 Ultra
Rated Input	5.0V --- 1.0A
Charging Method	Wireless charging
Maximum Transmission Power	-9.09dBuA/m@10m
Frequency Band	110.76-115.65kHz

This device has been tested and meets applicable limits for radio frequency (RF) exposure.



Power Adapter (Specific Models Only)

Power Adapter		
Input: 100-240 V~, 50/60 Hz, 0.2 A		Output: 5.0 V ---1.0 A, 5.0 W
Information Published	Value and Precision	Unit
Manufacturer's Name	Jiangxi Jian Aohai Technology Co., Ltd.	—
Model Identifier	A520-050100U-EU1	—
Input Voltage	100-240	V~
Input AC Frequency	50/60	Hz
Output Voltage	5.0	V ---
Output Current	1.0	A
Output Power	5.0	W
Average Active Efficiency	73.77	%
Efficiency at Low Load (10%)	64.588	%
No-load Power Consumption	0.075	W

Using Your Soocas NEOS II Ultra Electric Toothbrush



Step 1: Install the Brush Head

- Ensure the bristle direction is aligned with the front of the toothbrush handle, then gently push the brush head onto the motor shaft.
- Note: A small gap between the bristles and the handle is intentional, helping to reduce friction caused by motor vibration.

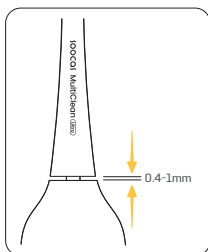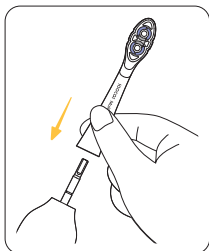Step 2: Fill the Water Tank


- Open the water tank seal and fill with water.
- Close the flip-top water inlet securely after filling to enjoy the simultaneous cleaning experience of brushing and flossing.

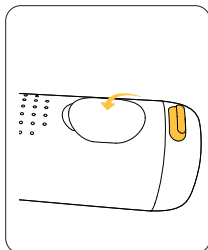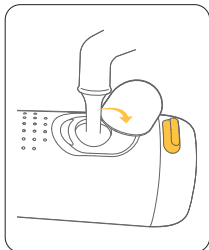Notes


- Use warm water below 40°C (104°F).
- Do not add toothpaste, salt, solid mouthwash tablets, or any chemical substances containing isopropyl myristate into the water reservoir to avoid clogging or damaging the device.

Step 3: Wet the brush head and apply an appropriate amount of toothpaste.

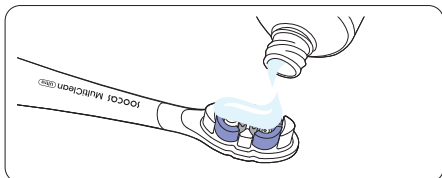
Step 4: Begin Oral Cleaning

1. Select Brushing Intensity

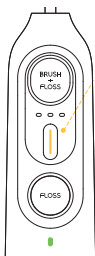Short press the Brushing Intensity Adjustment Button to select your desired setting: Low/Medium/High
Note: If this is your first time using the product, we recommend starting with Gentle mode.


2. Select Cleaning Mode

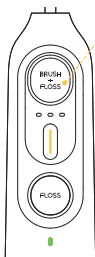
Deep Clean Mode


- Short press: Start / Pause Deep Clean Mode
- Long press: Enter standby mode
- When Deep Clean Mode is activated, the device will perform 2 minutes of brushing followed by 80 seconds of flushing, thoroughly cleaning tooth surfaces and gaps.

Brush your teeth first

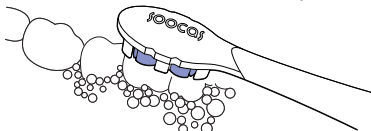Then brush your teeth while flossing

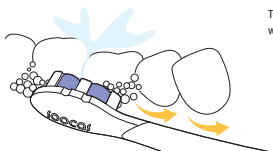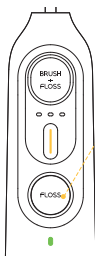
Quick Floss Mode

- Short press: Start / Pause
- Long press: Enter standby mode

When Quick Floss Mode is activated, the device will perform 80 seconds of combined flushing and brushing. Ensure the water tank is filled before activating Quick Floss Mode.

Oral Cleaning Tips

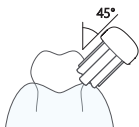
Position the bristles at a 45° angle against the gum line and move slowly. The Soocas NEOS II Ultra Electric Toothbrush will sweep up and down at a 20° angle to ensure thorough cleaning of the entire tooth and gum line.

Brush in short distances horizontally.


Avoid applying excessive pressure, as this may damage your gums and teeth. We recommend moving the brush head horizontally over short distances, allowing the toothbrush's high-frequency vibration and 20° sweeping motion to complete the cleaning.

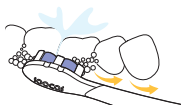When in flossing mode, position the toothbrush vertically toward the gaps between teeth and let the water flow slowly along the gum line.

Daily Maintenance and Cleaning


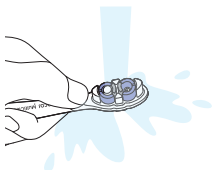
After using the Soocas NEOS II Ultra Electric Toothbrush, rinse the brush head thoroughly with clean water.


Tip: Replace the brush head every 3 months or sooner.

The water tank cover is located at the bottom of the handle. After the Quick Floss Mode ends, pull the yellow latch to detach the water tank cover and pour out any remaining liquid inside.

Reinstall the water tank bottom cover only after the inside of the water tank is completely dry, to reduce the risk of bacterial growth caused by moisture retention.

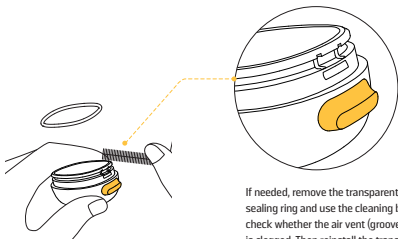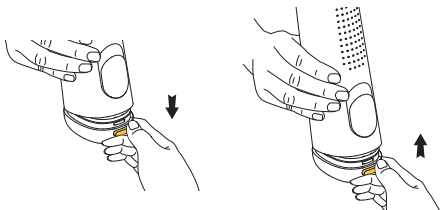
If needed, remove the transparent rubber sealing ring and use the cleaning brush to check whether the air vent (groove area) is clogged. Then reinstall the transparent rubber sealing ring securely.

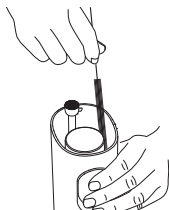Use the cleaning brush regularly to thoroughly clean the inside of the water tank.

Charging

- When the power indicator turns red, the battery level is low. Charge the device promptly to ensure normal operation. A full charge takes approximately 8 hours and provides up to 30 days of use (based on medium mode intensity setting).
- If the power indicator flashes red rapidly and the product cannot be powered on, the battery level is critically low. Charge immediately to restore functionality. A full charge takes approximately 8 hours; actual charging time may vary depending on the remaining battery level.

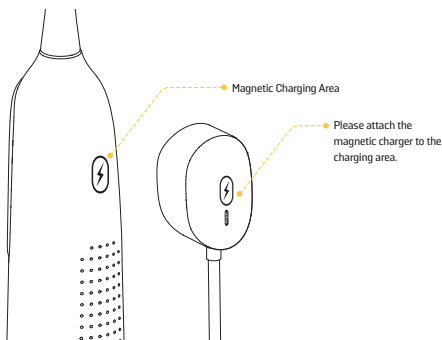Notes

- Do not immerse the magnetic charger in water or any other liquid.
- Use a 5V power adapter when charging.
- Ensure the charging port and power cord are completely dry before charging to avoid electric shock or short circuit.
- Do not place metal objects on the charger surface while it is powered on.
- Place the product on a stable surface during charging.
- Charge in environments between 0°C and 45°C. The product cannot charge in environments with extremely low or high temperatures.
- This product contains a built-in rechargeable battery. Keep it away from fire and avoid charging in high-temperature conditions.

Charging Indicator Lights

	Status	Battery Level
On (In Use/Paused)	Solid Green Light	30%~100%
	Solid Red Light	10%~30%
	Flashing Red Light	<10%
Charging	Solid Green Light	100%
	Breathing Green Light	30%~100%
	Breathing Red Light	<30%

Environmental Protection Policy

When this product reaches the end of its service life, please strictly follow local regulations for disposal:

First, remove the internal battery and recycle it in compliance with applicable standards. After the battery has been properly processed, the remaining components should be disposed of according to local waste classification guidelines.





If you encounter difficulty removing the battery yourself, please dispose of the product in accordance with local waste management regulations.

WARNINGS

- Do not throw the battery into fire to avoid explosion or the release of toxic substances.
- Do not short-circuit the battery to prevent burning and potential fire hazards.
- Do not disassemble any product components unless intentionally disposing of the product, as this may cause damage to the device.
- All procedures must strictly follow safety standards. Disconnect the power supply before removing the battery.
- The battery must be removed before disposal and handled separately.
- Dispose of used batteries only at designated recycling facilities.

How to Remove the Rechargeable Battery

- 1 Remove the bottom cover of the water tank and locate the round screw protection cap at the base of the motor. Use a flat-head screwdriver to pry it open.
- 2 Use a screwdriver to remove the fixing screw at the bottom of the motor.
- 3 Pry open the white plastic housing with a flat-head screwdriver, then detach the entire motor assembly.
- 4 Carefully remove the battery from inside the motor and place it in a designated recycling container.

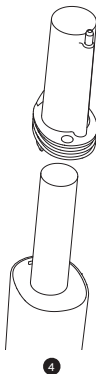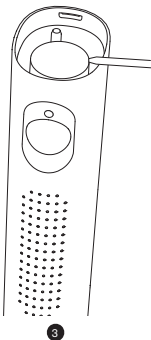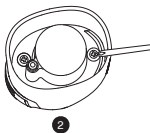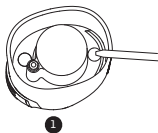Troubleshooting

Issue	Possible Cause	Solution
The electric toothbrush cannot be powered on.	The electric toothbrush is still connected to the charger.	Disconnect the charger.
	Battery level is too low.	Charge the battery until sufficient power is restored.
The brush head cannot be installed.	There are foreign objects or debris in the brush head socket.	Remove any debris from the brush head socket.
	There are foreign objects stuck in the drive shaft.	Clean the drive shaft.
	A non-genuine/incorrect brush head is used.	Use a genuine, compatible brush head.
Brush head vibration is weak.	Brush head is installed incorrectly or tilted.	Reinstall the brush head correctly.

Support E-mail: support@soocas.com
Soocas Official Website: soocas.com

Scan the QR code to
know more about Soocas

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