

A Closer Look


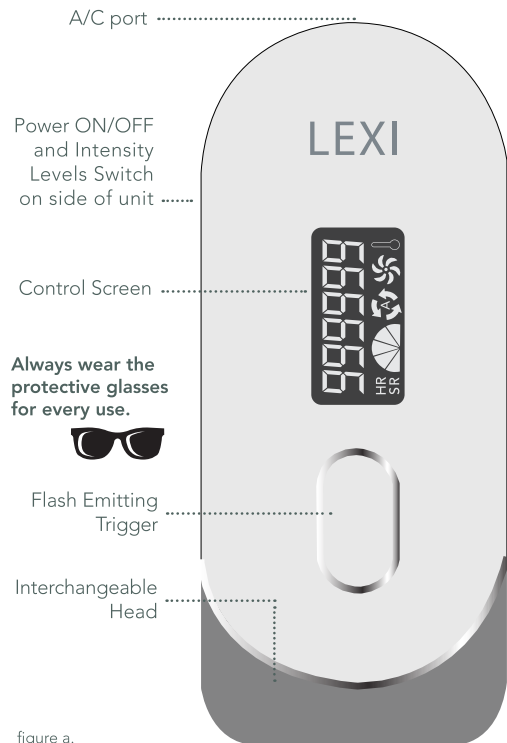


figure a.

Control Screen

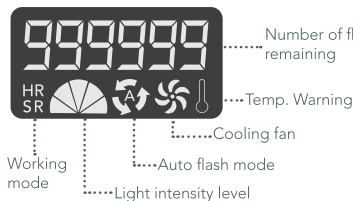

Light Exit Window

Tip:

Long Press the Flash Emitting Trigger to switch to Automatic Continuous Light mode.

Before Using Your Device

IMPORTANT BEFORE FIRST USE

IPL Treatment is designed to help break the cycle of hair growth. Light energy is transferred through the skin's surface and is converted to heat energy (below the surface of the skin), disabling the hair follicle. Treated hairs naturally fall out over the course of a few days to 1-2 weeks. Hair growth is cyclical and follows 3 distinct and consecutive phases. IPL treatment is only effective when hairs are in a particular phase (the Anagen-growth phase). Please follow all recommendations for best results.

Applicable User Reference Chart

		HAIR COLOR							
SKIN TONE									
		✓	✓	✓	✓	✓	✗	✗	✗
		✓	✓	✓	✓	✓	✗	✗	✗
		✓	✓	✓	✓	✓	✗	✗	✗
		✗	✗	✗	✗	✗	✗	✗	✗
		✗	✗	✗	✗	✗	✗	✗	✗

✓ = SUITABLE ✗ = NOT SUITABLE

THERAPEUTIC CYCLE REFERENCE TABLE:

1-4 weeks: Use 2-3 times a week. More than 75% of the hair begins to fall off naturally.

5-6 weeks: Use twice a week. New hair falls off and inhibits growth.

7-8 weeks: Use once a week. The hair becomes thinner and less dense.

More than 8 weeks: Use 1-2 times a week. Use it selectively according to the growth status of the hair.

Usage Guide

PREPARATION:

- First, choose an adequate location to conduct your IPL treatment: do NOT use in the bath or shower.
 - Shave the areas you wish to treat if your hair is longer than 1 mm. Do not use depilatory creams or waxing. Treatment of unshaved or improperly shaved areas can lead to a reduction in the efficacy of your device, skin reactions like redness or swelling and hair can burn resulting in unpleasant smell or pain.
- Note: Once the hair stops growing back, which usually occurs after 4-5 treatments, you no longer have to shave before your treatment.
- Clean your skin and make sure it is hair-free, entirely dry and free from oily substances (perfumes, creams, lotions, deodorants, sunscreens etc.)
 - Clean your device prior to use of trapped hair, dust and fibers.
 - Connect the power adapter to the a/c port of the device and plug into a wall outlet.

READY TO USE:

- Long press the Power switch until the display is turned on [figure a.] The fan will start as soon as the device is turned ON.
- At this point, you may use the main Power button to switch between the 5 Light intensity levels available. You may also turn on the Automatic Continuous Light Mode by Long pressing the Flash Emitting Trigger.

Note: Level 1 = SR mode for skin rejuvenation and hair removal. Levels 2 through 5 = SR and HR modes displayed, for hair removal only.

- Always start with the lowest light intensity level and choose one of the five light intensities in accordance to your comfort level. Reduce the light intensity if you experience any pain or discomfort. Persons with darker skin run a higher risk of skin irritation due to increased light absorption by their skin. Therefore lower intensities are recommended for persons with darker skin. Dark brown to brownish black/darker people should NOT use IPL treatment at all.

- TEST patch:** before your first treatment on each new body area, we recommend testing your skin for potential reactions on a small area. Wait 24 hours following the patch test, to ensure your skin is suitable for IPL treatment.

- Always wear the included safety eye gear to protect your vision from a series of bright flashlights.
- Press the device vertically on the skin at an angle of 90° to ensure that the skin is in full contact [figure b.]
- Press the Flash Emitting Trigger, and the device will send out a flash, accompanied by a slight "crackle sound", you will feel the warm feeling produced by the flash which is normal (if the temperature is too high and/or the display screen Temperature Warning comes on, it means that you selected to high of a light intensity level - it is recommended to always start from the lowest level.)
- Move the handset to the next area to be treated and repeat. Note that your device will need a few seconds to recharge between flashes. Ensure to neither leave gaps between treated areas nor to overlap treated areas [figure c.]
- For optimal results, repeat this treatment every week for 12 weeks to ensure all hair follicles, at all stages of their growth have been treated.

10. ALWAYS UNPLUG YOUR DEVICE AFTER USE.

figure b.


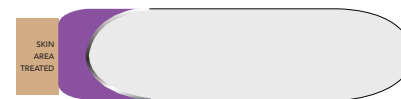


figure c.

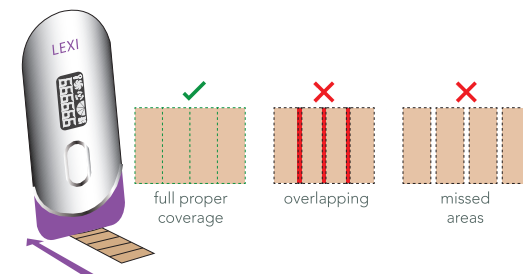

Contraindications & Precautions

DO NOT use this product on dark or very dark skin, as burns could occur;

DO NOT use if subject to bleeding disorder, and those who use isotretinoin, fruit acid (AHA), salicylic acid (BHA) or azalea acid are prohibited to use;

DO NOT use if you are pregnant, lactating, during menstrual periods, if you have body scars, diabetes and/or light-sensitive skin;

DO NOT use if you have severe heart disease, moderate to severe hypertension, malignant tumors, folliculitis, infectious skin disease, burned skin and anticoagulant disorders prohibit usage of this device;

DO NOT use if you have abnormal collagen metabolism, including history of scars and poor wound healing;

DO NOT use if you have a history of vascular diseases, such as varicose veins and vasodilation in the treatment area;

DO NOT use if your skin is sensitive to light and prone to rashes or allergies or psoriasis;

DO NOT use if your treatment area is infected, conditions such as herpes simplex, skin lesions and hematoma prohibit usage;

DO NOT use if subject to suppression of immune response diseases (including HIV infection, lupus erythematosus, porphyria, etc.);

Please use this product with caution for those who used chemical hair removal within 30 days, those who have used skin whitening products in a short time, or those who have used vitamin A ingredients. Do not use with essential oil products at the same time. Please ensure that the skin in the treated area is clean and dry, with no lotion, lubricant and other cosmetics.

DO NOT use while taking aspirin drugs, using immunosuppressive drugs, isotretinoin and other drugs. The recent use of light-sensitive drugs is prohibited (such as tetracycline and sulfa)

DO NOT use for the following body parts or the following conditions: 1/ Intra-oral, genitals, mucosa 2/ Plastic surgery 3/ Throat knot 4/ Eyeball, eyelid, around eyes 5/ Head 6/ Acne area 7/ Tattoo, body painting 8/ Eyebrows, hair 9/ Ear 10/ Nipple, Areola 11/ Inside the navel 12/ On large veins

(Note: This list is not exhaustive if you are not sure whether you can use this product we recommend that you consult a doctor)

During & After Usage Recommendations

DO NOT use antiperspirants immediately after treatment of your underarms;

Please regularly check and clean the dirt on the Light Exit Window so to avoid scalding and failure during use;

DO NOT place the Light Exit Window close to the items that are easy to be demagnetized so as to avoid malfunction and damage;

DO NOT block the fan vents while using, as it may cause overheating of the unit;

DO NOT point the Light Exit Window toward black objects so to avoid the unit running abnormally;

DO NOT use on clothes, long hair, and other items to avoid scorching, discoloration and damage;

DO NOT use on pets so to avoid injury;

DO NOT use this device around your eyes to avoid eye injury;

ALWAYS WEAR PROTECTIVE EYE GEAR (included) while using this device;

After use, please wipe off the hair and skin care products attached to hair removal to void malfunction and skin damage;

DO NOT use near children. Hair removal will emit strong light and easily hurt the child's eyes.

Please set the light intensity according to your skin condition. Be sure to use the first level for the first use;

Skin may become hot and red immediately after use;

DO NOT leave this device unattended after it is powered on.

Please shut off and unplug immediately after use;

After turning off the device, let it stand out for about 2 hours to cool down before storing;

DO NOT use near flammable and explosive materials (including sprayers and gas stations.etc.), which may cause fire and/or electric shock,

DO NOT use this device in an environment with an ambient temperature below 0° C (32° F) or above 45° C (110° F);

DO NOT sunbathe or expose your skin to direct ultraviolet rays within 24 hours after treatment;

During & After Usage Recommendations

DO NOT expose the product to direct sunlight or ultraviolet light for a long time which will seriously damage the product;

If the built-in fan does not work (no fan noise to be heard) after turning on the power, do not use it to avoid malfunctions such as overheating, smoke and fire.

WARNING: Do not use this product near water-filled bathtubs, shower basins, washbasin or other water-containing utensils and do not touch any internal parts to avoid electric shock.

Intended Use

LEXI IPL Hair Removal is an over the-counter device intended for removal of unwanted body and/or facial hair.

Cleaning & Maintenance

In order to ensure the best results and the longest service life of the device, please clean the product before and after each epilation to avoid long-term accumulation of hair and dust, especially in and around the Light Exit Window area.

Before cleaning, the device must be switched off and unplugged to allow the machine to cool.

Use a soft, clean and dry cloth to wipe and clean the mirror surface of the light outlet.

Store the device in a dry, dust-free environment at any temperature between -10 and 60°C.

Never wash the product and accessories in faucets or dishwashers.

Do not use steel balls, abrasive cleaners (toothpaste, etc.) or corrosive liquids (gasoline, motor oil, acetone, etc.) to clean the product.

Note: If the water does not clean well, you can use a few drops of high-concentration alcohol (not more than 96%) to clean the above parts.

Warnings & Warranty

BASIC SAFETY PRECAUTIONS.

When using electrical products, especially when children are present, basic safety precautions should always be taken. IF YOU THINK THE POWER IS OFF WHEN THE SWITCH IS OFF, YOU ARE WRONG. If the plug is in, the power is on. Unplug electrical appliances when you are not using them.

DANGERS.

To reduce the risk of electrocution: DO NOT get the Product or Charging Cable wet or immersed in water. Keep this product away from open flame or flammable aerosol products. The Device is not suitable for use in the presence of a flammable anesthetic mixture with air or oxygen or nitrous oxide.

WARNINGS.

To reduce the risk of burns, electrocution, fire or injury to persons: Close supervision is necessary when this product is used by, on, or near children or physically challenged persons. Use this product only for its intended use as described in this User Guide. DO NOT use this product where oxygen is being administered. The product has no serviceable parts. Never attempt to open the Device Handle for any reason. DO NOT disassemble, heat above 600°C (1400°F) or incinerate.

IMPORTANT CAUTIONS.

This is a medical device. If you are under a doctor's care for any skin condition. if you think you may have a skin condition or if you have very sensitive skin, consult your physician before using this product. Keep away from internal medical devices. This Device is electrical and contains small parts. Keep out of reach of children. THIS PRODUCT IS NOT DESIGNED OR INTENDED PRIMARILY FOR CHILDREN 12 YEARS OF AGE OR YOUNGER.

If you have any medical concerns about using LEXI, please consult your dermatologist before use.

If any of the components are damaged (LEXI Device Handle, or Charging Cord) DO NOT use. DO NOT attempt to open or service components. Email Spa Sciences customer service at: customerservice@spasciences.com

Warnings & Warranty

Warranty:

Your warranty is not effective unless you register the Product within 90 days of purchase. To register go online to www.warrantyspasciences.com and complete the information required for registration.

For a period of ONE (1) year from the date of purchase, Spa Sciences LP (the "Company") warrants the LEXI IPL Permanent Hair Removal System (the "Product") against damage, defects or failure resulting from ordinary consumer use under normal home conditions according to the instructions included with the Product (a "Covered Defect"). This warranty DOES NOT include damage, defects or failure caused by misuse or improper care of the Product; impact or insertion of foreign objects into or on the Product; application of harmful chemicals or abrasives to the Product; failure to follow manufacturer's instructions; normal wear and tear; alteration or overuse of the Product; or any professional or commercial use of the Product. If there is a Covered Defect during the one (1) year warranty period the Company agrees to exchange the Product with a new equivalent Product. This is the sole remedy available for breach of this limited warranty. Any claim for breach of this limited warranty must be commenced within one (1) year following the original purchase. This limited warranty applies only to a purchase of the Product from or an authorized reseller, is valid for the original end user of the Product only and cannot be transferred under any circumstances.

THE PROVISIONS OF THE FOREGOING WARRANTIES ARE IN LIEU OF ANY OTHER WARRANTY, WHETHER EXPRESS OR IMPLIED, WRITTEN OR ORAL (INCLUDING ANY WARRANTY OF MERCHANTABILITY OR FITNESS FOR A PARTICULAR PURPOSE). IN NO EVENT SHALL THE COMPANY BE LIABLE FOR ANY SPECIAL, INCIDENTAL, CONSEQUENTIAL, PUNITIVE, OR INDIRECT DAMAGES BASED UPON ANY LEGAL THEORY. SOME STATES DO NOT ALLOW LIMITS ON WARRANTIES, OR ON REMEDIES FOR BREACH IN CERTAIN TRANSACTIONS. IF APPLICABLE LAW DOES NOT PERMIT THE DISCLAIMER OF IMPLIED WARRANTIES, THEN THE DURATION OF SUCH IMPLIED WARRANTIES, ARE EXPRESSLY LIMITED TO THE DURATION OF THE EXPRESS LIMITED WARRANTY SET FORTH HEREIN.

Additional Info

Entire operating instructions including safety, care and warranty information is available at www.spasciences.com/products/lexi or hover camera over the QR CODE. NO APP NEEDED

