tria[®]

Tria FRX

ADVANCED LASER TECHNOLOGY

User guide & manual.

Contents

Introduction	4
What's in the box?.....	5
Indications for use.....	6
Contraindications	7
Warnings	8
Charging.....	12
How to Use	14
What to Expect	15
Device Care	16
Troubleshooting.....	17
Warranty.....	20
Technical Specifications	21
Français	24
Deutsch	42
Español.....	60
Italiano	78
繁體中文	96
简体中文	112
日本語.....	128
한국어.....	146

Introduction

Pioneers in laser-based skincare for over 20 years, Tria was founded by scientists who set the 'gold standard' for in-clinic laser hair removal. Their mission was to empower you to achieve clinical-grade results with handheld laser devices from the comfort and privacy of your own home.


Our priority has always been to give you excellent results, utilising the best technology. Our multi-award-winning laser devices utilise the same diode laser technology that is preferred by professionals around the world, so you can be confident, knowing that your Tria device is as close to in-clinic technology as possible in a handheld device.

Get clinically proven, skincare results at home with your Tria FRX. 89% of subjects reported smoother skin* and 78% of subjects reported a reduction in fine lines and wrinkles.

* Your Tria FRX uses patented non-ablative, fractional laser technology, working below the skin's surface to treat the multiple signs of aging, leaving you with visibly smoother, younger looking skin.

* Based on subject satisfaction questionnaires completed 4 weeks after finishing a 12-week treatment protocol (5x/week at highest power setting) in an independent clinical trial of 87 participants.

What's in the box?

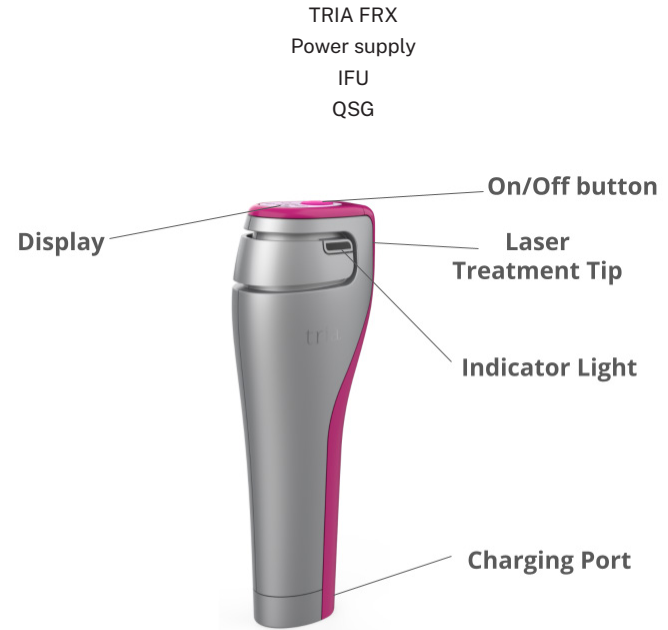Indications for use

The TRIA FRX is indicated for the treatment of periorbital wrinkles (crow's feet), which may result in smoother appearing skin in the treated area. If you are not sure whether the FRX Laser is right for you do not use your device until you have contacted Customer Care or a health care professional for more information.

Contraindications

DO NOT use the device on skin that is infected, burned, or cut as this could cause discomfort or irritation. Wait until the infection, burn, or cut has healed before using the laser.

Anyone with the following medical conditions should not use this device:

- Pimples: the laser may not be effective on pimples
- Veins: the laser may not be effective on veins
- Dyschromia: the laser may not be effective on dyschromia. Consult a doctor before use.
- Freckles: the laser may not be effective on freckles

Important safety information

Read the following important safety information before using the device.

If you have any questions contact The Beauty Tech Group Customer Service or a health care professional.

Warnings

DO NOT use to treat your eye or eyelids or aim the device directly into your eyes. Doing this could cause serious eye injury such as a permanent or temporary decrease in vision or permanent blindness. Treat around the eyes only where you can feel bone beneath the skin.

DO NOT use on tattooed skin or areas with permanent makeup. Doing this could result in a damaged tattoo or increased risk of skin injury.

DO NOT treat a zone on your face more than once per day (once every 24 hours). Doing this could result in increased risk of skin irritation.

After treatment is completed and the device automatically turns off, wait at least 24 hours to begin a new treatment session. Failure to do so may increase skin irritation.

DO NOT use the device if you are under the age of 18. The use of the device on or by anyone under the age of 18 has not been tested and could result in serious skin or eye injury.

DO NOT use if you are pregnant. The device has not been tested for safety on pregnant women and the risks are unknown.

DO NOT use the device if it is wet, or if it or the charger could fall into water while it is plugged in.

If the device or charger has fallen into water, DO NOT touch until you first unplug it. Touching the device while it is wet could result in possibly fatal electric shock.

DO NOT use the device if you experience skin color darkening. Failure to follow these instructions may cause further skin color darkening that may take more than 3 months to diminish.

DO NOT use the device to treat skin which is painful to touch, raised (higher than surrounding skin), crusting or bleeding. This may be a sign of a more serious condition and should be checked by your doctor.

DO NOT treat if your skin is significantly red or swollen from a previous treatment. Wait until the redness and swelling have disappeared before treating again. If you have skin redness or swelling that lasts for more than a week, stop using the device and contact a health care professional.

Invisible laser light. DO NOT view directly with optical instruments. Class 1 Laser Product.

Precautions

- DO NOT treat body areas other than the face. Using this device on other areas of the body has not been tested and the risks are unknown.
- The device has only been clinically tested on facial wrinkles. Using the device on other wrinkle areas may not improve their appearance.
- Applying treatment applications in the exact same spot without moving the device will increase redness and swelling and could potentially lead to blisters or scabs.
- DO NOT use the device if you have a skin condition or are recovering from professional cosmetic treatments, such as a skin peel, skin resurfacing, injectable toxins or fillers.
- If the area of skin for which you are using the device to treat wrinkles appears to be worsening, you should cease use and check with your doctor to make sure you do not have a more serious condition.
- Only store the device in a dry room at a temperature between -20°C and +50°C. (-4 °F and +122°F). Use the device only in a normal room temperature environment.
- DO NOT drop the device on hard surfaces, as this could damage it. Using a damaged device could result in serious skin or eye injury.
- Only use the provided charging cradle and power supply to recharge the device. Using different equipment could damage the device.
- DO NOT attempt to open or repair the device, charging cradle, or power supply.

Charging your device


1. The device may come with some charge, but we always recommend fully charging before first use.
2. Insert the cord into the port on the device. Plug the power supply into an outlet.
3. When your device is charging, the battery indicator will flash.
4. When fully charged, the indicator will turn off.

Sensitivity test


If using the TRIA FRX for the first time, we recommend you test the device on your outer forearm at the treatment level you wish to use. This will help you determine what the treatment feels like. Most users report becoming used to the sensation after several days of treatment. If a reaction such as blistering or scabbing occurs, do not use the TRIA FRX. This type of reaction is very rare.

1. Ensure your skin is clean and dry in the area you wish to treat.
2. Turn on your device by single pressing the power button.
3. Place the treatment window flat against your forearm and slowly glide side-to-side for several seconds.
4. Wait 24 hours and examine the tested area. If the skin appears normal, proceed with treatment.


How to use

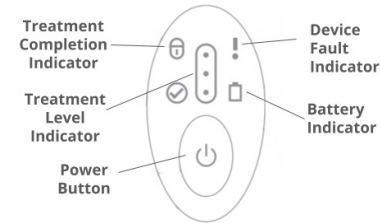1. Wash and dry your face before use.
2. Single press the power button on top of the device to turn it on.
3. Select your treatment level:
 - Training Level 1
 - Treatment Level 2
 - Treatment Level 3

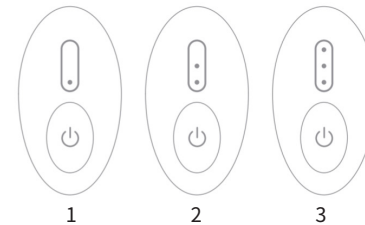4. Hold the treatment window flat against your skin. The indicator light will light up when in good contact with your skin and the device is moving correctly.

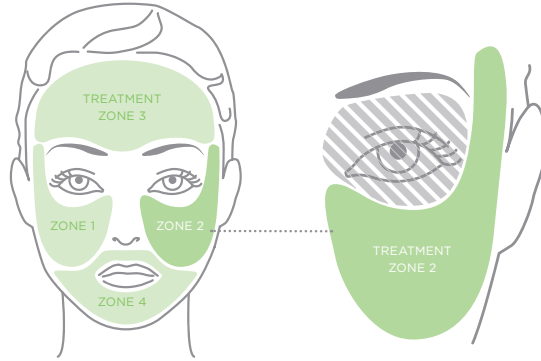
5. Slowly glide the laser side-to-side across your skin. Keep the device moving and in good contact with the skin to keep the indicator light on. See below user guidance image.


CORRECT
SKIN CONTACT


CORRECT
MOTION

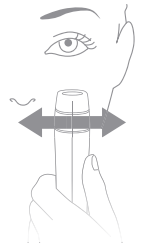X INCORRECT
SKIN CONTACT

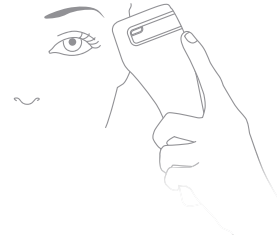X INCORRECT
MOTION

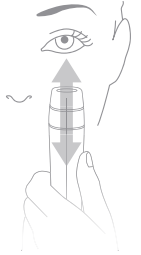6. Treat each zone until your device beeps.

- Level 1: 30 second intervals.
- Level 2: 1 minute intervals.
- Level 3: 2.5 minute intervals.

Your TRIA FRX will beep once after the first zone is treated, twice after the second zone is treated, and three times after the third zone is treated.

7. At the end of the treatment, the device will beep and turn off automatically. The treatment completion indicator will illuminate once complete.

WARNING: DO NOT use to treat your eye or eyelids or aim the device directly into your eyes. Doing this could cause serious eye injury such as a permanent or temporary decrease in vision or permanent blindness. Treat around the eyes only where you can feel bone beneath the skin.

WARNING: Do not treat each zone more than once per day (once every 24 hours). Doing this could result in increased risk of skin irritation.

Locking your device

When your device is locked, it will not power on. The icon on top of the device will illuminate when you try to power on.

To unlock your device:

- Press and hold the power button down until the device beeps 4 times.
- Release the power button.

To lock your device:

- Turn on your device.
- Press and hold the power button until the icon illuminates on top of the device.

Use your TRIA FRX 5 days a week for 12 weeks. Do not use for longer than a 12 week course of treatment.

What to expect

- We recommend making the TRIA FRX treatment part of your evening skincare routine. This allows the skin renewal to begin and side effects to diminish overnight.
- Users typically notice an improvement in their skin in as few as 2 weeks, with the best results occurring a few weeks after the 12 week course of treatment.
- During or immediately after treatment, you may feel mild discomfort such as stinging, itching, or prickling. There can also be a warming or heat sensation similar to that caused by sunburn. These sensations are normal and generally subside within minutes or hours and lessen with continued use.
- During or immediately after treatment, you may also develop temporary redness on your skin. This is normal and typically resolves overnight or within 3 days and lessens with continued use.
- Other side effects that may occur in normal use include:
 - Skin swelling. This generally diminishes overnight.
 - Flaking or dryness.
 - Temporary small brown dots on the skin.

NOTE: If you have more significant or persistent side effects, discontinue treatment and contact Customer Care or a health care professional.

Taking care of your Tria FRX

Cleaning

- For best results, keep the laser treatment tip clean. Oils, lotions, makeup, or other debris can block the laser light or keep the device from functioning properly.
- If the laser window or device tip does become visibly dirty, clean it with a lint-free cloth.
- If other parts of the device become dirty, clean as needed with a soft cloth moistened with rubbing alcohol.
- Avoid getting water or other liquid inside any part of the device.

Maintenance

- The device does not need regular or special maintenance.
- If you think your device is not working properly or is damaged, do not use the device.
- Contact Customer Care for assistance.

Storage

- Keep the device away from water, dampness, or debris that may get inside the device.

Recycling/disposal

- This device contains lithium batteries.
- Please follow the local regulations for recycling or disposing of electrical devices.

Travel

- The charging cradle will operate properly at all common household voltages (any supply between 100 and 240 VAC at 50/60 Hz).

Troubleshooting

My device won't turn on.

- a. Your device may not be charged.
- b. Your device may be locked. See above instructions.

My device won't turn off.

- a. Make sure you have pressed the power button properly.
- b. Your device may be faulty. Please contact Customer Care.

The device fault indicator (!) is lit.

- a. Turn your device off and back on.
- b. If the indicator remains lit, contact Customer Care.

My device won't charge.

- a. Check your cable is properly connected to the device and at the wall outlet.
- b. Try a different wall outlet.
- c. If the battery indicator remains off, your device is already charged and ready for use.

The indicator light won't light up or stops during treatment.

- a. Your device isn't charged.
- b. The treatment window isn't in complete contact with your skin. Ensure it is flat against your skin.
- c. Restart your treatment. If the problem persists, contact Customer Care.

The contact light illuminates intermittently.

- a. The treatment window isn't in complete contact with your skin. Ensure it is flat against your skin.
- b. Clean the laser tip with a dry lint-free cloth.

Warranty

Under normal use this product is guaranteed for 2 years after its original purchase date against manufacturing defects in materials or workmanship. This warranty extends only to the original purchaser. If the product should become faulty within the warranty period, contact The Beauty Tech Group.

- The Beauty Tech Group reserves the right to repair a defective product.
- The Beauty Tech Group reserves the right to replace a defective product with the most compatible product currently available.

This warranty does not cover damage caused by the following:

- Accident, misuse, abuse or alteration
- Servicing by unauthorised persons
- Connection to incorrect power sources
- Use with unauthorised accessories

The Beauty Tech Group shall not be responsible for any incidental, special or consequential damages resulting from the use of this product. All implied warranties, including but not limited to implied warranties of fitness and merchantability are limited in duration to two years from date of original purchase. This warranty gives you specific legal rights, and you may also have other rights that exclusion or limitation of incidental or consequential damages, or limitations on how long an implied warranty lasts, so some of the above limitations and exclusions may not apply to you. This warranty is given in lieu of all other warranties, written or oral, whether expressed by affirmation, promise, description, drawing, model or sample. Any and all warranties other than this one, whether expressed or implied, including implied warranties of merchantability and fitness for a particular purpose are hereby disclaimed. If the product is defective within the warranty period please contact The Beauty Tech Group.

Technical specifications

Model:	SRL
Laser type:	Diode laser
Laser class:	1 According with IEC 60825-1:2014
Note: This Class 1 laser emits invisible laser radiation. Do not treat the eye, eyelids, or point directly at the eye. This product is safe to eyes and skin if used in accordance with these Instructions for Use.	
Output wavelength:	1440 nm
Output laser energy:	5-12 mJ/pulse
Battery charger input:	100-240 V ~ , 50/60Hz, 0.15A
Battery charger output:	5.0 V DC, 1.0 A

Physical

Physical operating temperature:	 5°C to 30°C
Operating relative humidity:	Non-condensing 15-95%
Operating atmospheric pressure:	700 hPa to 1060 hPa
Transport and storage temperatures:	 -20°C to 50°C
Transport relative humidity:	Non-condensing 0-95%

Symbols

Manufacturer

“WEEE (Waste Electrical and Electronic Equipment)”. The waste products should be hand RGB legally.

Refer to instruction manual/booklet.

Protection against ingress of solids, it means the volume of the solids should less than 12mm. Protection against ingress of water, it means the device could protected against vertical falling water drop 15°

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The Beauty Tech Group Ltd.