


Waterpik™ Water Flosser WF-02 Series

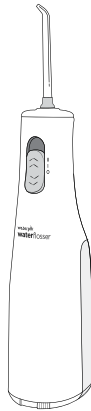

LET US HELP YOU GET STARTED!

We have many helpful tools to provide you with the best possible experience. Scan the QR code or visit: waterpik.com/express.

IMPORTANT SAFEGUARDS

READ ALL INSTRUCTIONS BEFORE USING.

- This product is for household use.

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

DANGER:

To reduce the risk of electrocution:

- Do not place or store product where it can fall or be pulled into a tub or sink.
- Do not place in or drop into water or other liquid.

WARNING:

To reduce the risk of burns, electrocution, fire, or injury to persons:

- This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety.
- Children should be supervised to ensure that they do not play with the appliance.
- Use this product only as indicated in these instructions for use or as recommended by your dental professional.
- Only use tips and accessories that are recommended by Water Pik, Inc.
- Never operate this product if it is not working properly, if it has been dropped or damaged, or dropped into water. Contact Water Pik, Inc. at 1-800-525-2774 or visit www.waterpik.com.
- Fill reservoir with warm water or other dental professional recommended solution only.

- Do not use iodine, saline (salt solution), or water insoluble concentrated essential oils in this product. Use of these can reduce product performance and will shorten the life of the product. Use of these products may void warranty.
- Do not expose the battery or appliance to excessive temperatures.
- Do not use lithium batteries in this product. These batteries are capable of producing pressures that may be uncomfortable to the user and may also result in damage to the appliance. Use of lithium batteries will void the warranty.
- Never drop or insert any object into any opening or hose.
- Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
- Do not direct water under the tongue, into the ear, nose or other delicate areas. This product is capable of producing pressures that may cause serious damage in these areas.
- Do not direct water into the nose or ear. The potentially deadly amoeba, *Naegleria fowleri*, may be present in some tap water or unchlorinated well water and may be fatal if directed into the nose or ear.

SAVE THESE INSTRUCTIONS

INTENDED USE

Waterpik™ water flossers with PRECISIONPULSE TECHNOLOGY are proven to deliver an effective combination of water pressure and pulsations to remove plaque bacteria that can cause gingivitis, cavities, and bad breath.

Symbols Glossary (Not all symbols apply to all products)	
Explanation of symbols used on Water Pik, Inc. product labels and packaging. Applicable symbols for a specific product are found on that product's labels and packaging.	
Symbol	Definition
	Refer to instruction manual
	Manufacturer
IPX7 or IPX4	Ingress Protection Code
	IEC Class II Equipment
	Detachable Power Supply
	Direct Current (DC) Voltage
	Alternating Current (AC) Voltage
	Do not immerse this device in any liquid.
	Do not use this device in a shower, bathtub and do not immerse in any liquid.
	Product contains a lithium ion battery
	CE marking of conformity
	Regulatory Compliance Mark
	Separate collection for waste electrical and electronic equipment.
	Contains a recyclable Li-ion battery
	Contains a recyclable Ni-MH battery

Limited One Year Warranty

Water Pik, Inc. warrants to the original consumer of this new product that it is free from defects in materials and workmanship for 1 year from the date of purchase. Consumer will be required to submit the original purchase receipt as proof of purchase date and if requested, the entire product, to support a warranty claim. Water Pik, Inc. will replace at its discretion the product, provided the product has not been abused, misused, altered or damaged after purchase, was used according to instructions, and was used only with accessories or consumable parts approved by Water Pik, Inc. This limited warranty excludes accessories or consumable parts such as tips, etc. Use with improper voltage will damage product and VOID warranty.

This warranty gives you specific legal rights. You may also have other rights that vary, depending on the law in your place of residence.

WATER PIK, INC. DISCLAIMS ALL OTHER EXPRESS OR IMPLIED WARRANTIES TO THE EXTENT PERMITTED BY LAW INCLUDING, BUT NOT LIMITED TO, ALL WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE. ALL IMPLIED WARRANTIES ARE LIMITED TO THE DURATION OF THIS LIMITED WARRANTY. SOME STATES OR PROVINCES MAY NOT ALLOW LIMITATIONS ON HOW LONG AN IMPLIED WARRANTY LASTS, SO THE ABOVE MAY NOT APPLY TO YOU.

Limitation of Liability

This limited warranty is your sole remedy against Water Pik, Inc. and it shall not be liable for any damages, whether direct, indirect, incidental, special, consequential, exemplary or otherwise, including but not limited to property damage, arising out of any theory of recovery, including statutory, contract or tort.

Arbitration

BINDING ARBITRATION AND WAIVER OF RIGHT TO ASSERT MATTERS IN COURT AND RIGHT TO JURY TRIAL

By purchasing, installing, or using this product you agree that any controversy or claim arising out of or relating to your purchase, use, modification, installation, or removal of this product shall be settled by binding arbitration. The arbitrator's decision will be final and binding. Judgment on the award rendered by the arbitrator may be entered in any court having jurisdiction over the matter. Further, you agree that you will not have the right to litigate any claim arising out of or in relation to your purchase, use, modification, installation, or removal of this product in any court or to have a jury trial on that claim.

Arbitration of any matters shall be administered by the American Arbitration Association and conducted by one arbitrator in accordance with the consumer arbitration rules of the American Arbitration Association in effect at the time of filing of the demand. The venue for all such arbitrations shall be Colorado and any such arbitration shall be subject to the Federal Arbitration Act and Colorado law.

WAIVER OF RIGHT TO PARTICIPATE IN CLASS ACTIONS OR PURSUE CLAIMS IN A REPRESENTATIVE CAPACITY.

Further, by purchasing, installing, or using this product you agree that arbitration must be on an individual basis. This means neither you, nor we, may join or consolidate claims in arbitration by or against other purchasers or users of Water Pik, Inc. products or the general public.

Additionally, you may not litigate in court or arbitrate any claims as a representative or member of a class or in a purported representative capacity on behalf of the general public, other purchasers or users of this product, other persons or entities similarly situated, or in a private attorney general capacity.

FCC STATEMENT (USA)

This device complies with part 15 of the FCC Rules. Operation is subject to the following two conditions:

- (1) This device may not cause harmful interference, and
- (2) this device must accept any interference received, including interference that may cause undesired operation.

Warning: Changes or modifications not expressly approved by Water Pik, Inc. could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation.

This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

waterpik™

Congratulations!

You are now on your way to better oral health!

✓ REGISTER AND SAVE

Register your product and receive over \$50 in savings. [waterpik.com/registration](https://www.waterpik.com/registration)

✓ WE'RE HERE TO HELP

We hope you are loving your new Water Flosser. If for any reason you are not satisfied or have any product questions, please contact our Customer Service Team based in Fort Collins, Colorado.




[waterpik.com/contact](https://www.waterpik.com/contact)


Instructions for Use

GETTING STARTED

1 Remove reservoir.

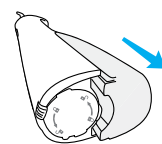


2 To remove battery door, twist, and pull.

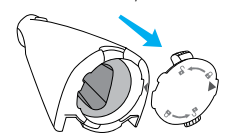


3 Insert 3 AA batteries.




4 To replace battery door, line-up arrows, push and twist.

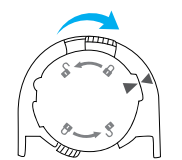

5 Insert tip firmly until it clicks into place.

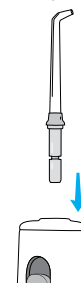


6 If the reservoir was removed, insert reservoir back on the unit. Fill reservoir with warm water.

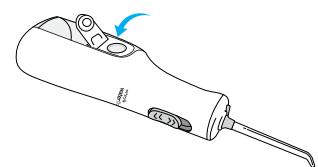

WATER FLOSS


1 Lean over sink and close lips enough to prevent splashing while still allowing water to flow from mouth into the sink.

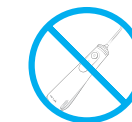

Keep unit upright during use to ensure continuous water flow.

2 Start on low pressure, by sliding up once to I.

Note: For high pressure, slide up to II.

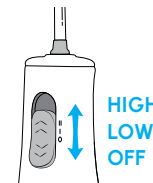


TECHNIQUE

3 Aim water at the gumline at a 90 degree angle. Follow the gumline and pause briefly between teeth.

Tips:
• **LOOK INTO SINK, not mirror, to reduce splashing.**

TO AVOID MESS:
KEEP TIP IN MOUTH WHILE UNIT IS ON.

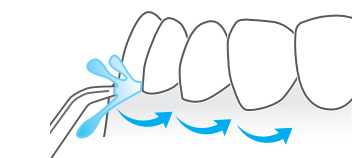

Replacement Parts and Accessories

Orthodontic Tip (OD-100E)
Replace every 3 months

Pik Pocket™ Tip (PP-100E)
Replace every 3 months

Plaque Seeker™ Tip (PS-100E)
Replace every 3 months

Handheld Cordless Precision Tip (JT-450E)
Replace every 6 months

Tongue Cleaner (TC-100E)
Replace every 6 months

Implant Denture Tip (DT-100E)
Replace every 6 months

Full line of replacement tips/attachments may be purchased online at www.waterpik.com, or by phone from Water Pik, Inc. at 1-800-525-2774.

Device Specifications

Rated Input: 4.5 V=
(3X 1.5V AA Alkaline Batteries)

Have Questions or Need More Information?

• **CLEANING**

• **TIPS**

• **TROUBLESHOOTING**

• **COMMON QUESTIONS**

Let Us Help You Get Started!

We have many helpful tools to provide you with the best possible experience. Scan the QR code or visit [waterpik.com/express](https://www.waterpik.com/express).



Service Maintenance

Waterpik™ water flossers have no consumer serviceable electrical items and do not require routine service maintenance. For all parts and accessories go to www.waterpik.com. You may also call our toll-free customer service number at 1-800-525-2774. Refer to the serial and model numbers in all correspondence. These numbers are located on the bottom of your unit.

Still have questions? Go to www.waterpik.com.

Cleaning

Never immerse unit in water. Before cleaning, unplug from the electrical outlet. Clean the product when necessary by using a soft cloth and mild non-abrasive cleanser to wipe exterior.

The water reservoir is removable for easy cleaning and is top rack dishwasher safe. We recommend air drying instead of a heated dry cycle.

It is recommended that the reservoir is removed and cleaned weekly.

Removing Hard Water Deposits / General Cleaning

Hard water deposits may build up in your unit, depending on the mineral content of your water. If left unattended, it can hinder performance.

Cleaning internal parts: add 2 tablespoons of white vinegar to a full reservoir of warm water. Point the handle and tip into sink. Turn unit ON and run until reservoir is empty. Rinse by repeating with a full reservoir of clean warm water.

This process should be done every 1 to 3 months to ensure optimal performance.



[waterpik.com/contact](https://www.waterpik.com/contact)



CONNECT WITH US @WATERPIK



Call: 1-800-525-2774

www.waterpik.com

waterpik™

Water Pik, Inc.,
a subsidiary of
Church & Dwight Co., Inc.
1730 East Prospect Road
Fort Collins, CO 80553-0001 USA
1-800-525-2774
www.waterpik.com