ENGLISH	Upper Arm Blood Pressure Monitor	
	Instructions for use	3
ESPAÑOL	Tensiómetro de brazo	
	Manual de instrucciones.....	17
FRANÇAIS	Tensiomètre au bras	
	Mode d'emploi	32

For Canada only/Solo para Canadá/
Pour le Canada seulement:
Manufacturer/ Fabricante/ Fabricant:
Beurer GmbH, Söflinger Str. 218
89077 Ulm Germany

Distributed by/Distribuido por/
Distribué par :
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020, USA
www.shop-beurer.com

-  **Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.**
-  **¿Preguntas o comentarios? Llame a nuestro servicio al cliente en Estados Unidos gratis al teléfono 1-800-536-0366.**
-  **Questions ou commentaires? Appelez gratuitement notre service clientèle situé aux États-Unis au 1-800-536-0366.**

**READ THIS MANUAL COMPLETELY AND CAREFULLY BEFORE
USING THIS PRODUCT**

Keep this manual in a safe location for future reference

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ENGLISH

Contents

1. Important safety notes	3	8. Displaying and Deleting Measurements	12
2. Intended Use	7	9. Care, Maintenance and Disposal	13
3. Package Contents	7	10. Troubleshooting Guide	13
4. Parts and Controls	7	11. Technical Specifications	14
5. Preparing for Measurement	8	12. Replacement Parts	14
6. Measuring Blood Pressure	10	13. FCC Compliance Information	14
7. Evaluating Results	11	14. Warranty	15

SAVE THESE INSTRUCTIONS FOR FUTURE REFERENCE

If the instruction manual is damaged or if you no longer have the instruction manual in your possession, please contact customer service. See warranty for service contact.

Dear Valued Customer,

Thank you for choosing one of our products. Our name stands for high-quality, thoroughly tested products for application in the areas of heating, weight management, blood pressure, body temperature, pulse, gentle therapy, massage, and beauty.

Register your product today to receive further benefits:

- added protection for your purchase
- health and wellness tips
- new product information
- continual outstanding customer service and support.

Visit www.registerbeurer.com or contact customer service at **1-800-536-0366**.

Please read these instructions for use carefully and keep them for later use.

Make them accessible to other users and observe the information they contain.

With kind regards,

Your Beurer team

1. Important safety notes

Signs and symbols

Whenever used, the following signs identify safety and property damage messages and designate a level of hazard or seriousness.

READ THIS ENTIRE MANUAL, THE SAFETY SECTION AND ALL INSTRUCTIONS AND WARNINGS COMPLETELY AND CAREFULLY BEFORE USING THIS PRODUCT. FOLLOW ALL SAFETY INSTRUCTIONS AND WARNINGS TO AVOID HAZARDOUS SITUATIONS AND TO MAKE CORRECT USE OF THIS PRODUCT.

	This is the safety alert symbol. It alerts you to potential personal injury hazards. Obey all safety messages that follow this symbol to avoid possible injury or death.
WARNING	Indicates a potentially hazardous situation, which, if not avoided, could result in death or serious injury.
CAUTION	Indicates a potentially hazardous situation, which, if not avoided, may result in minor/moderate injury.
NOTICE	Addresses practices not related to personal injury, such as product and/or property damage.
	Product information Note on important information
	Follow instructions for use
	Type BF, applied part
	Direct current
Storage/Transport 	Permissible storage temperature and humidity
Operating 	Permissible operating temperature and humidity
	Atmospheric pressure, limit Indicates the range of atmospheric pressures to which the medical device can be safely exposed
IP22	Device protected against foreign objects ≥ 12.5 mm and against water dripping at an angle
	Serial number

WARNING – To reduce the risk of fire, electric shock, or serious personal injury:

- For consistency, always measure your blood pressure at the same time each day.
- Before each measurement, relax for about five minutes.
- When performing multiple measurements on an individual, wait at least 1 minute between each measurement.
- Do not take a measurement within 30 minutes of eating, drinking, smoking or exercising.
- Repeat the measurement if you are unsure of the results.
- Measurements are for your information only – they are no substitute for a medical examination. Discuss the measurements with your doctor and never base any medical decisions on them (e.g. medicines and their administration).
- Do not use the blood pressure monitor on newborns or patients with pre-eclampsia. Consult your doctor before using the blood pressure monitor during pregnancy.
- This device is not intended for use by people (including children) with restricted physical, sensory or mental

skills or a lack of experience and/or a lack of knowledge, unless they are supervised by a person who has responsibility for their safety or they receive instructions from this person on how to use the device. Supervise children around the device to ensure they do not play with it.

- Cardiovascular diseases may lead to incorrect measurements or have a detrimental effect on measurement accuracy. The same also applies to very low blood pressure, diabetes, circulatory disorders, and arrhythmias as well as chills or shaking.
- The blood pressure monitor must not be used in connection with a high-frequency surgical unit.
- The blood pressure monitor is made from precision and electronic components. The accuracy of the measured values and service life of the device depend on its careful handling:
 - Protect the device from impacts, humidity, dirt, marked temperature fluctuations and direct sunlight.
 - Do not drop the device.
 - Do not use the device in the vicinity of strong electromagnetic fields and keep it away from radio systems or mobile telephones.
 - Only use the cuff included with the delivery or original replacement cuffs. Otherwise incorrect measured values will be recorded.
- We recommend that the batteries be removed if the device is not to be used for a prolonged period.
- Using the blood pressure monitor outside your home environment or whilst on the move (e.g. whilst travelling in a car, ambulance or helicopter, or whilst undertaking physical activity such as playing sport) can influence the measurement accuracy and cause incorrect measurements.
- This device is not intended for use as an ambulatory blood pressure monitor for carrying out long-term measurements at predefined time intervals.
- People with pacemakers or other electrical implants should consult their doctor before using the device.
- Do not use the cuff on people who have undergone a mastectomy.
- Do not place the cuff over wounds as this may cause further injury.
- Make sure that the cuff is not placed on an arm in which the arteries or veins are undergoing medical treatment, e.g. intravascular access or intravascular or therapy, or an arteriovenous (AV) shunt.
- Please note that when inflating, the functions of the limb in question may be impaired.
- During the blood pressure measurement, the blood circulation must not be stopped for an unnecessarily long time. If the device malfunctions remove the cuff from the arm.
- Avoid any mechanical restriction, compression or bending of the cuff line.
- Do not perform measurements more frequently than necessary. The resulting restriction of the blood flow may cause injury.
- Place the cuff on your upper arm only. Do not place the cuff on other parts of the body.
- The device is only intended for the purpose described in these instructions for use. The manufacturer is not liable for damage resulting from improper or careless use.
- The temperature of the cuff can reach up to 107.8 °F (42.1 °C) during normal operation at an ambient temperature of 104 °F (40 °C).

CAUTION – To reduce the risk of personal injury or product/ property damage:

- Remove any kind of upper arm jewelry or the like before taking a measurement. This could cause bruises.
- Do not use the device at the same time as other medical electrical devices (ME equipment). This could lead to a malfunction of the device and/or an inaccurate measurement.
- Do not use the device outside of the specified storage and operating conditions. This could lead to incorrect measurements.
- Small parts may present a choking hazard for small children if swallowed. They should therefore always be supervised.

NOTICE:

- Do not wash the unit in a washing machine, dish washer or tumble dryer.
- Never immerse and or spill water or any other liquid onto the monitor or any components, otherwise liquid will enter it and cause damage.
- Never attempt to repair, open and/or disassemble the unit (including arm cuff and optional accessory) or adjust it yourself. This may damage the unit and impair the functions. If you need to have the unit repaired, please contact our customer service. Please see warranty for service contact. Before submitting any complaint, first check the batteries and replace them if necessary.
- Changes or modifications to the device will nullify the user warranty.
- Do not press the buttons with excessive force or with pointed objects.
- Only use the device on people who have the specified upper arm measurement for the device.

Battery Handling Safety Precautions

- Use only the size and type of batteries specified.
- Be sure to follow the correct polarity when installing the batteries. Reversed batteries may cause damage to the device.
- Do not mix different types of batteries together (e.g. Alkaline and Carbon-zinc or rechargeable batteries) or old batteries with fresh ones. Always replace batteries as a simultaneous set.
- If the batteries in the device are depleted or the device will not be used for a long period of time, remove the batteries to prevent damage or injury from possible battery leakage.
- Do not try to recharge batteries not intended to be recharged; they can overheat and rupture (follow battery manufacturer's directions.)
- Do not dispose of batteries in fire, batteries may explode or leak.
- Clean the battery contacts and also those of the device prior to battery installation.
- Remove discharged batteries from the product and dispose/recycle in compliance with all applicable laws.
- Keep batteries away from children and pets. Batteries may be harmful if swallowed. Should a child or pet swallow a battery, seek medical assistance immediately.
- If your skin or eyes come into contact with battery fluid, rinse the affected areas with water and seek medical assistance.
- Do not disassemble, open or crush the batteries.
- If a battery has leaked, put on protective gloves and clean the battery compartment with a dry cloth.

Notes on electromagnetic compatibility

- The device is suitable for use in all environments listed in these instructions for use, including domestic environments.
- The use of the device may be limited in the presence of electromagnetic disturbances. This could result in issues such as error messages or the failure of the display/device.
- Avoid using this device directly next to other devices or stacked on top of other devices, as this could lead to faulty operation. If, however, it is necessary to use the device in the manner stated, this device as well as the other devices must be monitored to ensure they are working properly.
- The use of accessories other than those specified or provided by the manufacturer of this device can lead to an increase in electromagnetic emissions or a decrease in the device's electromagnetic immunity; this can result in faulty operation.
- Failure to comply with the above can impair the performance of the device.

2. Intended Use

This blood pressure monitor performs non-invasive measurement and monitoring of arterial blood pressure and heart rate on adult patients.

3. Package Contents


1 x Upper Arm Blood Pressure Monitor with integrated easyLock® universal cuff (9.4 in - 15.7 in / 24 cm - 40 cm)

4 x 1.5 V LR03 AAA batteries

1 x Quick Start Guide


1 x Instruction manual

4. Parts and Controls

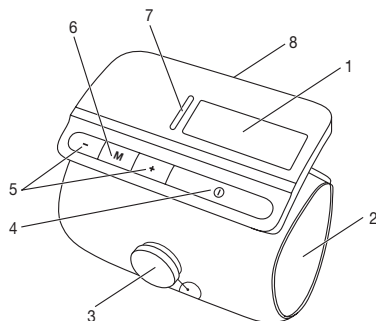1. LCD Display
2. Upper Arm Cuff
3. EasyLock® Button
4. Start/Stop Button
5. +/- Button
6. Memory **M** Button
7. LED Risk Indicator
8. Battery Compartment Cover

LCD Display Information:

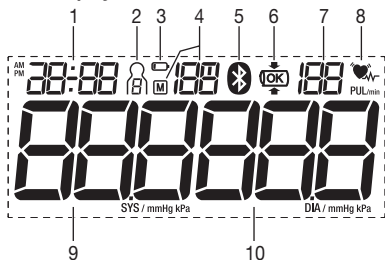
1. Time and Date
2. User Number
3. Low Battery Icon
4. Memory Space Number/Average ($\bar{R}$), Morning Reading (R_M), Evening Reading (P_M)
5. *Bluetooth*® Icon
6. Cuff Tightness Indicator
7. Pulse Rate
8. Pulse/Irregular Heartbeat Indicator
9. Systolic Pressure
10. Diastolic Pressure

5. Preparing for Measurement

Inserting/changing batteries:

1. Remove the Battery Compartment Cover from the rear of the device.
2. Insert four "AAA" alkaline batteries (included) according to the polarity markings inside the compartment. **NOTE:** Do not use rechargeable batteries.
3. Reattach the Battery Compartment Cover.
4. All LCD Display elements will briefly appear and $t2$ h will flash. Set the date and time as described below. Any saved measurements will be retained.

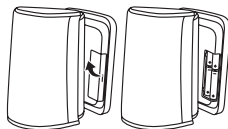If the battery replacement symbol  is permanently displayed, you can no longer perform any measurements and must replace all batteries. Once the batteries have been removed from the device, the date and time must be set again. Any saved measurements are retained.

Setting the hour format, date, and time:

If date and time are not set, stored measurements will not have a valid timestamp associated with them.

Hour format	If $t2$ h is not flashing in the LCD Display, press and hold the M Button for five seconds while the unit is switched off. Select an hour format (12 or 24 hour display) with the +/- Buttons and press the M Button to confirm.	 
Date	The year flashes on the LCD Display. Set the year with the +/- Buttons and press the M Button to confirm. The month flashes on the LCD Display. Set the month with the +/- Buttons and press the M Button to confirm. The day flashes on the LCD Display. Set the day with the +/- Buttons and press the M Button to confirm. NOTE: If the hour format is set to 12h, the day/month display sequence is reversed.	  
Time	The hour flashes on the LCD Display. Set the hour with the +/- Buttons and press the M Button to confirm. The minutes flash on the LCD Display. Set the minutes with the +/- Buttons and press the M Button to confirm.	 
Bluetooth® settings	The Bluetooth® icon flashes on the LCD Display. Activate or deactivate Bluetooth® data transfer with the +/- Buttons and press the M Button to confirm. Note that operation with Bluetooth® activated will increase battery drain.	 
User	The User Memory icon will appear in the LCD Display. Select a User Memory using the +/- and press the M Button to confirm and exit; the device will shut off	

Bluetooth® Setup

In addition to displaying and saving measurements locally on the blood pressure monitor itself, you can also transfer measurements to your mobile device via *Bluetooth*®. To do so, you will need to install the app "beurer HealthManager Pro". The app is freely available in either the Apple App Store or Google Play.

For a list of system requirements and compatible devices, scan the QR Code below:

When connecting for the first time, a randomly generated six-digit PIN code is displayed on the LCD Display of the blood pressure monitor while an input field appears on the smartphone to enter this six-digit PIN code. After you successfully enter the code, the device will be connected to your smartphone.

- Make sure the „beurer HealthManager Pro“ app on your smartphone is always activated and open before you start data transfer from the blood pressure monitor.
- Data transfer is in progress when the *Bluetooth*® symbol appears on the LCD Display.

To transfer the measured values, proceed as follows:

Step 1:

Activate the *Bluetooth*® function on your device as described in section "Setting the hour format, date, time and *Bluetooth*®".

Step 2:

Start the app. Select ELITE 950 in the app and follow the instructions.

Step 3:

Take a measurement.

Step 4:

(Data transfer immediately following the measurement): If the *Bluetooth*® function has been activated, the data transfer to the app "beurer HealthManager Pro" starts automatically after the measurement.

Step 4:





(Data transfer at a later point in time): Go into the memory recall mode for the desired user memory (see "Displaying and Deleting Measurements"). The data transfer starts automatically.

6. Measuring Blood Pressure

The blood pressure monitor needs to be at room temperature when measuring. Measurements can be performed on either arm.

Positioning the Cuff:

1. Place your arm inside the cuff and slide the cuff up to your upper arm about 1 in - 1½ in (2-3 cm) above the bend of the elbow. Align the cuff so that the "ARTERY"  mark is directly over the artery. Ensure blood circulation is not restricted by tight clothing or other objects.
2. Press the easyLock® Button until it clicks into place to activate the locking mechanism. You should be able to fit two fingers under the cuff after it is in place.
3. To adjust the cuff more precisely, turn the easyLock® Button to increase or decrease tightness.

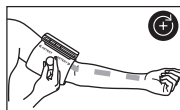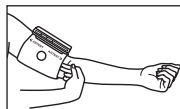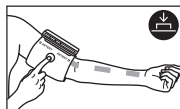
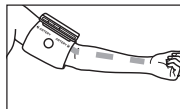
Cuff Position Indicator

This device features a cuff tightness indicator to check the correct positioning of the cuff on your upper arm during measurement. When the cuff is correctly positioned, the  symbol will appear in the LCD Display during measurement and will also display along with measurement results.

If the cuff is positioned incorrectly, the measurement is cancelled and an error message appears in the LCD Display. If this happens, refer to the "Troubleshooting Guide" section for solutions.

Correct Posture:

1. Rest for at least five minutes before each measurement.
2. Perform the measurement either sitting or lying down. Always make sure that the cuff is level with your heart.
3. To perform a measurement, sit comfortably with your arms and back supported. Do not cross your legs and place your feet flat on the floor.
4. Keep still during the measurement and do not talk.

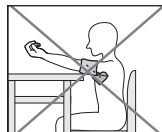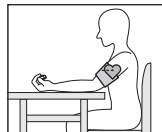
Performing a Measurement:

- To select a desired user memory, press the +/- Buttons until the desired user memory icon appears in the LCD Display. Confirm by pressing the Start/Stop Button; the device will then shut off. Each of the two user memories can store 120 sets of measurements. If you do not select a user memory, the measurement will be stored in the most recently-used user memory.
 - With the cuff attached properly and in a correct posture, press the **Start/Stop** Button. All LCD Display elements will briefly display and measurement will begin after about three seconds. To cancel a measurement at any time, press the Start/Stop Button.
 - As a measurement concludes, air pressure in the Cuff will release. The pulse symbol (♥) will display, along with systolic and diastolic pressure. The LED Risk Indicator will also illuminate, indicating the measurement's risk category.
- 
- 
- "Er_" (where the blank space is a number 0-7) will appear if the measurement was not performed properly. Refer to the "Troubleshooting Guide" section for details about specific error message numbers and then repeat the measurement.
 - If *Bluetooth*® data transfer is active, the unit will transmit measurement data to the app "beurer HealthManager Pro" installed on your paired mobile device. The *Bluetooth*® icon will display while the unit establishes communication with your mobile device and will shut off when a link is successful. After data is transferred, the unit will shut off. Also, if the unit cannot establish a *Bluetooth*® connection, the icon will shut off after about 60 seconds and the unit will shut off after about a minute.
 - Press the **Start/Stop** Button to shut off the blood pressure monitor if needed. The measurement is stored in the selected user memory. The device will also shut itself off after about one minute.
 - Wait at least five minutes before taking another measurement.

Unlocking the cuff and removing the device:

After the measurement has been taken successfully, you can take off the device again. To do so, proceed as follows:

- Pull the easyLock® Button upwards with your fingers to unlock the cuff.
- You can now slightly expand the tube-shaped cuff and then remove it from your upper arm.

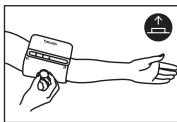7. Evaluating Results

Irregular Heartbeat:

This device can detect irregular heartbeat patterns during measurement and will indicate them by displaying the Irregular Heartbeat Icon . This could be an indicator for arrhythmia, a condition where the heart rhythm is abnormal. These symptoms (omitted or premature heartbeats, slow or excessively fast heart rate) may be caused, among other things, by heart disease, age, physical predisposition, excessive use of stimulants, stress, or lack of sleep. Only your doctor can determine if you do have arrhythmia.

Repeat the measurement if the Irregular Heartbeat Icon displays after a measurement. Rest for at least five minutes between measurements and do not talk or move during measurement. If the icon appears often, contact your doctor, since self-diagnosis and treatment based solely on the test results may be dangerous. It is vital to follow your doctor's instructions.

Risk Classification

Measurements can be classified and assessed by the table below. However, these standard values are only a general guideline, since individual blood pressure varies in different people and different age groups, etc.

Hypertension Category	Systole (in mmHg)	Diastole (in mmHg)	LED Risk Indicator
Grade 3: Severe hypertension	≥ 180	≥ 110	Red
Grade 2: Moderate hypertension	160–179	100–109	Orange
High blood pressure stage 1 (mild)	140–159	90–99	Yellow
High normal	130–139	85–89	Green
Normal	120–129	80–84	Green
Optimal	< 120	< 80	Green

Source: WHO, 1999 (World Health Organization)

It is important to consult your doctor regularly for advice. Your doctor will tell you your individual values for normal blood pressure as well as the value above which your blood pressure is classified as dangerous.

The classification on the display and the scale on the unit show which category the recorded blood pressure values fall into. If the systolic and diastolic values fall into two different categories, the Risk Classification Indicator on the device always shows the higher category.

8. Displaying and Deleting Measurements

User memory

Every successful measurement is stored along with its date and time. If there are more than 120 measurements in a particular user memory bank, the oldest measurement is deleted.

With the device turned off, select a desired user memory by pressing the +/- Buttons until the desired user memory icon appears in the LCD Display. Confirm by pressing the Start/Stop Button; the device will then shut off.


Average Values

With the device turned off, press the **M** Button; “**A**” flashes on the display. The average of all saved measurements in the selected user memory is displayed.


Press the **M** Button again. “**AM**” flashes on the display.

The average value of the morning measurements for the last seven days is displayed (morning: 5:00 AM – 9:00 AM).

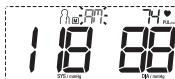Press the **M** Button a third time. “**PM**” flashes on the display.

The average value of the evening measurements for the last seven days is displayed (evening: 6:00 PM – 8:00 PM).

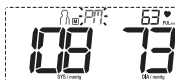Individual Measurements

- Press the **M** Button again and the last individual measurement is displayed.
- To shut the device off, press the **Start/Stop** Button.

You can exit the menu at any time by pressing the Start/Stop Button.

- To clear one of the user memories, you must first select it.
- Follow the instructions described above to select the desired user memory.
- With the device turned off, press the Memory Button once. You should see the Average measurement for the selected user. Press and hold the +/- Buttons for five seconds. All values in the selected user memory will be deleted.

9. Care, Maintenance and Disposal

- Clean the device and cuff with a slightly moistened cloth.
- Do not use detergents or solvents.
- Never immerse the device in water or other liquids.
- When storing the device, make sure no heavy objects are placed on top of it. Remove the batteries. Ensure the Cuff line does not have any sharp kinks.

Disposal

Observe the local regulations for material disposal. Dispose of the device in accordance with local regulations. If you have any questions, please contact the local authorities responsible for waste disposal.

10. Troubleshooting Guide

Error message	Possible cause	Solution
ER1	Pulse not detected.	Please wait one minute and repeat the measurement.
ER2	User moves or speaks during measurement.	Ensure that you do not speak or move during the measurement.
ER3 	Cuff is positioned incorrectly.	Refer to the "Positioning the Cuff" section for details.
ER4	An unknown error occurs during measurement or measured values are outside the measurement range.	Please wait one minute and repeat the measurement. Ensure that you do not speak or move during the measurement. If the error occurs repeatedly, consult a doctor to check you are healthy.
ER5	Cuff does not inflate properly.	Please take another measurement to check whether the cuff can be correctly inflated. Make sure that neither your arm nor other heavy objects are pressing on the line, and that the line is not bent.
ER6 	Batteries are almost depleted.	Replace batteries.
ER7 	Bluetooth® data transfer not completed.	Refer to the "Bluetooth® Setup" subsection for details.
In most of these cases, simply repeat the measurement. Do not move or speak.		

11. Technical Specifications

Measurement Method	Oscillometric, non-invasive blood pressure measurement on the upper arm
Measurement Range	Cuff pressure 0–300 mm Hg, systolic 50–250 mm Hg, diastolic 30–200 mm Hg, Pulse 40–180 beats per minute
Display Accuracy	Pressure ± 3 mmHg, Pulse ± 5 %
Measurement inaccuracy	Max. permissible standard deviation according to clinical testing: systolic 8 mmHg/diastolic 8 mmHg
Memory	2 x 120 memory spaces
Dimensions	5.6 in L x 3.6 in W x 0.9 in H (143 mm L x 91 mm W x 22 mm H)
Weight	9.9 oz (280 g) without batteries
Cuff size	9.4 in - 15.7 in (24 - 40 cm)
Operating Range	41° F - 104° F (5° C - 40° C), 15-90% relative humidity (non-condensing) Atmospheric pressure range of 700 hPa to 1060 hPa
Storage Range	-13° F - 158° F (-25° C - 70° C), $\leq 90\%$ relative humidity (non-condensing) Atmospheric pressure range of 700 hPa to 1060 hPa
Power Supply	4 AAA batteries
Battery Life	Approx. 300 measurements
Expected Service life	at least 10.000 measurements

The serial number is located on the device or in the battery compartment.

12. Replacement Parts

Description	Item number
easyLock® universal cuff (24-40 cm)	164.290

13. FCC Compliance Information

Upper Arm Blood Pressure Monitor ELITE 950

Responsible Party – U.S. Contact Information

Beurer North America LP, 1 Oakwood Boulevard, Suite 255, Hollywood, FL 33020. United States

1-800-536-0366, info@beurer.com

FCC Compliance Statement

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:

(1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

NOTE: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

14. Warranty

Your Beurer Upper Arm Blood Pressure Monitor, model ELITE 950, is warranted to be free from defects in materials and workmanship for the life of the product under normal conditions of intended use and service. This warranty extends only to the original retail purchaser and does not extend to retailers or subsequent owners.

We will, at our option, repair or replace the Beurer Upper Arm Blood Pressure Monitor, model ELITE 950, without additional charge, for any part or parts covered by these written warranties. No refunds will be given. Repair or replacement is our only responsibility and your only remedy under this written warranty. If replacement parts for defective materials are not available, Beurer reserves the right to make product substitutions in lieu of repair or replacement.

For warranty service contact our customer service department at 1-800-536-0366 or at info@beurer.com to provide a description of the problem. If the problem is deemed to be within the scope of the limited lifetime warranty, you will be asked to mail the product at your costs in its original package with proof of purchase, your name, address and phone number. If the problem is not deemed to be within the scope of the limited lifetime warranty, we will provide a quotation for repair respectively replacement and return shipping fee.

This warranty does not cover damage caused by misuse or abuse; accident; the attachment of unauthorized accessory; alteration to the product; improper installation; misapplication; lack of reasonable care with respect to the product; unauthorized repairs or modifications; improper use of electrical/power supply; old worn batteries; normal wear; loss of power; dropped product; malfunction or damage of an operating part as a result of failure to comply with instructions for use or to provide manufacturer's recommended maintenance; transit damage; theft; neglect; vandalism; or environmental conditions; loss of use during the period the product is at a repair facility or otherwise awaiting parts or repair; or any other conditions whatsoever that are beyond the control of Beurer. This warranty is void if the product is ever used in a commercial or business environment. The maximum liability of Beurer under this warranty is limited to the purchase price actually paid by the customer for the product covered by the warranty, as confirmed by proof of purchase, regardless of the amount of any other direct or indirect damage suffered by the customer.

This warranty is effective only if the product is purchased and operated in the country in which the product is purchased. A product that requires modifications or adaptation to enable it to operate in any other country than the country for which it was designed, manufactured, approved and/or authorized, or repair of products damaged by these modifications is not covered under this warranty.

THE WARRANTY PROVIDED HEREIN SHALL BE THE SOLE AND EXCLUSIVE WARRANTY. ANY IMPLIED WARRANTIES, OBLIGATIONS, OR LIABILITIES, INCLUDING BUT NOT LIMITED TO THE IMPLIED WARRANTY OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, ARE LIMITED IN DURATION TO THE DURATION OF THIS APPLICABLE WRITTEN WARRANTY. Some states do not allow limitations on how long an implied warranty lasts, so the above limitations may not apply to you.

IN NO EVENT SHALL BEURER BE LIABLE FOR ANY SPECIAL, INCIDENTAL, INDIRECT OR CONSEQUENTIAL DAMAGES FOR BREACH OF THIS OR ANY OTHER WARRANTY, EXPRESS, IMPLIED OR ANY OTHER THEORY OF LIABILITY, WHATSOEVER. Some states do not allow the exclusion or limitation of special, incidental, or consequential damages, so the above limitation may not apply to you. Beurer does not authorize anyone, including, but not limited to, retailers, the subsequent consumer purchaser of the product from a retailer or remote purchaser, to obligate Beurer in any way beyond the terms set forth herein.

This warranty does not extend to the purchase of opened, used, repaired, repackaged and/or resealed products, including but not limited to sale of such products on Internet auction sites and/or products by surplus or bulk resellers. Any and all warranties or guarantees shall immediately cease and terminate in connection with any products or parts thereof which are repaired, replaced, altered, or modified, without the prior explicitly written consent of Beurer.

This warranty gives you specific legal rights, and you may also have other rights which may vary from state to state.

For more information regarding our product line in the USA, please visit: www.beurer.com

For Canada only:
Manufacturer:
Beurer GmbH
Soefflinger Str. 218
89077 Ulm, Germany

Distributed by:
Beurer North America LP
1 Oakwood Boulevard, Suite 255
Hollywood, FL 33020, USA
www.shop-beurer.com

Questions or comments? Call our US-based customer service toll free at 1-800-536-0366.

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Assemblé en Chine

The *Bluetooth*® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by Beurer North America LP is under license. Other trademarks and trade names are those of their respective owners. Apple and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. App Store is a service mark of Apple LLC., registered in the U.S. and other countries.

Google Play and the Google Play logo are trademarks of Google LLC. Android is a trademark of Google LLC.