Digital Graphic Printer

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EU: Sony Europe B.V.
Da Vincilaan 7-D1, 1930 Zaventem, Belgium
UK: Sony Europe B.V.
The Heights, Brooklands, Weybridge,
Surrey KT13 0XW, United Kingdom
CH: Sony Europe B.V., Hoofddorp,
Schlieren/Switzerland Branch
Wiesenstrasse 5, 8952 Schlieren, Switzerland

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Product name: Digital Graphic Printer

Before operating the unit, please read this manual thoroughly and retain it for future reference.

Indications for Use/ Intended Use

The Sony UP-D899MD printer is a compact, medical grade black and white digital printer. It is designed to be connected to medical imaging systems such as mobile c-arm, ultrasound, cardiac catheterization laboratory and other compatible imaging systems. The printer produces hard copy prints of still images captured by these systems for the patient record or for referrals.

Regulatory Information relevant under the European Medical Device Regulation 2017/745 as amended (MDR)

The printer is intended to produce hard copy prints of still images captured by medical imaging systems for review and to support further clinical decisions.

Notes

- Output images from this equipment cannot be used for diagnostic use.
- This equipment is for medical professionals.
- This equipment is intended for use in medical environments, such as clinics, examination rooms, and operating rooms.

Contraindications

There are no known specific situations that contraindicate the use of this device.

Patient population

Patient population is not subject to any restrictions.

WARNING

To reduce the risk of fire or electric shock, do not expose this equipment to rain or moisture.

To avoid electrical shock, do not open the cabinet. Refer servicing to qualified personnel only.

No modification of this equipment is allowed.

WARNING

To avoid the risk of electric shock, this equipment must only be connected to a supply mains with protective earth.

WARNING

To disconnect the main power, unplug the power plug. When installing the unit, incorporate a readily accessible disconnect device in the fixed wiring, or connect the power plug to an easily accessible socket-outlet near the unit. Do not position the ME equipment where it is difficult to unplug the power plug. If a fault should occur during operation of the unit, operate the disconnect device to switch the power supply off, or disconnect the power plug.

Symbols on the product

Consult the instructions for use

Follow the directions in the instructions for use for parts of the unit on which this symbol appears.

This symbol indicates the manufacturer, and appears next to the manufacturer's name and address.

This symbol indicates the Importer, and appears next to the Importer's name and registered office address.

This symbol indicates the European Community representative, and appears next to the European Community representative's name and address.

This symbol indicates the UK Responsible Person, and appears next to the UK Responsible Person's name and address.

This symbol indicates the Swiss authorized representative, and appears next to the Swiss authorized representative's name and address.

This symbol indicates the medical device in the European Community.

This symbol indicates the date of manufacture.

This symbol indicates the serial number.

This symbol indicates the Unique Device Identifier (UDI), and appears next to the bar code representation of the Unique Device Identification.

This symbol indicates the equipotential terminal which brings the various parts of a system to the same potential.

Storage and transport temperature

This symbol indicates the acceptable temperature range for storage and transport environments.

Storage and transport humidity

This symbol indicates the acceptable humidity range for storage and transport environments.

Storage and transport pressure

This symbol indicates the acceptable atmospheric pressure range for storage and transport environments.

This symbol indicates sharp edges.

For accessories

This symbol indicates the batch code.

This symbol indicates "Do not reuse".

For customers in the U.S.A.

This equipment has been tested and found to comply with the limits for a Class A digital device, pursuant to part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference when the equipment is operated in a commercial environment. This equipment generates, uses and can radiate radio frequency energy and, if not installed and used in accordance with the instruction manual, may cause harmful interference to radio communications. Operation of this equipment in a residential area is likely to cause harmful interference in which case the user will be required to correct the interference at his own expense.

You are cautioned that any changes or modifications not expressly approved in this manual could void your authority to operate this equipment.

All interface cables used to connect peripherals must be shielded in order to comply with the limits for a digital device pursuant to Subpart B of part 15 of FCC Rules.

If you have any questions about this product, you may call;
Sony Customer Information Service
Center 1-800-222-7669 or <http://www.sony.com/>

Supplier's Declaration of Conformity

Trade Name	: SONY
Model	: UP-D899MD
Responsible party	: Sony Electronics Inc.
Address	: 16535 Via Esprillo, San Diego, CA 92127 U.S.A.
Telephone Number	: 858-942-2230

For customers in the U.S.A.

Caution

Federal law (United States of America) restricts this device to sale by or on the order of a licensed healthcare practitioner.

Rx
ONLY

For customers in the U.S.A.

For the glossary of symbols, visit:
<https://rd1.sony.net/help/med/symbols/en/>

For customers in Canada

This unit has been certified according to Standard CAN/CSA-C22.2 No. 60601-1.

For customers in Europe

Any serious incident that has occurred in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and/or patient is established.

Important safeguards and notices for use in the medical environments

1. All devices connected to the unit must be certified or compliant according to IEC 60601-1, IEC 60950-1, IEC 62368-1 and IEC 60065 standards and other IEC/ISO standards applicable to the devices.
2. Furthermore, the system as a whole must comply with IEC 60601-1 standards. All peripheral devices connected to the signal input/output sections of the unit constitute the medical-use system, and therefore, the user is responsible for ensuring that the system as a whole complies with IEC 60601-1 standards. If in doubt, consult qualified Sony service personnel.
3. Connecting the unit to other devices may increase the leakage current.
4. For all peripheral devices connected to the unit that operate on commercial power supplies and do not comply with IEC 60601-1 standards, incorporate an isolation transformer that complies with IEC 60601-1 standards and connect to the commercial power supply via the transformer.
5. The unit generates, uses, and may radiate radio frequency energy. If it is not installed and used in accordance with the instruction manual, it may cause interference on other devices. If the unit causes interference (which can be determined by disconnecting the power cord from the unit), try the following.

- Relocate the unit with respect to the affected devices.
- Connect the unit and the affected devices to different branch circuits.

For more information, consult qualified Sony service personnel. (Applicable standard: IEC 60601-1-2)

Important EMC notices for use in medical environments

- The UP-D899MD needs special precautions regarding EMC and needs to be installed and put into service according to the EMC information provided in the instructions for use.
- The UP-D899MD is intended for use in a professional healthcare facility environment.
- Portable and mobile RF communications equipment, such as cellular phones, can affect the UP-D899MD.

Warning

- Portable RF communications equipment should be used no closer than 30 cm (12 inches) to any part of the UP-D899MD. Otherwise, degradation of the performance of this equipment could result.
- If the UP-D899MD will be used adjacent to or stacked with other equipment, normal operation of the UP-D899MD under such configurations should be verified via observation.
- The use of accessories and cables other than those specified, with the exception of replacement parts sold by Sony Corporation, may result in increased emissions or decreased immunity of the UP-D899MD.

Guidance and manufacturer’s declaration – electromagnetic emissions		
The UP-D899MD is intended for use in the electromagnetic environment specified below. The customer or the user of the UP-D899MD should assure that it is used in such an environment.		
Emission test	Compliance	Electromagnetic environment – guidance
RF emissions CISPR 11	 Group 1	The UP-D899MD uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11 CISPR 32	 Class B	The UP-D899MD is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Harmonic emissions IEC 61000-3-2	 Class A	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	 Complies	

Guidance and manufacturer's declaration – electromagnetic immunity			
The UP-D899MD is intended for use in the electromagnetic environment specified below. The customer or the user of the UP-D899MD should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Electrostatic discharge (ESD)	±8 kV contact	±8 kV contact	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, a relative humidity of at least 30% is recommended.
IEC 61000-4-2	±15 kV air	±15 kV air	
Electrical fast transient/burst	±2 kV for power supply lines	±2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
IEC 61000-4-4	±1 kV for input/output lines	±1 kV for input/output lines	
Surge	±1 kV line(s) to line(s)	±1 kV differential mode	Mains power quality should be that of a typical commercial or hospital environment.
IEC 61000-4-5	±2 kV line(s) to earth	±2 kV common mode	
Voltage dips, short interruptions and voltage variations on power supply input lines	0% U_T (100% dip in U_T) for 0.5/1 cycles ^a	0% U_T (100% dip in U_T) for 0.5/1 cycles ^a	Mains power quality should be that of a typical commercial or hospital environment. If the user of the UP-D899MD requires continued operation during power mains interruptions, it is recommended that the UP-D899MD be powered from an uninterruptible power supply or a battery.
IEC 61000-4-11	40% U_T (60% dip in U_T) for 5 cycles	40% U_T (60% dip in U_T) for 5 cycles	
	70% U_T (30% dip in U_T) for 25/30 cycles ^a (for 0.5 sec)	70% U_T (30% dip in U_T) for 25/30 cycles ^a (for 0.5 sec)	
	0% U_T (100% dip in U_T) for 250/300 cycles ^a (for 5 sec)	0% U_T (100% dip in U_T) for 250/300 cycles ^a (for 5 sec)	
Power frequency (50/60 Hz) magnetic field	30 A/m	30 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
IEC 61000-4-8			

Proximity magnetic fields IEC 61000-4-39	65 A/m 134.2 kHz Pulse modulation 2.1 kHz 7.5 A/m 13.56 MHz Pulse modulation 50 kHz	65 A/m 134.2 kHz Pulse modulation 2.1 kHz 7.5 A/m 13.56 MHz Pulse modulation 50 kHz	Proximity magnetic fields of the UP-D899MD should be at levels characteristic of a typical location in a typical commercial or hospital environment.
NOTE: $\mathcal{U}_T$ is the a.c. mains voltage prior to application of the test level.			
a For example, 10/12 means 10 cycles at 50 Hz or 12 cycles at 60 Hz.			

Guidance and manufacturer’s declaration – electromagnetic immunity			
The UP-D899MD is intended for use in the electromagnetic environment specified below. The customer or the user of the UP-D899MD should assure that it is used in such an environment.			
Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment – guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz outside ISM bands ^c	3 Vrms	Portable and mobile RF communications equipment should be used no closer to any part of the UP-D899MD, including cables, than the recommended separation distance calculated from the equation appliance to the frequency of the transmitter. Recommended separation distance $d = 1.2 \sqrt{P}$
	6 Vrms 150 kHz to 80 MHz in ISM bands ^c	6 Vrms	

Radiated RF	3 V/m	3 V/m	IEC 60601-1-2: 2007
IEC 61000-4-3	80 MHz to 2.7 GHz		$d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.5 GHz IEC 60601-1-2: 2014 + A1: 2020 $d = 2.0 \sqrt{P}$ 80 MHz to 2.7 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, ^a should be less than the compliance level in each frequency range. ^b Interference may occur in the vicinity of equipment marked with the following symbol: 
NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.			
NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.			

- a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the UP-D899MD is used exceeds the applicable RF compliance level above, the UP-D899MD should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating the UP-D899MD.
- b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.
- c The ISM (industrial, scientific and medical) bands between 150 kHz and 80 MHz are 6.765 MHz to 6.795 MHz; 13.553 MHz to 13.567 MHz; 26.957 MHz to 27.283 MHz; and 40.66 MHz to 40.70 MHz.

Recommended separation distances between portable and mobile RF communications equipment and the UP-D899MD

The UP-D899MD is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the UP-D899MD can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the UP-D899MD as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m				
	IEC 60601-1-2 : 2007			IEC 60601-1-2 : 2014 + A1: 2020	
	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 1.2 \sqrt{P}$	800 MHz to 2.5 GHz $d = 2.3 \sqrt{P}$	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 2.7 GHz $d = 2.0 \sqrt{P}$
0.01	0.12	0.12	0.23	0.12	0.20
0.1	0.38	0.38	0.73	0.38	0.63
1	1.2	1.2	2.3	1.2	2.0
10	3.8	3.8	7.3	3.8	6.3
100	12	12	23	12	20

For transmitters rated a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Guidance and manufacturer's declaration – electromagnetic immunity					
The UP-D899MD is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. Portable RF communications equipment should be used no closer than 30 cm (12 inches) to any part of the UP-D899MD. Otherwise, degradation of the performance of this equipment could result.					
Immunity test	Band ^a	Service ^a	Modulation	IEC 60601 test level	Compliance level
Proximity fields from RF wireless communications equipment IEC 61000-4-3	380 – 390 MHz	TETRA 400	Pulse modulation 18 Hz	27 V/m	27 V/m
	430 – 470 MHz	GMRS 460 FRS 460	FM ±5 kHz deviation 1 kHz sine	28 V/m	28 V/m
	704 – 787 MHz	LTE Band 13, 17	Pulse modulation 217 Hz	9 V/m	9 V/m
	800 – 960 MHz	GSM 800/ 900 TETRA 800 iDEN 820 CDMA 850 LTE Band 5	Pulse modulation 18 Hz	28 V/m	28 V/m
	1,700 – 1,990 MHz	GSM 1800 CDMA 1900 GSM 1900 DECT LTE Band 1, 3, 4, 25 UMTS	Pulse modulation 217 Hz	28 V/m	28 V/m
	2,400 – 2,570 MHz	Bluetooth WLAN 802.11 b/g/n RFID 2450 LTE Band 7	Pulse modulation 217 Hz	28 V/m	28 V/m
	5,100 – 5,800 MHz	WLAN 802.11 a/n	Pulse modulation 217 Hz	9 V/m	9 V/m
NOTE: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.					
^a For some services, only the uplink frequencies are included.					

Caution

When you dispose of the unit or accessories, you must obey the laws in the relative area or country and the regulations in the relative hospital regarding environmental pollution.

Warning on power connections

Use a proper power cord for your local power supply.

1. Use the approved Power Cord (3-core mains lead) / Appliance Connector / Plug with earthing-contacts that conforms to the safety regulations of each country if applicable.
2. Use the Power Cord (3-core mains lead) / Appliance Connector / Plug conforming to the proper ratings (Voltage, Ampere).

If you have questions on the use of the above Power Cord / Appliance Connector / Plug, please consult a qualified service personnel.

Warning on power connections for medical use

Customers in the U.S.A. and Canada should use the following type of power cord.

Customers in other countries or regions should use the power cord prescribed by their country or region.

	U.S.A. and Canada
Plug type	HOSPITAL GRADE*

	U.S.A. and Canada
Cord type	Min. Type SJT
	Min. 18 AWG
Minimum rating for plug and appliance couplers	10 A / 125 V
Safety approval	UL Listed and CSA

* Note: Grounding reliability can only be achieved when the equipment is connected to an equivalent receptacle marked "Hospital Only" or "Hospital Grade".

Important

The nameplate is located on the bottom.

Warning Using this unit for medical purposes

The connectors on this unit are not isolated.

Do not connect any device other than one which conforms to IEC 60601-1 standards.

When an information technology device or AV device that uses an alternating current is connected, current leakage may result in an electric shock to the patient or operator.


If use of such a device is unavoidable, isolate its power supply by connecting an isolation transformer, or by connecting an isolator between the connecting cables.

After implementing these measures, confirm that the reduced risk now conforms to IEC 60601-1 standards.

Caution

Do not use the device in a MR (Magnetic Resonance) environment. It may cause a malfunction, fire, and unwanted movement.

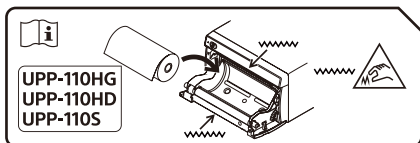Caution

Paper cutters are mounted in the unit. When loading paper, be careful not to touch the paper cutters. Touching the paper cutters may cause injury.

Caution

Do not touch the inside of the door and the patient at the same time. If the unit malfunctions, it may generate voltage that could be harmful to the patient.

Caution

To prevent internal heat built-up, maintain sufficient space around the printer (at least 10 cm (4 in.)). Do not obstruct the ventilation holes on the right and left sides of the printer and on the rear panel.

This equipment is not suitable for use in locations where children are likely to be present.

For customers in Canada (including in the province of Quebec)

- Use the unit with the specified mains voltage. If you use the unit with a different mains voltage from the unit specification, there is a risk of fire or electric shock.
- Damage to the power cord may cause a fire or electric shock.
 - When installing, do not wedge the unit between a wall, rack, shelf, etc.
 - Do not modify or damage the power cord.
 - Do not place heavy objects on the unit or pull the unit.
 - Do not place the unit close to heating appliances or otherwise expose it to heat.
 - When unplugging the power cord, always hold the plug.
- If the unit becomes wet, it may cause an electric shock or fire due to electrical leakage.
- Not using the specified power cord may cause a fire or electric shock.
- Installing the unit on an unstable or inclined surface may cause the unit to fall and result in injury. Make sure that the installation/mounting location is sufficiently strong.
- If you connect/disconnect the power plug with wet hands, it may cause an electric shock.
- Do not step on the unit or place heavy objects on the unit. The unit may tip over, fall, or break, resulting in injury.
- Tripping over the power cord or connection cables may cause the unit to fall or tip over and result in injury. Exercise caution when connecting and installing the unit.
- Failing to insert the power cord plug/connector all the way in straight may cause a fire or electric shock.
- Always remove the power cord and connection cables when moving

the unit. If any connection cable becomes snagged, it may cause the unit to tip over or fall. Moving the unit with cord still connected may damage the cords and result in a fire or electric shock.

FOR CUSTOMERS IN CANADA (INCLUDING IN THE PROVINCE OF QUEBEC)

ALL INSTRUCTIONS AND STATEMENTS WHICH ARE NECESSARY FOR CANADIAN CUSTOMERS ARE PROVIDED IN ENGLISH AND FRENCH. OTHER INSTRUCTIONS AND STATEMENTS NOT PROVIDED IN ENGLISH AND FRENCH ARE NOT FOR CANADIAN CUSTOMERS (INCLUDING IN THE PROVINCE OF QUEBEC).

Reduction in the Use of Hazardous Substances in Electrical & Electronic Equipment (Applicable in Republic of India)

This product and its components, consumables, parts or spares comply with the hazardous substances restriction of India's E-Waste (Management) Rules. The maximum allowable concentrations of the restricted substances are 0.1% by weight in homogenous materials for Lead, Mercury, Hexavalent Chromium, Polybrominated Biphenyls (PBB) and Polybrominated Diphenyl Ethers (PBDE), and 0.01% by weight in homogenous materials for Cadmium, except for the exemptions specified in Schedule II of the aforesaid Rules.

For details on operation, refer to the manuals stored on the supplied CD-ROM.

Using the CD-ROM manuals

The manuals can be viewed on a computer with Adobe Acrobat Reader installed.

You can download Adobe Acrobat Reader for free from the Adobe website.

1. Open the "Instructions For Use" folder on the CD-ROM.
2. Select the language of the manual you want to view.

Note

If you damage or lose the CD-ROM, contact your dealer or Sony service representative.

For the customers in the U.S.A.**SONY LIMITED WARRANTY** -

Please visit www.sony.com/psa/warranty for important information and complete terms and conditions of Sony's limited warranty applicable to this product.

For the customers in Canada**SONY LIMITED WARRANTY** -

Please visit www.sony.com/psa/warranty for important information and complete terms and conditions of Sony's limited warranty applicable to this product.

For the customers in Europe

Sony Europe B.V. - Standard Warranty and Exceptions on Standard Warranty. Please visit <https://pro.sony/support-services/primessupport/support-professional-solutions-europe-standard-product-warranty> for important information and complete terms and conditions.

For the customers in Korea**SONY LIMITED WARRANTY** -

Please visit https://pro.sony/ko_KR/support-services for important information and complete terms and conditions of Sony's limited warranty applicable to this product.

Precautions

On condensation

If the unit is suddenly taken from a cold to a warm location, or if ambient temperature suddenly rises, moisture may form on the outer surface of the unit and/or inside of the unit. This is known as condensation. If condensation occurs, turn off the unit and wait until the condensation clears before operating the unit. Operating the unit while condensation is present may damage the unit.

On consumable parts


The life expectancy of the electrolytic capacitor is about 5 years under normal operating temperatures and normal usage (8 hours per day; 25 days per month).

If usage exceeds the above normal usage frequency, the life expectancy may be reduced correspondingly.

Use with electrosurgical knives and similar devices

If this unit is used together with an electrosurgical knife, etc., the picture may be disturbed, warped or otherwise abnormal as a result of strong radio waves or voltages from the device. This is not a malfunction. When you use this unit simultaneously with a device from which strong radio waves or voltages are emitted, confirm the effect of this before using such devices, and install this unit in a way that minimizes the effect of radio wave interference.

Cleaning the Cabinet

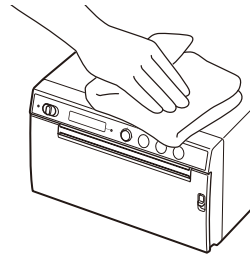Note

Before cleaning the printer, be sure to set the standby/ON switch to standby and then **disconnect the power cord**.

When solvents such as benzene or thinner, or acid, alkaline or abrasive detergent, or chemical cleaning cloths are used on the printer surface, the surface finish may be damaged. Take care with respect to the following:

- Clean the printer surface with a 50 to 70 v/v% concentration of isopropyl alcohol or a 76.9 to 81.4 v/v% concentration of ethanol.
- Stubborn stains may be removed with a soft cloth such as a cleaning cloth lightly dampened with mild detergent solution and then cleaned using the above chemical.
- Do not use unnecessary force to rub the printer surface with a stained cloth. The printer surface may be scratched.
- Do not keep the printer surface in contact with a rubber or vinyl resin product for a long period of time. The surface finish may deteriorate or the coating may come off.

Specifications

Power Requirements	100 V to 240 V ~, 50/60 Hz
Input Current	0.75 A to 0.5 A
Operating Temperature	10 °C to 40 °C (50 °F to 104 °F)
Operating Humidity	20% to 80% (no condensation allowed)
Operating Pressure	700 hPa to 1,060 hPa
Storage and Transport Temperature	-20 °C to +60 °C (-4 °F to +140 °F)
Storage and Transport Humidity	20% to 90% (no condensation allowed)
Storage and Transport Pressure	700 hPa to 1,060 hPa
Dimensions	154 × 88 × 208 mm (w/h/d) (6 1/8 × 3 1/2 × 8 1/4 inches)
Mass	2.0 kg (4 lb 6.5 oz) (printer only)
Thermal head	Thin-film thermal head, 1280 dots
Picture size	(Standard print size) 100 × 75 mm
Printing speed	About 2 sec./image (960 × 1280 dots) (when [HI] is selected from [-SPEED] of [CONFIG.] menu) About 3.8 sec./image (960 × 1280 dots) (when [LOW] is selected from [-SPEED] of [CONFIG.] menu)
Interface	Hi-Speed USB (USB 2.0)

Supplied accessories

Thermal head cleaning
sheet (1)
CD-ROM (including multi-
lingual Instructions for
Use, Readme and
printer driver) (1)
Before Using this Printer (1)
Service Contact List (1)
Information for Customers
in Europe (1)

Optional accessories

Thermal Print Media
UPP-110HG (High
Glossy)
UPP-110HD (High
Density)
UPP-110S (High Quality)

The design and these specifications
are subject to change without prior
notification.

Medical Specifications

Protection against electric
shock:
Class I
Protection against harmful
ingress of water:
Ordinary
Degree of safety in the
presence of a
flammable anesthetic
mixture with air or with
oxygen or nitrous oxide:
Not suitable for use in
the presence of a
flammable anesthetic
mixture with air or with
oxygen or nitrous oxide
Mode of operation:
Continuous

Notes

- Always verify that the unit is operating properly before use. SONY WILL NOT BE LIABLE FOR DAMAGES OF ANY KIND INCLUDING, BUT NOT LIMITED TO, COMPENSATION OR REIMBURSEMENT ON ACCOUNT OF FAILURE TO PRINT CONTENT OF ANY TYPE OR LOSS OF DATA DUE TO FAILURE OF THIS UNIT OR ITS PRINTING MEDIA, ASSOCIATED SOFTWARE, EXTERNAL STORAGE, OR OTHER EXTERNAL DEVICE.
- SONY WILL NOT BE LIABLE FOR DAMAGES OF ANY KIND INCLUDING, BUT NOT LIMITED TO, COMPENSATION OR REIMBURSEMENT ON ACCOUNT OF THE LOSS OF PRESENT OR PROSPECTIVE PROFITS DUE TO FAILURE OF THIS UNIT, EITHER DURING THE WARRANTY PERIOD OR AFTER EXPIRATION OF THE WARRANTY, OR FOR ANY OTHER REASON WHATSOEVER.
- SONY WILL NOT BE LIABLE FOR CLAIMS OF ANY KIND MADE BY USERS OF THIS UNIT OR MADE BY THIRD PARTIES.
- SONY WILL NOT BE LIABLE FOR THE TERMINATION OR DISCONTINUATION OF ANY SERVICES RELATED TO THIS UNIT THAT MAY RESULT DUE TO CIRCUMSTANCES OF ANY KIND.

Trademarks

- Adobe and Adobe Acrobat Reader are trademarks of Adobe Systems Incorporated in the United States and/or other countries.
- The products or system names appearing in this document are trademarks or registered trademarks of their respective owners. Further, the ® or ™ symbols are not used in the text.

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© 2025 Sony Corporation