


CYBER-RELAX

PERFORMANCE

Massage Chair : **OH-8000**

FUJIIRYŌKI

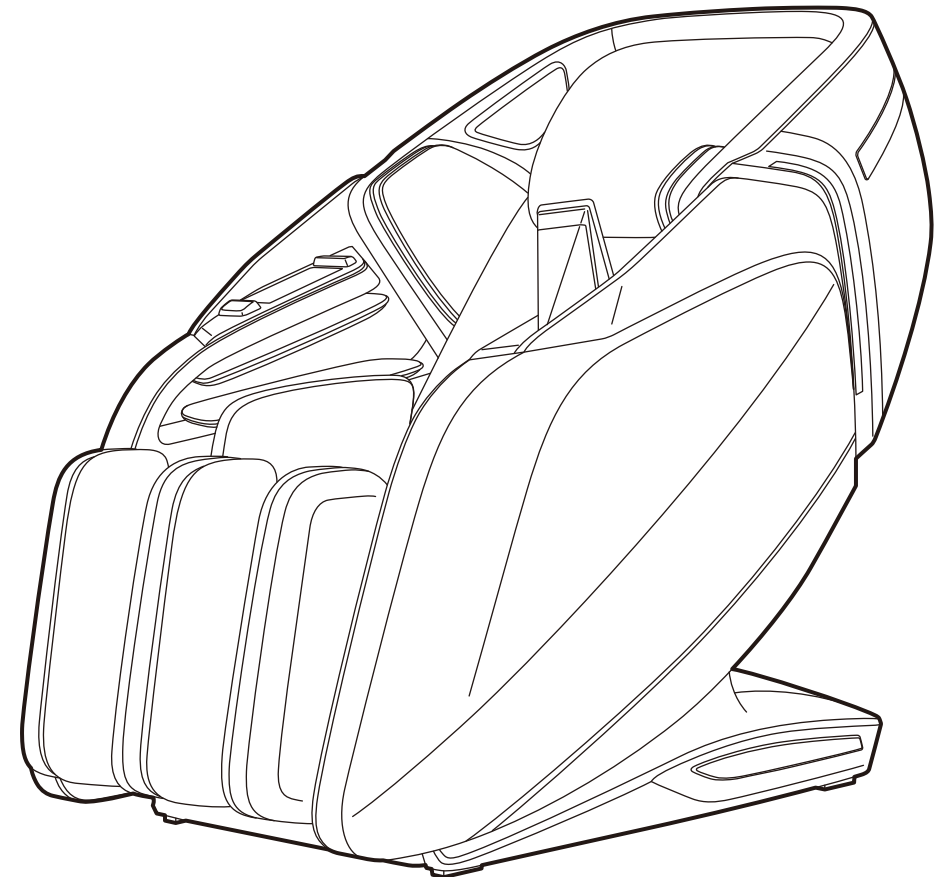

Thank you for purchasing our massage chair. Please read this manual thoroughly before use, paying special attention to the "IMPORTANT SAFETY INSTRUCTIONS" section for proper operation. Keep this manual accessible for future reference.

Contents

- IMPORTANT SAFETY INSTRUCTIONS 4
- Place Massage Chair 12
- Assembly Instruction 14
- Before Use 19
- Names and Functions of Each Part..... 20
- How to Use Auto Course..... 24
- After completing the Massage..... 27
- Manager mode 28
- Language 30
- Troubleshooting..... 31
- Cleaning and Maintenance..... 32
- Product Specifications..... 34

IMPORTANT SAFETY INSTRUCTIONS

- Read these "Important Safety Instructions" before use and use this product properly
- These important safety instructions are prescribed to ensure correct usage of the product, and to prevent injury and damage to the user and other persons.
- Followings are the display and contents.
- Please keep to the regulations as these are important instructions for your safety.

 DANGER	Improper usage may result in serious injury or death.
 WARNING	Improper usage may result in injury or damage to property.*

*Property refers to houses, household belongings, domestic animals and pets.

	 Prohibited actions The prohibited action is illustrated inside the illustration (for example, the illustration on the left refers to "Prohibiting disassembling")
	 Directives, actions which are compulsory The prohibited action is illustated inside the illustration (for example, the illustration on the left refers to "unplugging the power plug")

- Always keep this operation manual where it will be readily available for referral.
- Attention: When disposing of this product and/or its accessories, do so in accordance with local laws and regulations. For details, please refer to your local authority.

 **DANGER**

	The appliance has a heated surface. Persons insensitive to heat must be careful when using the appliance. Pay full attention when it falls under the following. •elderly, ailing person, those who cannot operate by themselves •persons who took medicine that cause drowsiness(sleeping pills, cold medicine) •persons with heavy drinking, heavy tiredness •persons with hematogenous disorder (diabetes, etc.) •persons who recieve a note from the doctor •people with weak temperature sense. There is a risk of low temperature burns, dehydration.
	The following persons must consult a medical physician before using this product. •Persons with implanted medical devices (e.g., pacemakers) that may be affected by electromagnetic interference. •persons with malignant tumors. •persons with heart diseases. •persons who are pregnancy or those immediately after childbirth. •persons who suffer from paresthesia related to peripheral circulatory disturbance due to diabetes etc. •persons with wounds on their skin. •persons in need of repose. •persons with Temperature sensory loss. •persons with body temperatures exceeding 38°C. (when strong acute inflammation symptoms, malaise, hills, blood pressure fluctuations are observed, when debilitated) •persons with herniated discs. •a person with tendovaginitis. •persons who suffer from osteoporosis, spinal column fractures, or acute pains. •persons with abnormalities in the spine, or those whose spines are bent. •Persons with thrombosis, phlebitis, or suspected blood clots, Persons at risk of deep vein thrombosis (DVT). •persons who feel an abnormality on the body not prescribed above, persons in treatment by a doctor.

 **DANGER**

	If the chair stops moving or any abnormalities occur, unplug the power plug and contact the point of purchase for an inspection/repairs. Electrification, short circuits and fire may occur.
	Take care of the movements of the knead balls when massaging around neck. Avoid excessive massaging or massaging the front of the neck. Accidents or injuries may occur.
	When reclining or lowering the footrest, always check underneath and surroundings the footrest for people, pets, and other items. Accidents or injuries may occur.
	When reclining or lowering the footrest, be careful not to get hands, fingers, feet, head caught underneth the footrest or backrest. •underneath the chair itself •backrest, seat, clearance of armrest •clearance of footrest, seat, armrest •underneath the footrest, the back of footrest •clearance between arm and armrest •behind the backrest Accidents or injuries may occur.
	Before use, check underneath the back pad to make sure that there are no rips or tears in the cover cloth of the seat back. (Should there be any rips or tears, no matter how small, immediately discontinue use, disconnect the power cord, and ask for repairs to be conducted). Continued use with rips or tears in the cloth of the seat back may cause electric shocks or injuries.
	Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
	Not long term use of heater and same spot on the body If heat is used on the same spot for a long time, low temperature burns may occur even if the temperature is comparatively low in 40-60°C. About low temperature burn •low temperature burn effect the deep part of the skin and show red specks or blisters. •In case red specks or blisters are found, consult a medical physician immediately.
	The following persons are strongly recommended not to use this product. Certain prohibited from receiving a massage by a medical physician. (for example those suffering from thrombosis, severe aneurysm, acute vericose veins, various kinds of dermatitis and skin infections (including inflammations of subsutaneous tissue), etc.)
	Do not remove, rip or cut the cover of the chair(backrest, footrest, armrest) on purpose. Never use the chair with the cover removed or torn as clothes and hair may get caught in the knead balls.
	Do not use the chair other than on the designated voltage. Fire or electrification may occur.
	Do not use the chair when the power cord or its plug is damaged, or when the receptacle is loose. Contact your point of purchase for repairs if the power cord or power plug are damaged in anyway. Continued usage may result in electric shocks, shortage, or even fire.
	Surely plug the power plug to the root and keep away from the dust. Doing so may cause damage, fire and electrification.
	Do not allow unattended children or persons unable to declare their own intentions to use this chair. Do not allow infants near the chair. Electrification or injuries may occur.
	Do not allow children to play or get on the chair. Accidents or injuries may occur.
Do not connect or disconnect the power plug with wet hands. Electrification may occur.	

IMPORTANT SAFETY INSTRUCTIONS


 DANGER	
	Never use pins or other metallic fasteners with this appliance.
	Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.
	Do not use or store in the bathrooms or places of high humidity or moisture. Electrification, fire, accidents, and moulding may occur.
	Never remodel or decompose/repair by yourself. Fire or injury due to abnormal action may occur.

 WARNING	
	Usage of this product should be limited to 30 minutes for each use. Continual use on the same part of the body should be limited to 5 minutes. Prolonged use may damage muscles and nerves.(Only use for 30 minutes per day.)
	The following persons should consult a medical physician before using this product. ●persons whose muscles have become weak due to increasing age, and those who are of slim build. ●persons suffering from lower back pain caused by bones or internal organs. ●persons easily bruised, susceptible to sprains. ●persons who suffer from intense motion sickness. ●persons who have undergone internal organ or heart operations.
	If any rashes, reddening of the skin, itches or similar symptoms appear after using this product, immediately stop using and consult a physician. Accident and ill health may occur.
	An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts. To disconnect, turn all controls to the off position, then remove plug from outlet. This may result in accidents due to children's mischief.
	Use the chair on a level surface. Do not use outdoors. Damage or accidents may occur.
	Make sure operation buttons and timers are operating normally before use. Accidents or injury may occur.
	Read the operation manual once again and confirm the product is operating normally before using the chair after a prolonged period of non-use. Accidents or injury may occur.
	Confirm no items get caught in the chair (backrest, footrest, armrest) before use. Never drop or insert any object into any opening. Accidents or injury may occur.
	Do not place feet or hands on the footrest when the footrest is retracting. Accidents or injury may occur.
	When any physical abnormalities appear or are felt during operation, stop using the product immediately and consult a physician.
	Sit down slowly after confirming the knead balls are located at the storage position. Accidents or injury may occur.
	Use weak stimulation massage when you start to use. Injury may occur.

 WARNING	
	Make sure the ground terminal is installed. Damages or electrification may occur when ground terminal is not installed.
	Turn the power switch OFF and unplug the power cord immediately in times of power shortages. When power is restored accidents may occur.
	Unplug the power plug by holding the actual power plug and not the cord. This may cause electrification, short circuits, or fire.
	Unplug the power plug when not using the chair. This may cause a fire due to insulation degradation by dust or humidity.
	When maintaining the chair always unplug the power cord. This may cause electrification or injury.
	Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
	Never drop or insert any object into any opening.
	To disconnect, turn all controls to the off position, than remove plug from outlet.
	The appliance is intended for commercial use.
	Temperatures high enough to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
	Do Not Crush - Avoid sharp folds.
	Do not damage the power cord or power plug. ●Do not hurt ●Do not damage ●Do not process ●unreasonably bend or pull. ●Do not twist ●Do not tie ●Do not place heavy items on the power cord. ●Do not get the power cord caught, especially when moving the chair. Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair. Continued usage may result in electric shocks, shortage, or even fire.
	Do not use immediately after eating or consuming alcohol. Nausea, illness, accidents or injury may occur.
	Do not close to stoves or other sources of fire, nor smoke whilst using the chair. Do not use on heated carpets or other geating equipment. Do not operate under blanket or pillow. Excessive heating can occur and cause fire, Electric shock, or injury to persons. Keep the cord away from heated surfaces. This may cause fire.
	Do not stand on the seat, or get on the seat back, footrest, arm or shoulder portion of the chair. The user or chair may rollover and injury or accidents may result.
	Do not stand or sit down on upper part of the backrest, the front of the seat or on the footrest when the chair is reclining. The user or chair may rollover and injury or accidents may result.
	Do not try to sit or get up from the chair when the chair is reclining or the footrest is elevated. The user or chair may rollover and injury or accidents may result.

IMPORTANT SAFETY INSTRUCTIONS

!WARNING	
	Do not use accessories other than those which accompany the chair. Doing so may cause product malfunction.
	Do not use the chair whilst wearing earrings, necklaces, scarves, or other ornaments. This may cause injury.
	Do not use the chair whilst wearing hard objects on the head, such as hair ornaments. This may cause injury.
	Do not use thigh or hip massage with hard objects in pockets. This may cause accidents or injury.
	Do not use while naked. Direct massages on the naked body may damage the skin.
	Do not use on the elbow, knee, head, abdomen areas. Do not get hands or feet caught in the knead balls. Injuries may occur.
	Do not pull out the power cord, or turn off the power switch whilst using the chair. This may cause accidents or injury.
	Do not sleep whilst using the chair. This may cause ill health or injury.
	Do not move the chair while people or items are on the chair. Damage or injuries may occur.
	Make sure no more than 1 person is on the chair at any one time. Doing so may result in damage to the chair or injury.
	Do not wipe with benzine, paint thinner or alcohol, nor spray with any pharmaceutical sprays. This may cause fire or electrification.
	Not to be used for purposes other than massage purposes. This may cause damage or accidents.
	Do not use other therapeutic equipment whilst using this product. This may cause ill health or accidents.
	Do not haul the chair itself or move with the casters on the floor where it easily gets scratches.
	Do not overturn the chair or subject the chair to any strong shocks. Damage or injury may occur.
	Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
	Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
	Do not carry this appliance by supply cord or use cord as a handle.
	Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
	Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.

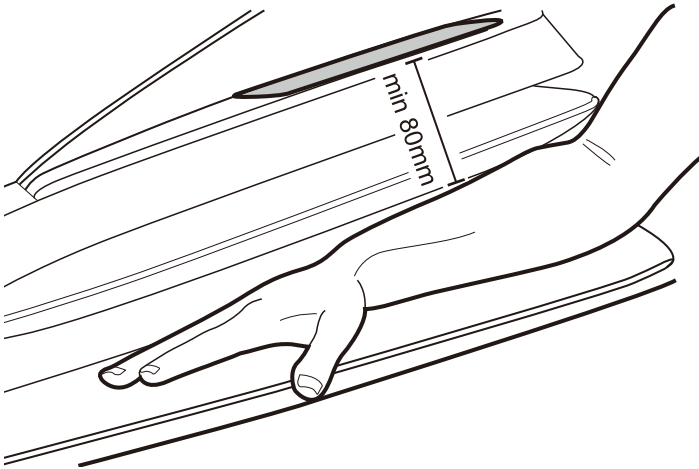!WARNING	
	Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
	Keep the cord away from heated surfaces.
	Use heated surfaces carefully. May cause serious burns. Do not use over insensitive skin areas or in the presence of poor circulation. The unattended use of heat by children or incapacitated persons may be dangerous.
	Keep children away from extended foot support, back support (or other similar parts).
	Connect this appliance to a properly grounded outlet only. See Grounding Instructions.
	This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the appliance by a person responsible for their safety. Children should be supervised to ensure that they do not play with the appliance.
	This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
	Cleaning and user maintenance shall not be made by children without supervision
	Please do not listen at high volume.
	When using the massage chair, the arm should not be placed on the wireless charger. The minimum distance between the arm and the wireless charger is 80mm.
	

IMPORTANT SAFETY INSTRUCTIONS

When using an electrical appliance, basic precautions should always be followed, including the following:

Read all instructions before using this appliance.

DANGER – To reduce the risk of electric shock:

- 1) Always unplug this appliance from the electrical outlet immediately after using and before cleaning.
- 2) Never use pins or other metallic fasteners with this appliance.
- 3) Carefully examine the covering before each use. Discard the appliance if the covering shows any sign of deterioration, such as checking, blistering, or cracking.
- 4) Keep Dry – Do not operate in a wet or moist condition.

WARNING – To reduce the risk of burns, fire, electric shock, or injury to persons:

- 1) An appliance should never be left unattended when plugged in. Unplug from outlet when not in use, and before putting on or taking off parts.
- 2) Do not operate under blanket or pillow. Excessive heating can occur and cause fire, electric shock, or injury to persons.
- 3) Close supervision is necessary when this appliance is used by, on, or near children, invalids, or disabled persons.
- 4) Use this appliance only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
- 5) Never operate this appliance if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the appliance to a service center for examination and repair.
- 6) Do not carry this appliance by supply cord or use cord as a handle.
- 7) Keep the cord away from heated surfaces.
- 8) Never operate the appliance with the air openings blocked. Keep the air openings free of lint, hair, and the like.
- 9) Never drop or insert any object into any opening.
- 10) Do not use outdoors.
- 11) Do not operate where aerosol (spray) products are being used or where oxygen is being administered.
- 12) To disconnect, turn all controls to the off position, then remove plug from outlet.
- 13) Temperatures sufficiently high to cause burns may occur regardless of the control setting. Do not use on an infant or invalid or on a sleeping or unconscious person. Do not use on insensitive skin or on a person with poor blood circulation. Check the skin in contact with the heated area of the appliance frequently to reduce the risk of blistering.
- 14) Do Not Crush – Avoid sharp folds.
- 15) Do not stand on or in appliance. Use only while seated.

16) Keep children away from extended foot support (or other similar parts).

17) Connect this appliance to a properly grounded outlet only. See Grounding Instructions.

18) Do not use massager in close proximity to loose clothing or jewelry.

19) Keep long hair away from massager while in use.

20) This appliance is intended for commercial use.

SAVE THESE INSTRUCTIONS


GROUNDING INSTRUCTIONS

This product must be grounded. If it should malfunction or breakdown, grounding provides a path of least resistance for electric current to reduce the risk of electric shock. This product is equipped with a cord having an equipment-grounding conductor and a grounding plug. The plug must be plugged into an appropriate outlet that is properly installed and grounded in accordance with all local codes and ordinances.

DANGER – Improper connection of the equipment-grounding conductor can result in a risk of electric shock. Check with a qualified electrician or serviceman if you are in doubt as to whether the product is properly grounded. Do not modify the plug provided with the product

– if it will not fit the outlet, have a proper outlet installed by a qualified electrician.


This product is for use on a nominal 120-V circuit and has a grounding plug that looks like the plug illustrated in sketch A in below Figure. Make sure that the product is connected to an outlet having the same configuration as the plug. No adapter should be used with this product.

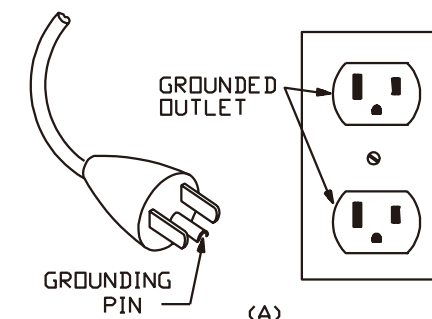

Place Massage Chair

Required space for installation

Ensure there is enough space for the massage chair to recline. The back of the chair should be at least 4 inches (10 cm) from the wall.
Do not expose the chair to direct sunlight or high temperatures, such as near heating equipment.
For long-term use in a fixed location, place a protective mat underneath. The mat should be large enough to cover both the chair's contact area and where the user's feet rest.

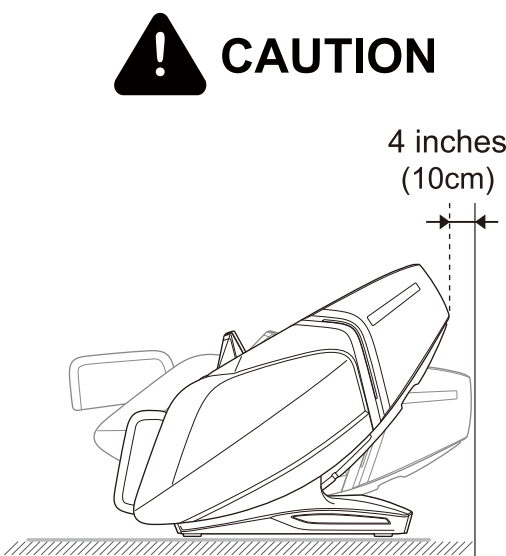

How to lift the massage chair

Grip the bottom sides of the side panels and the rear handle, then lift and move. The chair is heavy—use caution to avoid back strain. (Lifting this machine requires Three or more people.)

1. Please lift the massage chair completely when moving it on fragile surfaces such as wooden floors.
2. When moving on wooden floors, lift the chair to prevent the casters from damaging the surface. Keep your hands and feet clear while lifting, moving, and lowering the chair to avoid serious injury. Do not release the chair until it is completely flat and stable. It is recommended to place protective mats or other protectors on the floor before lowering the chair into position.




How to move the massage chair using casters

Use the rear casters to move the massage chair. According to the method shown in the picture, tilt the chair downwards and backwards about 45° to move it.

Note: Do not press down too much to avoid product damage or injury to the carrying person by falling down the chair. (Carrying this machine using the casters must be carried by two or more persons.)


1. Before moving, you need to shut down the machine correctly, completely restore the backrest and legrest position, and then remove the power supply from the product.
2. The product must be left empty when moving. Do not allow people, pets or debris to stay on the product.
3. When using the casters to move the chair, avoid obstacles exceeding 1 inch (2 cm) in height or 1 inch (2 cm) in width.
4. The casters are not designed for long-distance movement. For distances exceeding 166 feet, move the product slowly and steadily in short segments.

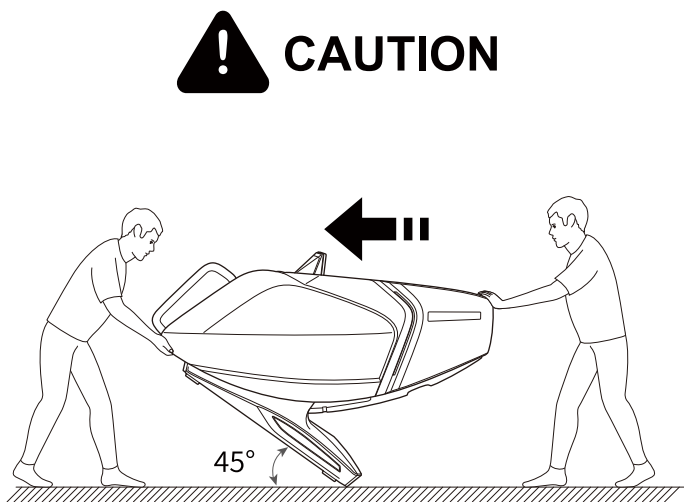


Assembly Instruction

■Accessories List

- Armrest retaining screws x 4pcs (2 for each side)

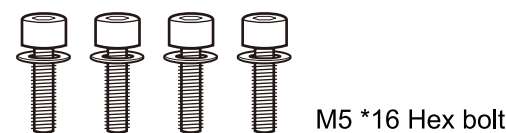

- Legrest retaining screws x 2pcs
- *The screws are attached to the main chair unit

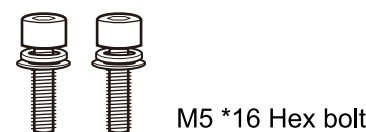


- Hex wrench x 1pcs



- Power cord x 1pc


- Zero-Wall Frame x 1pc

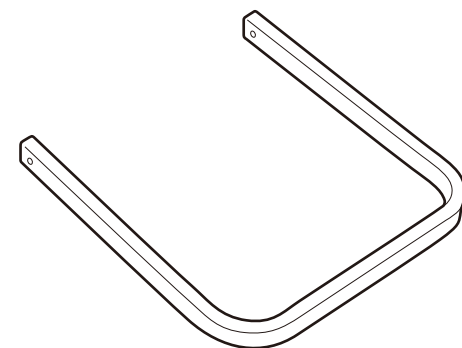

■Assembly instruction for armrest


Step 1

Connect the wiring harness assemblies and air hose from the steel frame to their corresponding slots on the armrest. Ensure all connections are secure.

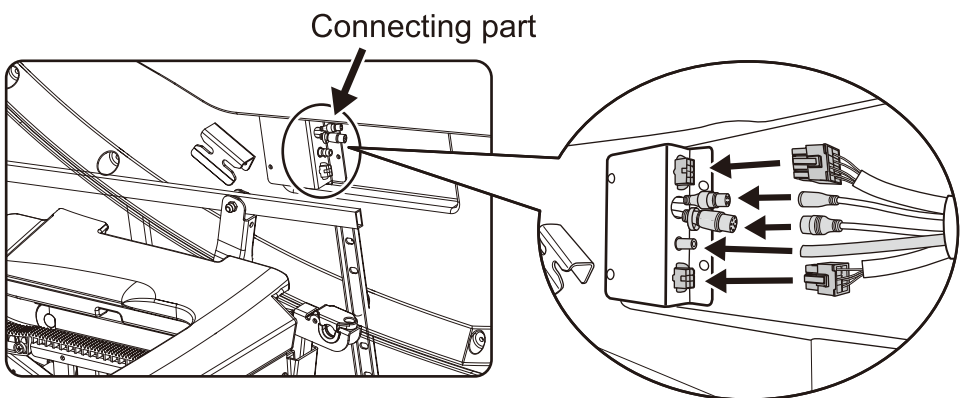

Step 2


Lift and tilt the front of the armrest. Align the rear end (where the roller is located) with the slot on the main body. Push the armrest backward and upward until the roller reaches the end of the slot.

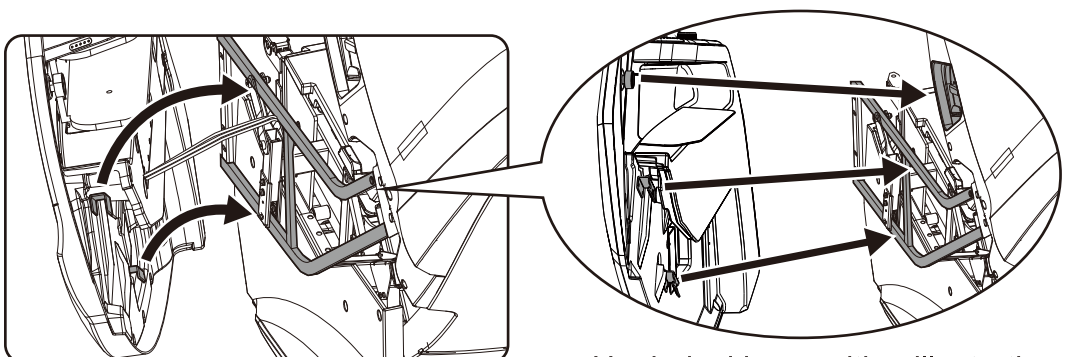

Hook docking position illustration

Step 3



Press the front end of the armrest down so the hook on the armrest engages with the tubing of the steel frame assembly. Align the screw position on the front side of the armrest with the steel frame assembly to confirm it is properly assembled. Use 4 armrest retaining screws to lock it in place.

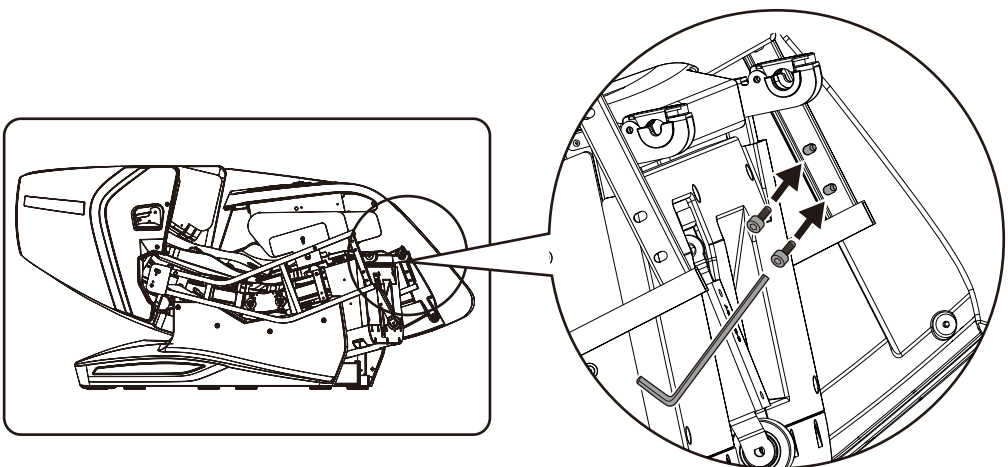


Assembly Instruction

■Assembly instruction for Legrest

Step 1 Turn on the power switch—the chair will automatically rise.

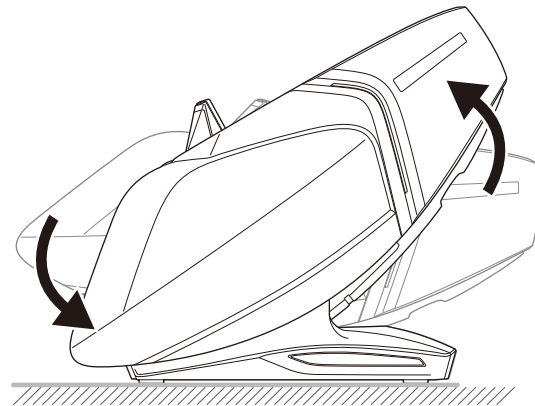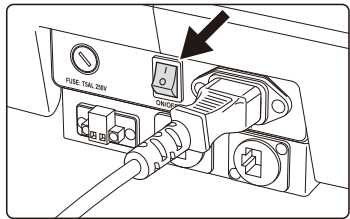


Step 2 Unzip the seat pad and lift it up.

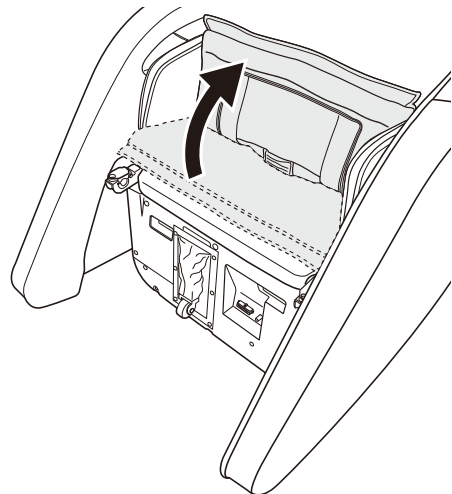


*It will be easier to work with if you lift it up.

Step 3 Connect the wire assembly and air hose from the legrest to their corresponding slots on the main chair body.

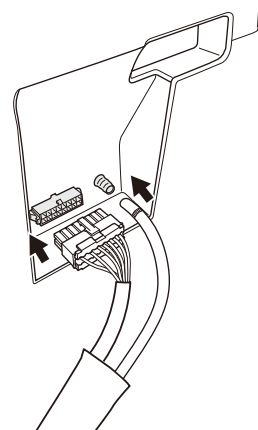

*When connecting the wiring harness assembly, hold the connector—not the wires directly—and insert it firmly all the way into the slot.

Step 4 After completing the wire connection, open the fixing covers of the legrest on both sides of the steel frame assembly. Lift the legrest and place it onto the slot. Then flip the fixing cover down and fasten it with the 2 legrest retaining screws.

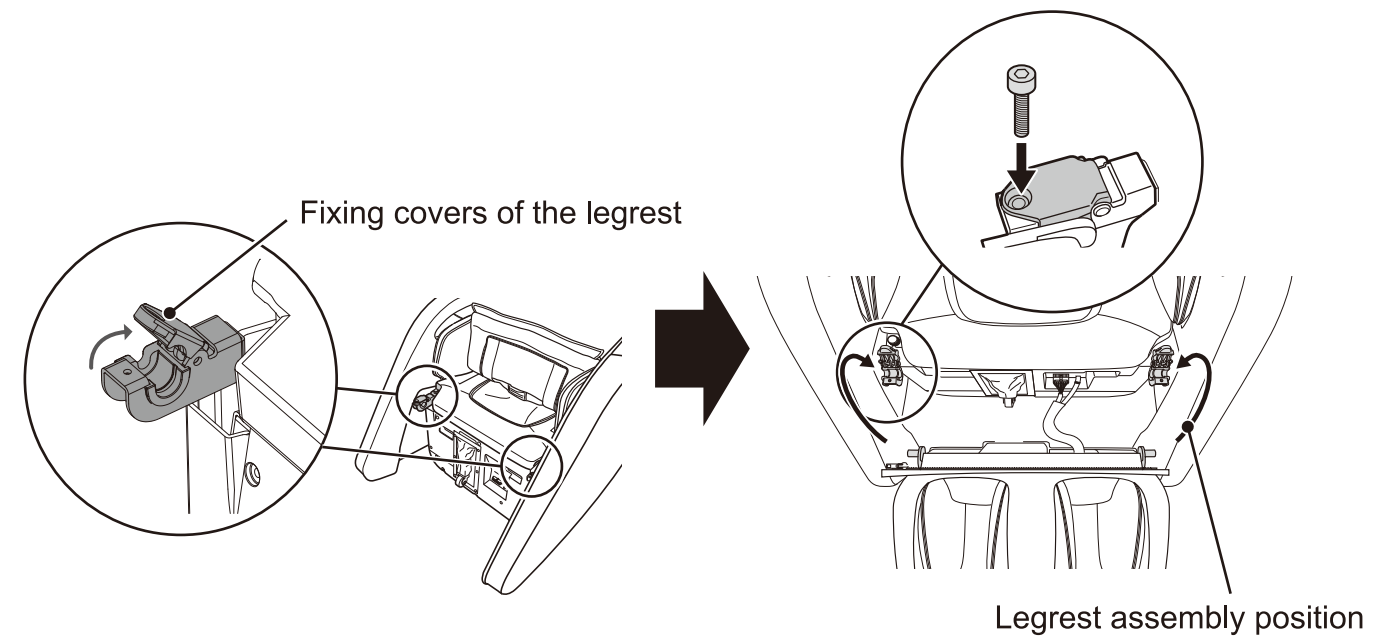


Assembly Instruction


■Installing the Zero-Wall Frame

After installing the chair, remove the casters and attach the frame. This prevents the chair from hitting the wall when reclining.


Step 1 Remove the 2 screws (①) and remove the back cover (②).

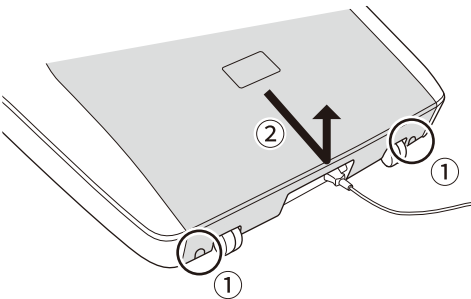


Step 2 Remove the pin securing the casters, then remove the casters.

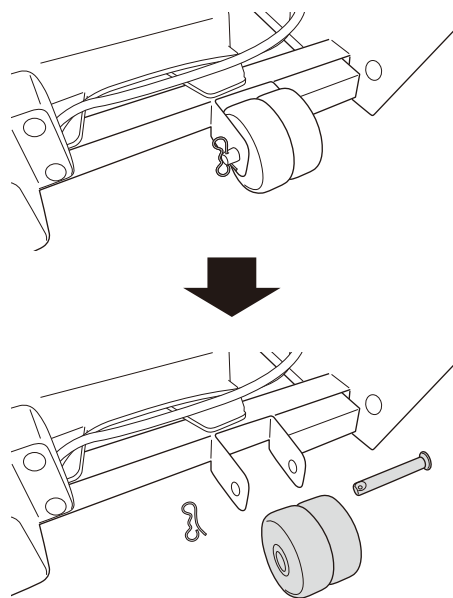

Step 3 Align the frame and secure it with the pin.

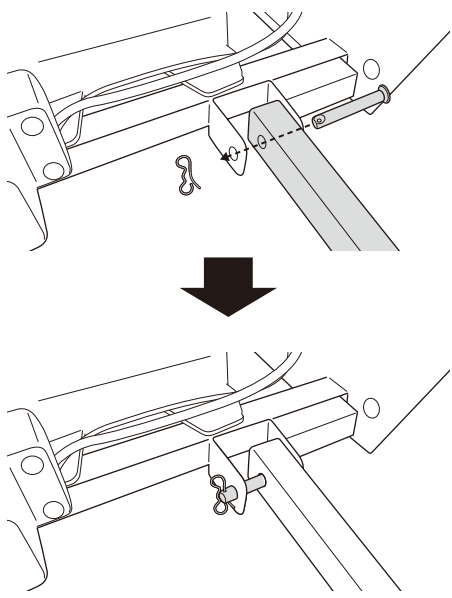

Step 4 Attach the back cover and secure it with 2 screws. Apply the warning label in the desired language.

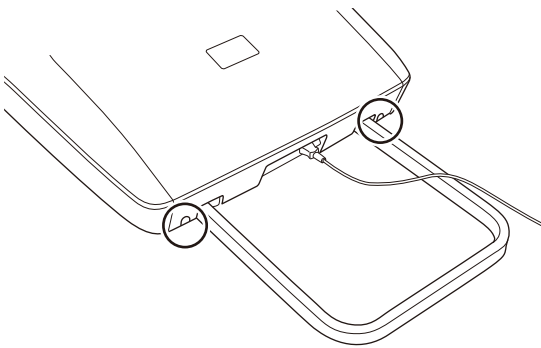


Before Use

Check around


Make sure there are no people, pets, or other objects around the unit.

! CAUTION

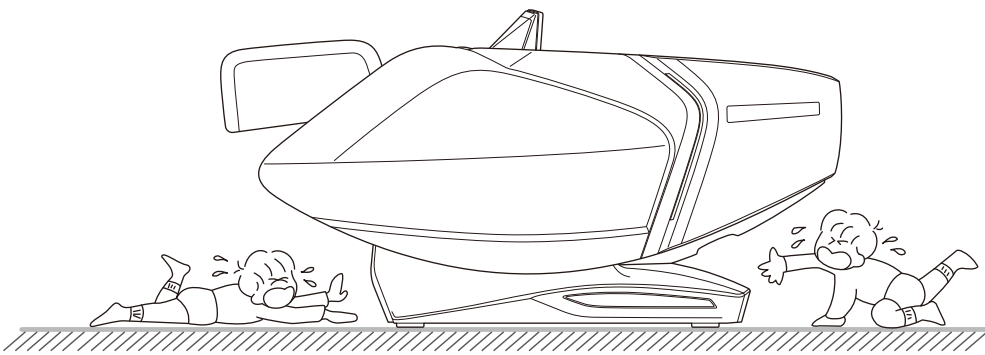


Turn on the power


! CAUTION

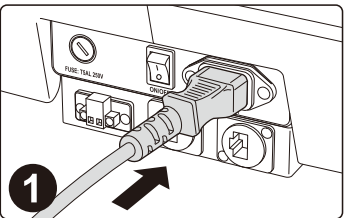

1 Connect the power cord on the chair

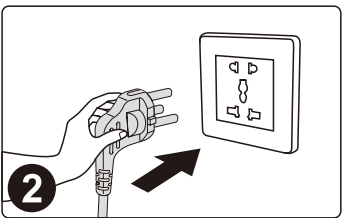

2 Connect to power source

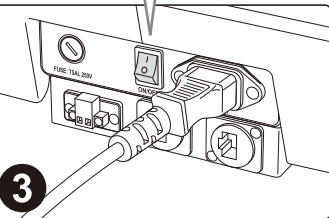

3 Turn on the power switch ("I" position)


! DANGER

Please check whether the power cord and power plug are damaged before use.
Please supply power to this product in accordance with the calibrated working power requirements of this product.
It is strictly prohibited to use power supply other than the rated power supply of this product.
It is strictly prohibited to connect various power conversion devices to power this product without authorization.

Names and Functions of Each Part


Main unit

Front side

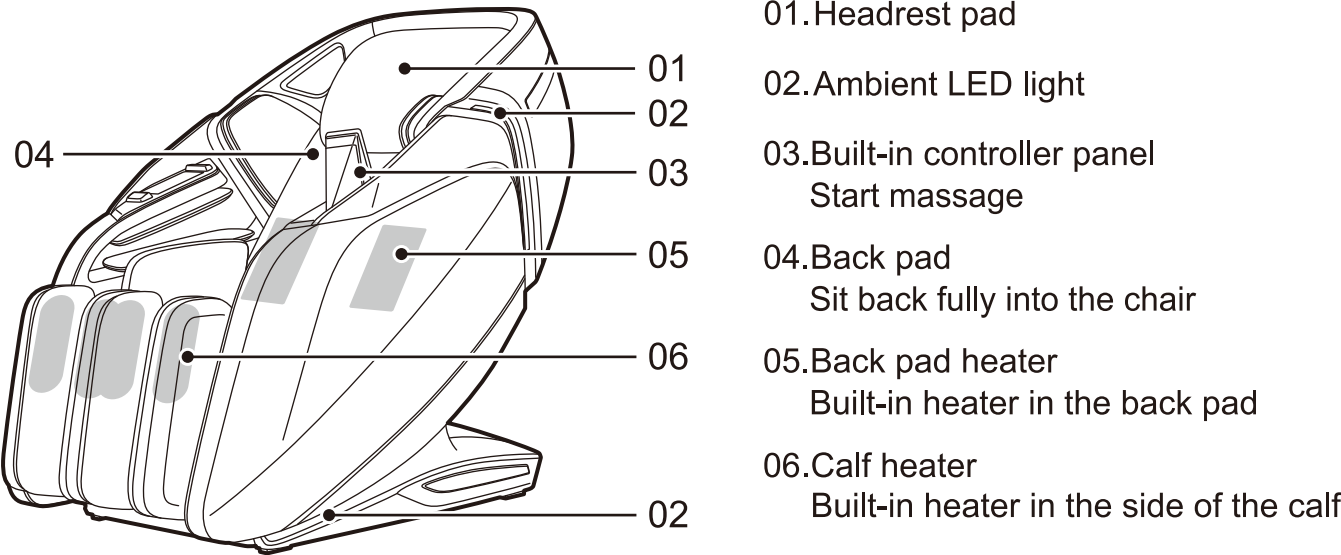


Main unit


Back

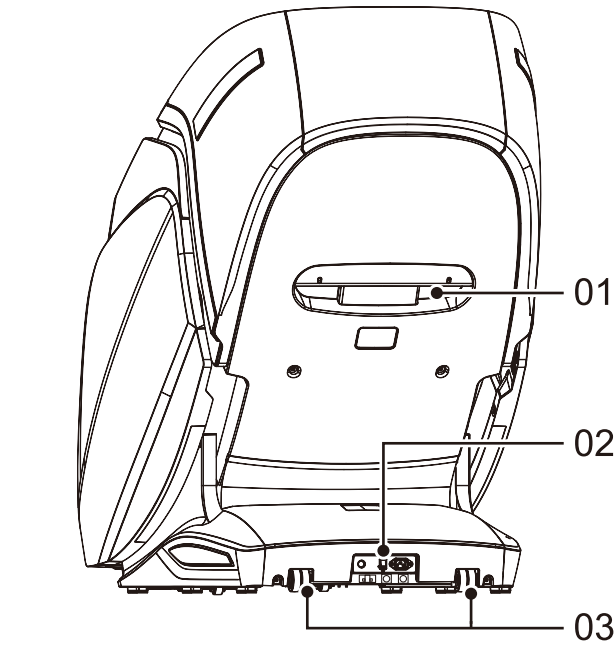

- 01.Handle
 Lifting massage chair
- 02.Power switch, Power cord input,
 RJ45 port, RJ22/RJ9 port,
 Screw terminal port

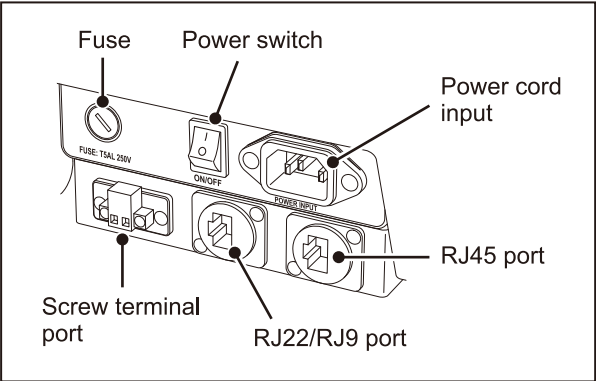


Massage area

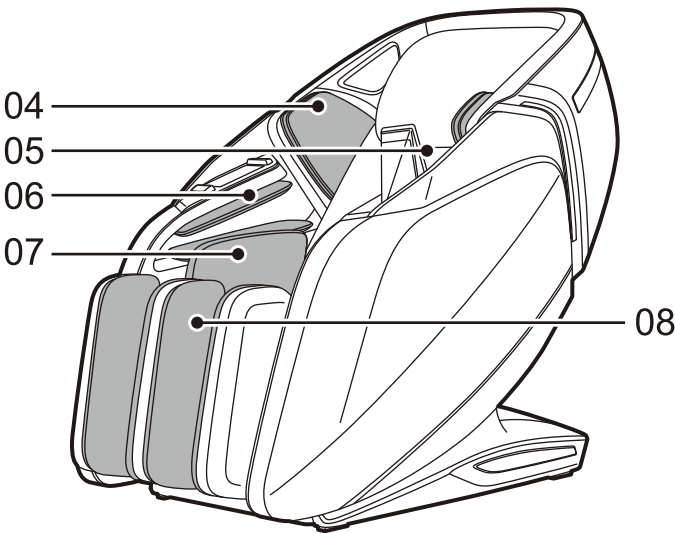

— Air compression
 massage areas


Names and Functions of Each Part

Remote control

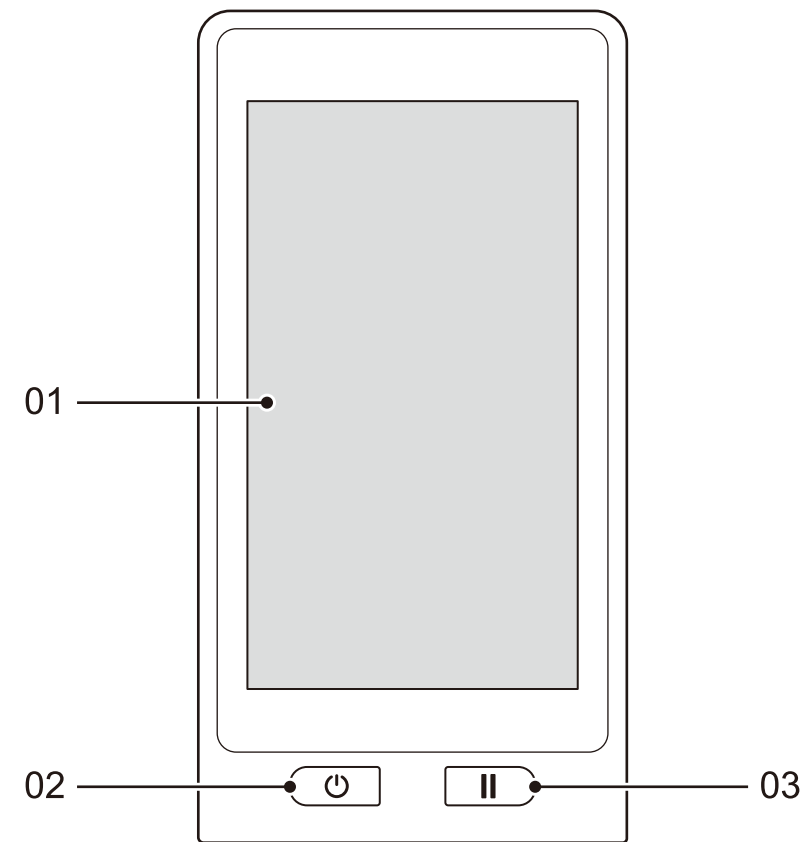

- 01. Touch display : Controls the user interface
- 02. Power button : Turns the power on and off
- 03. Pause button : Pause or resume the operation

Charging port

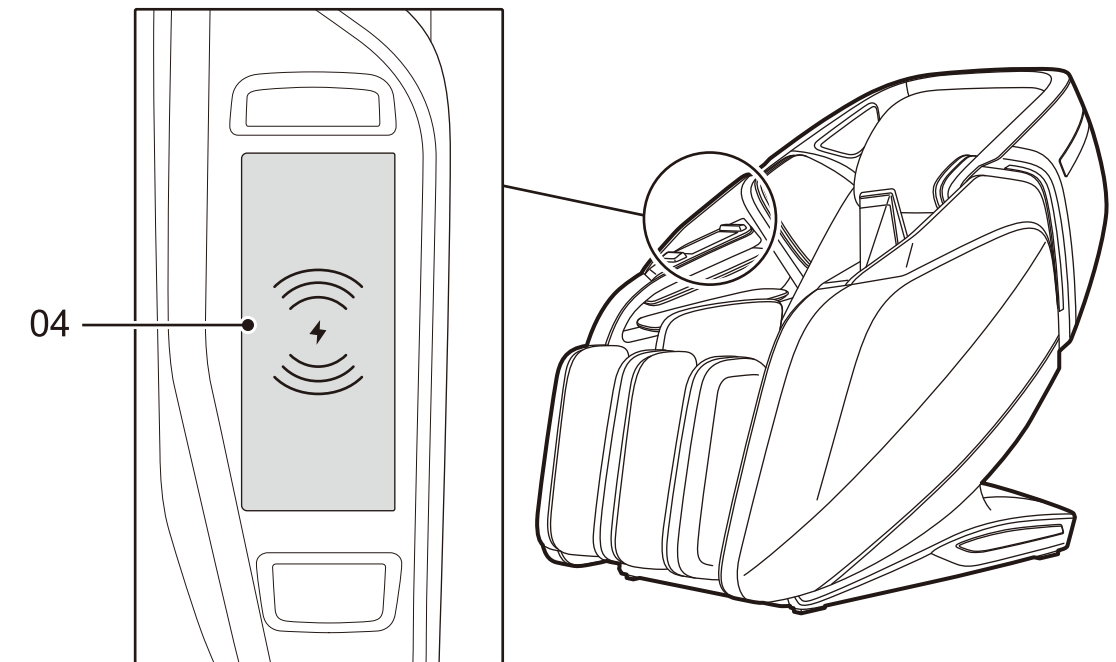

04. Mobile phone wireless charger

*If wireless charging does not work, please remove the smartphone case and try again.

How to Use Auto Course

■Auto course




①	RECOVERY	A standard massage course combining kneading and tapping techniques
②	RELAXATION	A slow-paced massage course focusing on soothing kneading techniques
③	VIGOROUS	A fast-paced massage course focusing on tapping techniques
④	STRETCH	A mixture of various stretch movements to relieve the entire body. * During the stretch movement, recline and air intensity are automated.
⑤	UPPER BODY	An upper body-focused massage course with varied techniques
⑥	LOWER BODY	A lower body-focused massage course with varied techniques

■Start an auto course


1. Press the power button on the controller to activate the touch screen.

* If you have the scanner option installed on the right armrest, scan the QR code/barcode or touch the RFID tag on the scanning area.

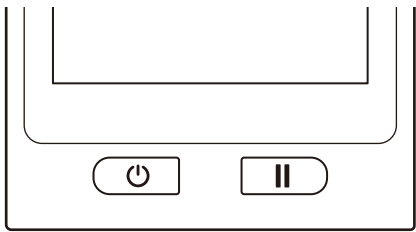

2. Choose your preferred massage course from the 6 available options.
3. When you select a course, the chair will automatically perform a body scan and recline.
 - Lean back fully into the chair during the body scan.
 - To skip the body scan, tap the skip icon on the screen.

4. Once the body scan is complete, the massage will begin.
 - After the auto body scan, you can adjust the shoulder position by pressing the ^ v icon and (✓) button to confirm.

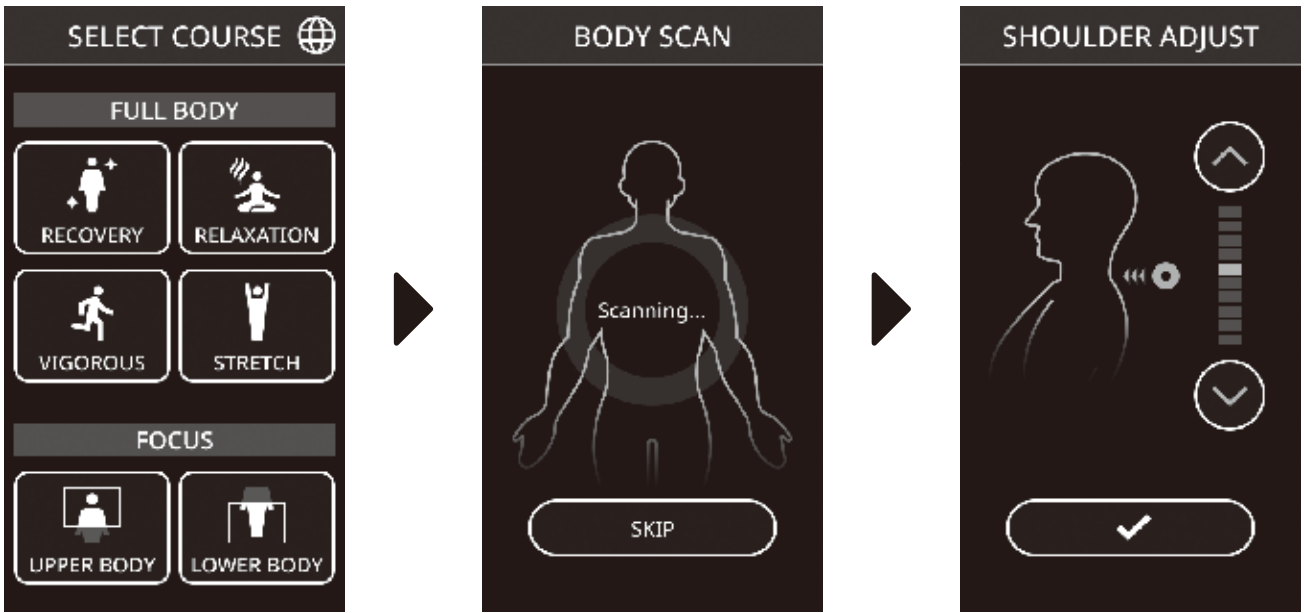

5. The selected massage course and operation screen will appear.

How to Use Auto Course



Items Adjustable During Massage

①		Recline level	Adjust the recline angle (5 levels)
②		Hand intensity	Adjust the massage roller 3D protrusion (3 levels)
③		Pressure level	Adjust the air compression intensity (3 levels)
④		Home	Return to the massage selection screen. * The massage will continue until you select a new course.
⑤		Adjustment	• Backrest and legrest angles can be adjusted independently. • Shoulder position can be re-adjusted if necessary.
⑥		Heater	Turn the heater on or off.

6. When the timer elapses, the massage will end.

After completing the Massage

■Massage Timer

The massage will end after the set timer elapses.

- The backrest, legrest, and massage roller will automatically return to their starting positions after the massage.
- Massage time can be configured in manager mode.

■Pause and Resume

Pressing the button during massage operation temporarily pauses the massage.

- Press the button again or tap the icon on the screen to resume.




■To end the massage early


When you press the button during the massage, the massage will end.

- All massage movements will stop, and the backrest, legrest, and massage roller will be restored.
- Timer will be canceled.
- Wait until the chair returns completely to its original position before exiting to prevent accidents.


■Turn off the power to the entire machine

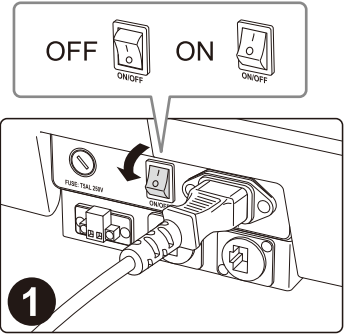

Turn off the power switch

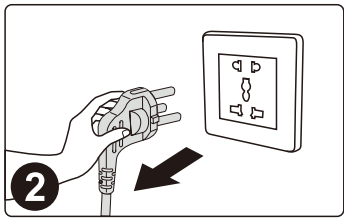

Remove the power plug from the outlet

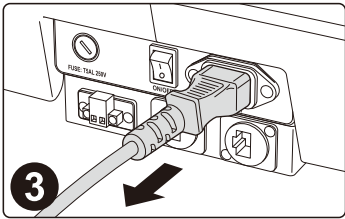


Unplug the power cord from the product

Manager mode

■How to Enter Manager mode

On the language selection screen, tap the top-right area more than 6 times consecutively to enter manager mode.

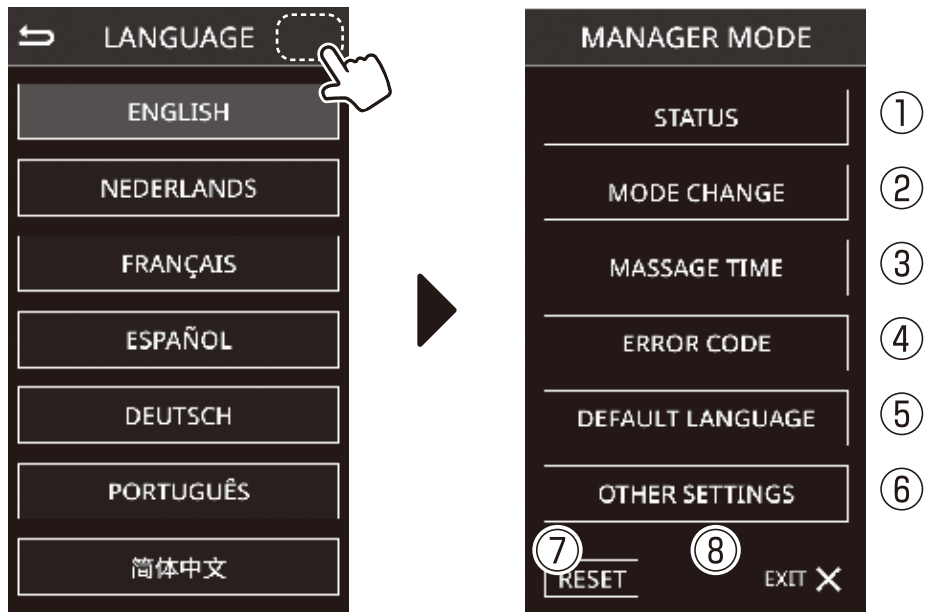

■Manager mode Description

①. STATUS

- You can check the current status.
 - Chair status: Displays chair information, usage data, and error status.
 - * Reset will not erase this information.
 - Setting status: Displays the current chair configuration.

②. MODE CHANGE

- You can switch the mode.
 - Open mode: Anyone can start a massage by pressing the power button on the controller.

- Member mode: Massage starts only for authorized members who scan a QR code or RFID tag on the chair management module or other external devices.
 - * When using the chair management module or devices such as T-Max or Gantner, set the chair mode to member mode.
- You can turn power saving mode (sleep mode) on or off and set the auto timeout duration to trigger the mode.

③. MESSAGE TIME

- You can configure the message time in open mode. (10-40 minutes)
 - * In member mode, the message time is dependent on the chair management module or other external devices.
- You can set the auto-start time in member mode. When enabled, a designated message course will start automatically after the specified time. This lets new users begin the massage immediately without needing to learn the interface. (OFF / 1–99 seconds)

④. ERROR CODE

- You can view current error codes and error history in sequence.
 - * If error codes persist in the Current Error tab, please contact your local service provider with this information.

⑤. DEFAULT LANGUAGE

- You can set the default language.
 - * When a session restarts or the power is recycled, the language will always return to the default.

⑥. OTHER SETTINGS



- You can turn the ambient LED lights on or off (side and bottom LEDs).
- You can adjust the brightness of the controller screen.

⑦. RESET

- You can reset the chair settings.
 - * Chair usage history (chair status data) will not be erased.

⑧. EXIT

- Exit manager mode.

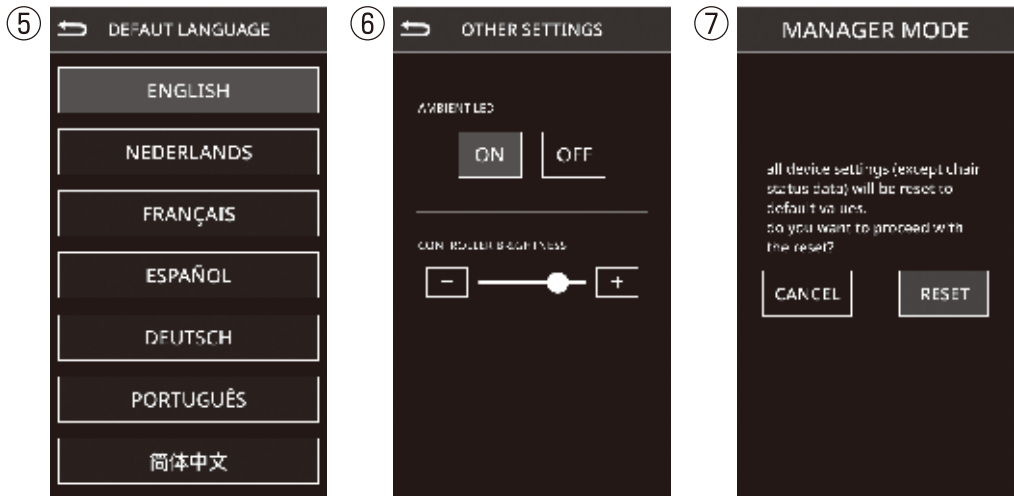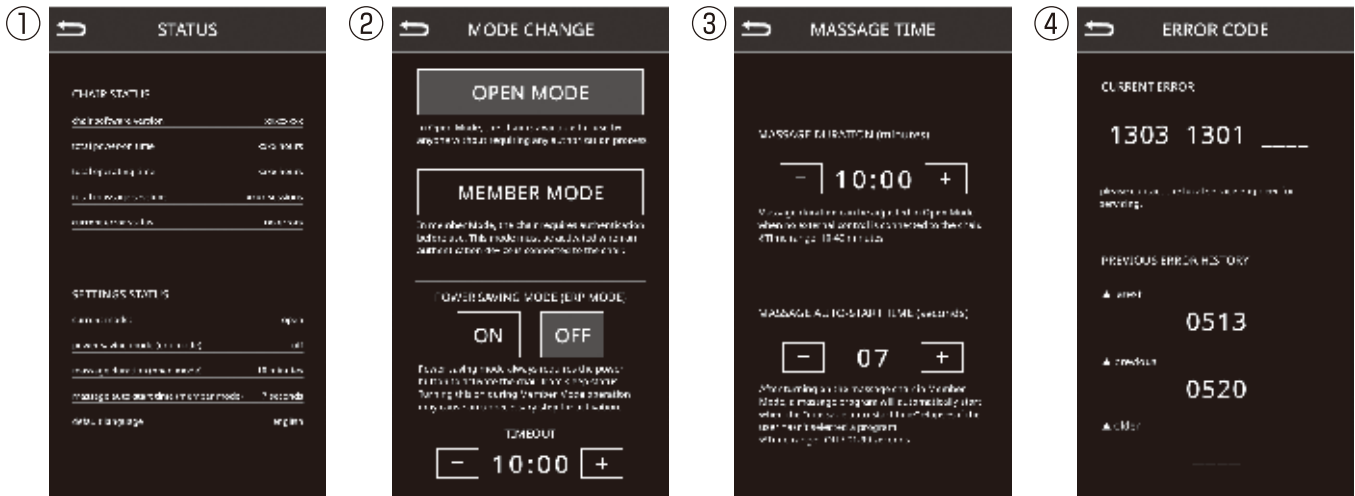


Language

■Language

On the selection screen, tap the  to select your desired language.

- You can change the user interface language on the controller.
- When the massage session is complete and you restart, the language returns to the default setting.
- The default language can be configured in manager mode.

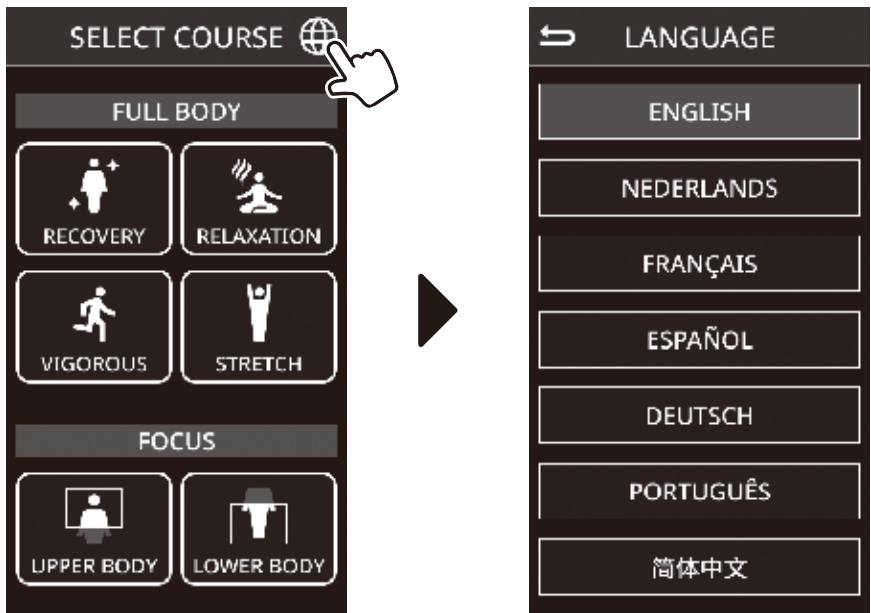

Troubleshooting

Issue	When the calf mechanism or the backrest reaches a certain position,press the corresponding adjustment button and it will no longer change,and a continuous "beep, beep" sound will be emitted.
Solution	Due to product structure design and ergonomic requirements, this product has a limited protection design for the limit state of attitude adjustment. At the same time, a reminder sound will sound when the limit is reached.

Issue	The controller is not turned on.
Solution	Check if the mode is set correctly. In open mode, the controller activates when you push the power button. In member mode, the controller requires the chair management module or other external devices to activate. Enter Manager mode to check the current status.

Issue	The machine makes a leather rubbing sound when adjusting its posture.
Solution	This product adopts high-grade anti-corrosion leather that meets environmental protection requirements. Due to the inherent characteristics of leather materials, there will be a normal slight noise when rubbing at the relative movement, which is a normal phenomenon. Note: If you hear fabric or leather tearing sound when the machine is running, please stop using it immediately and seek professional maintenance.

Issue	The massage roller does not reach your shoulder or neck properly.
Solution	The shoulder position may be set lower than your actual position after the initial body scan. Sit back fully against the backrest and tap the Shoulder Adjust icon in the Adjustment screen.

Issue	The machine is damaged. The power cord or power plug is abnormally heating.
Solution	To prevent accidents, please contact the authorized service center immediately.

For any questions about this product, please contact our service provider.

Cleaning and Maintenance

■Upholstery Care

Proper maintenance ensures hygiene, durability, and optimal appearance.

Daily Cleaning (Before & After Each Use) :

- Wipe surfaces with a soft microfiber cloth moistened with mild soap water or 70% or less isopropyl alcohol.
- Focus on high-contact areas: headrest, backrest, armrests, seat, legrest
- Allow 1-2 minutes to air dry

Key Guidelines :


- Use only alcohol or approved disinfectants safe for synthetic leather
- Avoid bleach, acetone, ammonia, and abrasive cleaners
- Use soft microfiber cloths only
- Do not oversaturate the cloth; use only a light dampening.

Weekly : Deep clean all surfaces and inspect for wear, cracks or damage

Monthly : Comprehensive inspection for damage, Use commercial-grade synthetic leather conditioner (alcohol-compatible)


Avoid : Sharp objects, direct sunlight, excessive moisture, and rough scrubbing

*Note : Regular alcohol-based disinfection is essential in commercial environments. This microfabric synthetic leather is designed to withstand frequent alcohol cleaning with minimal deterioration when proper protocols are followed.

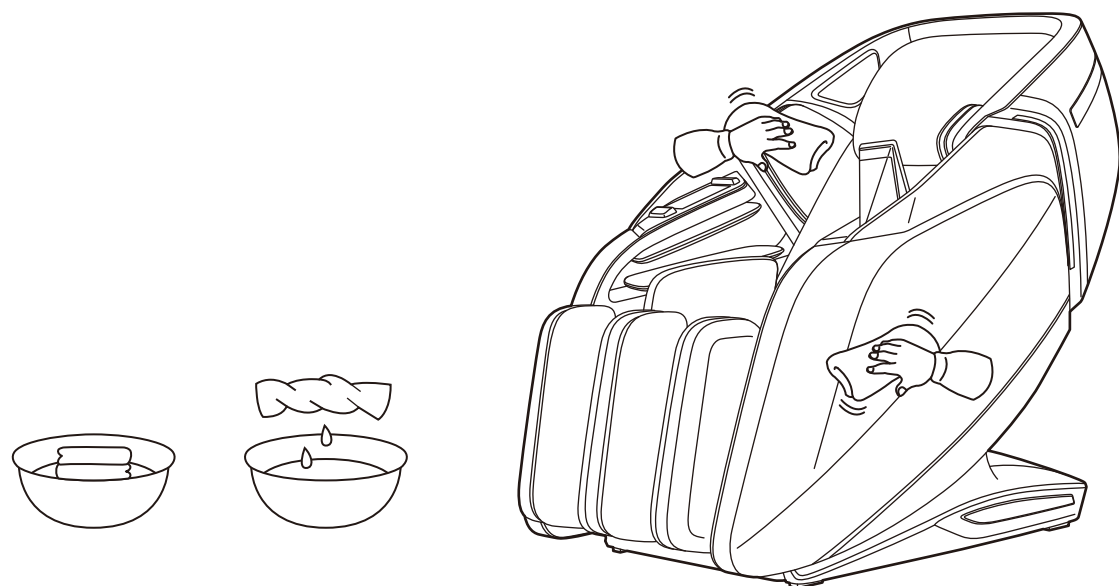

■Replaceable Parts

All upholstery components that contact skin use zipper-type construction for easy replacement. Contact an authorized supplier for replacement parts when worn.

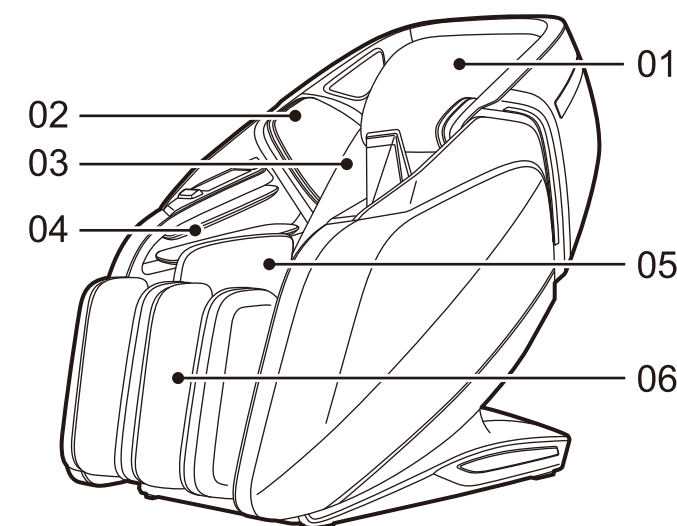

- 01. Headrest pad
- 02. Shoulder cover set
- 03. Back and seat cover set
- 04. Armrest cover set
- 05. Pelvis cover set
- 06. Calf cover set

Product Specifications

Model name	OH-8000
Power supply	AC 100-240V, 50/60Hz
Rated time	30min.
Power consumption	100W
Safety structure	Class I
Max. User Weight / Height (Approx.)	150kg / 195cm (330 lbs / 6'5")
Country of origin	Made in China
Max. User Shoulder/Seat/Calf width	56 / 48 / 17 cm (22 / 19 / 6.7 inch)
Frame Length / Material	130.5cm (51 inch) / Steel, POM
Power Recline Angle	129.5°~ 164°
Massage Mecha Travel Distance	90cm (35.4 inch)
Upholstery Material	Synthetic Leather (Microfiber PU)
Controller	Embedded 5" touchscreen control panel
Massage Strength level	3 levels (Protrusion / Compression)
Airbags	Shoulder - Arms - Hip - Calf (36cells)
External Port	RJ45, RJ22/RJ9, Screw Terminal Block
Upright Chair Size (LxWxH)	153.5 x 83 x 123 cm (60.4 x 32.7 x 48.4 inch)
Reclined Chair Size (LxWxH)	188 x 83 x 92.5 cm (74 x 32.7 x 36.4 inch)
Product Net Weight	121 kg (267 lbs)

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions:
(1) this device may not cause harmful interference, and
(2) this device must accept any interference received, including interference that may cause undesired operation.
Changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

Note: This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses, and can radiate radio frequency energy, and if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation. If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.